User Guide

iE

Revision 0 (Sep 2026)

Document No.: UG-00022

CONTENTS

02	English
16	Български
30	中文(繁體)
44	Hrvatski
58	čeština
72	Dansk
86	Dutch
100	Eesti keel
114	Suomi
128	Français
142	Deutsch
156	Ελληνικά
170	magyar
184	íslenska
198	Italiano
212	한국어
226	Latviešu
240	Lietuvis
254	Norsk
268	Português
282	Română
296	русский
310	Slovenčina
324	Slovenščina
338	Español
352	Svenska
366	Türk
380	Srpski
394	עברית
408	língua portuguesa
422	Shqip

English

About This Guide	4
1 Introduction and Overview	4
1.1 Intended Use	4
1.2 Indication for Use	4
1.3 Contraindications	4
1.4 Qualifications of the Operator	4
1.5 Symbols	4
1.6 iE Components Overview	5
1.6.1 Additional Components (Sold Separately)	6
1.7 Setting Up the iE System	6
1.7.1 Basic Settings of iE (Medit Plug & Scan)	6
2 Medit Scan for Clinics Overview	7
2.1 Introduction	7
2.2 Feature Availability	7
2.3 Installation	7
2.3.1 System Requirements	7
2.3.2 Software Installation Guide	8
2.3.3 Medit Scan for Clinics User Guide	8
2.4 Error Messages	8
3 Maintenance	8
3.1 Calibration	8
3.1.1 How to Calibrate iE	8
3.2 Cleaning, Disinfection, Sterilization Procedure	9
3.2.1 Reusable Tip	9
3.2.2 Mirror	9
3.2.3 Handpiece	9
3.2.4 Other Components	10
3.3 Disposal	10
3.4 Updates in Medit Scan for Clinics	10
4 Safety Guide	10
4.1 System Basics	10
4.2 Proper Training	10
4.3 In Case of Equipment Failure	11
4.4 Hygiene	11
4.5 Electrical Safety	11
4.6 Eye Safety	11
4.7 Explosion Hazards	11
4.8 Pacemaker and ICD Interference Risk	11
4.9 Cybersecurity Safety	11
5 Electro-Magnetic Compatibility Information	12
5.1 Electromagnetic Emissions	12
5.2 Electromagnetic Immunity	12
6 Specifications	15

About This Guide

Convention in This Guide

This guide uses various symbols to highlight important information to ensure correct usage, prevent injury to the user and others, and prevent property damage. The meanings of the symbols used are described below.

 **WARNING**

The WARNING symbol indicates information that, if ignored, could result in a medium risk of personal injury.

 **CAUTION**

The CAUTION symbol indicates safety information that, if ignored, could result in a slight risk of personal injury, property damage, or damage to the system.

 **TIPS**

The TIPS symbol indicates hints, tips, and additional information for optimal operation of the system.

1. Introduction and Overview

1.1 Intended Use

The iE system is an intraoral 3D scanner intended to record topographical characteristics of teeth and surrounding tissues digitally. The iE system produces 3D scans for use in computer-assisted design and manufacturing of dental restorations.

1.2 Indication for Use

The iE system is for scanning the patient’s intraoral features. Various factors (intraoral environment, operator’s expertise, and laboratory workflow) may affect the final scan results when using the iE system. The iE system can also be used in full-arch scans, but various factors (intraoral environment, operator’s expertise, and laboratory workflow) may affect the final results.

1.3 Contraindications

The iE system is not intended to be used to create images of the internal structure of teeth or the supporting skeletal structure.

1.4 Qualifications of the Operator

 **CAUTION**

- The iE system is designed for use by individuals with professional knowledge in dentistry and dental laboratory technology.
- The user of the iE system is solely responsible for determining whether or not this device is suitable for a particular patient’s case and circumstances.
- The user is solely responsible for the accuracy, completeness, and adequacy of all data entered the iE system and the provided software.
- The user should check the accuracy of the results and assess each individual case.
- The iE system must be used in accordance with its accompanying User Guide.
- Improper use or handling of the iE system will void its warranty. If you require additional information on the proper use of the iE system, please contact your local distributor.
- The user is not allowed to modify the iE system.

1.5 Symbols

No.	Symbol	Description
1		Serial number
2		Medical device
3		Date of manufacture
4		Manufacturer
5		Caution
6		Warning
7		Read the user guide
8		European Conformity marking
9		Authorized representative in the European Union
10		BF type of applied part
11		WEEE mark
12		Prescription use (U.S.A)
13		MET mark
14		AC
15		DC
16		Temperature limitation: -10 – 50°C (14 – 122°F)
17		Humidity limitation

18		Atmospheric pressure limitation
19		Fragile
20		Keep dry
21		This way up
22		Stacking limit: 14 layers
23		Consult instructions for use
24		UK Conformity Assessed marking (UKCA)
25		Authorized representative in Switzerland
26		Authorized representative in the United Kingdom
27		Model Number
28		Quantity
29		Unique device identifier

1.6 iE Components Overview

No.	Component (Model No.)	Qty	Illustration
1	iE Handpiece (MO8-iE)	1ea	
2	iE Handpiece Cover (MO1-HC1)	1ea	
3	Reusable Tip (Large) (MO1-RTL)	2ea	
4	Reusable Tip (Medium) (MO1-RTM)	2ea	
5	Calibration Tool (MO1-CT1)	1ea	
6	Practice Model	1ea	
7	Desktop Cradle (MO1-DC)	1ea	
8	Wall Mount Holder (MO1-WH1)	1ea	
9	Power Delivery Cable (2.5 m)	1ea	
10	Power Delivery Cable (2 m)	1ea	
11	USB Flash Drive (Installer of Medit Scan for Clinics included)	1ea	
12	User Guide	1ea	

1.6.1 Additional Components (Sold Separately)

No.	Component (Model No.)	Qty	Illustration
1	Reusable Tip (Large) (MO1-RTL)	4ea	
2	Reusable Tip (Medium) (MO1-RTM)	4ea	

- All components listed in 1.6 & 1.6.1 can be purchased separately.
- Availability of items for sale may vary depending on the medical device registration status in each country or region. Please contact Medit or your local distributor for the availability of specific items.


 CAUTION

- Keep the practice model in a cool place away from direct sunlight. A discolored practice model may affect the results of the practice mode.
- Medit Scan for Clinics is included in the USB drive. This product is optimized for PC, and using other devices is not recommended. Do not use anything other than a USB port. It may cause malfunction or fire.
- The small reusable tip (MO1-RTS) is not compatible with the iE system.

1.7 Setting Up the iE System


1.7.1 Basic Settings of iE (Medit Plug & Scan)

You can connect iE directly to a PC with the power delivery cable.

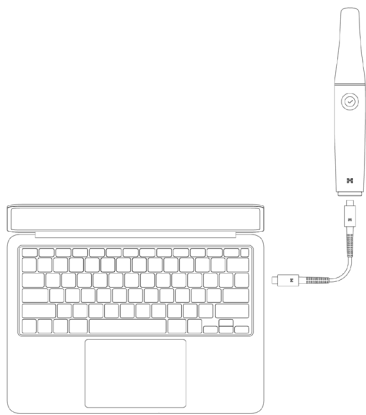Turning On the iE

- ① Connect the iE to your PC with the power delivery cable, and the scanner will automatically power on.
- ② When power is applied, the rear LED lights up.

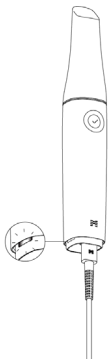
Turning Off the iE


When you disconnect the power delivery cable, the scanner switches off.

Desktop Cradle

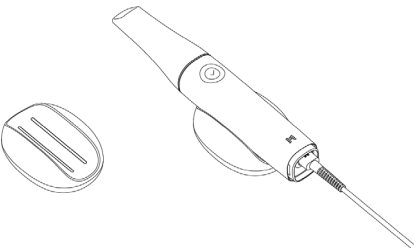Wall Mount Holder

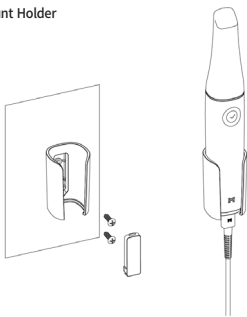2. Medit Scan for Clinics Overview

2.1 Introduction

Medit Scan for Clidfriendly working interface to digitally record topographical characteristics of teeth and surrounding tissues using the iE system.

2.2 Feature Availability

The following features and services are not supported on this scanner model:

- Medit Link Express and Medit Link Express for Mac
- Air Mouse

2.3 Installation

2.3.1 System Requirements

Recommended System Requirements

Windows			macOS
Laptop		Desktop	Laptop/Desktop
Processor	Intel Core i7 – 14700HX Intel Core i7 – 13700H Intel Core Ultra 7 – 155H	Intel Core i7 – 14700 Intel Core i7 – 13700	M1 Pro (10-core CPU, 16-core GPU) M2 (8-core CPU, 10-core GPU) M2 Pro (10-core CPU, 16-core GPU) M3 (8-core CPU, 10-core GPU) M3 Pro (8-core CPU, 10-core GPU)
	AMD Ryzen 7 8745H AMD Ryzen 7 7735H	AMD Ryzen 7 8700F AMD Ryzen 7 7700X	
RAM	32 GB		24 GB
Graphic	NVIDIA GeForce RTX 5060 (VRAM 8 GB or higher) *Compatible with Medit Scan for Clinics v1.13.3 or later. NVIDIA GeForce RTX 4060 (VRAM 8 GB or higher) NVIDIA GeForce RTX 3070 (VRAM 8 GB or higher) NVIDIA RTX A3000 (VRAM 8 GB or higher) *AMD Radeon is not supported.		-
OS	Windows 10 64-bit Windows 11 (recommended for Intel Core 12th Gen or later)		Sonoma 14, Sequoia 15 or higher

Minimum System Requirements

Windows			macOS
Laptop		Desktop	Laptop/Desktop
Processor	Intel Core i5 – 14500HX Intel Core i5 – 13500H Intel Core Ultra 5 – 125H	Intel Core i5 – 14400 Intel Core i5 – 13400	M3 (8-core CPU, 8-core GPU) M4 (10-core CPU, 8-core GPU)
	AMD Ryzen 5 8640U AMD Ryzen 5 7640H	AMD Ryzen 5 8600G AMD Ryzen 5 7600	
RAM	16 GB		16 GB
Graphics	NVIDIA GeForce RTX 4050 (VRAM 6 GB or higher) NVIDIA GeForce RTX 3060 (VRAM 6 GB or higher) NVIDIA RTX A2000 (VRAM 6 GB or higher) *AMD Radeon is not supported.		-
OS	Windows 10 64-bit Windows 11 (recommended for Intel Core 12th Gen or later)		Sonoma 14, Sequoia 15 or higher

For accurate and up-to-date system requirements, please visit www.meditlink.com.

- Use a PC and monitor certified to IEC 62368-1 (or IEC 60950-1), IEC 55032, and IEC 55024.
- To prevent electric shock, connect the PC only to a properly grounded power source.
 - Ensure that the PC's power plug remains easily accessible at all times.

The device may not work when using cables other than the USB 3.0 cable provided by Medit. Medit is not responsible for any problems caused by cables other than the USB 3.0 cable provided by Medit. Be sure to use only the USB 3.0 cable included in the package.

2.3.2 Software Installation Guide

- ① Connect the enclosed USB flash drive to a PC.
- ② Run the installation file.
- ③ Select the setup language and click "Next."
- ④ Choose the installation path.
- ⑤ Read the "License Agreement" carefully, check "I agree to the License terms and conditions," and then click "Install."
- ⑥ The installation process may take several minutes. Please do not shut down the PC until the installation is complete.
- ⑦ After the installation is complete, restart the PC to ensure optimal program operation.

The installation will not be processed if the iE system is connected to a PC. Please turn off the scanner before starting the installation.

2.3.3 Medit Scan for Clinics User Guide

Please refer to the User Guide of Medit Scan for Clinics: Medit Scan for Clinics > Menu > User Guide.

2.4 Error Messages

The Medit iE scanner system presents error messages, enabling users to intuitively identify hardware and system conditions. Some messages include solutions to aid users in self-diagnosis and problem resolution. Error messages are presented in plain text and supplemented with images when needed to enhance comprehension.

Connection failed.

Failed to connect to the cameras. Please follow the instructions below and try again.

- Press the "Retry" button. Device will automatically reboot and try to connect again.
- Turn off the device and disconnect all cables from the PC. Connect cables again and restart the device.
- If connection problem persists, please install the program again.
- If connection problem persists, please contact your service provider.

CancelRetry

Turn Off the External Light

External light was detected while scanning.

This scan might have proceeded with an external light source, which would increase the chances of noisy and inaccurate scanned data.

We recommend turning the external light off before the next scan.

IgnoreClose

CAUTION

If the instructions provided in the error message do not address the issue, please contact the local distributor from whom you acquired your scanner system or support@medit.com.

3. Maintenance

CAUTION

- Equipment maintenance should only be carried out by a Medit employee or Medit-Certified company or personnel.
- In general, users are not required to perform maintenance work on the iE system besides calibration, cleaning, and sterilization. Preventive inspections and other regular maintenance are not required.

3.1 Calibration

Periodic calibration is required to produce precise 3D models. You should perform calibration when:

- The quality of the 3D model is not reliable or accurate when compared to previous results.
- Environmental conditions such as temperature have changed.
- The calibration period has expired. You can set the calibration period in the Menu > Settings > Calibration Period (Days).

The calibration panel is a delicate component.

Do not touch the panel directly. Check the calibration panel if the calibration process is not performed properly. If the calibration panel is contaminated, please contact your service provider.

If the calibration tool is long exposed to low temperatures below Medit's recommendations mentioned in the use and storage guidance, the rotation motion of the device may be compromised. In that case, forcibly turning the device may cause damage. To avoid it, keep the calibration device in the recommended temperature environment before use.

We recommend performing a calibration periodically. You can set the calibration period via Menu > Settings > Calibration Period (Days). The default calibration period is 90 days.

3.1.1 How to Calibrate iE

- ① Turn on the iE and launch the Medit Scan for Clinics.
- ② Run the Calibration Wizard at the bottom of the main toolbar panel in Medit Scan for Clinics.
- ③ Prepare the calibration tool and the iE handpiece.
- ④ Turn the dial of the calibration tool to the starting position.
- ⑤ Put the iE handpiece into the calibration tool.
- ⑥ Click "Next" to start the calibration process.
- ⑦ If the calibration tool is mounted properly in the correct position, the system automatically acquires data.
- ⑧ When data acquisition is complete at the starting position, turn the dial to the next position.
- ⑨ Repeat the steps to the last position.
- ⑩ When data acquisition is complete at the last position, the system automatically calculates and shows the calibration results.

3.2 Cleaning, Disinfection, Sterilization Procedure

3.2.1 Reusable Tip

The reusable tip is used intraorally during scanning and is intended for limited reuse. Prior to reuse between patients, the tip must be cleaned and either disinfected or sterilized to prevent cross-contamination.

Cleaning & Disinfection

- Prepare a cleaning solution.
 - » Dilute a general neutral detergent at a 1:100 ratio before use.
- Clean the reusable tips with a cleaning solution and a brush.
 - » Ensure the mirror of the tip is completely clean and stain-free. If the mirror appears stained or foggy, repeat the cleaning process.

CAUTION

- » The reusable tip has a complex structure, and automatic cleaning may not achieve complete cleaning; therefore, do not clean the reusable tip in an automatic washer.
- Rinse the reusable tips three times with purified water.
- Remove moisture with a paper towel and allow them to air-dry completely at room temperature for at least 80 minutes.
- Disinfect the reusable tips using a disinfectant containing 15% or less Isopropyl Alcohol (IPA) for 1 minute. Then, ensure they are thoroughly dried at room temperature for at least 5 minutes.
 - » Before using a disinfectant product, please refer to the product user manual for proper use.
 - » You can find the list of recommended disinfectants in the Medit Help Center at <http://support.medit.com/hc>.

Cleaning & Sterilization

- Complete the cleaning procedure described in the section above.
- Insert the tip into a paper sterilization pouch and seal it, making sure that it is airtight. Use either a self-adhesive or heat-sealed pouch.
- Sterilize the wrapped tip in an autoclave with the following conditions:
 - » Sterilize in a gravity autoclave at 135°C (275°F) for 10 minutes and dry for 30 minutes.
 - » Sterilize in a pre-vacuum autoclave at 134°C (273.2°F) for 4 minutes and dry for 20 minutes.
 - » Sterilize in a gravity autoclave at 121°C (249.8°F) for 30 minutes and dry for 15 minutes.
- Use an autoclave program that dries the wrapped tip before opening the autoclave.
- Scanner tips can be re-sterilized up to 150 times. After reaching this limit, they must be disposed of according to the guidelines in the disposal section.
- Autoclave times and temperatures may vary depending on the autoclave type and manufacturer. For this reason, it may not be able to meet the maximum number of times. Please refer to the user's manual of the autoclave manufacturer you are using to determine whether the required conditions are met.

CAUTION

- The mirror of the scanner tip is a delicate optical component that should be handled with care to ensure optimal scan quality. Be careful not to scratch or smudge it as any damage or blemishes may affect the data acquired.
- Wrap the scanner tip in cotton gauze before autoclaving to prevent permanent mirror stains. Refer to your autoclave's operating manual for specific sterilization parameters.
- Tips that have been cleaned, disinfected, and sterilized must remain sterile until they are used on the patient.
- Medit is not responsible for any damage, including tip distortion, caused by cleaning, disinfection, or sterilization procedures that do not comply with the guidelines outlined above.

3.2.2 Mirror

The presence of impurities or smudges on the tip mirror may lead to poor scan quality and an overall poor scanning experience. In such situations, clean the mirror following the steps below:

- ① Disconnect the scanner tip from the iE handpiece.
- ② Pour alcohol on a clean cloth or cotton-tipped swab and wipe the mirror. Make sure to use alcohol that is free of impurities or it may stain the mirror. You can use either ethanol or propanol (ethyl-/propyl alcohol).
- ③ Wipe the mirror dry using a dry, lint-free cloth.
- ④ Make sure the mirror is free of dust and fibers. Repeat the cleaning process as necessary.

3.2.3 Handpiece

After treatment, clean and disinfect all other surfaces of the iE handpiece except for the scanner front (optical window) and end (air vent hole). Cleaning and disinfecting must be done with the device turned off. Use the device only after it is completely dry.

The recommended cleaning and disinfecting solution is denatured alcohol (ethyl alcohol or ethanol) – typically 60 – 70% Alc/Vol.

The general cleaning and disinfecting procedures are as follows:

- ① Turn off the power by unplugging the power delivery cable from the device.
- ② Clean the filter on the front end of the iE handpiece.
 - » If alcohol is poured directly into the filter, it may seep inside the iE handpiece and cause a malfunction.
 - » Do not clean the filter by pouring alcohol or cleaning solution directly into the filter. The filter must be gently wiped with a cotton or soft cloth moistened with alcohol.
 - » Do not wipe by hand or apply excessive force.
 - » Medit is not responsible for any damage or malfunction that occurs during cleaning that does not follow the guidelines above.
- ③ After cleaning the filter, put the cover on the front of the iE handpiece.
- ④ Pour the disinfectant onto a soft, lint-free, and non-abrasive cloth.
- ⑤ Wipe the scanner surface with the cloth.
- ⑥ Dry the surface with a clean, dry, lint-free, and non-abrasive cloth.

CAUTION

- Do not clean the iE handpiece when the device is turned on as the fluid may enter the scanner and cause malfunction.
- Use the device after it is completely dry.
- Chemical cracks may appear if improper cleaning and disinfecting solutions are used during cleaning.

3.2.4 Other Components

- Pour the cleaning and disinfecting solution onto a soft, lint-free, and non-abrasive cloth.
- Wipe the component surface with the cloth.
- Dry the surface with a clean, dry, lint-free, and non-abrasive cloth.

CAUTION

- Chemical cracks may appear if improper cleaning and disinfecting solutions are used during cleaning.

3.3 Disposal

CAUTION

- The scanner tip must be sterilized before disposal. Sterilize the tip as described in the section "3.2.1 Reusable Tip."
- Dispose of the scanner tip as you would any other clinical waste.
- Other components are designed to conform with the following directives:
 - RoHS, Restriction of the Use of Certain Hazardous Substances in Electrical and Electronic Equipment. (2011/65/EU)
 - WEEE, Waste Electrical and Electronic Equipment Directive. (2012/19/EU)

3.4 Updates in Medit Scan for Clinics

Medit Scan for Clinics automatically checks for updates when the software is in operation. If a new version of the software is released, the system will automatically download it.

4. Safety Guide

Please adhere to all the safety procedures as detailed in this User Guide to prevent human injury and equipment damage. This document uses the words WARNING and CAUTION when highlighting precautionary messages.

Carefully read and understand the guidelines, including all preventive messages as prefaced by the words WARNING and CAUTION. To avoid bodily injury or equipment damage, make sure to adhere strictly to the safety guidelines. All instructions and precautions as specified in the Safety Guide must be observed to ensure the system's proper functionality and personal safety.

The iE system should only be operated by dental professionals and technicians who are trained to use the system. Using the iE system for any purpose other than its intended usage as outlined in the section "1.1 Intended Use" may result in injury or damage to the equipment. Please handle the iE system according to the guidelines in the safety guide.

Any serious incident involving the device should be reported to the manufacturer and the competent authority of the member state in which the user and patients are located.

4.1 System Basics

The iE system is a high-precision optical medical device. Acquaint all the following safety and operating instructions before the installation, usage, and operation of iE.

CAUTION

- If the product has been stored in a cold environment, give it time to adjust to the temperature of the environment before use.
- If used immediately, condensation may occur, which may damage the electronic parts inside the unit.
- Ensure that all components provided are free from physical damage. Safety cannot be guaranteed if there is any physical damage to the unit.
- Before using the system, check that there are no issues such as physical damage or loose parts. If there is any visible damage, do not use the product and contact the manufacturer or your local representative.
- Check the iE handpiece and its accessories for any sharp edges.
- When not in use, the iE system should be kept mounted on a desk stand or a wall mount stand.
- Do not install the desk stand on an inclined surface.
- Do not place any object on the iE system.
- Do not place the iE system on any heated or wet surface.
- Do not block the air vents located at the rear of the iE system. If the equipment overheats, the iE system may malfunction or stop working.
- Do not spill any liquid on the iE system.
- The iE handpiece and other included components are made of electronic components. Do not allow any kind of liquid or foreign objects to enter.
- Do not pull or bend the cable connected to the iE system.
- Carefully arrange all the cables so that you or your patient do not trip or get caught in the cables. Any pulling tension on the cables may cause damage to the iE system.
- Always keep an eye on the product and your patient while using the product to check for abnormalities.
- Proceed with calibration, cleaning, disinfection, and sterilization in accordance with the contents of the user guide.
- If you drop the tips on the floor, do not attempt to reuse them. Discard the tip immediately as there is a risk that the mirror attached to the tip may have been dislodged.
- Due to its fragile nature, the tips should be handled with care. To prevent damage to the tip and its internal mirror, be careful to avoid contact with the patient's teeth or restorations.
- If the iE system is dropped on the floor or if the unit is impacted, it must be calibrated before use. If the instrument is unable to connect to the software, consult the manufacturer or authorized resellers.
- If the equipment fails to operate normally, such as having issues with accuracy, stop using the product, and contact the manufacturer or authorized resellers.
- Install and use only approved programs to ensure the proper functionality of the iE system.
- In the event of a severe accident involving the iE system, notify the manufacturer and report it to the competent national authority of the country or region where the user and patient reside.
- If the PC with the software installed does not have security software or if there is a risk of malicious code intrusion into the network, the PC may be breached with malware (malicious software such as viruses or worms that damage your computer).
- The software for this product must be used in compliance with medical and personal information protection laws.
- To provide electrical insulation and maintain electrical safety, a coating is applied to insulate the device except on the areas where the USB ports are located.
- Check the product's markings, such as certification symbols or manufacturer details, under lighting of about 500 lux for clear visibility.

4.2 Proper Training

WARNING

Before using your iE system on patients:

- You should have been trained to use the system or read and fully understand this User Guide.
- You should be familiar with the safe use of the iE system, as detailed in this User Guide.
- Before use or after changing any settings, the user should check that the live image is displayed properly in the camera preview window of the program.

4.3 In Case of Equipment Failure

WARNING

If your iE system is not working properly or if you suspect that there is a problem with the equipment:

- Remove the device from the patient's mouth and discontinue use immediately.
- Disconnect the device from the PC and check for errors.
- Contact the manufacturer or authorized resellers.
- Modifications to the iE system are prohibited by law as they may compromise the safety of the user, patient, or a third party.

4.4 Hygiene

WARNING

For clean working conditions and patient safety, ALWAYS wear clean surgical gloves when:

- Making any direct contact with the tip.
- Using the iE system on patients.
- Touching the iE system.

WARNING

The iE system and its optical window should be always kept clean. Before using the iE system on a patient, be sure to:

- Sterilize the iE system as described in the section "3.2 Cleaning, Disinfection, Sterilization Procedure."
- Use a sterilized tip.

4.5 Electrical Safety

WARNING

- The iE system is a Class II device.
- The iE system only uses RF energy internally. The amount of RF radiation is low and does not interfere with surrounding electromagnetic radiation.
- There is a risk of electric shock if you attempt to access the inside of the iE system. Only qualified service personnel should access the system.
- Do not place liquids such as beverages near the iE system and avoid spilling any liquid on the system.
- Never spill liquid of any kind on the iE system.
- Condensation due to changes in temperature or humidity can cause moisture buildup inside the iE system, which may damage the system. Before connecting the iE system to a power supply, be sure to keep the iE system at room temperature for at least two hours to prevent condensation. If condensation is visible on the product surface, the iE should be left at room temperature for more than 8 hours.
- The iE system doesn't have a physical power button. Turn off the scanner by unplugging the power delivery cable.
- The EMISSIONS characteristics of this equipment make it suitable for use in industrial areas and hospitals (CISPR 11 Class A). If used in a residential environment (for which CISPR 11 Class B is normally required), this equipment might not offer adequate protection to radio-frequency communication services.
- Avoid pulling on the power delivery cables used with the iE system.
- Do not touch the connectors of the device and the patient simultaneously.

4.6 Eye Safety

WARNING

- The iE system projects a bright light from its tip during scanning.
- The bright light projected from the tip of the iE is not harmful to the eyes. However, you should not look directly at the bright light nor aim the light beam into the eyes of others. Generally, intense light sources can cause eyes to become brittle and the likelihood of secondary exposure is high. As with other intense light source exposure, you may experience a temporary reduction in visual acuity, pain, discomfort, or visual impairment, all of which increase the risk of secondary accidents.
- There is an LED that emits UV-C wavelengths inside the iE handpiece. It is irradiated only inside the iE handpiece and does not go outside. The blue light visible inside the iE handpiece is for guidance, not UV-C light. It is harmless to the human body.
- The UV-C LED operates with a wavelength of 270 – 285 nm.
- Disclaimer for Risks Involving Patients with Epilepsy The Medit iE should not be used on patients that have been diagnosed with epilepsy due to the risk of seizures and injury. For the same reason, dental staff who have been diagnosed with epilepsy should not operate the Medit iE.

4.7 Explosion Hazards

WARNING

- The iE system is not designed to be used near flammable liquids, gases, or in environments with high oxygen concentrations.
- There is a risk of explosion if you use the iE system near flammable anesthetics.

4.8 Pacemaker and ICD Interference Risk

WARNING

- Implantable Cardioverter Defibrillators (ICDs) and pacemakers may have interference due to some devices.
- Maintain a moderate distance from the patient's ICD or pacemaker when using the iE system.
- For more information on peripherals used with iE, check the respective manufacturer's manuals.

4.9 Cybersecurity Safety

- If a cybersecurity incident occurs, stop using the scanner and software immediately. Power off the scanner and log out of the software.
- Promptly report the incident to our support team through email, phone, or other available means of contact. Please refer to the last page of the User Guide for contact information.
- When reporting an incident, please provide as much information as possible, including the time of occurrence and any unusual behavior you noticed. This information will assist us in resolving the issue quickly.

5. Electro-Magnetic Compatibility Information

5.1 Electromagnetic Emissions

The iE system is intended for use in the electromagnetic environment as specified below. The customer or the user of the iE system should ensure that it is used in such an environment.

Guidance and Manufacturer's Declaration – Electromagnetic Emission		
Emission Test	Compliance	Electromagnetic Environment – Guidance
RF Emissions CISPR 11	Group 1	The iE uses RF energy only for its internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment. The iE system is suitable for use in all establishments other than domestic establishments and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes.
RF Emissions CISPR 11	Class A	
Harmonic Emissions IEC 61000-3-2	Class A	
Voltage Fluctuations / Flicker Emissions	Complies	

WARNING
This iE system is intended for use by healthcare professionals only. This equipment/system may cause radio interference or may disrupt the operation of nearby equipment. It may be necessary to take mitigation measures, such as re-orienting or relocating the iE or shielding the location.

5.2 Electromagnetic Immunity

Guidance 1

The iE system is intended for use in the electromagnetic environment as specified below. The customer or the user of the iE system should ensure that it is used in such an environment.

Guidance and Manufacturer's Declaration – Electromagnetic Immunity			
Immunity Test	IEC 60601 Test Level	Compliance Level	Electromagnetic Environment – Guidance
Electrostatic Discharge (ESD) IEC 61000-4-2	± 8 kV contact ± 2 kV, ± 4 kV, ± 8 kV, ± 15 kV air	± 8 kV contact ± 2 kV, ± 4 kV, ± 8 kV, ± 15 kV air	Floors should be made of wood, concrete, or ceramic tiles. If floors are covered with a synthetic material, relative humidity of at least 30% is recommended.
Electrical Fast Transient / Burst IEC 61000-4-4	± 2 kV for power supply lines ± 1 kV for input/output lines	± 2 kV for power supply lines ± 1 kV for input/output lines	The mains power quality should be that of a typical commercial or hospital environment.
Surge IEC 61000-4-5	± 0.5 kV, ± 1 kV differential mode	± 0.5 kV, ± 1 kV differential mode	The mains power quality should be that of a typical commercial or hospital environment.
	± 0.5 kV, ± 1 kV, ± 2 kV common mode	± 0.5 kV, ± 1 kV, ± 2 kV common mode	
Voltage Dips, Short Interruptions, and Voltage Variations on Power Supply Input Lines IEC 61000-4-11	0% Ur (100% dip in Ur) for 0.5/1 cycles 70% Ur (30% dip in Ur) for 25/30 cycles 0% Ur (100% dip in Ur) for 250/300 cycles	0% Ur (100% dip in Ur) for 0.5/1 cycles 70% Ur (30% dip in Ur) for 25/30 cycles 0% Ur (100% dip in Ur) for 250/300 cycles	The mains power quality should be that of a typical commercial or hospital environment. If the user of the iE system requires continued operation during power mains interruptions, it is recommended that the iE system be powered from an uninterruptible power supply or a battery.
Power Frequency Magnetic Fields (50/60Hz) IEC 61000-4-8	30 A/m	30 A/m	Power frequency magnetic fields should be at levels characteristic of a location in a typical commercial or hospital environment.
Proximity Magnetic Fields in the Frequency Range 9 kHz to 13.56 MHz Immunity IEC 61000-4-39	8 A/m 30 kHz CW modulation	8 A/m 30 kHz CW modulation	Resistance to magnetic fields was tested and applied only to surfaces of enclosures or accessories accessible during intended use.
	65 A/m 134.2 kHz PM 2.1 kHz	65 A/m 134.2 kHz PM 2.1 kHz	
	7.5 A/m 13.56 MHz PM 50 kHz	7.5 A/m 13.56 MHz PM 50 kHz	

NOTE: Ur is the main voltage (AC) prior to the application of the test level.

• **Guidance 2**

Recommended Separation Distances Between Portable and Mobile Communication Equipment and the iE		
Rated Maximum Output Power of the Transmitter [W]	Separation Distance According to the Frequency of Transmitter [M]	
	IEC 60601-1-2:2014	
	150 kHz to 80 MHz $d = 1.2 \sqrt{P}$	80 MHz to 2.7 GHz $d = 2.0 \sqrt{P}$
0.01	0.12	0.20
0.1	0.38	0.63
1	1.2	2.0
10	3.8	6.3
100	12	20

For transmitters rated at a maximum output power not listed above, the recommended separation distance (d) in meters (m) can be estimated using the equation applicable to the frequency of the transmitter, where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.

NOTE 1: At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies.

NOTE 2: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects, and people.

• **Guidance 3**

The iE system is intended for use in the electromagnetic environment specified below. The customer or the user of the iE system should ensure that it is used in such an environment.

Guidance and Manufacturer's Declaration – Electromagnetic Immunity			
Immunity Test	IEC 60601 Test Level	Compliance Level	Electromagnetic Environment – Guidance
Conducted RF IEC 61000-4-6	3 Vrms 150 kHz to 80 MHz Outside ISM Bands amateur	3 Vrms	Portable and mobile RF communications equipment should not be used closer to any part of the Medit iE system, including cables, than the recommended separation distance. This is calculated using the equation applicable to the frequency of the transmitter. Recommended Separation Distance (d): $d = 1.2 \sqrt{P}$ IEC 60601-1-2:2007 $d = 1.2 \sqrt{P}$ 80 MHz to 800 MHz $d = 2.3 \sqrt{P}$ 80 MHz to 2.5 GHz IEC 60601-1-2:2014 $d = 2.0 \sqrt{P}$ 80 MHz to 2.7 GHz Where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer, d is the recommended separation distance in meters (m). Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey should be less than the compliance level in each frequency range. Interference may occur in the vicinity of equipment marked with the following symbol:
	6 Vrms 150 kHz to 80 MHz In ISM Bands amateur	6 Vrms	
Radiated RF IEC 61000-4-3	3 V/m 80 MHz to 2.7 GHz	3 V/m	

NOTE 1: At 80 MHz and 800 MHz, the higher frequency range applies.

NOTE 2: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects, and people.

NOTE 3: The ISM (Industrial, Scientific, and Medical) bands between 150 kHz and 80 MHz are 6.765 MHz to 6.795 MHz; 13.553 MHz to 13.567 MHz; 26.957 MHz to 27.283 MHz; and 40.66 MHz to 40.70 MHz.

- Guidance 4**
The iE system is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled.
Portable RF communications equipment should be used no closer than 30 cm (12 inches) to any part of the iE system. Otherwise, degradation of the performance of this equipment could result.

Guidance and Manufacturer's Declaration – Electromagnetic Immunity					
Immunity Test	Band ¹	Service ¹	Modulation	IEC 60601 Test Level	Compliance Level
Proximity Fields from RF Wireless Communications IEC 61000-4-3	380 – 390 MHz	TETRA 400	Pulse Modulation 18 Hz	27 V/m	27 V/m
	430 – 470 MHz	GMRS 460; FRS 460	FM ±5 kHz Deviation 1 kHz sine	28 V/m	28 V/m
	704 – 787 MHz	LTE Band 13, 17	Pulse Modulation 217 Hz	9 V/m	9 V/m
	800 – 960 MHz	GSM 800/900; TETRA 800; iDEN 820; CDMA 850; LTE Band 5	Pulse Modulation 18 Hz	28 V/m	28 V/m
	1700 – 1990 MHz	GSM 1800; CDMA 1900; GSM 1900; DECT; LTE Band 1, 3, 4, 25; UMTS	Pulse Modulation 217 Hz	28 V/m	28 V/m
	2400 – 2570 MHz	Bluetooth; WLAN 802.11b/g/n; RFID 2450; LTE Band 7	Pulse Modulation 217 Hz	28 V/m	28 V/m
	5100 – 5800 MHz	WLAN 802.11a/n	Pulse Modulation 217 Hz	9 V/m	9 V/m

NOTE: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects, and people.

 WARNING

- Use of the iE adjacent to or on other equipment must be avoided as it may result in improper operation.
If this use is necessary, it is advisable that this and the other equipment be observed to verify that they are operating normally.
- The use of accessories, transducers, and cables other than those specified or provided by the Medit of the iE could result in high electromagnetic emissions or reduced electromagnetic immunity from this equipment and result in improper operation.

¹ For some services, only the uplink frequencies are included.

6. Specifications

Model Name		MO8-iE	
Trade Name	iE		
Packing Unit	1 set		
Rating	5 V ~~, 3 A		
Classifications for Protection Against Electric Shock	Class II, Type BF Applied Parts (Reusable Tip)		
* This product is a medical device.			
Handpiece (including medium-size tip)			
Dimension	223.4 x 36.7 x 35.3 mm (L x W x H)		
Weight	165 g		
Reusable Tip			
Dimension - Full Tip	Large	36.1 x 34.1 x 90.8 mm (W x H x L)	
	Medium	36.1 x 34.1 x 90.4 mm (W x H x L)	
Dimension - Tip Head	Large	26.9 x 19.7 mm (W x H)	
	Medium	22.4 x 16.3 mm (W x H)	
Calibration Tool			
Dimension	160 x 48.5 mm (H x Ø)		
Weight	205 g		
Operating, Storage, and Transport Conditions			
Operating Condition	Temperature	18 – 28°C (64.4 – 82.4°F)	
	Humidity	20 – 75% relative humidity (non-condensing)	
	Air Pressure	800 – 1,100 hPa	
Storage Condition	Temperature	-10 – 50°C (14 – 122°F)	
	Humidity	20 – 80% relative humidity (non-condensing)	
	Air Pressure	800 – 1,100 hPa	
Transport Condition	Temperature	-10 – 50°C (14 – 122°F)	
	Humidity	20 – 80% relative humidity (non-condensing)	
	Air Pressure	620 – 1,200 hPa	
Emission limits per environment			
Environment	Hospital environment		
Conducted and radiated RF Emissions	CISPR 11		

EU REP

Meditrial Srl
Via Po 9 00198, Rome Italy

CH REP

Meditrial Europe Ltd
Bahnhofstrasse 23 6300 Zug, Switzerland

UK REP

Meditrial UK Ltd
19 The Circle, Queen Elizabeth Street, London, SE1 2JE, United Kingdom

Medit Corp.
9F, 10F, 13F, 14F, 16F, 8, Yangpyeong-ro 25-gil, Yeongdeungpo-gu,
Seoul, 07207, Republic of Korea Tel: +82-02-2193-9600

Contact for Product Support
Email: support@medit.com
Tel: +82-02-2193-9600

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- Press the “Retry” button. Device will automatically reboot and try to connect again.
- Turn off the device and disconnect all cables from the PC. Connect cables again and restart the device.
- If connection problem persists, please install the program again.
- If connection problem persists, please contact your service provider.

Cancel Retry

Turn Off the External Light

External light was detected while scanning

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