MEDIT Link

ME-UG-702H
Revision 13 (29.09.2026)
SW version 3.5.0

Table of contents

Medit Link

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Administrator

Dentist, Manager

Have the authority to create a new dental clinic account and have full access to Medit Link.

- ✓ Scan/Order/Manage orders
- ✓ Manage members

Register

Member

Dental Hygienist, Dental Staff

Belong to the dental clinic their admin created and have limited access to Medit Link after the admin's approval.

- ✓ Scan/Order/Manage orders

Register

Lab

Administrator

Laboratory Director, Manager

Have the authority to create a new dental lab account and have full access to Medit Link.

- ✓ Scan/Manage ordered cases
- ✓ Manage members

Register

Member

Dental Technician, Laboratory Staff

Belong to the dental lab their admin created and have limited access to Medit Link after the admin's approval.

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Email Address

Check

 Please ensure your email address is active. Medit Link will send a verification email.

 Registering with a disposable and/or fake e-mail address can make you vulnerable to security issues and make it difficult to receive information from Medit Link.

Password

 Password must be 8 to 16 characters long, including a combination of at least three letters, numbers, and special characters. (Passwords are accepted in English language only.)

Confirm Password

 Please type in the same password as above.

Select Country/Region

☐ I agree to all of the Terms and Conditions.

☐ Personal Information Collection Agreement (Mandatory)

View

☐ Terms & Conditions (Mandatory)

View

☐ Agree to the Marketing and Promotion Agreement (optional)

View

 Check this option to receive marketing emails, which include product updates, software news, and promotional events.

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
Welcome, stage-test-wh!

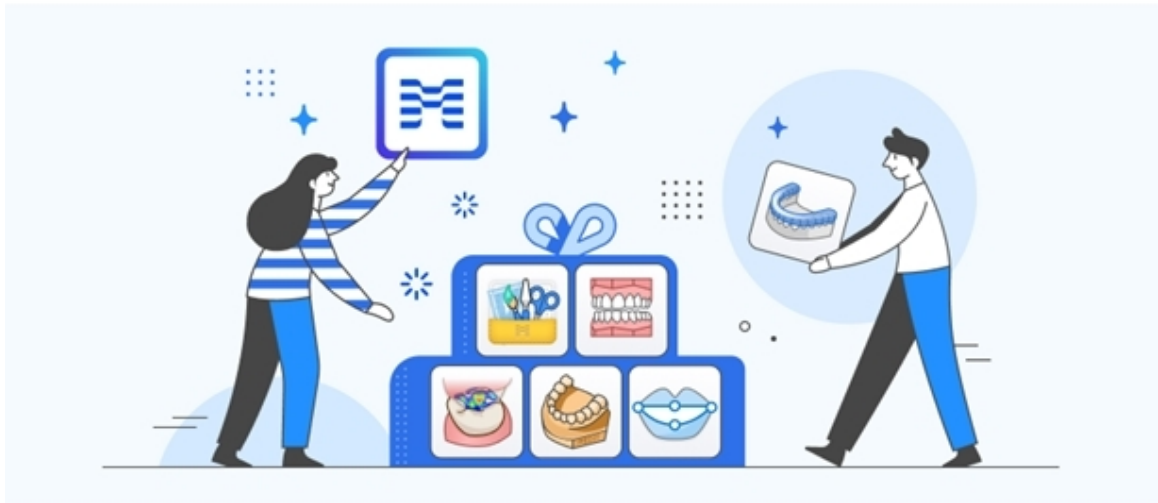Utilize Various Apps

Download Medit Link to your desktop to utilize scanning programs and various apps.
All data generated is synced to the cloud and is also available on the web.

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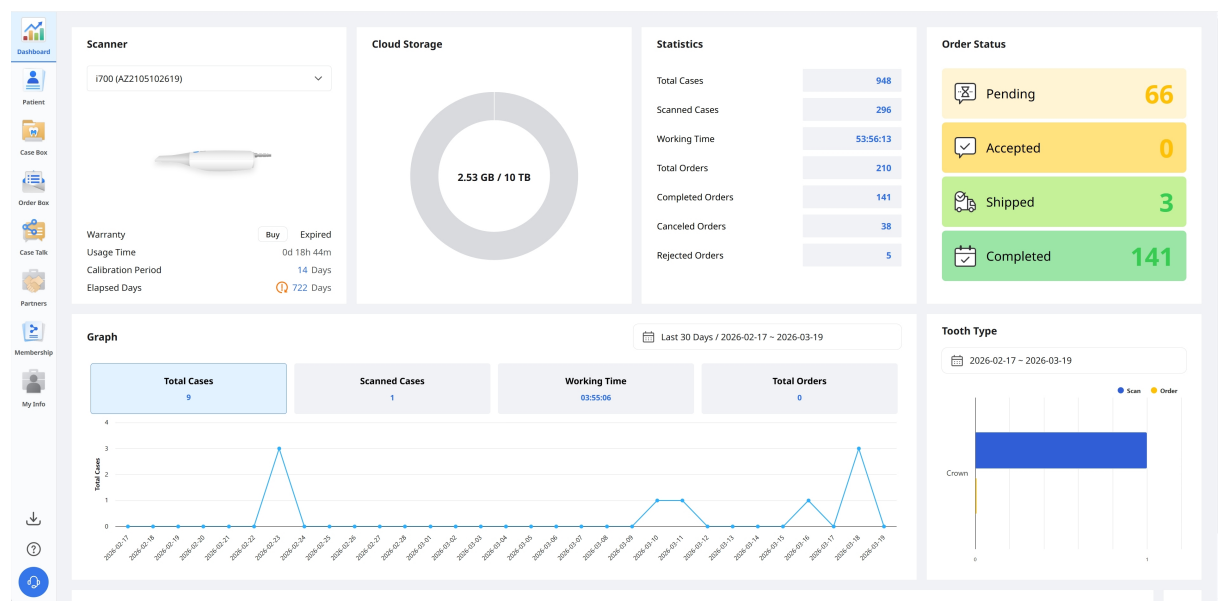
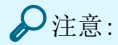
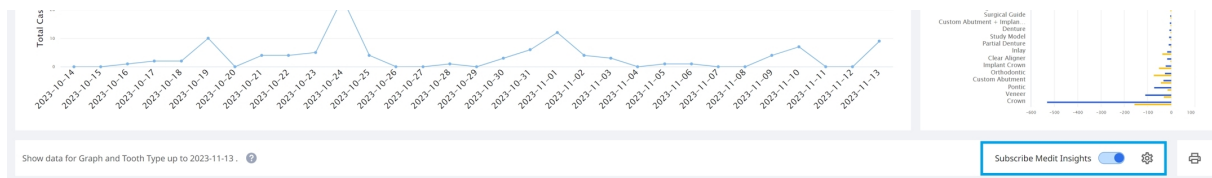
☒ Medit Link account


☐ Another email address

The report will be sent on the 2nd day of every month.

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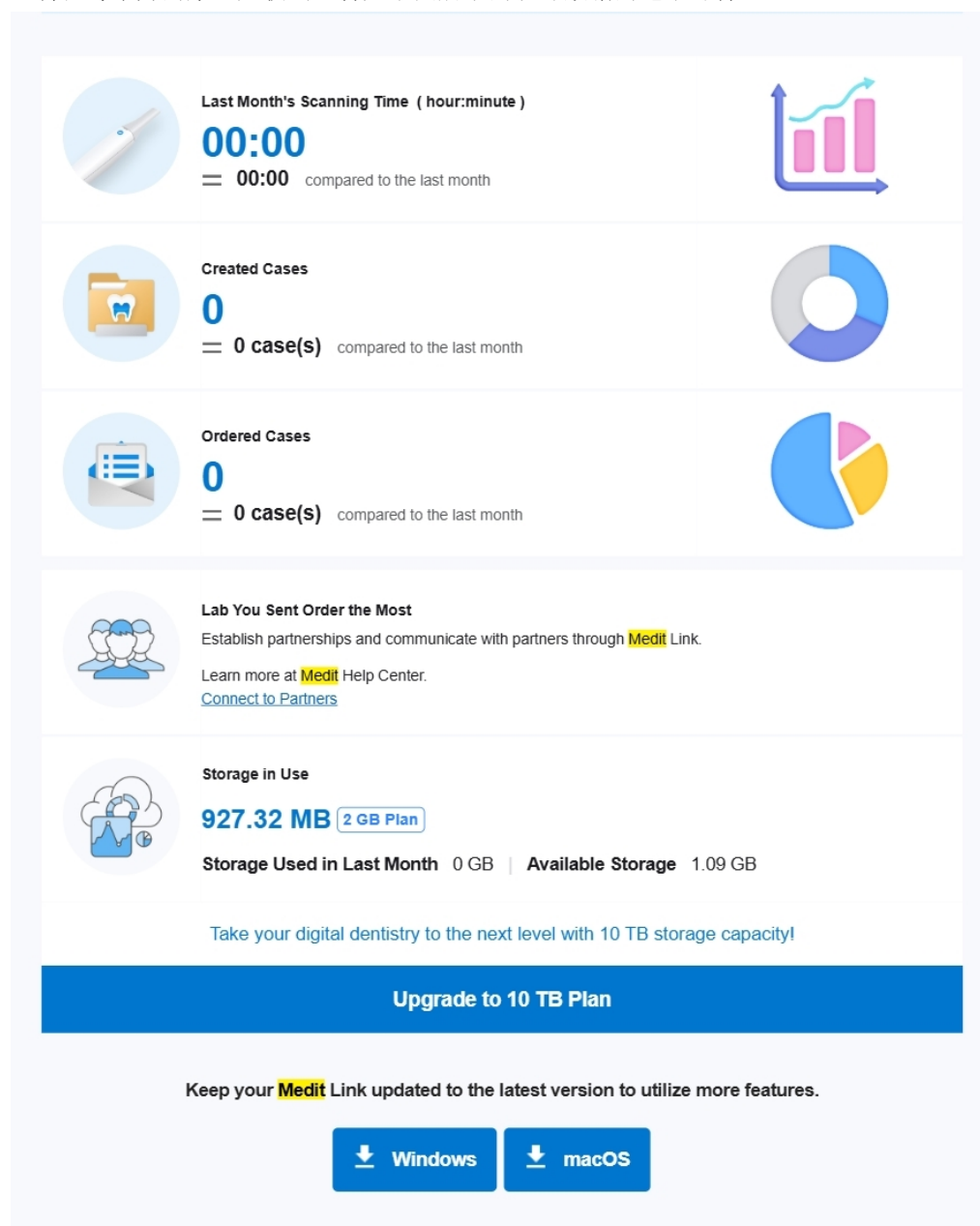This is an automatically generated email. Please do not reply to this email.
 Please contact us at insights@medit.com for more information.
[Unsubscribe](#) If you no longer wish to receive this email, click "Unsubscribe." You can change the email subscription settings in **Medit Link Dashboard**.
 Copyright © **Medit** Corp. All rights reserved.

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
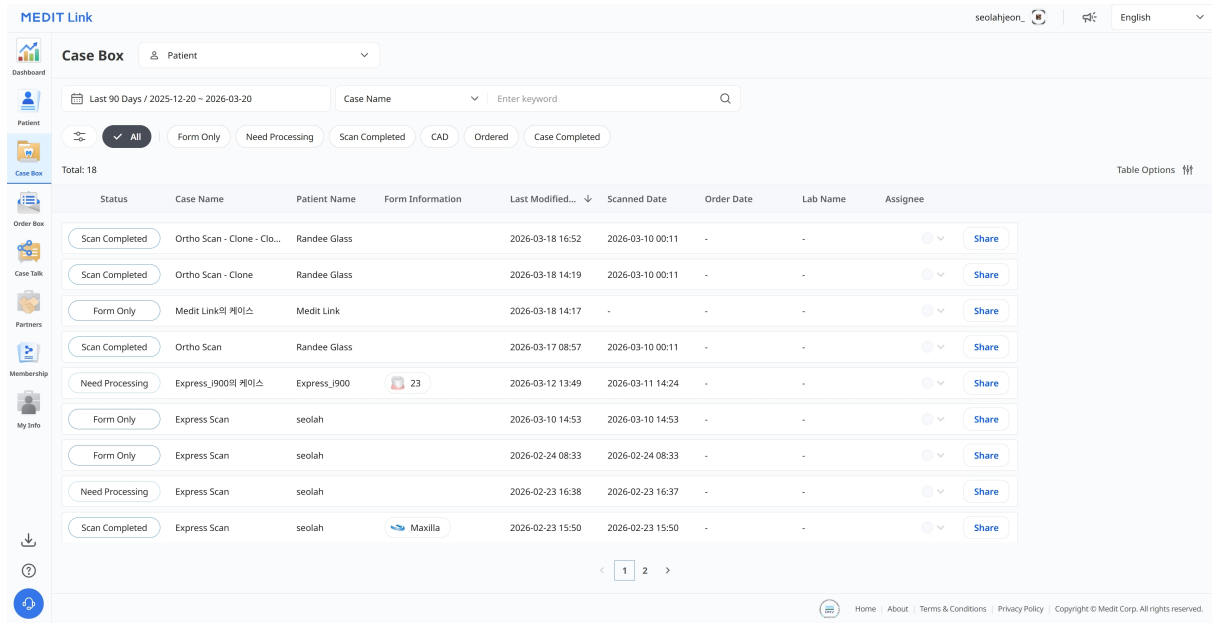
The screenshot shows the Medit Link Case Box interface. At the top, there's a header with the Medit Link logo and user information. Below it, a sidebar on the left contains navigation icons for Dashboard, Patient, Case Box, Order Box, Case Task, Partners, Membership, and My Info. The main area is titled 'Case Box' and features a search bar, a date range selector (Last 90 Days / 2025-12-20 ~ 2026-03-20), and a dropdown for 'Case Name'. Below these are filter buttons: All (selected), Form Only, Need Processing, Scan Completed, CAD, Ordered, and Case Completed. A table lists 18 cases with columns for Status, Case Name, Patient Name, Form Information, Last Modified, Scanned Date, Order Date, Lab Name, and Assignee. Each row includes a 'Share' button. The bottom of the interface has a footer with links to Home, About, Terms & Conditions, Privacy Policy, and Copyright information.

Status	Case Name	Patient Name	Form Information	Last Modified...	Scanned Date	Order Date	Lab Name	Assignee
Scan Completed	Ortho Scan - Clone - Clo...	Randee Glass		2026-03-18 16:52	2026-03-10 00:11	-	-	Share
Scan Completed	Ortho Scan - Clone	Randee Glass		2026-03-18 14:19	2026-03-10 00:11	-	-	Share
Form Only	Medit Link의 케이스	Medit Link		2026-03-18 14:17	-	-	-	Share
Scan Completed	Ortho Scan	Randee Glass		2026-03-17 08:57	2026-03-10 00:11	-	-	Share
Need Processing	Express_i900의 케이스	Express_i900	23	2026-03-12 13:49	2026-03-11 14:24	-	-	Share
Form Only	Express Scan	seolah		2026-03-10 14:53	2026-03-10 14:53	-	-	Share
Form Only	Express Scan	seolah		2026-02-24 08:33	2026-02-24 08:33	-	-	Share
Need Processing	Express Scan	seolah		2026-02-23 16:38	2026-02-23 16:37	-	-	Share
Scan Completed	Express Scan	seolah	Maxilla	2026-02-23 15:50	2026-02-23 15:50	-	-	Share

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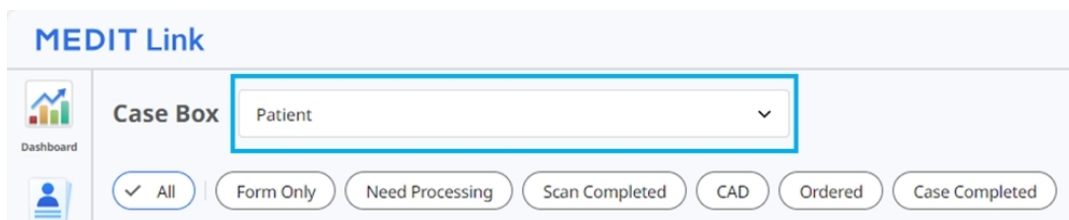The screenshot shows the Medit Link Case Box interface with a focus on the patient selection dropdown menu. The dropdown is open, showing 'Patient' as the selected option. Below the dropdown are filter buttons: All (selected), Form Only, Need Processing, Scan Completed, CAD, Ordered, and Case Completed.

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
The screenshot shows the Medit Link Case Box interface. The top navigation bar includes 'Dashboard', 'Patient', 'Case Box', 'Order Box', 'Case Talk', 'Partners', 'Membership', and 'My Info'. The 'Case Box' section is active, displaying a list of cases. The table has columns for Status, Case Name, Patient Name, Form Information, Last Modified Date, Scanned Date, Order Date, and Lab Name. Each row has a 'Share' button on the right.

Status	Case Name	Patient Name	Form Information	Last Modified Date	Scanned Date	Order Date	Lab Name	Share
Form Only	test's Case	test		2023-11-03 12:50	-	-	-	Share
Scan Completed	wonhoonchoi-test1's Case - Clone - Clone	wonhoonchoi-test1	23, 24, 25, 26	2023-11-02 19:37	2023-11-02 09:01	-	-	Share
Scan Completed	wonhoonchoi-test1's Case - Clone	wonhoonchoi-test1	23, 24, 25, 26	2023-11-02 19:21	2023-11-02 09:01	-	-	Share
Ordered	wonhoonchoi-test1's Case	wonhoonchoi-test1	23, 24, 25, 26	2023-11-02 19:21	2023-11-02 09:01	2023-11-02 ...	Wonhoonchoi_lab2	Share
Scan Completed	G11598492311021004395	test test (12345)	21	2023-11-02 19:09	-	-	-	Share
Scan Completed	test's Case	test		2023-11-01 13:32	-	-	-	Share
Scan Completed	G11598492311010319123	test test (test)	21	2023-11-01 12:20	-	-	-	Share
Scan Completed	test's Case	test		2023-10-26 17:25	2023-10-26 17:24	-	-	Share
Scan Completed	wonhoonchoi-test1's Case	wonhoonchoi-test1		2023-10-26 17:22	2023-10-26 17:21	-	-	Share
Ordered	test's Case	test	14, 15, 16, 17	2023-10-24 15:45	-	2023-10-24 ...	Wonhoonchoi_lab2	Share

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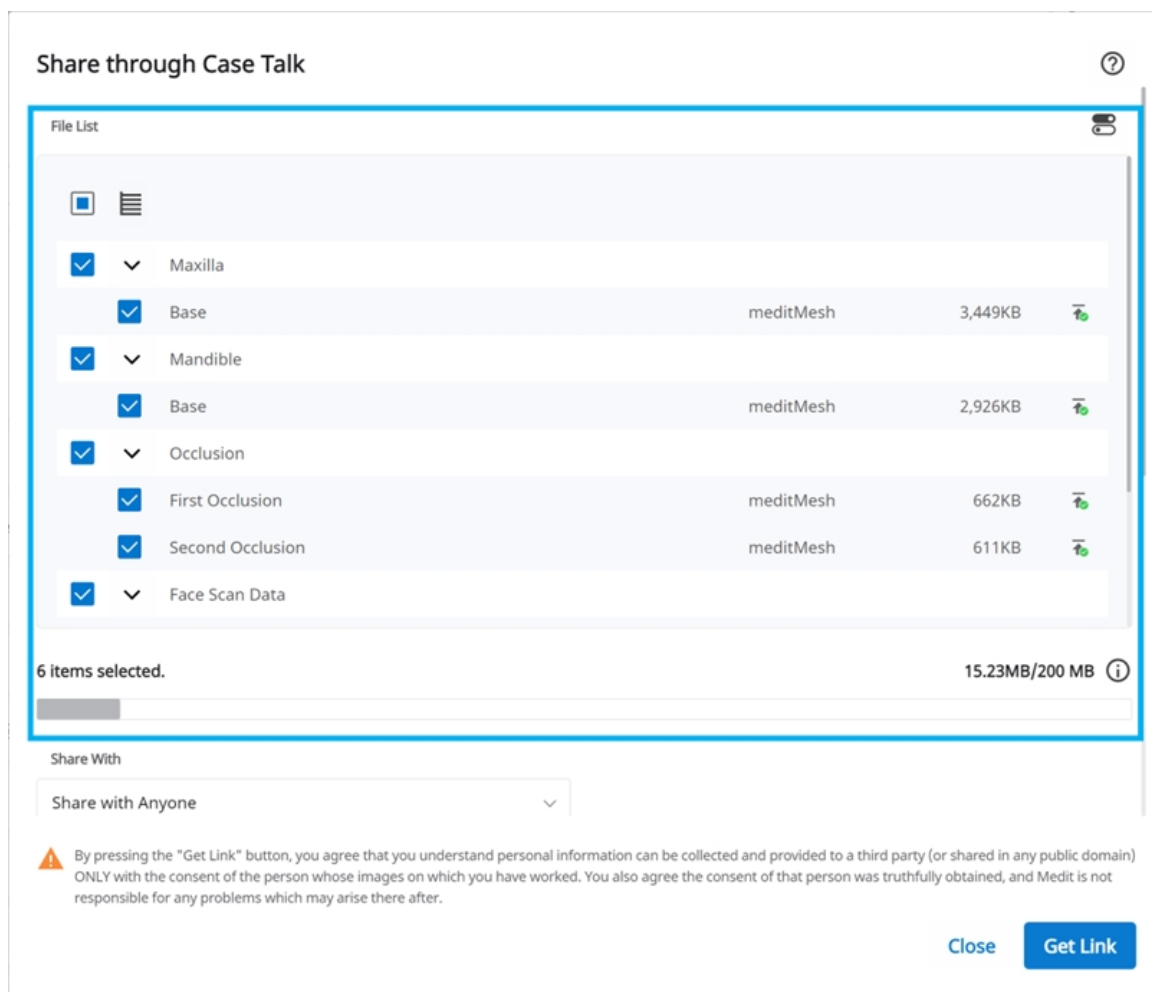Invite Partner

Invite your desired partner to join Medit Link. When the lab accepts your invitation, a partnership is automatically formed.

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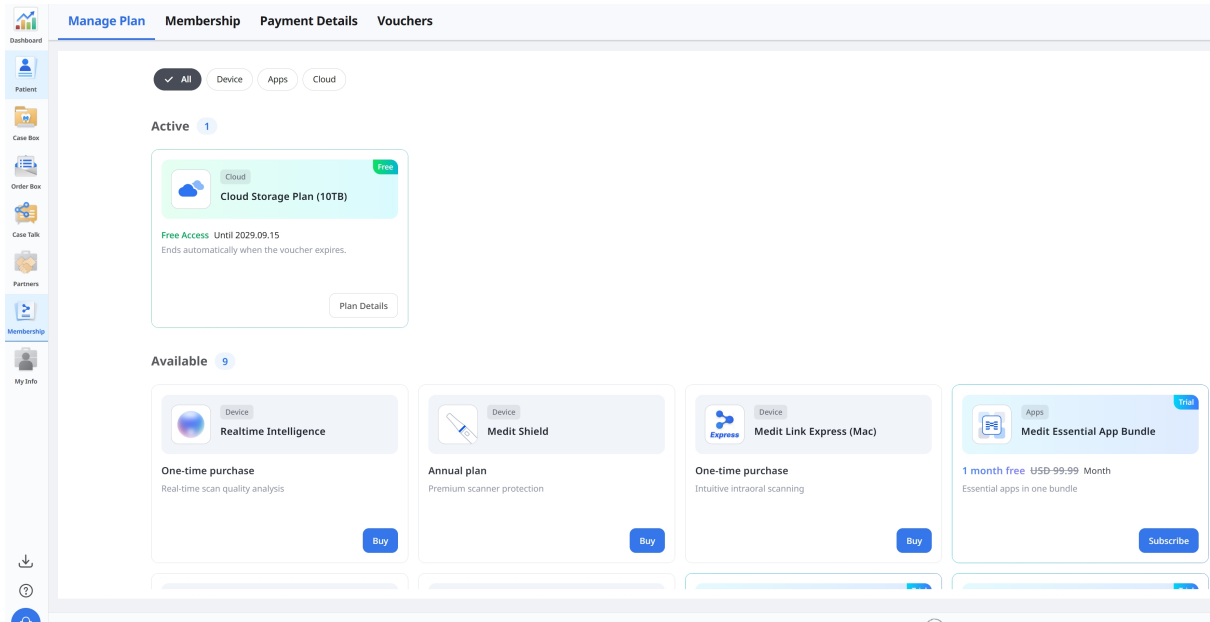
[Manage Plan](#) **Membership** **Payment Details** **Vouchers**

✓ All Device Apps Cloud

Active 1

Cloud

Free

Cloud Storage Plan (10TB)

Free Access Until 2029.09.15
Ends automatically when the voucher expires.

Plan Details

Available 9

Device

Realtime Intelligence

Device

Me

Plan Details

×

Free Access

Your free access is active.

Product	Cloud Storage Plan (10TB)
Status	Free
Expiration Date	2029.09.15
Regular Price	USD 9.99/Month

⚠ Important


- You can't subscribe to this product while the voucher is active.
- Free access ends automatically when the voucher expires. A paid subscription is required to continue.

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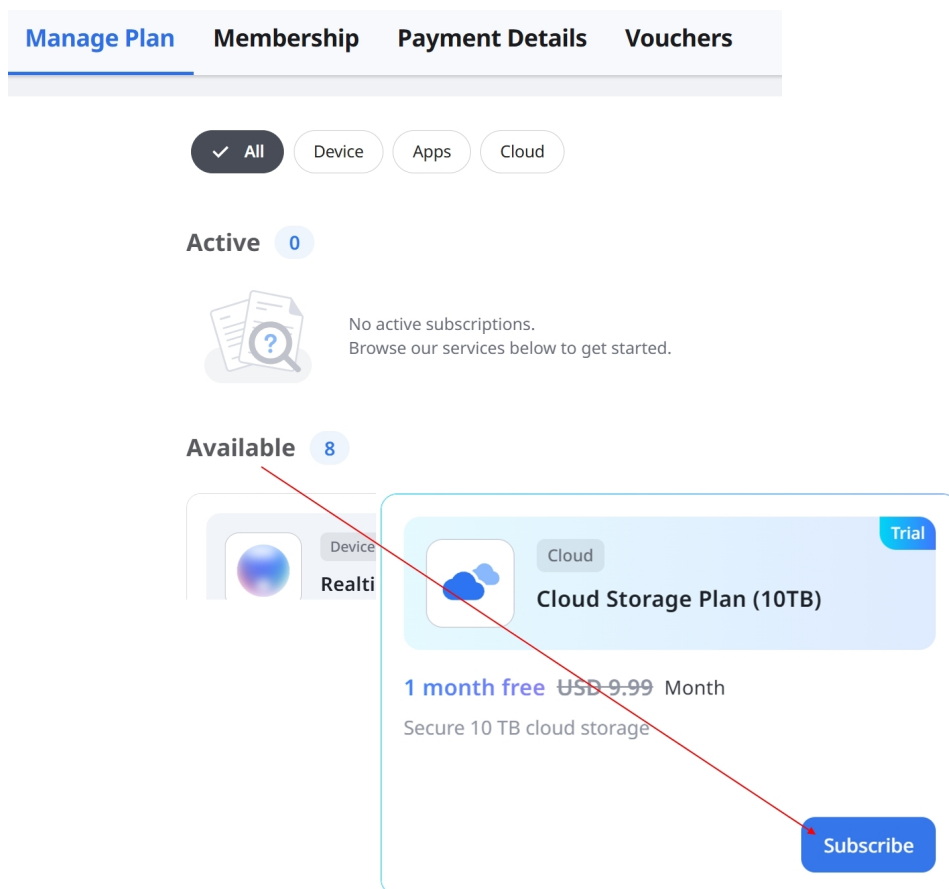~~USD 9.99~~ Month

✓ Includes 10 TB of cloud storage

✓ Securely store large cases and other data

✓ Streamline your Medit Link data management in the cloud

✓ Store design files without separate local backups

• Actual charges may vary depending on your bank or card issuer's fees and policies.

• For more information about Cloud Storage Plan (10TB), see [here](#).

Trial

Cloud Storage: First Month Free

1 month free >

☐ I agree to automatic renewal and recurring charges.

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Subscribe to Medit Cloud Storage 10TB Plan

\$0.00

Then \$9.99 per day after coupon expires

Medit Cloud Storage 10TB Plan for Clinics

Medit Cloud Storage Plan

Billed daily

\$9.99

Subtotal

\$9.99

1 Month Free for New Subscribers

100% off for a day

-\$9.99

Tax

\$0.00

Total due today

\$0.00

Pay with card

Email

yhgam0829@gmail.com

Card information

1234 1234 1234 1234

MM / YY

CVC

Cardholder name

Full name on card

Country or region

South Korea

Subscribe

By confirming your subscription, you allow Medit USA Inc to charge you for future payments in accordance with their terms. You can always cancel your subscription.

Powered by [Stripe](#) | [Terms](#) | [Privacy](#)

Subscribe

Your Cloud Storage Plan (10TB) subscription is now active.

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Manage Plan

Promotion Applied

Regular pricing will apply after the promotion ends.

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
Payment Date	Product	Payment Type	Duration	Details	Amount of Payment	Status
2026-09-21	Cloud Storage Plan (10TB)	Subscription	2026-09-21 - 2026-10-21	Monthly payment	USD 0.00	Payment Completed
2026-09-21	Cloud Storage Plan (10TB)	Subscription	2026-09-21 - 2026-10-21	Monthly payment	USD 0.00	Payment Completed
2026-09-21	Cloud Storage Plan (10TB)	Subscription	2026-09-21 - 2026-10-21	Monthly payment	USD 0.00	Payment Completed
2026-09-21	Cloud Storage Plan (10TB)	Subscription	2026-09-21 - 2026-10-21	Monthly payment	USD 0.00	Payment Completed
2026-09-21	Cloud Storage Plan (10TB)	Subscription	2026-09-21 - 2026-10-21	Monthly payment	USD 0.00	Payment Completed

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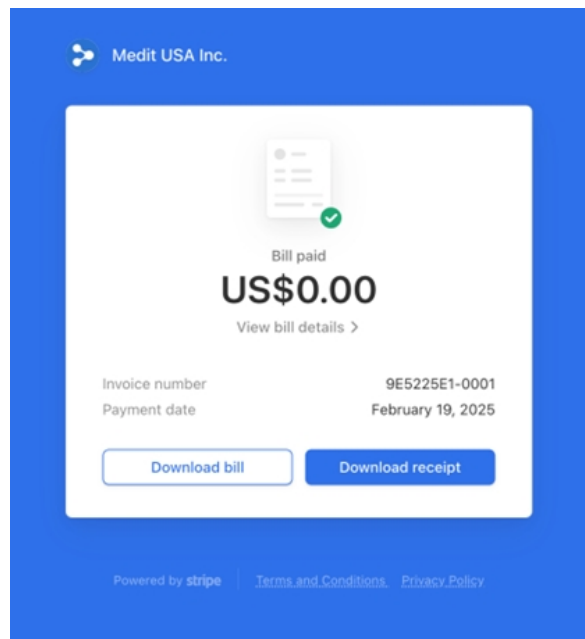ITEM

FEBRUARY 19 - FEBRUARY 20, 2025

Medit Cloud Storage 10TB Plan for Clinics
Quantity 1 piece

US\$0.60

subtotal

US\$0.60

1 Month Free for New Subscribers(100% 할인)

-US\$0.60

Payment amount

US\$0.00

Amount paid

US\$0.00

Remaining amount

US\$0.00

Do you have any questions?Contact Medit USA Inc.

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Apply the voucher at checkout.

Redeem

Activate the voucher for immediate access.

Order

Apply the voucher when checking out in the app.

Vouchers

13

All

Product	Promotion	Code Number	Registration Date	Valid Until	From	To	Status
Cloud Storage Plan	Basic EJ_Basic_Cloud Storage Plan (10TB)	2QJWR0FP	2026-09-21	2026-12-10	2026-10-21	2026-11-21	Used
Cloud Storage Plan	Basic EJ_Basic_Cloud Storage Plan (10TB)	QQ30XJ3T	2026-09-21	2026-12-10	-	-	Registered
Cloud Storage Plan	Basic EJ_Basic_Cloud Storage Plan (10TB)	RSP4GZOI	2026-09-21	2026-12-10	-	-	Registered
Cloud Storage Plan	Basic EJ_Basic_Cloud Storage Plan (10TB)	JYAYH1XH	2026-09-21	2026-12-10	-	-	Registered
Cloud Storage Plan	Basic EJ_Basic_Cloud Storage Plan (10TB)	RSMP5SA4	2026-09-21	2026-12-10	-	-	Registered
Cloud Storage Plan	Redeem EJ_Redeem_Cloud Storage Plan (10TB)	ZGYUDFRY	2026-09-21	2100-01-01	-	-	Registered
Cloud Storage Plan	Redeem EJ_Redeem_Cloud Storage Plan (10TB)	1MGYSVHF	2026-09-21	2100-01-01	2026-09-21	2026-10-21	Used
Cloud Storage Plan	Redeem EJ_Redeem_Expired1	3FWWHOAA	2026-09-21	2026-09-21	-	-	Expired
Cloud Storage Plan	Redeem EJ_Redeem_Expired1	COUR8QKR	2026-09-21	2026-09-21	-	-	Expired
Cloud Storage Plan	Redeem EJ_Redeem_Terminated1	PZPJXKH8	2026-09-21	2026-12-21	2026-09-21	2026-10-21	Used
Cloud Storage Plan	Redeem SQA_Redeem_Cloud Storage Plan_36M	GFZKBWJS	2026-09-15	2100-01-01	2026-09-15	2029-09-15	Activated
Cloud Storage Plan	Redeem SQA_Redeem_Cloud Storage Plan_36M	VE00ZB46	2026-09-14	2100-01-01	2026-09-14	2029-09-14	Used

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Vouchers 2

All

Product	Promotion	Code Number
Cloud Storage Plan	Redeem SQA_Redeem_Cloud Storage Plan_36M	GFZl

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30 Days

⚠ Important

- Once activated, the voucher cannot be canceled.
- You can't subscribe to this product while the voucher is active.

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Language
English ▼

Country/Region
Korea, Republic of ▼

ⓘ Your chosen language is applied to the email sent by MEDIT Link and the notification.

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Edit

Cover and Logo

Looking for digitalized dental laboratories?

Find the most ordered Labs around your clinic with 'Featured Part

1

2

3

Services

Select

Specialties (choose up to 3)

Select

Select the services your business specializes in from the list.

About Us

e.g., Welcome to OOO Dental Clinic.

We are a licensed dental clinic situated in California specializing mainly in aesthetic restorations and implantology.

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The scanner information is automatically registered once your scanner is connected to a Medit scan program.

Home

About

Terms & Conditions

Privacy Policy

Copyright © Medit Corp. All rights reserved.

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MEDIT Link

My Account My Lab My Scanner **Order Forms**

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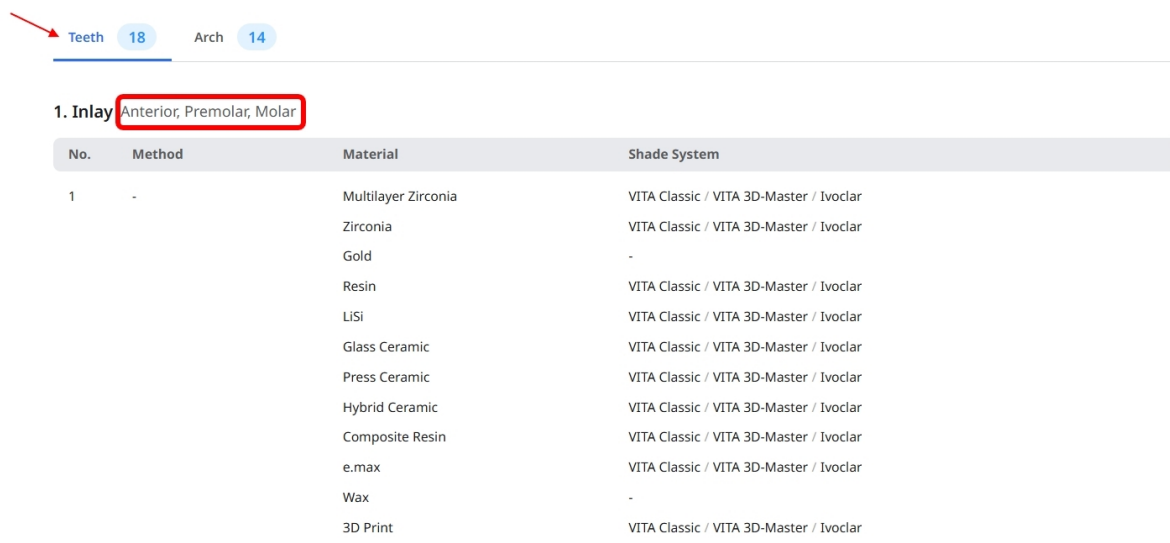
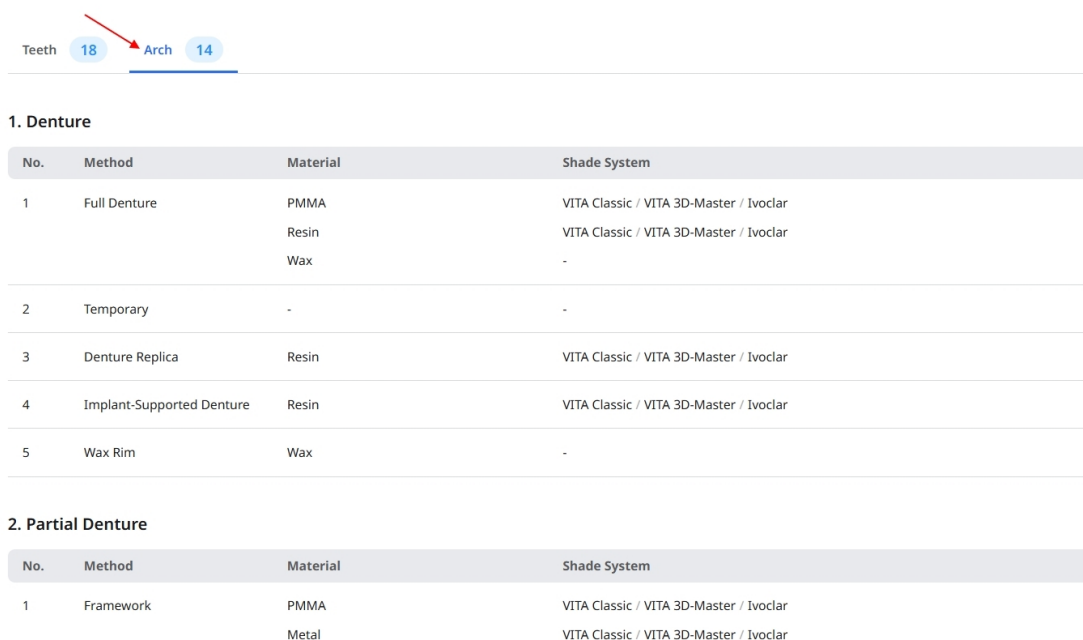
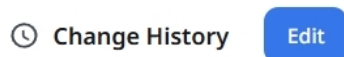
Methods ¹

Search

Screw-Retained **Full-Contour** Offset Cement-Retained

SCRIP Temporary Cut-back Pontic Press Preform Screw-Retained Layered

Screw-Retained PFZ Screw-Retained PFG Screw-Retained PFM Cement-Retained Layered

[More](#)

☐ Do not use

Full-Contour

Materials ²

Search

Multilayer Zirconia Zirconia PMMA Gold Resin

LiSi Glass Ceramic Press Ceramic Hybrid Ceramic Composite Resin e.max

Titanium Metal Stainless Steel Co-Cr (Cobalt Chromium) Ni-Cr (Nickel-Chromium) Wax

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The screenshot shows a software interface for configuring a 'Crown'. It has two main sections: 'Allowed Teeth' and 'Methods'.

Allowed Teeth: Three checkboxes are present: 'Anterior' (checked), 'Premolar' (checked), and 'Molar' (checked).


Methods: A search bar is followed by a grid of buttons: 'Screw-Retained', 'Full-Contour' (highlighted in dark grey), 'Offset', 'Cement-Retained', 'SCRIP', 'Temporary', 'Cut-back', 'Pontic', 'Press', 'Preform', 'Screw-Retained Layered', 'Screw-Retained PFZ', 'Screw-Retained PFG', 'Screw-Retained PFM', and 'Cement-Retained Layered'. A 'More' link is below the grid. A 'Do not use' checkbox is at the bottom left.

A dropdown menu is open for 'Full-Contour', showing a list of materials. The 'Materials' section includes a search bar and buttons for: 'Multilayer Zirconia', 'Zirconia', 'PMMA', 'Gold', 'Resin', 'LiSi', 'Glass Ceramic', 'Press Ceramic', 'Hybrid Ceramic', 'Composite Resin', 'e.max', 'Titanium', 'Metal', 'Stainless Steel', 'Co-Cr (Cobalt Chromium)', 'Ni-Cr (Nickel-Chromium)', 'Wax', 'PEEK', '3D Print', and 'Zirconia HT' (highlighted in dark grey with a red arrow pointing to it). A '+ Add Material' link is at the bottom right of the list. A 'Hide' link is at the bottom right of the dropdown. A 'Do not use' checkbox is at the bottom left of the dropdown.

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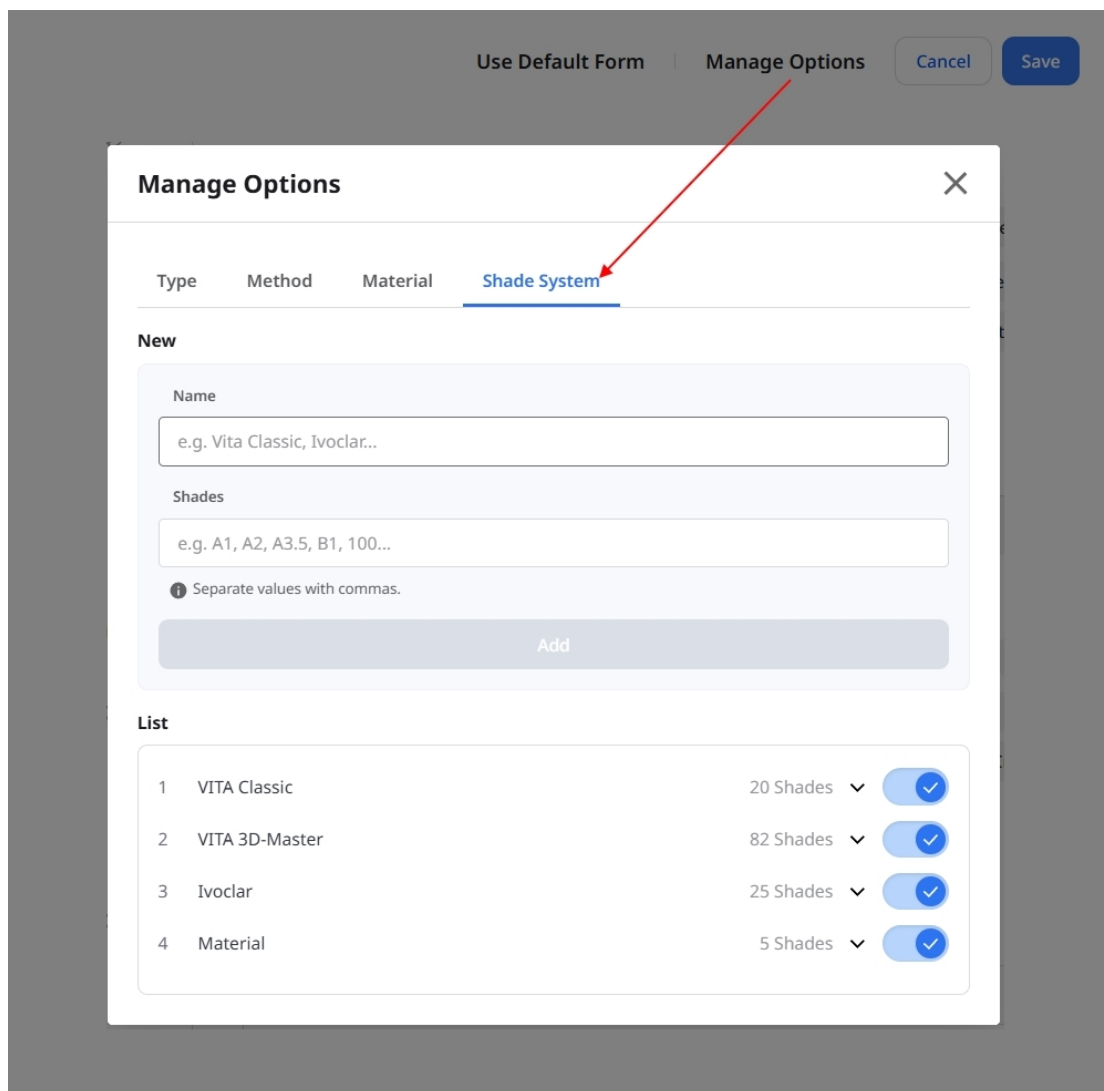
The screenshot shows a 'Manage Options' dialog box with a dark grey background. At the top, there are two tabs: 'Use Default Form' and 'Manage Options', with the latter being active. To the right of the tabs are 'Cancel' and 'Save' buttons. The dialog box itself has a title bar 'Manage Options' and a close button 'X'. Inside, there are four tabs: 'Type', 'Method', 'Material', and 'Shade System', with 'Shade System' being the active tab. Below the tabs, there are two sections: 'New' and 'List'. The 'New' section has two text input fields: 'Name' (with placeholder text 'e.g. Vita Classic, Ivoclar...') and 'Shades' (with placeholder text 'e.g. A1, A2, A3.5, B1, 100...'). Below the 'Shades' field is a small icon and text 'Separate values with commas.' and an 'Add' button. The 'List' section shows a table with four rows: 1 VITA Classic (20 Shades), 2 VITA 3D-Master (82 Shades), 3 Ivoclar (25 Shades), and 4 Material (5 Shades). Each row has a dropdown arrow and a toggle switch.

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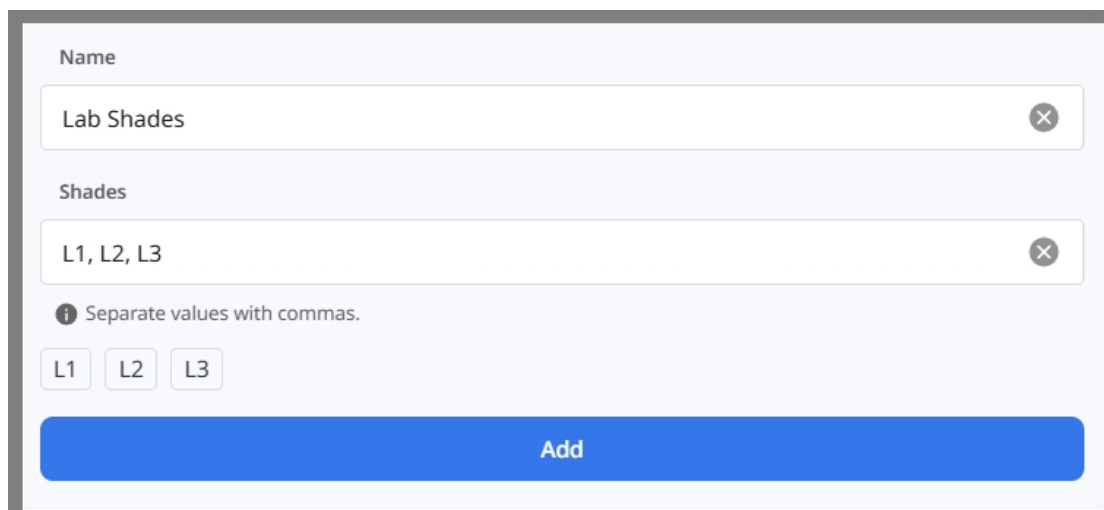The screenshot shows the 'Manage Options' dialog box with the 'Shade System' tab selected. The 'Name' field is filled with 'Lab Shades' and the 'Shades' field is filled with 'L1, L2, L3'. Below the 'Shades' field, there are three buttons labeled 'L1', 'L2', and 'L3'. At the bottom, there is a large blue 'Add' button.

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The screenshot shows a 'Crown' configuration window. Under 'Allowed Teeth', 'Anterior', 'Premolar', and 'Molar' are all checked. In the 'Methods' section, 'Full-Contour' is selected among various options like 'Screw-Retained', 'Offset', and 'Cement-Retained'. Below this, a 'Do not use' checkbox is present. A section titled 'Full-Contour 1' is expanded, showing a 'Materials' list where 'Zirconia HT' is selected. Below the materials, the 'Zirconia HT' section is expanded, showing 'Shade Systems'. In this list, 'Lab Shades' is selected with a checkmark, and a red arrow points to the '+' icon next to it. Other shade systems like 'VITA Classic' and 'VITA 3D-Master' are also listed with their respective shade counts.

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Edit

2026.09.01 | 09:19

> Teeth Custom Abutment + Implant Crown Added

2026.09.01 | 09:19

> Teeth Custom Abutment + Implant Crown Removed

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2026.09.01 | 09:19

▼ Teeth Custom Abutment + Implant Crown Added

Anterior Added

2026.09.01 | 09:19

> Teeth Custom Abutment + Implant Crown Removed

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Notice
NFT Service Termination Notice
2023-04-19

To our users,

We want to thank all Medit Link users who have supported and used our NFT service.

Unfortunately, we must inform you that this service will be terminated following the upcoming program update.

Please see the details below for more information:

- **Feature in Question**
 - Share via NFT (in Medit Link Case Talk)
- **Effective Date**
 - Starting from the update in the first half of 2023 (scheduled for April 27)

※ The update schedule is subject to change.

We'll continuously work to provide better service.

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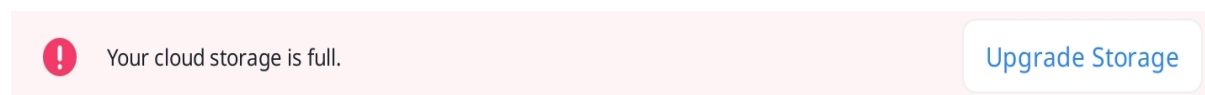
Register Offline Data

Register

Finds and registers offline data on this device that is not registered in Medit Link.

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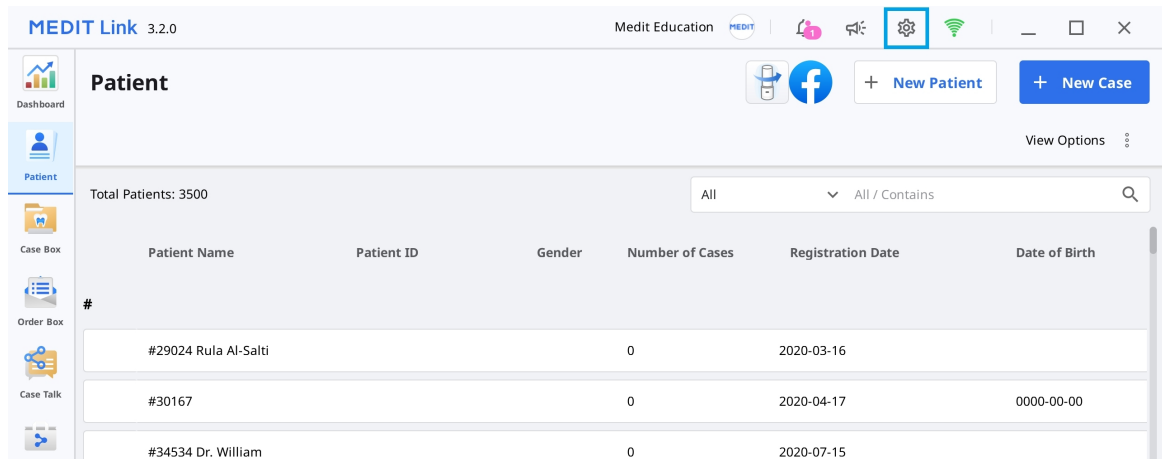
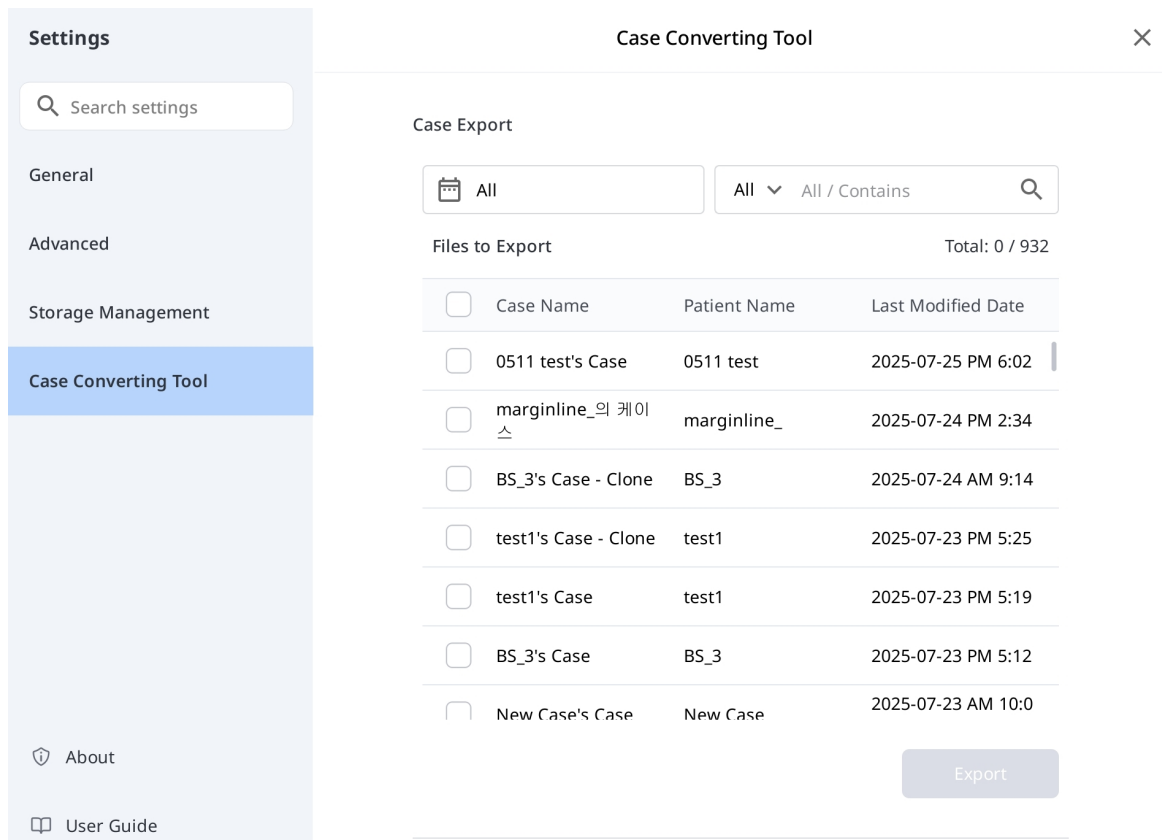
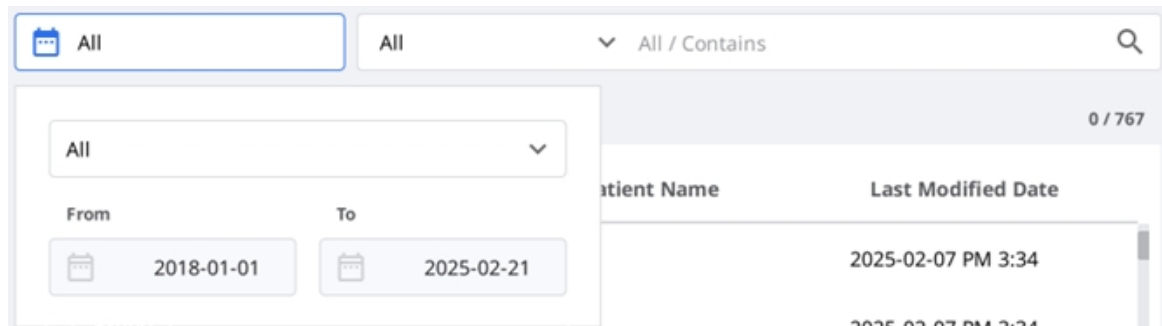
Question

Exporting successfully completed.

Do you want to open the folder containing the exported file?

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
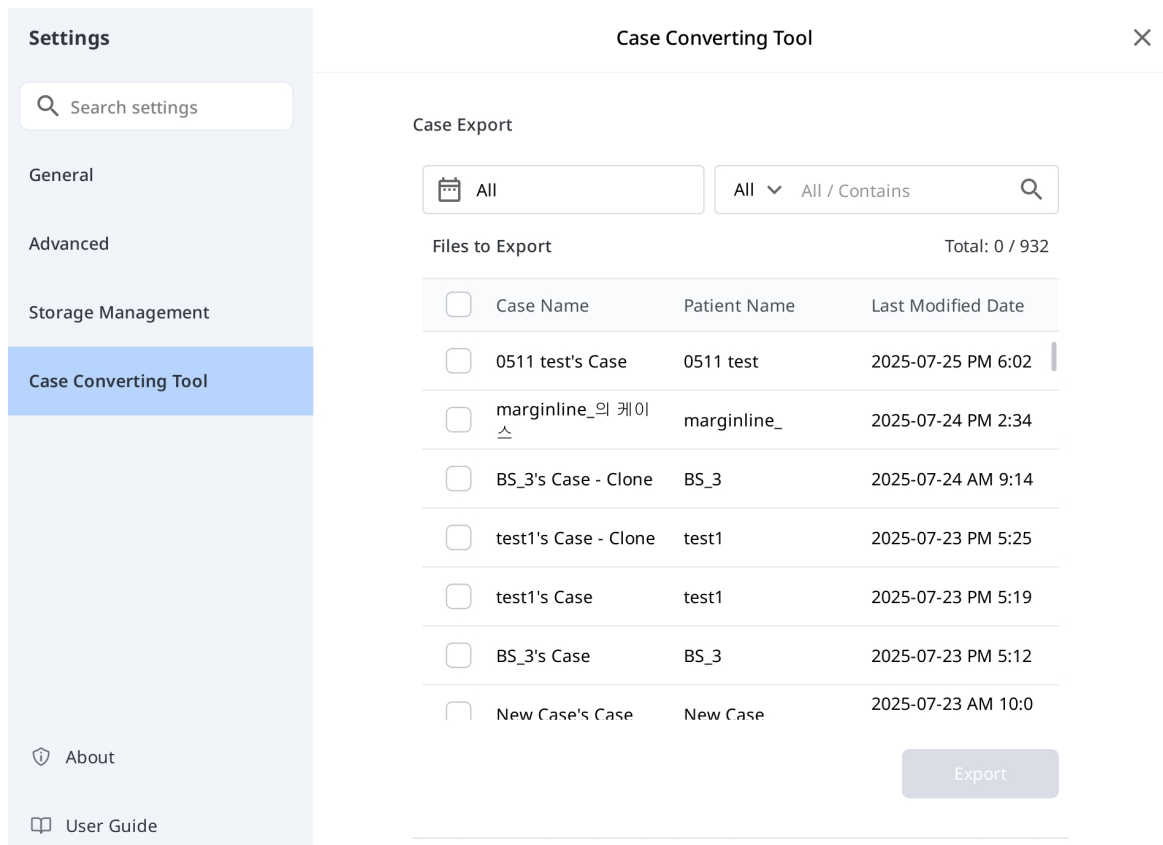
Case Import

A new case is created for every imported case. It contains case data, patient information, and scan file. However, it does not contain order information.

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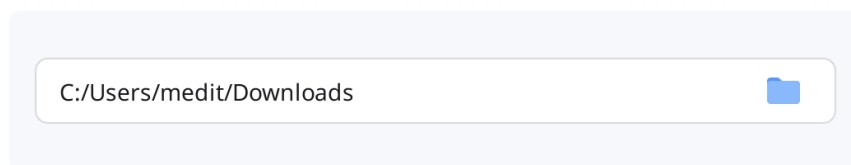[Terms & Conditions](#)

[Open Source License](#)

EC REP Meditrial Srl
Via Po 9 00198, Rome Italy

CH REP Meditrial Europe Ltd
Bahnhofstrasse 23 6300 Zug, Switzerland

UK REP Meditrial UK Ltd
19 The Circle, Queen Elizabeth Street, London, SE1 2JE, United Kingdom

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Order Box												
All / 2018-01-01 ~ 2026-04-08		All	All / Contains									
☒ All ☐ Pending ☐ Canceled ☐ Rejected ☐ Accepted ☐ Shipped ☒ Case Completed												
Total: 254		Table Options										
Status	Case Name	Patient Name	Order Date	Requested Delivery Date	Lab Name	Order ID						
☐ Accepted	texture on's Case - Clone	texture on	2026-04-08 PM 5:40		MEDIT_LAB_STAGE	6012449						
☒ Pending	Express Scan	Seol register	2026-04-02 PM 9:38	2026-04-09 AM 9:00	MEDIT_LAB_STAGE	6012242						
☒ Case Completed	전설아의 케이스	전설아	2025-07-16 PM 5:46	2025-07-17 PM 3:00	MEDIT_LAB_STAGE	6004719						
☒ Case Completed	27	1	2025-03-16 PM 3:11	2025-03-17 PM 3:00	MEDIT_LAB_STAGE	6002656						
☐ Case Completed	3.3.0.136 TEST's Case	3.3.0.136 TEST	2024-04-12 PM 3:25	2024-04-13 PM 2:00	MEDIT_LAB_STAGE	5998987						
☐ Case Completed	test의 케이스	Patient TEST	2024-04-09 PM 3:06	2024-04-10 PM 2:00	MEDIT_LAB_STAGE	5998955						
☐ Case Completed	Mandibular Movement Case - Clone - Clone	MEDIT	2023-06-12 PM 12:12	2023-06-13 PM 2:00	MEDIT_LAB_STAGE	5996680						
☐ Case Completed	Release New_2's Case	Release New_2	2023-04-25 PM 6:47	2023-04-25 PM 2:00	MEDIT_LAB_STAGE	5996180						
☐ Case Completed	cX 문의's Case - Clone	cX 문의	2023-05-16 PM 5:26	2023-05-17 PM 2:00	MEDIT_LAB_STAGE	5996393						
☐ Case Completed	cX 문의's Case	cX 문의	2023-05-15 PM 2:39	2023-05-16 PM 2:00	MEDIT_LAB_STAGE	5996362						
☐ Case Completed	310.187's Case	310.187	2023-05-15 PM 1:33	2023-05-16 PM 2:00	MEDIT_LAB_STAGE	5996360						
☐ Case Completed	310.187 test's Case - Clone			2023-05-03 PM 2:00	MEDIT_LAB_STAGE	5996306						

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Pending ⁴

Order List ⁶

No.	Type	Method	Material	Shade
23	Crown	Anatomic	Multi-Layered ZL...	None

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Watch the Video >

The Doctors' Choice

Medit Intraoral Scanner Trusted by the Majority

Click to find out more >

Featured Apps (Top 5)

Diagnosis/Consultation

Order Placement

Medit Orthodontic Suite

Install

Diagnosis/Consultation

CAM

CAD

Medit ClinicCAD

CAD solution for creating a wide range of dental prostheses

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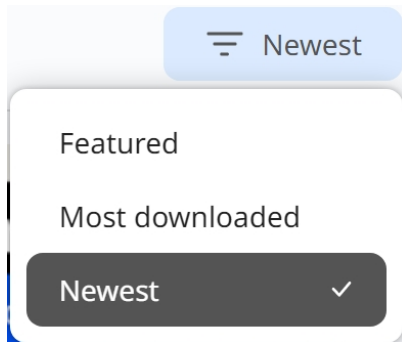
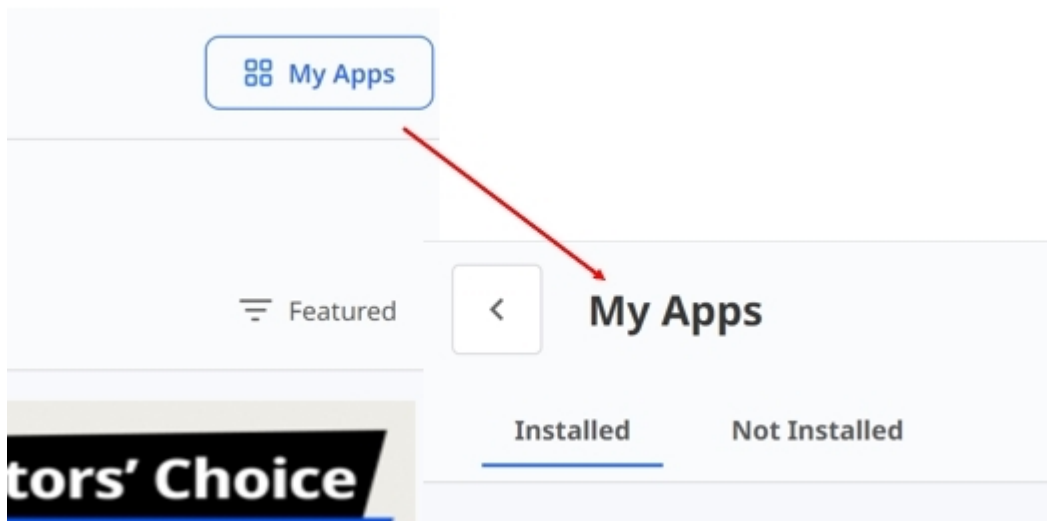
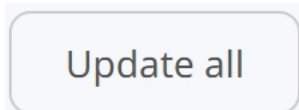
MEDIT Link 3.2.0 wonhoon choi

Trash Box

☐ Restore ☐ Selected Cases: 0 All / Contains

	Status	Case Name	Patient Name	Lab Name	Deletion Date ↓	Deleted By	Size
☐	Scan Completed	Jira, 임플란트 46, 47's Case - Clone - Clone	임플란트 46, 47		2023-04-26 오전 9:12	wonhoon choi	66,654 KB
☐	Form Only	G11598492304180810429	test test		2023-04-18 오후 5:15	wonhoon choi	0 KB
☐	Form Only	G11598492304180806313	test test		2023-04-18 오후 5:15	wonhoon choi	0 KB
☐	Form Only	G11598492304180811349	tes test		2023-04-18 오후 5:15	wonhoon choi	0 KB
☐	Form Only	G11598492304180754083	Glidewell_fastdesign		2023-04-18 오후 5:00	wonhoon choi	0 KB
☐	Form Only	G11598492304180733332	Glidewell_fastdesign		2023-04-18 오후 5:00	wonhoon choi	0 KB
☐	Form Only	G11598492304180648146	test test		2023-04-18 오후 5:00	wonhoon choi	0 KB
☐	Case Completed	Multi Occlusion Normal Case's Case - Clone	Multi Occlusion Normal Case		2023-04-18 오전 10:50	wonhoon choi	23,652 KB
☐	CAD	Multi Occlusion Normal Case's Case - Clone - Clone	Multi Occlusion Normal Case		2023-04-18 오전 10:50	wonhoon choi	23,652 KB
☐	Form Only	Occlusion Test's Case	Occlusion Test		2023-04-17 오후 6:31	wonhoon choi	0 KB
☐	Form Only	Test Occlusion 's Case	Test Occlusion		2023-04-17 오후 6:31	wonhoon choi	0 KB
☐	Form Only	Occlusion Test's Case	Occlusion Test		2023-04-17 오후 6:31	wonhoon choi	0 KB
☐	Form Only	Occlusion Test's Case	Occlusion Test		2023-04-17 오후 6:26	wonhoon choi	0 KB

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Data Type

All

☐	Case Name	Patient Name	Scanned Date↓	Size	Status
☐	BLATRIX Charline's Case	BLATRIX Charline	2025-01-06 PM 6:28	729.30 MB	Upload Required
☐	BLATRIX Charline's Case	BLATRIX Charline	2025-01-06 PM 6:28	729.30 MB	Upload Required
☐	BLATRIX Charline's Case - 복제	BLATRIX Charline	2025-01-06 PM 6:28	733.35 MB	Upload Required
☐	BLATRIX Charline's Case - Clone	BLATRIX Charline	2025-01-06 PM 6:28	734.01 MB	Upload Required
☐	test의 케이스	Patient TEST	2023-04-26 PM 7:58	1.25 GB	Upload Required
☐	i700 HD's Case	i700 HD	2023-04-24 PM 5:23	911.24 MB	Upload Required
☐	Mandibular Movement	MEDIT	2023-04-24 PM 1:59	1.88 GB	Upload Required
☐	Mandibular Movement Case	MEDIT	2023-04-24 PM 1:59	1.75 GB	Upload Required
☐	Mandibular Movement Case - Clo...	MEDIT	2023-04-24 PM 1:59	1.75 GB	Upload Required
☐	서라_Mandibular Movement_TEST	Mandibular Movement_export_test	2023-04-10 PM 3:39	34.73 MB	Upload Required
☐	서라_Mandibular Movement_TEST...	Mandibular Movement_export_test	2023-04-10 PM 3:39	731.43 MB	Upload Required

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
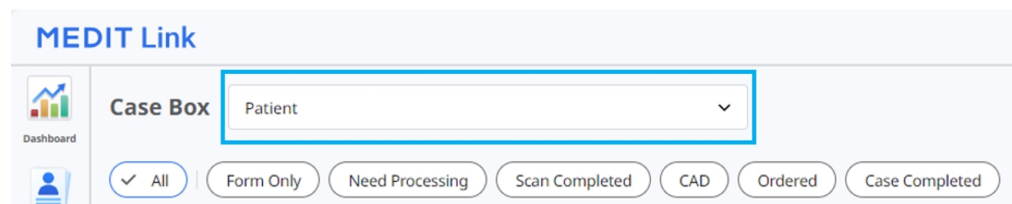
Case Box Patient

Dashboard

✓ All Form Only Need Processing Scan Completed CAD Ordered Case Completed

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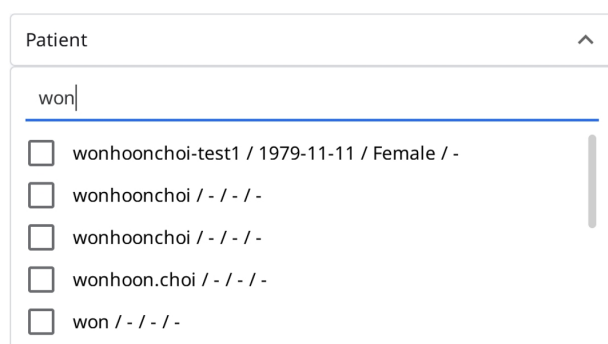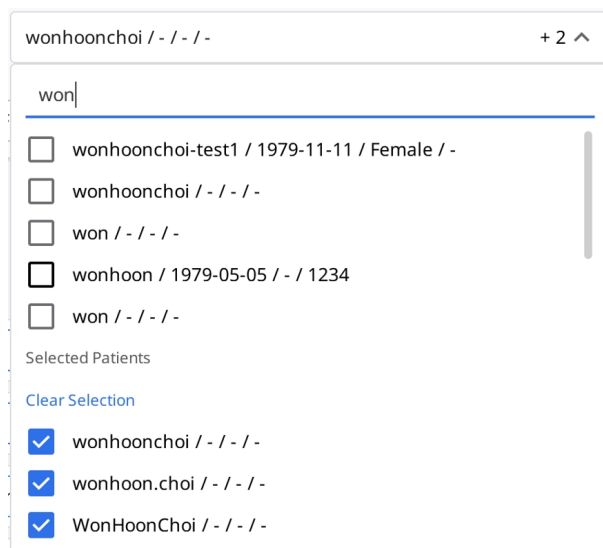
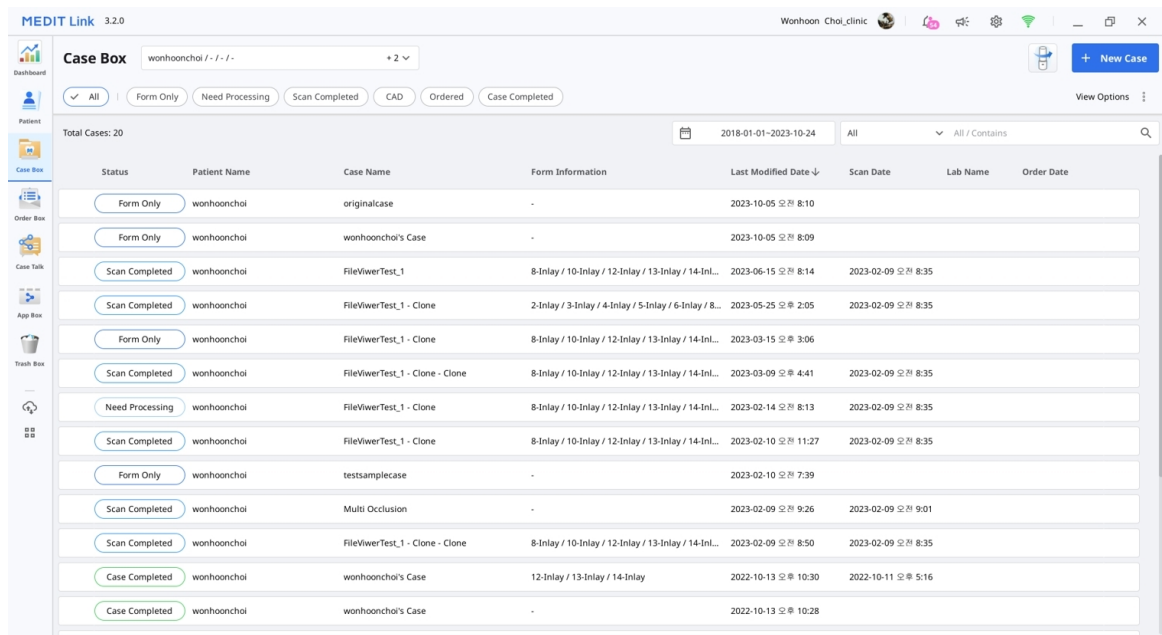
Status	Patient Name	Case Name	Form Information	Last Modified Date	Scan Date	Lab Name	Order Date
Form Only	wonhoonchoi	originalcase	-	2023-10-05 오전 8:10			
Form Only	wonhoonchoi	wonhoonchoi's Case	-	2023-10-05 오전 8:09			
Scan Completed	wonhoonchoi	FileVierTest_1	8-Inlay / 10-Inlay / 12-Inlay / 13-Inlay / 14-Inl...	2023-06-15 오전 8:14	2023-02-09 오전 8:35		
Scan Completed	wonhoonchoi	FileVierTest_1 - Clone	2-Inlay / 3-Inlay / 4-Inlay / 5-Inlay / 6-Inlay / 8...	2023-05-25 오후 2:05	2023-02-09 오전 8:35		
Form Only	wonhoonchoi	FileVierTest_1 - Clone	8-Inlay / 10-Inlay / 12-Inlay / 13-Inlay / 14-Inl...	2023-03-15 오후 3:06			
Scan Completed	wonhoonchoi	FileVierTest_1 - Clone - Clone	8-Inlay / 10-Inlay / 12-Inlay / 13-Inlay / 14-Inl...	2023-03-09 오후 4:41	2023-02-09 오전 8:35		
Need Processing	wonhoonchoi	FileVierTest_1 - Clone	8-Inlay / 10-Inlay / 12-Inlay / 13-Inlay / 14-Inl...	2023-02-14 오전 8:13	2023-02-09 오전 8:35		
Scan Completed	wonhoonchoi	FileVierTest_1 - Clone	8-Inlay / 10-Inlay / 12-Inlay / 13-Inlay / 14-Inl...	2023-02-10 오전 11:27	2023-02-09 오전 8:35		
Form Only	wonhoonchoi	testsamplecase	-	2023-02-10 오전 7:39			
Scan Completed	wonhoonchoi	Multi Occlusion	-	2023-02-09 오전 9:26	2023-02-09 오전 9:01		
Scan Completed	wonhoonchoi	FileVierTest_1 - Clone - Clone	8-Inlay / 10-Inlay / 12-Inlay / 13-Inlay / 14-Inl...	2023-02-09 오전 8:50	2023-02-09 오전 8:35		
Case Completed	wonhoonchoi	wonhoonchoi's Case	12-Inlay / 13-Inlay / 14-Inlay	2022-10-13 오후 10:30	2022-10-11 오후 5:16		
Case Completed	wonhoonchoi	wonhoonchoi's Case	-	2022-10-13 오후 10:28			

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Patient Name*

Case Name*

Creating patient information involves consent from the patient.
Patient information can be shared with 3rd parties while placing an order.

Partner

Select partner

Assignee

Select assignees

You can select up to 5 assignees.

Tags

You can add the entered tag with Enter or Spacebar.

You can add up to 10 tags of 100 characters or less.

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Patient List

Select a patient in the list to set the patient of the case. You can add or manage the patient information in Medit Link.

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
qa_lab_01 Uses Lab Form

Full Services Dental Esthetics Crown & Bridge

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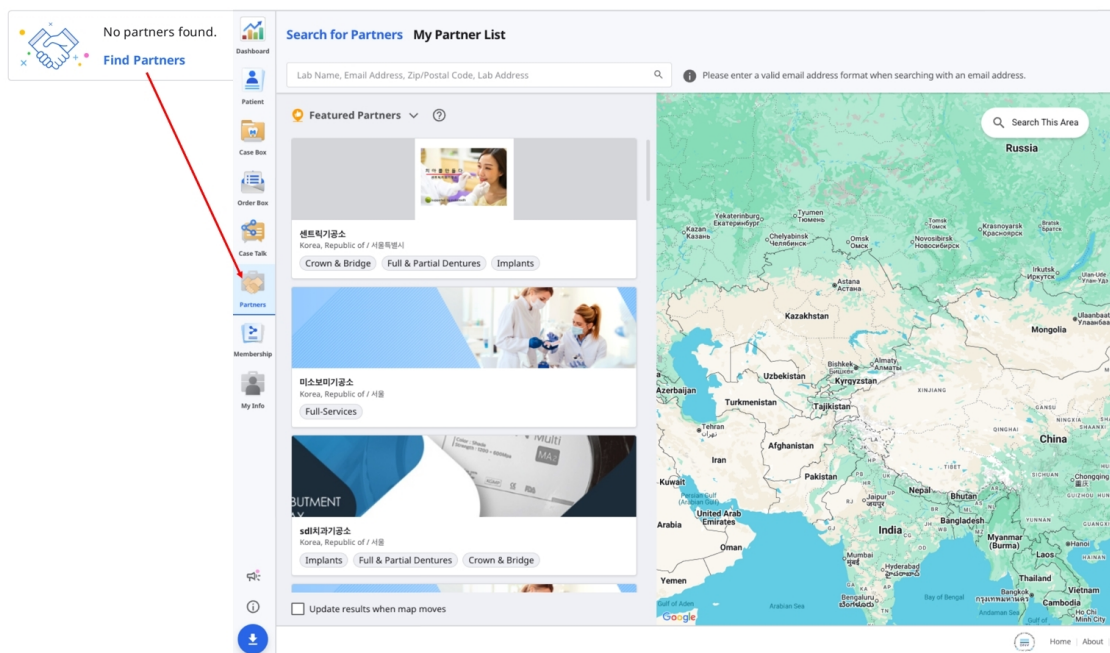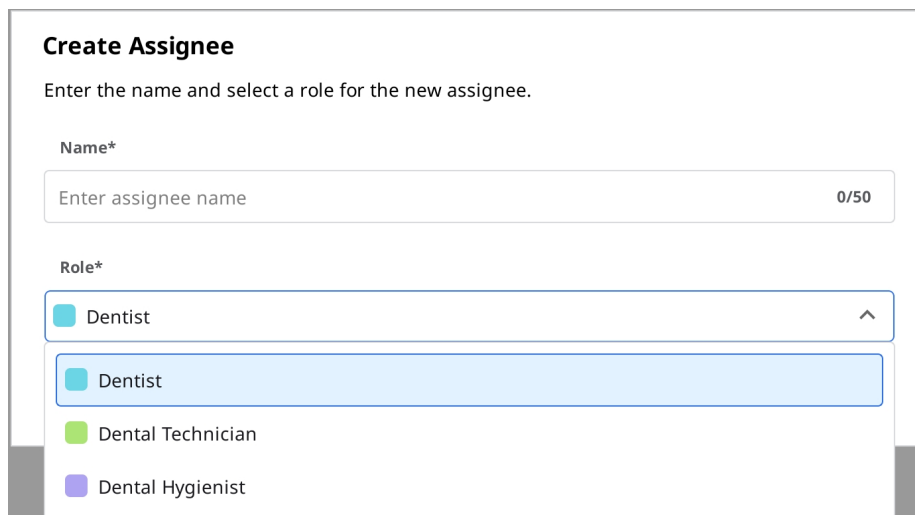
Create Assignee

Enter the name and select a role for the new assignee.

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MEDIT Link 3.2.0

Case Box wonhoonchoi-test1 / 1979-11-11 / Female / - + 2 ▾

☒ All |
 ☐ Form Only |
 ☐ Need Processing |
 ☐ Scan Completed |
 ☐ CAD |
 ☐ Ordered |
 ☐ Case Completed

☒ **Export** ▾ |
 Status ▾ |
 |
 Selected Cases: 4

	Status	Patient Name	Case Name	Form Information
☒	☐ Scan Completed	wonhoonchoi-test1	wonhoonchoi-test1's Case	-
☒	☐ Form Only	wonhoonchoi-test1	wonhoonchoi-test1's Case	-
☒	☐ Form Only	wonhoonchoi	originalcase	-
☒	☐ Scan Completed	wonhoonchoi-test1	wonhoonchoi-test1's Case	-

Dashboard
 Patient
Case Box
 Order Box
 Case Talk
 App Box

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Question

Exporting successfully completed.
Do you want to open the folder containing the exported file?

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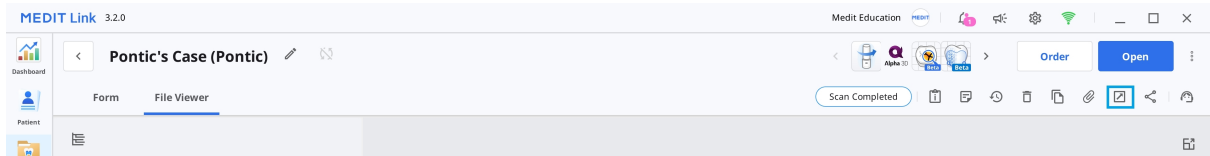
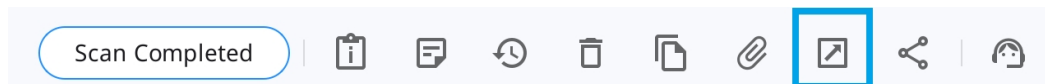
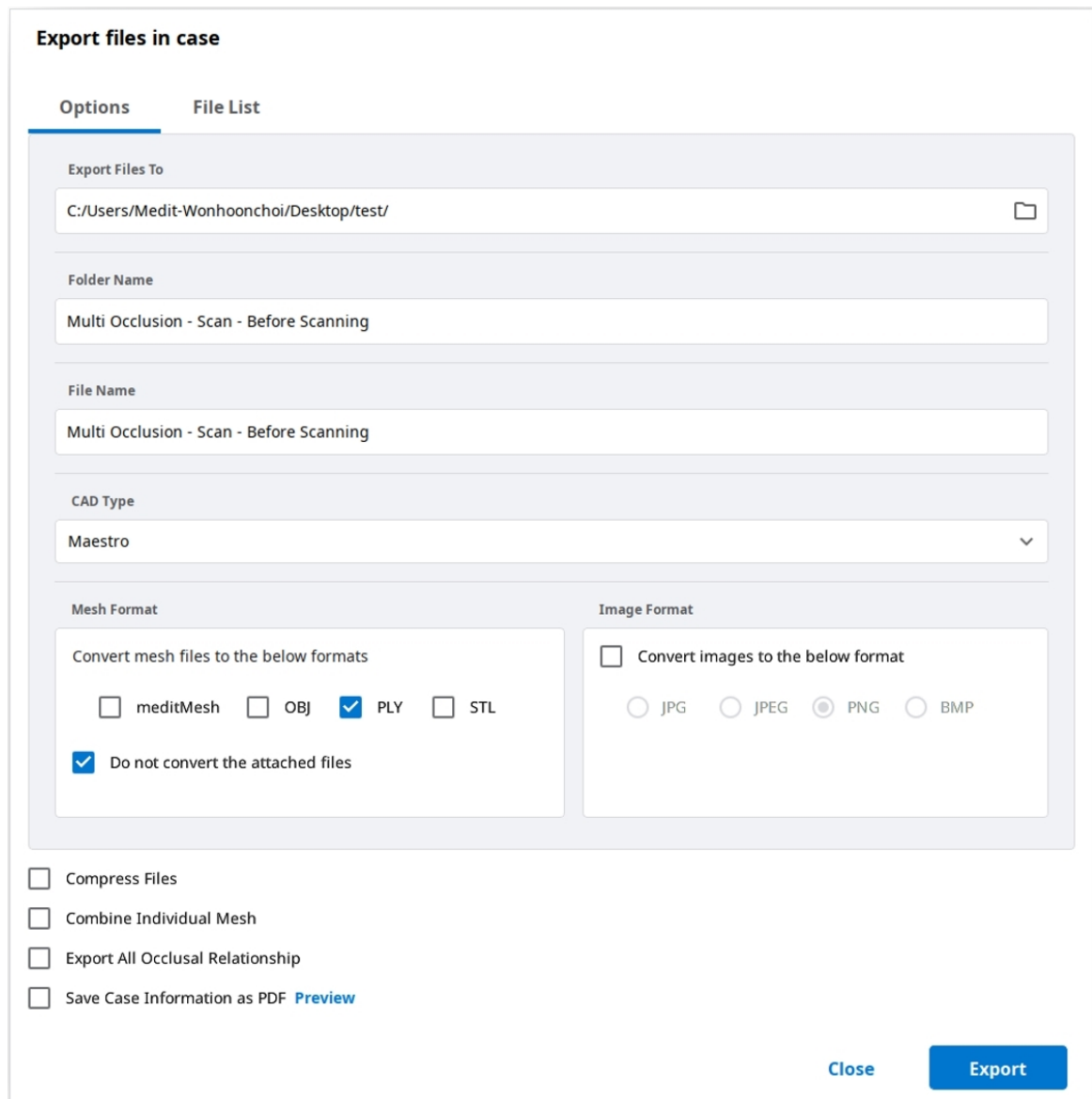
Options **File List**

Export Files To
C:/Users/Medit-Wonhoonchoi/Desktop/test/

Folder Name
Multi Occlusion - Scan - Before Scanning

File Name
Multi Occlusion - Scan - Before Scanning

CAD Type
Maestro


Mesh Format
Convert mesh files to the below formats
☐ meditMesh ☐ OBJ ☒ PLY ☐ STL
☒ Do not convert the attached files

Image Format
☐ Convert images to the below format
☐ JPG ☐ JPEG ☒ PNG ☐ BMP


☐ Compress Files
☐ Combine Individual Mesh
☐ Export All Occlusal Relationship
☐ Save Case Information as PDF [Preview](#)

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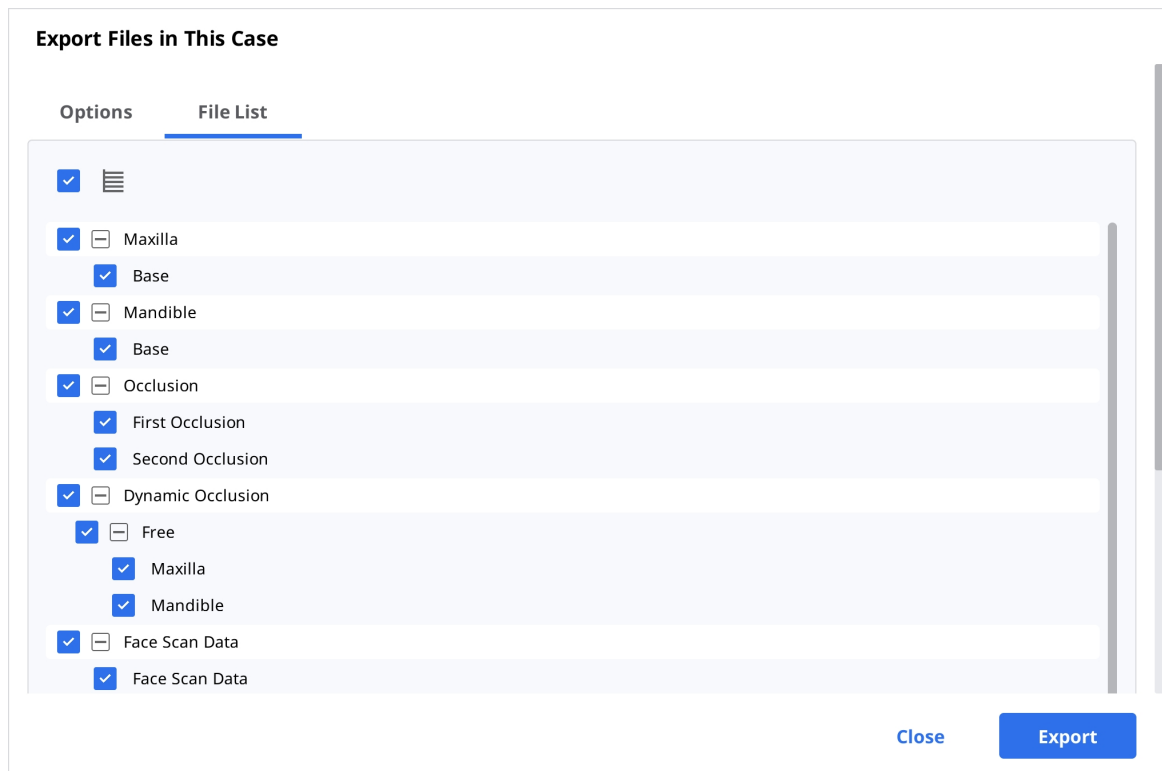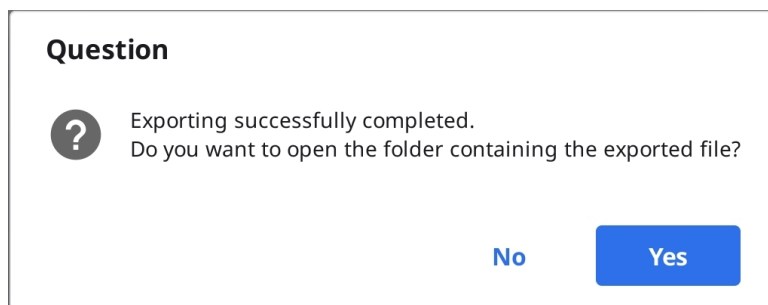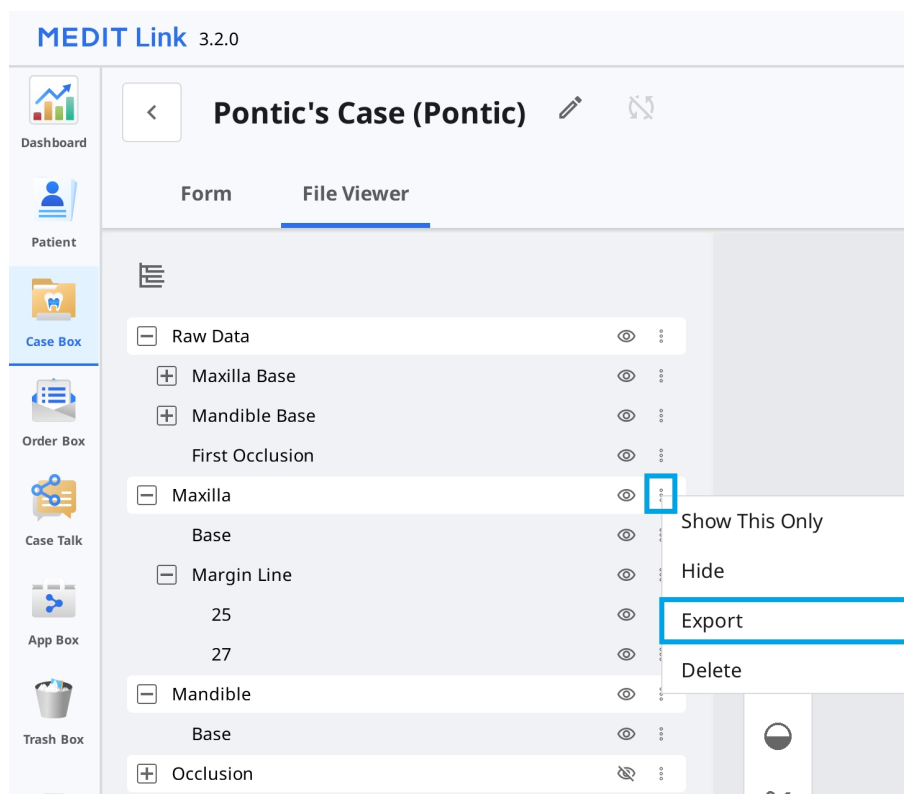
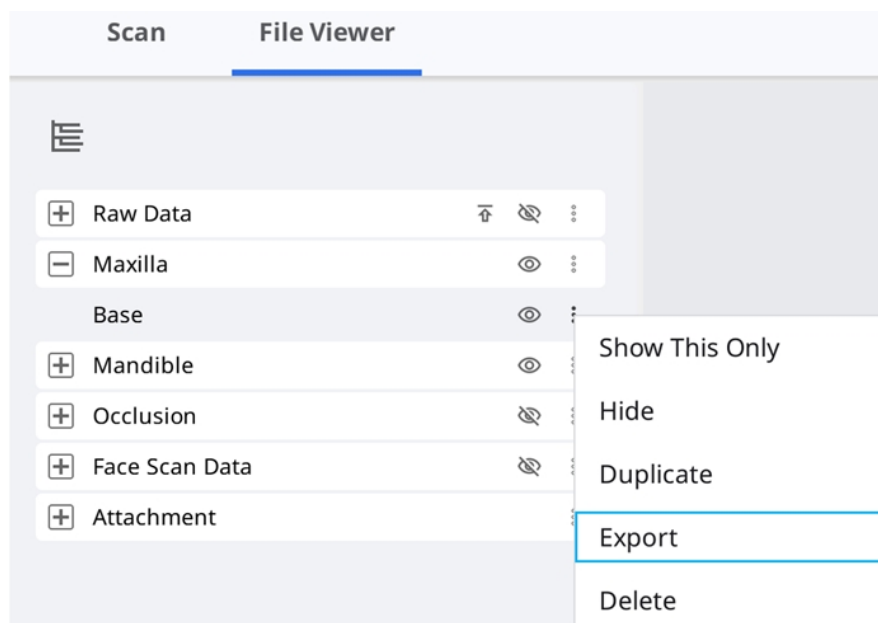
C:/Users/Medit-Wonhoonchoi/Desktop/test/

Folder Name

Multi Occlusion - Scan - Before Scanning

File Name

Multi Occlusion - Scan - Before Scanning

CAD Type

Maestro

Mesh Format

Convert mesh files to the below formats

☐ meditMesh

☐ OBJ

☒ PLY

☐ STL

☒ Do not convert the attached files

Image Format

☐ Convert images to the below format

☐ JPG

☐ JPEG

☒ PNG

☐ BMP

☐ Compress Files

☐ Combine Individual Mesh

☐ Export All Occlusal Relationship

☐ Save Case Information as PDF [Preview](#)

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Options File List

☐ Maxilla
 ☐ Base
 ☒ Mandible
 ☒ Base

☐ Occlusion
 ☐ First Occlusion
 ☐ Second Occlusion

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Question

? Exporting successfully completed.
Do you want to open the folder containing the exported file?

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Case Box Patient

All / 2018-01-01 ~ 2026-09-23 All / Contains

☒ All
 ☐ Form Only
 ☐ Need Processing
 ☐ Scan Completed
 ☐ CAD
 ☐ Ordered
 ☐ Case Completed

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Share through Case Talk

Share with Anyone

Create a shared link and share it with anyone.

Expiration Date

2022-04-27

You can extend the period up to 90 days from today.

PIN Code

☐ PIN Code

Enter PIN Code

Set the PIN code using six digits for privacy protection.

☒ Share the patient's name

☒ Grant Download Permission Premium Feature

If you check this option, all users who have access to the shared link can download the files.

The premium feature is provided to all users for a limited time.
Feature availability is subject to change depending on the membership plan.

☒ Allow comments

If you check this option, you will be able to communicate through the shared link.

Comment

Please enter your comment.

0/1024


By pressing the "Get Link" button, you agree that you understand personal information can be collected and provided to a third party (or shared in any public domain) ONLY with the consent of the person whose images on which you have worked. You also agree the consent of that person was truthfully obtained, and Medit is not responsible for any problems which may arise there after.

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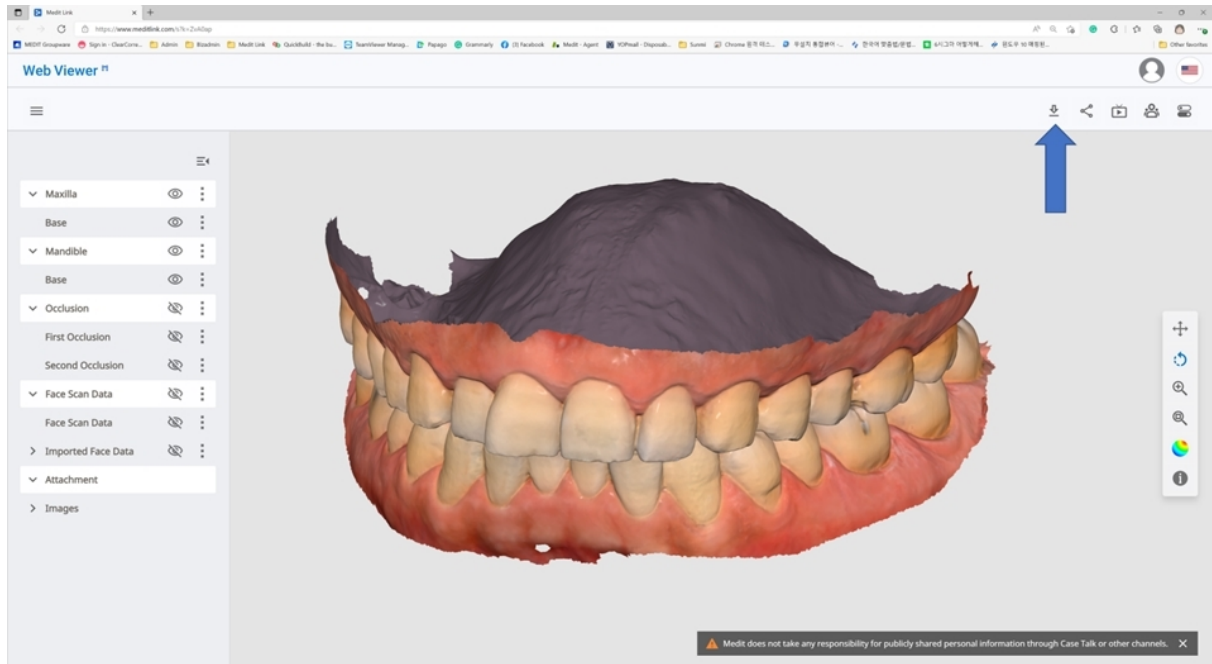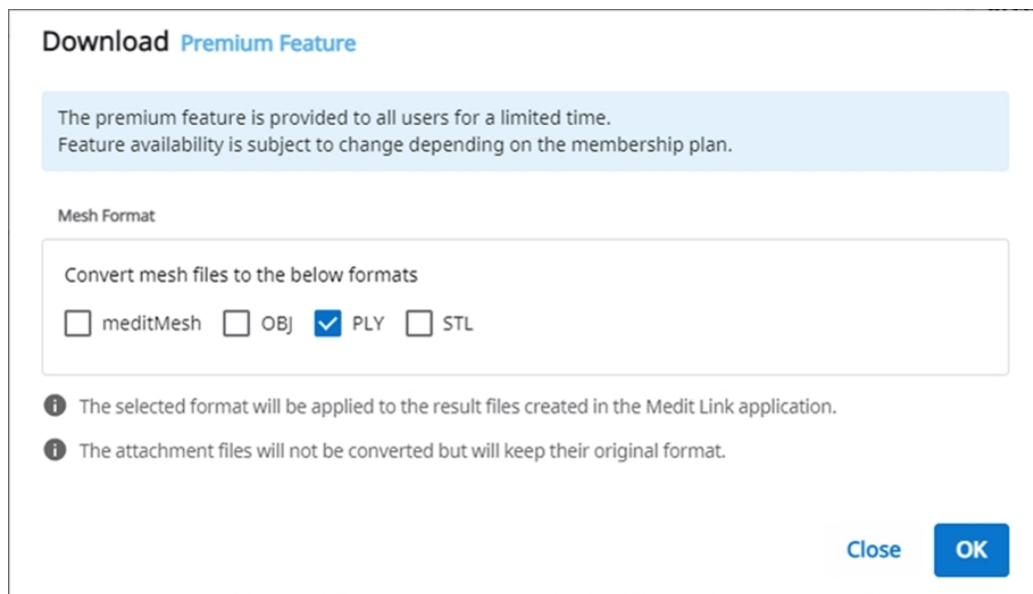
Uses Lab Form

San Diego, US

[92131] Residence Inn San Diego Rancho Bernardo

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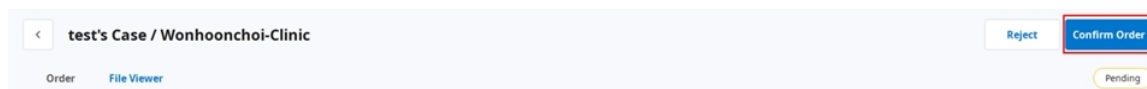
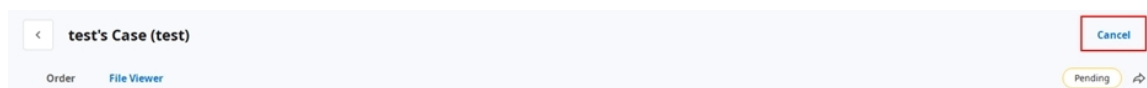
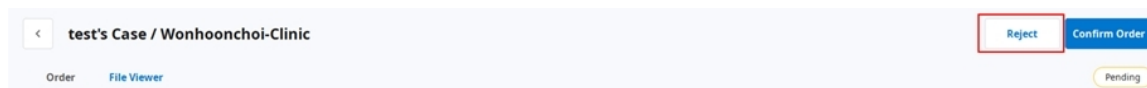
The screenshot shows a web interface for managing a case. The title is "test1's Case (test1)". On the left is a sidebar with icons for Dashboard, Patient, Case Box, Order Box, Case Talk, and Trash Box. The main area has two tabs: "Order" and "File Viewer". The "Order" tab is active, showing a form with the following fields:

- Order ID: 6004858
- Assignee: (dropdown menu)
- Patient Name: test1
- Partner Name: Elysee Dental Test
- Scanned Date: (empty)
- Order Date: 2025-07-23 PM 5:17
- Requested Delivery Date: 2025-07-30 PM 05:00
- Expected Delivery Date: 2025-07-30 PM 05:00
- Memo: (empty text area)
- Tags: test

On the right, there is a "Shipment" section with a table titled "Order List". The table has columns: No., Type, Method, Material, and Shad.

No.	Type	Method	Material	Shad
15	Crown	Anatomic	Zirconia	A4

At the top right of the main area, there is a "Cancel Order" button and a "Pending" status indicator.

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Shipment Registration

Select a delivery method and enter the required information to register.
The information you enter is saved in the Shipment tab and can be checked at any time.

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Case Name
mandibular movement - Clone - Clone

Patient Name
Test Occlusion

☒ Share patient information

Requested Delivery Date
2023-04-13

Delivery Time
오전 8:15 AM 8:00 PM 8:00

Memo

Tags
You can add the entered tag with Enter or Spacebar.

☒ Include all occlusion relationships in the order
☒ Include Tags In The Order

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Tags
tag

Order List Shipment

No.	Type	Method	Material	Shad
25	Crown	Anatomic	Multi-Layered Zirconia	A3.5
26	Crown	Anatomic	Multi-Layered Zirconia	A3.5

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