MEDIT Link

ME-UG-702H
Revision 13 (29.09.2026)
SW version 3.5.0

Table of contents

Medit Link

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Administrator

Dentist, Manager

Have the authority to create a new dental clinic account and have full access to Medit Link.

- ✓ Scan/Order/Manage orders
- ✓ Manage members

Register

Member

Dental Hygienist, Dental Staff

Belong to the dental clinic their admin created and have limited access to Medit Link after the admin's approval.

- ✓ Scan/Order/Manage orders

Register

Lab

Administrator

Laboratory Director, Manager

Have the authority to create a new dental lab account and have full access to Medit Link.

- ✓ Scan/Manage ordered cases
- ✓ Manage members

Register

Member

Dental Technician, Laboratory Staff

Belong to the dental lab their admin created and have limited access to Medit Link after the admin's approval.

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Email Address

Check

 Please ensure your email address is active. Medit Link will send a verification email.

 Registering with a disposable and/or fake e-mail address can make you vulnerable to security issues and make it difficult to receive information from Medit Link.

Password

 Password must be 8 to 16 characters long, including a combination of at least three letters, numbers, and special characters. (Passwords are accepted in English language only.)

Confirm Password

 Please type in the same password as above.

Select Country/Region

☐ I agree to all of the Terms and Conditions.

☐ Personal Information Collection Agreement (Mandatory)

View

☐ Terms & Conditions (Mandatory)

View

☐ Agree to the Marketing and Promotion Agreement (optional)

View

 Check this option to receive marketing emails, which include product updates, software news, and promotional events.

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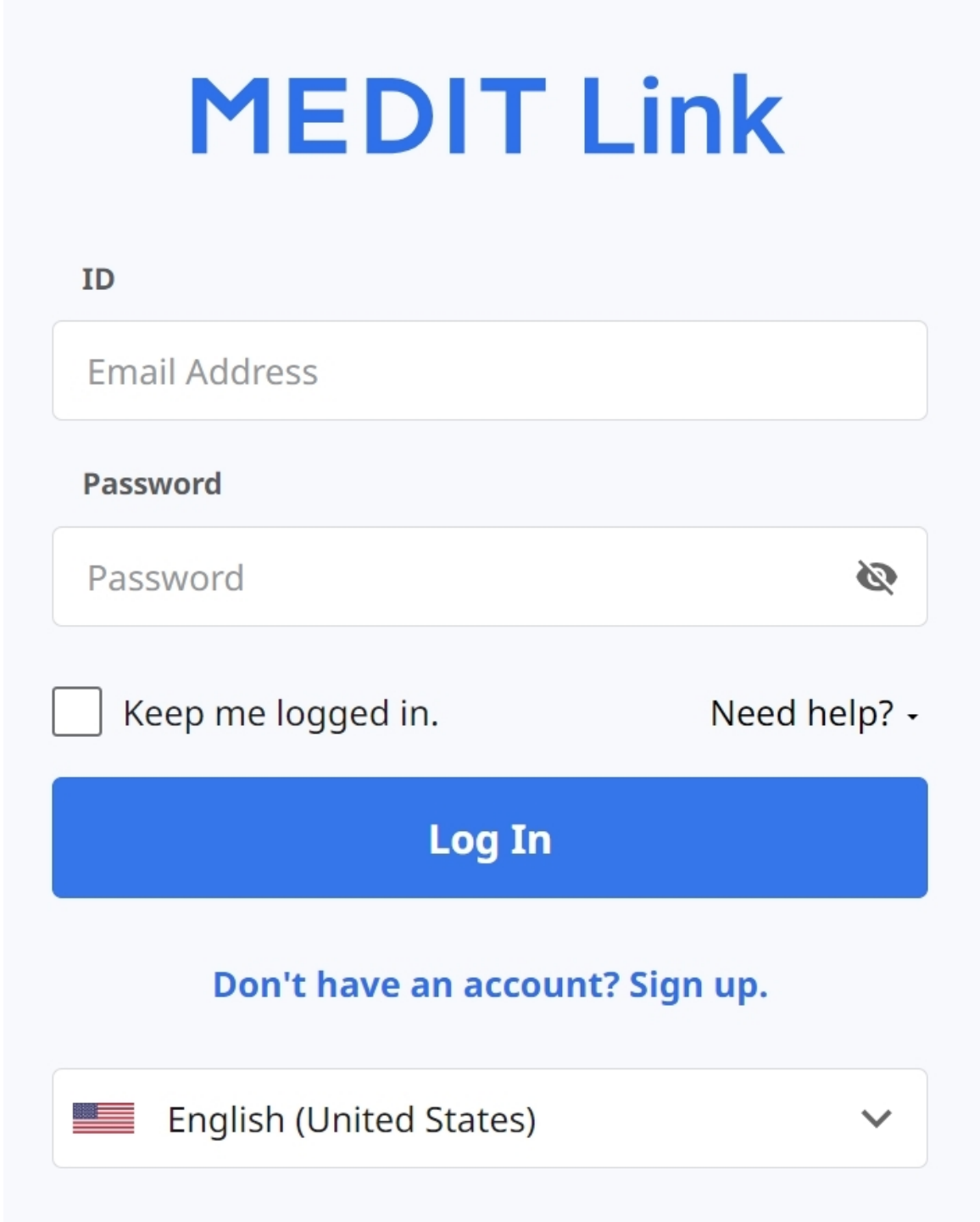
The image shows the login page for MEDIT Link. It features a light blue background with the title 'MEDIT Link' in large blue letters. Below the title, there are two input fields: 'ID' (Email Address) and 'Password'. The password field has a toggle icon for visibility. Below these fields, there is a checkbox for 'Keep me logged in.' and a link for 'Need help?'. A large blue 'Log In' button is centered below the checkbox. At the bottom, there is a link for 'Don't have an account? Sign up.' and a language selection dropdown menu showing 'English (United States)' with a US flag icon and a downward arrow.

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
Welcome, stage-test-wh!

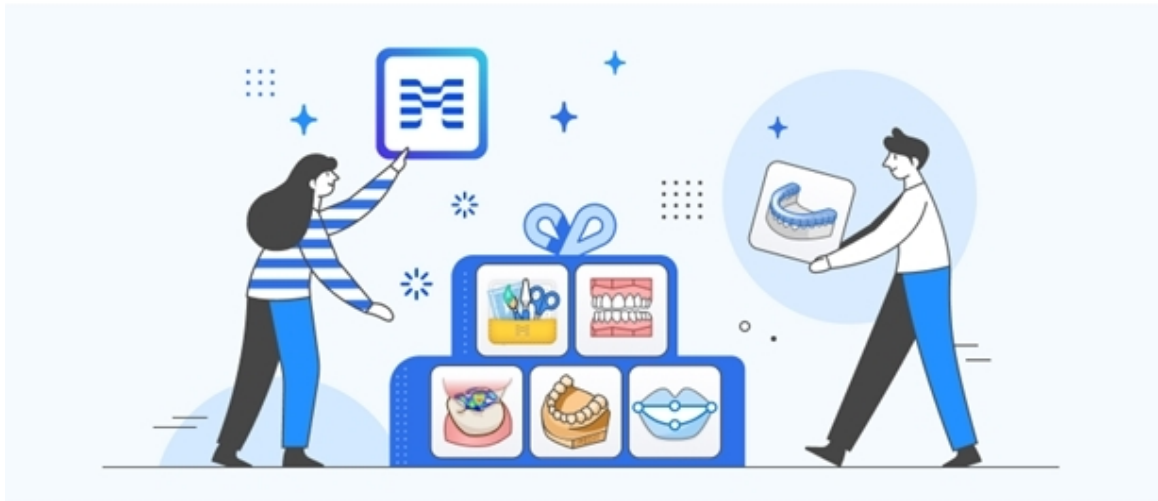Utilize Various Apps

Download Medit Link to your desktop to utilize scanning programs and various apps.
All data generated is synced to the cloud and is also available on the web.

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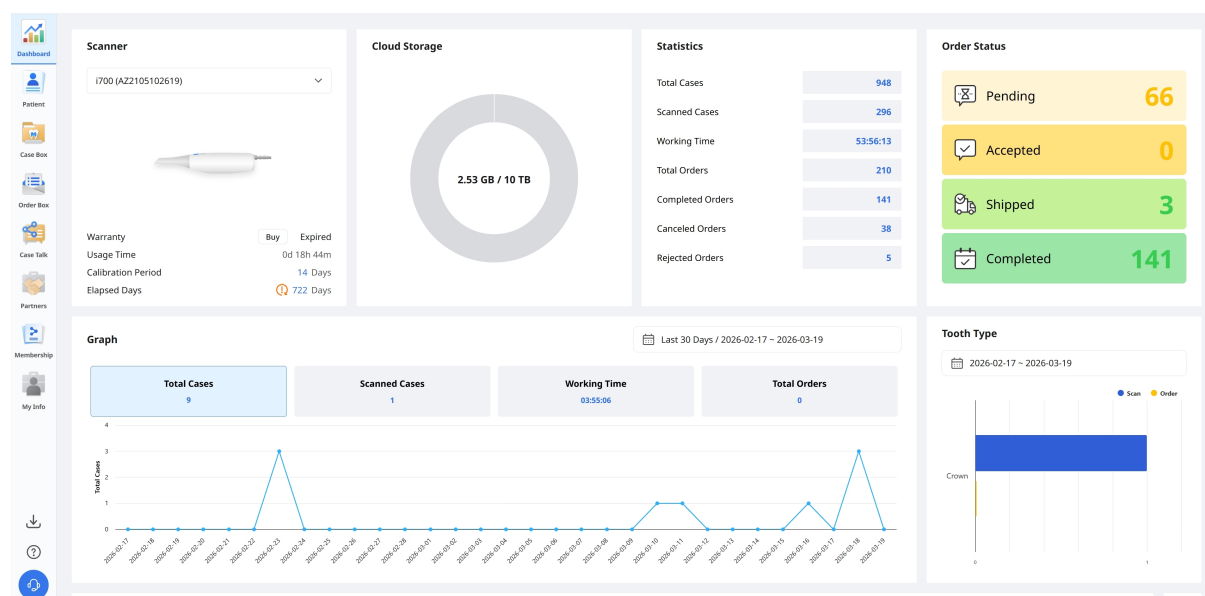
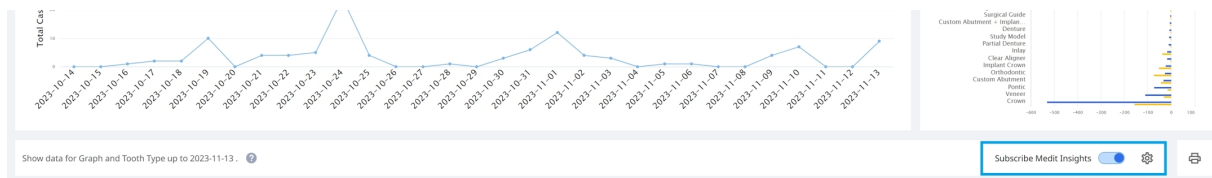
☒ Medit Link account


☐ Another email address

Enter another email address to receive the report email.

The report will be sent on the 2nd day of every month.

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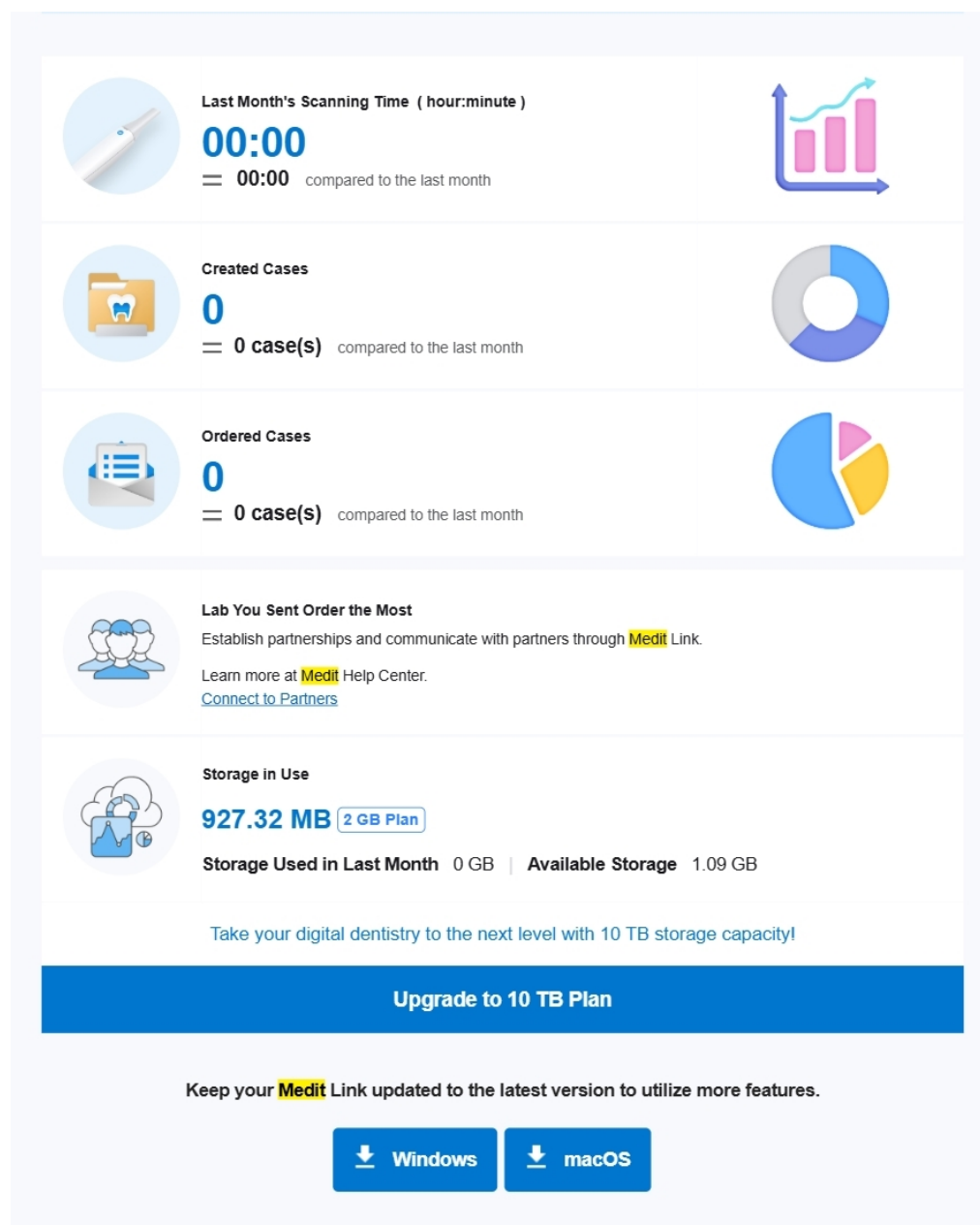This is an automatically generated email. Please do not reply to this email.
Please contact us at insights@medit.com for more information.
[Unsubscribe](#) If you no longer wish to receive this email, click "Unsubscribe." You can change the email subscription settings in [Medit Link Dashboard](#).
Copyright © Medit Corp. All rights reserved.

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
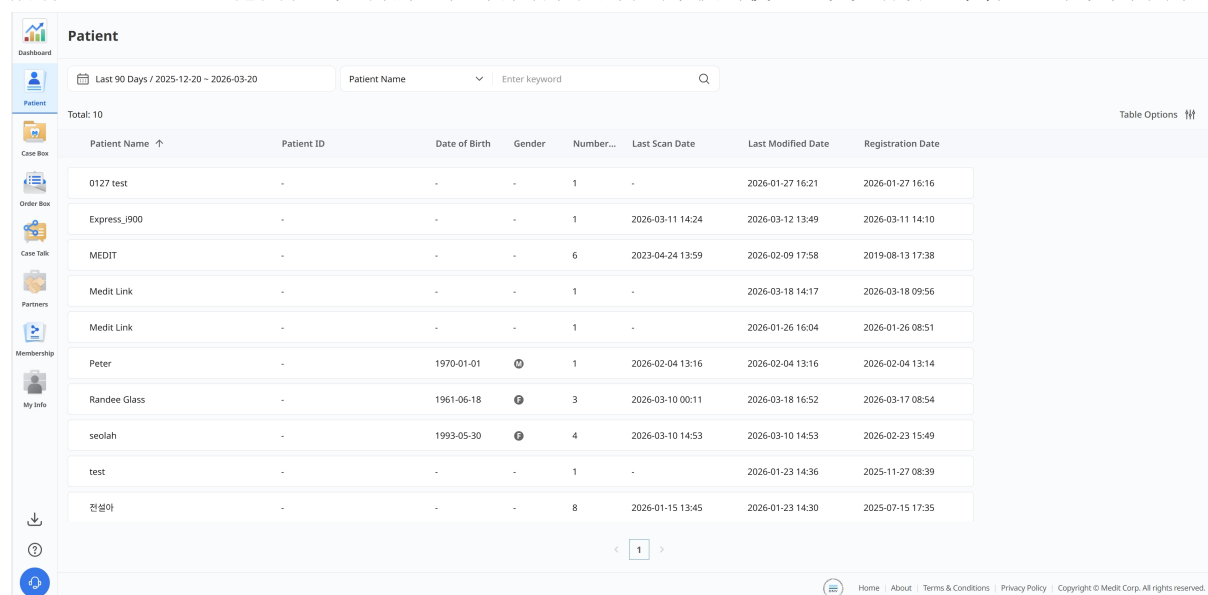
< 1 >

Home | About | Terms & Conditions | Privacy Policy | Copyright © Medit Corp. All rights reserved.

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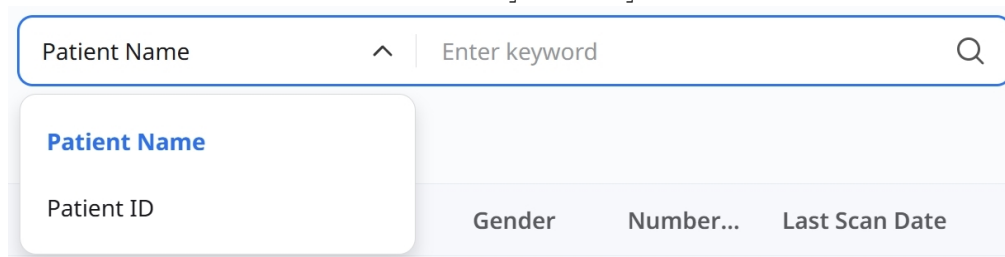The screenshot shows the Medit Link Case Box interface. At the top, there's a header with the Medit Link logo, a user profile (seolahjeon), and a language dropdown (English). Below the header, the 'Case Box' section is active, showing a filter for 'Patient'. A search bar is present with a date range 'Last 90 Days / 2025-12-20 ~ 2026-03-20' and a 'Case Name' dropdown. A filter bar includes 'All' (selected), 'Form Only', 'Need Processing', 'Scan Completed', 'CAD', 'Ordered', and 'Case Completed'. The main area displays a table of 18 cases. The table has columns for Status, Case Name, Patient Name, Form Information, Last Modified, Scanned Date, Order Date, Lab Name, and Assignee. Each row includes a status badge, a case name, a patient name, and a 'Share' button. The footer contains a navigation bar with links to Home, About, Terms & Conditions, Privacy Policy, and Copyright © Medit Corp. All rights reserved.

Status	Case Name	Patient Name	Form Information	Last Modified...	Scanned Date	Order Date	Lab Name	Assignee
Scan Completed	Ortho Scan - Clone - Clo...	Randee Glass		2026-03-18 16:52	2026-03-10 00:11	-	-	Share
Scan Completed	Ortho Scan - Clone	Randee Glass		2026-03-18 14:19	2026-03-10 00:11	-	-	Share
Form Only	Medit Link의 케이스	Medit Link		2026-03-18 14:17	-	-	-	Share
Scan Completed	Ortho Scan	Randee Glass		2026-03-17 08:57	2026-03-10 00:11	-	-	Share
Need Processing	Express_i900의 케이스	Express_i900	23	2026-03-12 13:49	2026-03-11 14:24	-	-	Share
Form Only	Express Scan	seolah		2026-03-10 14:53	2026-03-10 14:53	-	-	Share
Form Only	Express Scan	seolah		2026-02-24 08:33	2026-02-24 08:33	-	-	Share
Need Processing	Express Scan	seolah		2026-02-23 16:38	2026-02-23 16:37	-	-	Share
Scan Completed	Express Scan	seolah	Maxilla	2026-02-23 15:50	2026-02-23 15:50	-	-	Share

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The screenshot shows the Medit Link Case Box interface. The 'Case Box' section is active, and the 'Patient' filter is selected. A red box highlights the 'Patient' dropdown menu. Below the filter bar, there are buttons for 'All' (selected), 'Form Only', 'Need Processing', 'Scan Completed', 'CAD', 'Ordered', and 'Case Completed'.

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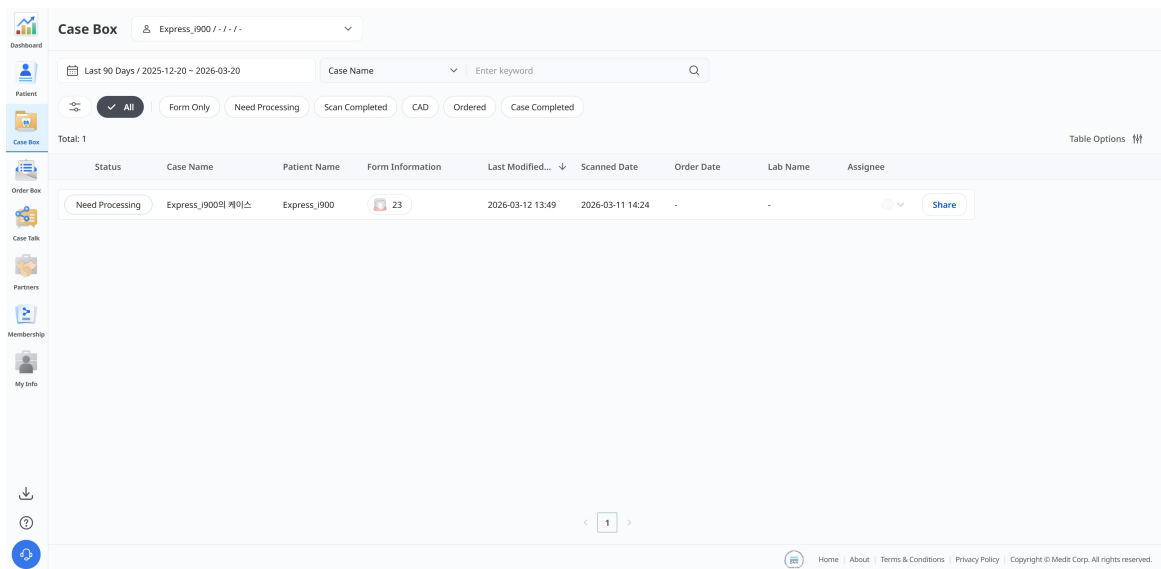
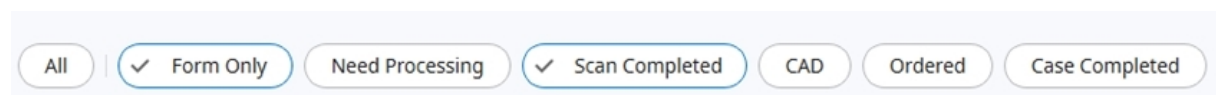
Case Box(診所)	所有 / Form Only / Need Processing / Scan Completed / CAD / Ordered / Case Completed
Work Box(技工所)	所有 / Form Only / Need Processing / Scan Completed / CAD / 完成

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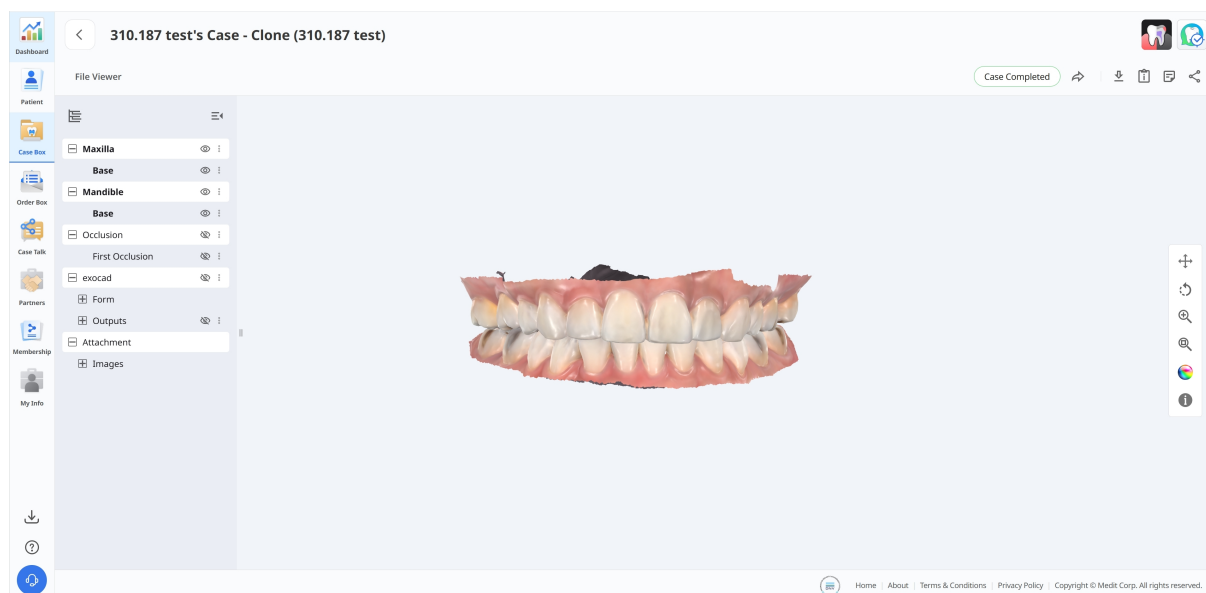
The screenshot displays the 'Order Box' interface in the Medit Link application. It features a sidebar with navigation icons for Dashboard, Patient, Case Box, Order Box, Case Task, Partners, Membership, and My Info. The main content area is titled 'Order Box' and includes a filter section with a date range (Last 90 Days / 2026-01-08 - 2026-04-08), a search bar for Case Name, and status filters (All, Pending, Canceled, Rejected, Accepted, Shipped, Case Completed). A table below shows one order with details: Status (Pending), Case Name (Express Scan), Patient Name (Seol register), Order Date (2026-04-02 21:38), Requested Deliv... (2026-04-09 09:00), Last Modified... (2026-04-02 21:38), Lab Name (MEDIT_LAB_STAGE), Order ID (6012242), and Assignee. The interface also includes a 'Total: 1' indicator and a 'Table Options' button.

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Invite Partner

Invite your desired partner to join Medit Link. When the lab accepts your invitation, a partnership is automatically formed.

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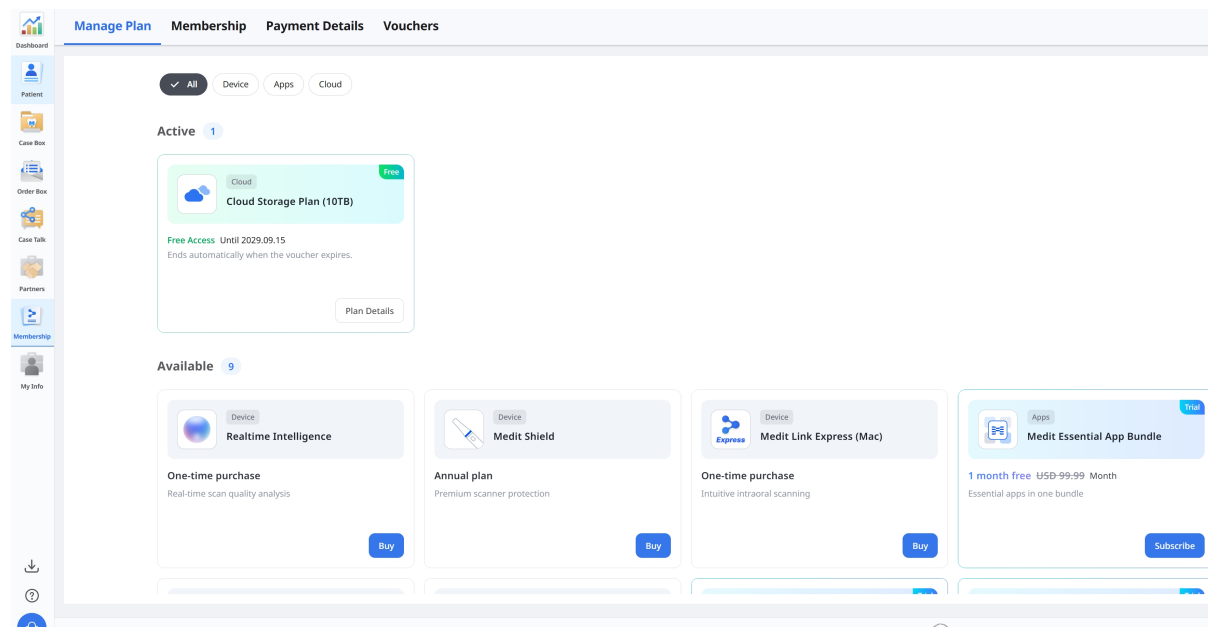
Manage Plan **Membership** **Payment Details** **Vouchers**

✓ All Device Apps Cloud

Active 1

Cloud

Free

 **Cloud Storage Plan (10TB)**

Free Access Until 2029.09.15
Ends automatically when the voucher expires.

Plan Details

Available 9

Device

 **Realtime Intelligence**

Device

 **Me**

Plan Details

×

Free Access
Your free access is active.

Product	Cloud Storage Plan (10TB)
Status	Free
Expiration Date	2029.09.15
Regular Price	USD 9.99/Month

 **Important**


- You can't subscribe to this product while the voucher is active.
- Free access ends automatically when the voucher expires. A paid subscription is required to continue.

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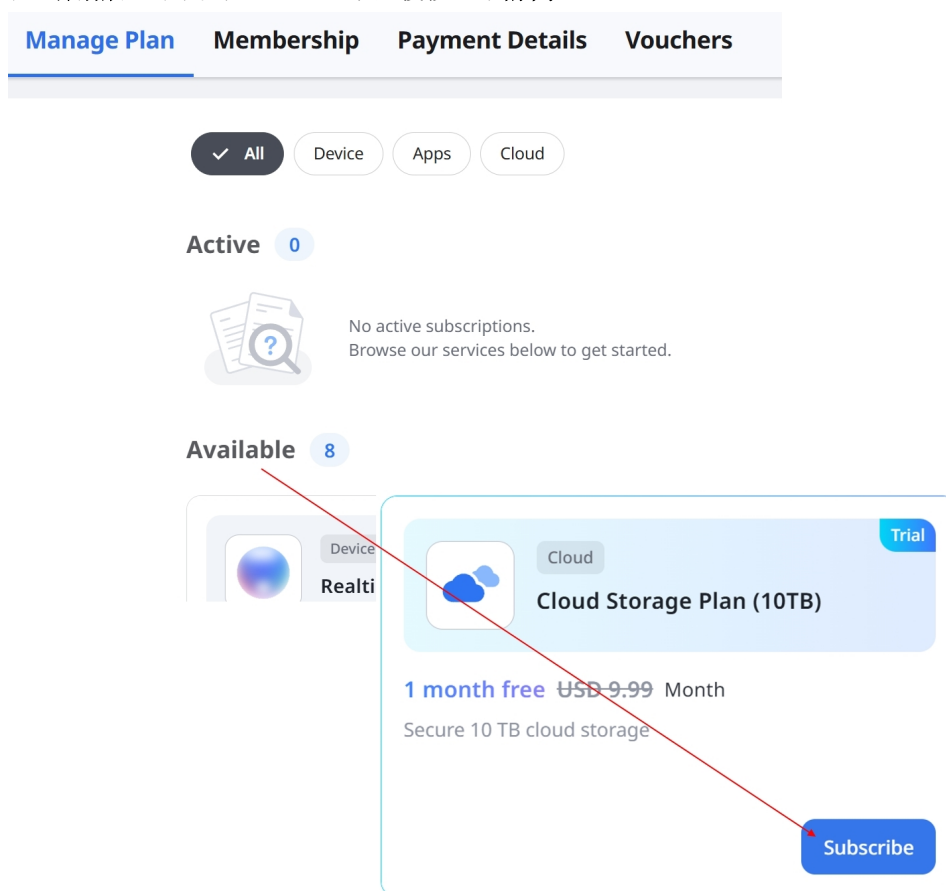1 month free

~~USD 9.99~~ / Month

- ✓ Includes 10 TB of cloud storage
- ✓ Securely store large cases and other data
- ✓ Streamline your Medit Link data management in the cloud
- ✓ Store design files without separate local backups

- Actual charges may vary depending on your bank or card issuer's fees and policies.
- For more information about Cloud Storage Plan (10TB), see [here](#).

Trial

Cloud Storage: First Month Free

1 month free >

☐ I agree to automatic renewal and recurring charges.

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The image shows a two-step subscription process. The first step is a payment summary screen with a blue header 'Subscribe to Medit Cloud Storage 10TB Plan'. It displays a price of \$0.00, noting that the rate will be \$9.99 per day after a coupon expires. A table lists the 'Medit Cloud Storage 10TB Plan for Clinics' at \$9.99, a '1 Month Free for New Subscribers' discount of -\$9.99, and a tax of \$0.00, resulting in a 'Total due today' of \$0.00. The second step is a 'Pay with card' form with fields for email (yhgam0829@gmail.com), card number (1234 1234 1234 1234), MM / YY, CVC, cardholder name, and country (South Korea). A red arrow points from the 'Subscribe' button in the payment form to the 'Confirm' button in a modal window.

Subscribe

Your Cloud Storage Plan (10TB) subscription is now active.
Enjoy full access to all features.

Amount of Payment	Free
Next Payment Date	2026.10.23
Last Payment Method	****

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The image shows a 'Manage Plan' modal window. It features a 'Promotion Applied' banner stating 'Regular pricing will apply after the promotion ends.' Below this, a table displays subscription details: 'Product' (Cloud Storage Plan (10TB)), 'Price' (Free USD 9.99/Month), 'Expiration Date' (2026.10.23), and 'Last Payment Method' (**** Edit). A red arrow points to the 'Unsubscribe' button, which is highlighted in red.


Manage Plan

Promotion Applied
Regular pricing will apply after the promotion ends.

Product	Cloud Storage Plan (10TB)
Price	Free USD 9.99/Month
Expiration Date	2026.10.23
Last Payment Method	**** Edit

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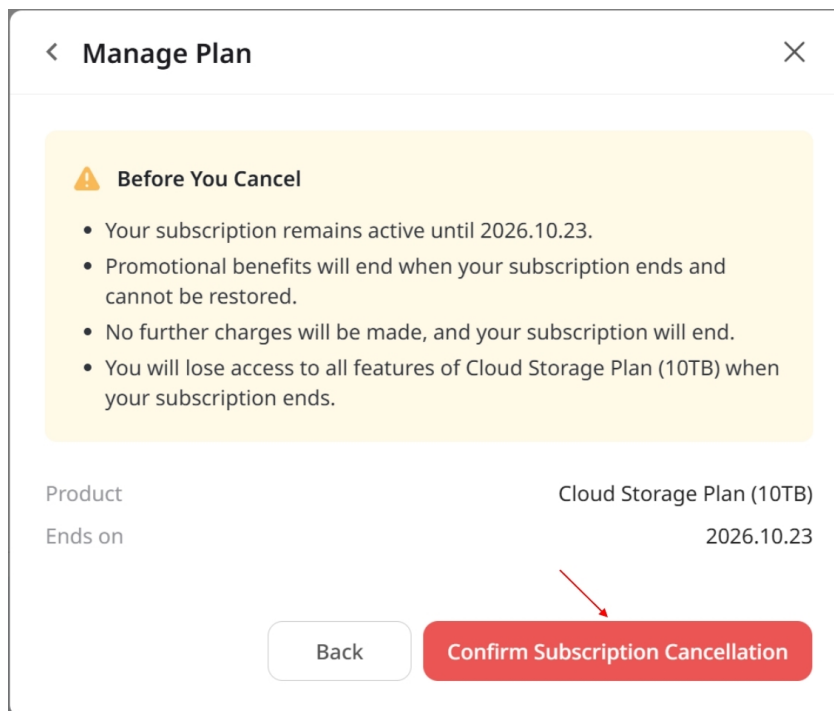The screenshot shows a 'Manage Plan' dialog box with a close button (X) in the top right corner. Inside, there is a yellow warning box titled 'Before You Cancel' with a warning icon. It contains four bullet points: 'Your subscription remains active until 2026.10.23.', 'Promotional benefits will end when your subscription ends and cannot be restored.', 'No further charges will be made, and your subscription will end.', and 'You will lose access to all features of Cloud Storage Plan (10TB) when your subscription ends.' Below this, the 'Product' is listed as 'Cloud Storage Plan (10TB)' and 'Ends on' is '2026.10.23'. At the bottom, there are two buttons: a white 'Back' button and a red 'Confirm Subscription Cancellation' button. A red arrow points to the red button.

< **Manage Plan** ×


⚠ Before You Cancel

- Your subscription remains active until 2026.10.23.
- Promotional benefits will end when your subscription ends and cannot be restored.
- No further charges will be made, and your subscription will end.
- You will lose access to all features of Cloud Storage Plan (10TB) when your subscription ends.

Product Cloud Storage Plan (10TB)
Ends on 2026.10.23

Back Confirm Subscription Cancellation

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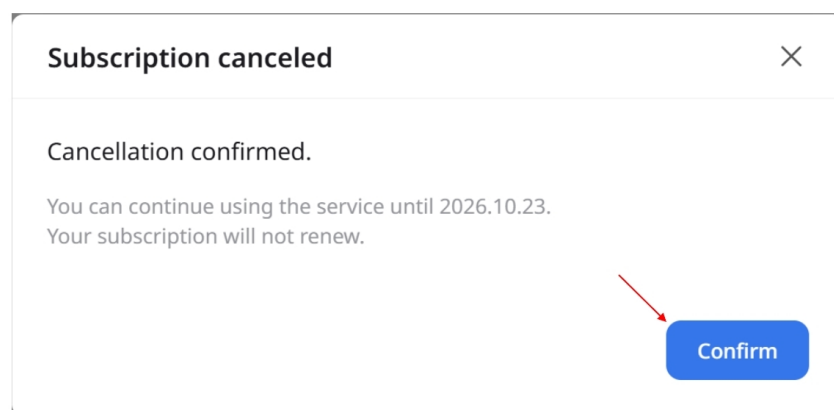The screenshot shows a 'Subscription canceled' dialog box with a close button (X) in the top right corner. It contains the text 'Cancellation confirmed.' followed by 'You can continue using the service until 2026.10.23.' and 'Your subscription will not renew.' At the bottom right, there is a blue 'Confirm' button. A red arrow points to the blue button.


Subscription canceled ×

Cancellation confirmed.

You can continue using the service until 2026.10.23.
Your subscription will not renew.

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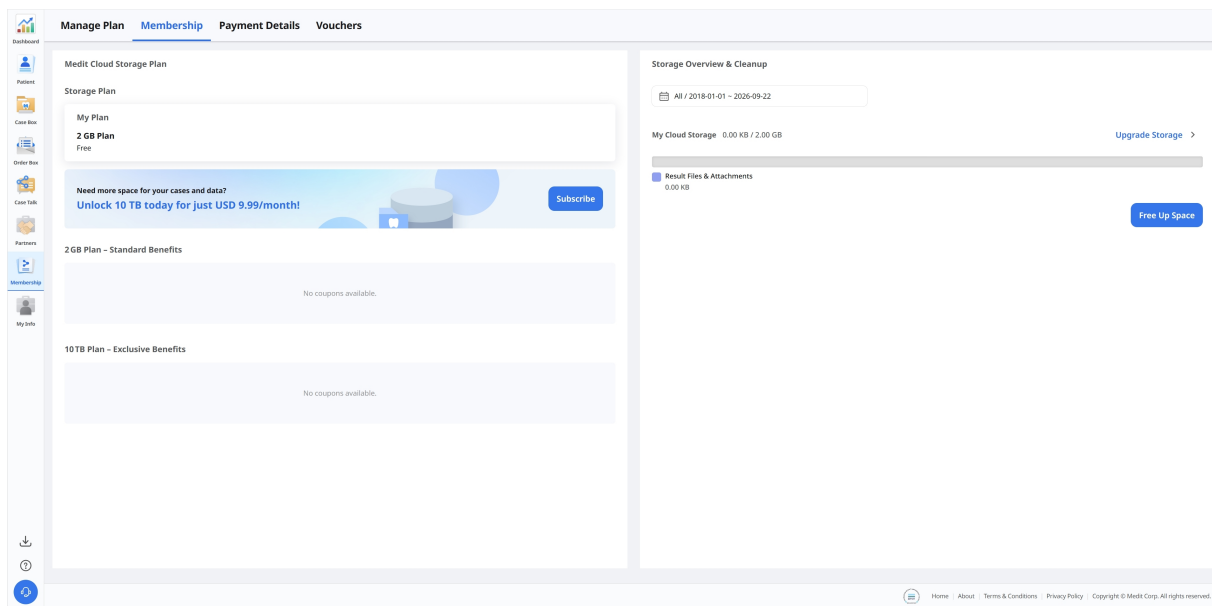
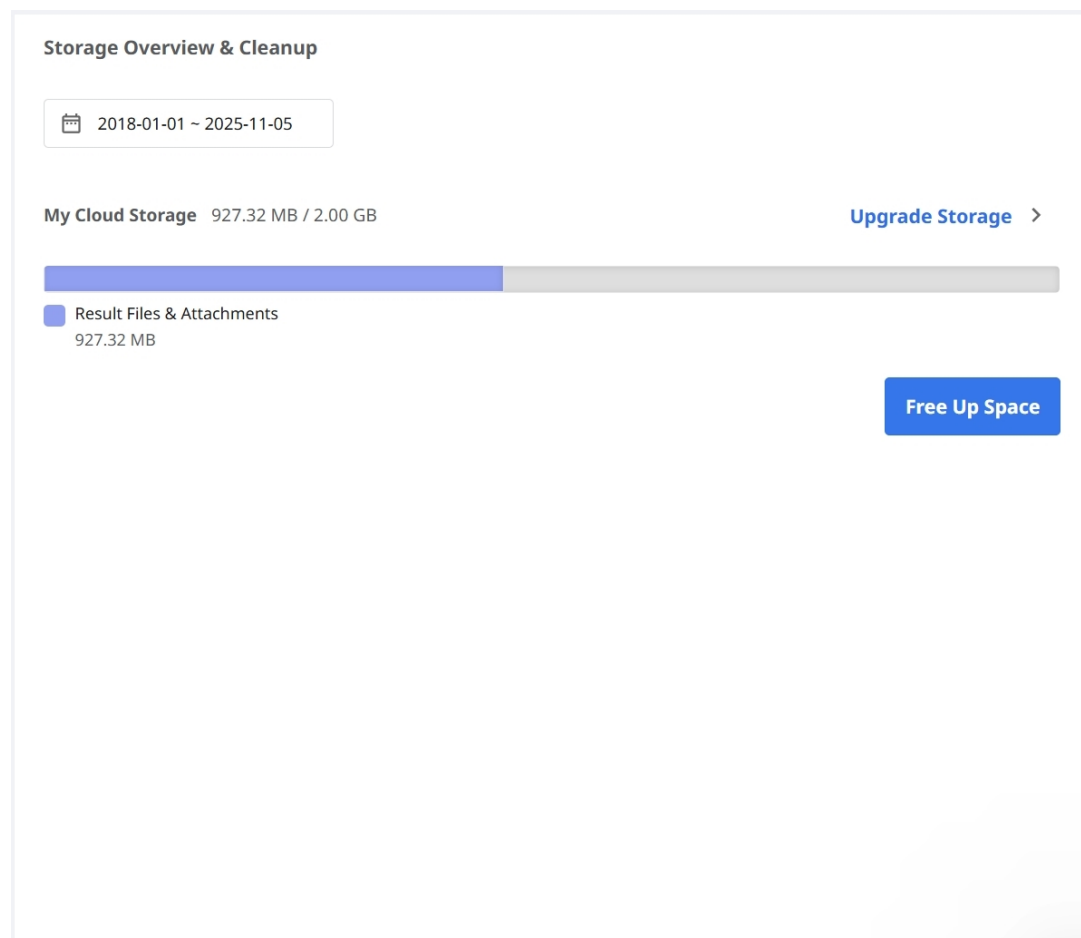
Payment Date	Product	Payment Type	Duration	Details	Amount of Payment	Status
2026-09-21	Cloud Storage Plan (10TB)	Subscription	2026-09-21 - 2026-10-21	Monthly payment	USD 0.00	Payment Completed
2026-09-21	Cloud Storage Plan (10TB)	Subscription	2026-09-21 - 2026-10-21	Monthly payment	USD 0.00	Payment Completed
2026-09-21	Cloud Storage Plan (10TB)	Subscription	2026-09-21 - 2026-10-21	Monthly payment	USD 0.00	Payment Completed
2026-09-21	Cloud Storage Plan (10TB)	Subscription	2026-09-21 - 2026-10-21	Monthly payment	USD 0.00	Payment Completed
2026-09-21	Cloud Storage Plan (10TB)	Subscription	2026-09-21 - 2026-10-21	Monthly payment	USD 0.00	Payment Completed

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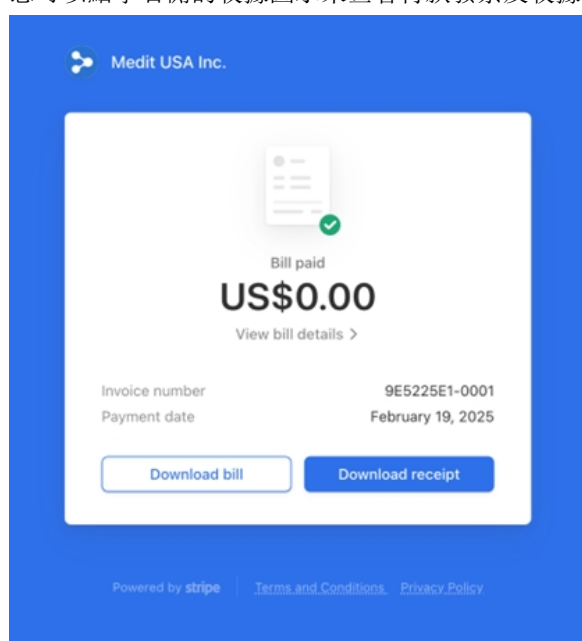ITEM

FEBRUARY 19 - FEBRUARY 20, 2025

Medit Cloud Storage 10TB Plan for Clinics

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Vouchers 2

Product	Promotion	Code Number
Cloud Storage Plan	Redeem SQA_Redeem_Cloud Storage Plan_36M	GFZ!

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Free Access30 Days

⚠ Important

- Once activated, the voucher cannot be canceled.
- You can't subscribe to this product while the voucher is active.

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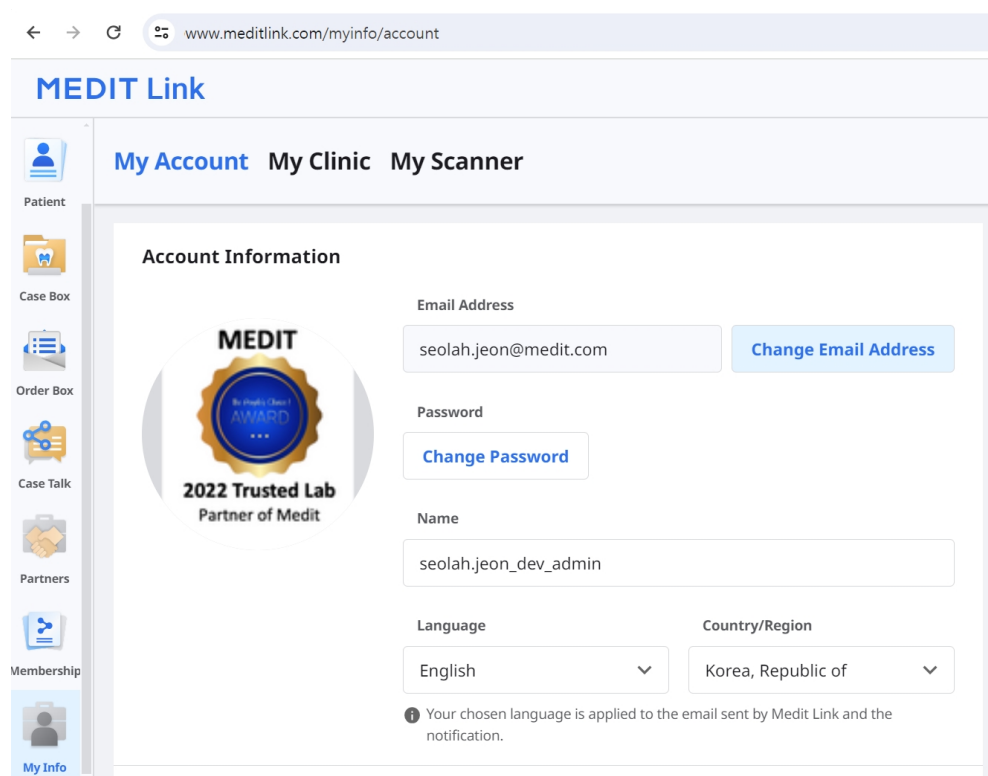
The screenshot shows the 'My Info' page of the MEDIT Link website. The browser address bar displays 'www.meditlink.com/myinfo/account'. The page has a light blue header with the 'MEDIT Link' logo. A sidebar on the left contains icons for 'Patient', 'Case Box', 'Order Box', 'Case Talk', 'Partners', 'Membership', and 'My Info' (which is highlighted). The main content area is titled 'Account Information' and contains the following fields:

- Email Address:** seolah.jeon@medit.com, with a 'Change Email Address' button.
- Password:** A 'Change Password' button.
- Name:** seolah.jeon_dev_admin
- Language:** English (dropdown menu).
- Country/Region:** Korea, Republic of (dropdown menu).

Below the dropdowns, a small information icon and text state: 'Your chosen language is applied to the email sent by Medit Link and the notification.'

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Edit

Cover and Logo

Looking for digitalized dental laboratories?

Find the most ordered Labs around your clinic with 'Featured Part

123

Services

Select

Specialties (choose up to 3)

Select

Select the services your business specializes in from the list.

About Us

e.g., Welcome to OOO Dental Clinic.

We are a licensed dental clinic situated in California specializing mainly in aesthetic restorations and implantology.

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The scanner information is automatically registered once your scanner is connected to a Medit scan program.

Home

About

Terms & Conditions

Privacy Policy

Copyright © Medit Corp. All rights reserved.

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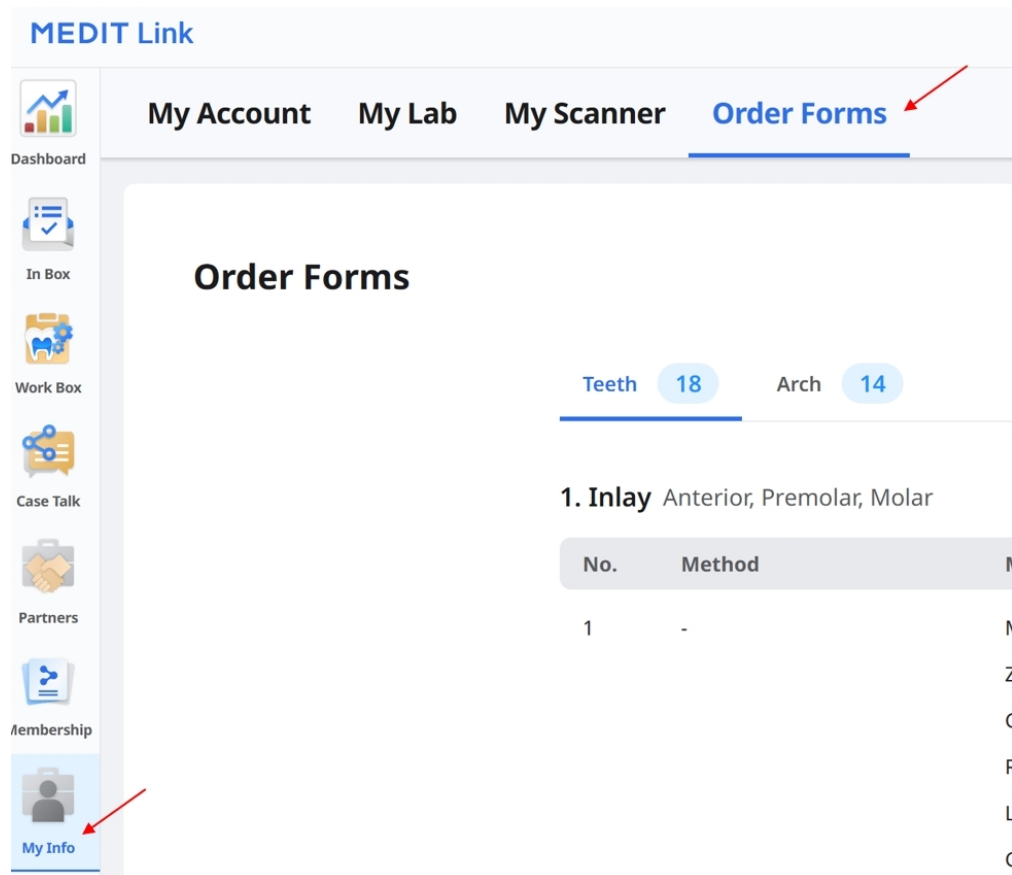
MEDIT Link

My Account My Lab My Scanner **Order Forms**

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Methods ¹

Search

Screw-Retained **Full-Contour** Offset Cement-Retained

SCRIP Temporary Cut-back Pontic Press Preform Screw-Retained Layered

Screw-Retained PFZ Screw-Retained PFG Screw-Retained PFM Cement-Retained Layered

[More](#)

☐ Do not use

Full-Contour

Materials ²

Search

Multilayer Zirconia Zirconia PMMA Gold Resin

LiSi Glass Ceramic Press Ceramic Hybrid Ceramic Composite Resin e.max

Titanium Metal Stainless Steel Co-Cr (Cobalt Chromium) Ni-Cr (Nickel-Chromium) Wax

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
The screenshot displays the 'Crown' configuration window. Under 'Allowed Teeth', 'Anterior', 'Premolar', and 'Molar' are all checked. In the 'Methods' section, 'Full-Contour' is the selected method among various options like 'Screw-Retained', 'Offset', and 'Cement-Retained'. Below the methods, a 'Do not use' checkbox is present. A dropdown menu for 'Full-Contour' is open, showing a list of materials. 'Zirconia HT' is highlighted with a blue selection bar and a red arrow points to it. Other materials visible include 'Multilayer Zirconia', 'Zirconia', 'PMMA', 'Gold', 'Resin', 'LiSi', 'Glass Ceramic', 'Press Ceramic', 'Hybrid Ceramic', 'Composite Resin', 'e.max', 'Titanium', 'Metal', 'Stainless Steel', 'Co-Cr (Cobalt Chromium)', 'Ni-Cr (Nickel-Chromium)', 'Wax', 'PEEK', and '3D Print'. A '+ Add Material' link is also visible. At the bottom of the materials list, there is another 'Do not use' checkbox and a 'Hide' link.

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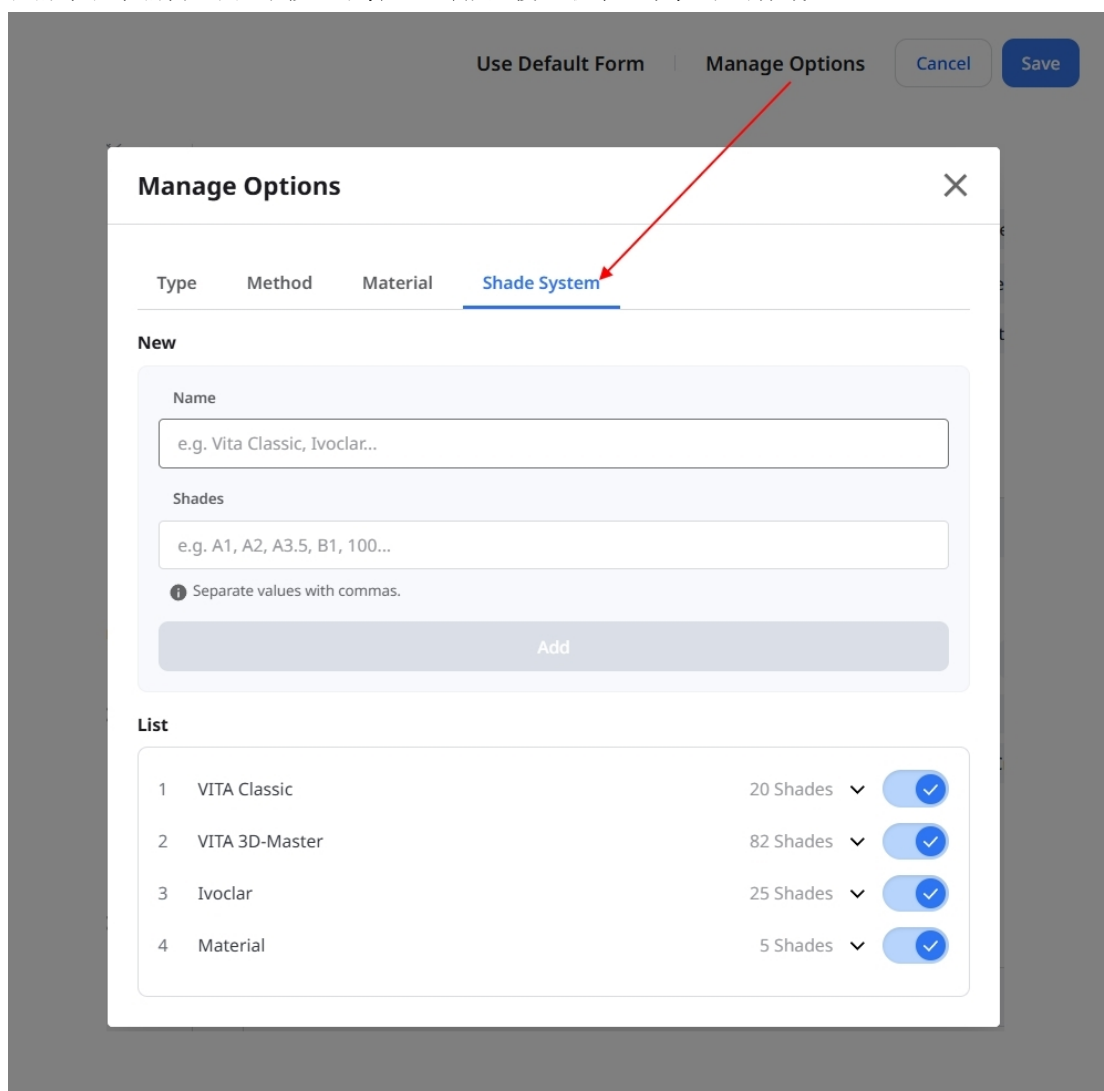
The screenshot shows the 'Manage Options' dialog box. At the top, there are buttons for 'Use Default Form', 'Manage Options' (highlighted with a red arrow), 'Cancel', and 'Save'. The dialog has four tabs: 'Type', 'Method', 'Material', and 'Shade System' (which is selected and underlined). Below the tabs, there are two sections: 'New' and 'List'. The 'New' section contains input fields for 'Name' (with placeholder 'e.g. Vita Classic, Ivoclar...') and 'Shades' (with placeholder 'e.g. A1, A2, A3.5, B1, 100...'). Below these fields is a note 'Separate values with commas.' and an 'Add' button. The 'List' section shows a table of existing shade systems:

ID	Name	Shades	Active
1	VITA Classic	20 Shades	☒
2	VITA 3D-Master	82 Shades	☒
3	Ivoclar	25 Shades	☒
4	Material	5 Shades	☒

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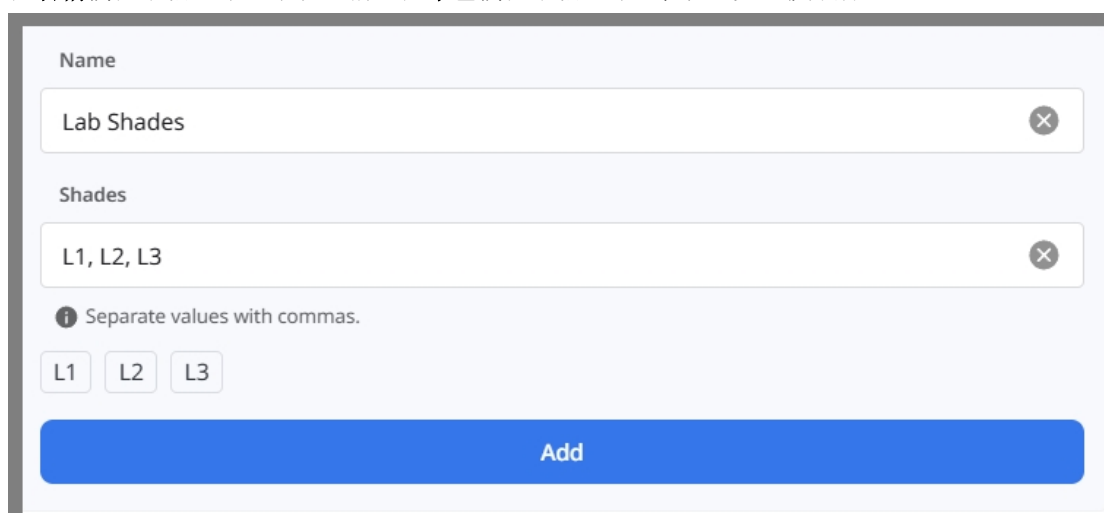This is a close-up of the 'New' section of the 'Manage Options' dialog box. It shows the 'Name' input field containing 'Lab Shades' and the 'Shades' input field containing 'L1, L2, L3'. Below the 'Shades' field is a note 'Separate values with commas.' and three buttons labeled 'L1', 'L2', and 'L3'. At the bottom is a large blue 'Add' button.

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The screenshot shows the 'Crown' configuration window. Under 'Allowed Teeth', 'Anterior', 'Premolar', and 'Molar' are all checked. In the 'Methods' section, 'Full-Contour' is selected among various options like 'Screw-Retained', 'Offset', 'Cement-Retained', etc. Below this, a 'Do not use' checkbox is present. A section titled 'Full-Contour 1' is expanded, showing a 'Materials' list where 'Zirconia HT' is selected. Below the materials, the 'Zirconia HT' section is expanded, showing 'Shade Systems'. In this list, 'Lab Shades' is selected with a blue checkmark, and a red arrow points to the '+' icon next to it. Other shade systems like 'VITA Classic', 'VITA 3D-Master', 'Ivoclar', and 'Material' are also listed with their respective shade counts.

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 Change History 

2026.09.01 | 09:19

> Teeth Custom Abutment + Implant Crown Added

2026.09.01 | 09:19

> Teeth Custom Abutment + Implant Crown Removed

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2026.09.01 | 09:19

 Teeth Custom Abutment + Implant Crown Added

Anterior Added

2026.09.01 | 09:19

> Teeth Custom Abutment + Implant Crown Removed

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Notice
NFT Service Termination Notice
2023-04-19

To our users,

We want to thank all Medit Link users who have supported and used our NFT service.

Unfortunately, we must inform you that this service will be terminated following the upcoming program update.

Please see the details below for more information:

- **Feature in Question**
 - Share via NFT (in Medit Link Case Talk)
- **Effective Date**
 - Starting from the update in the first half of 2023 (scheduled for April 27)

※ The update schedule is subject to change.

We'll continuously work to provide better service.

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Register Offline Data

[Register](#)

Finds and registers offline data on this device that is not registered in Medit Link.

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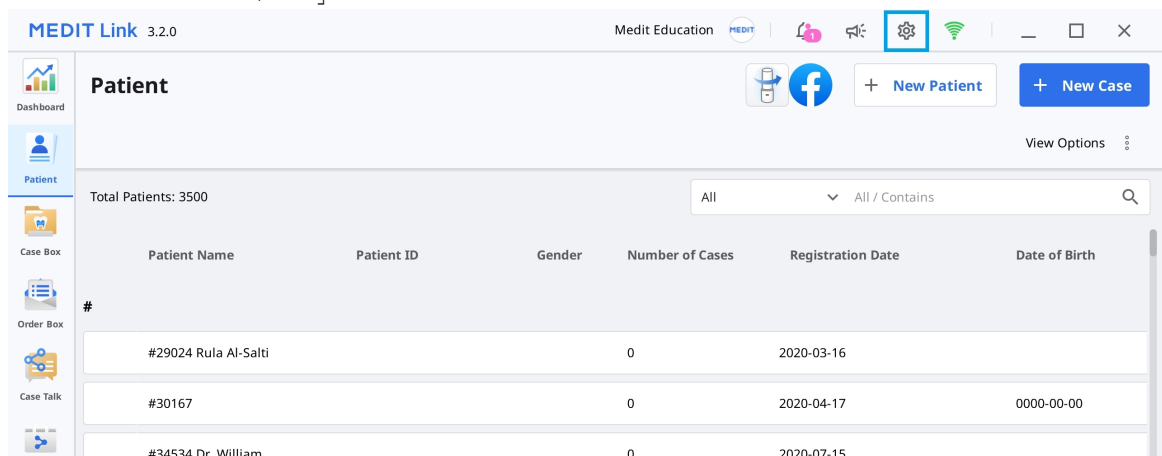
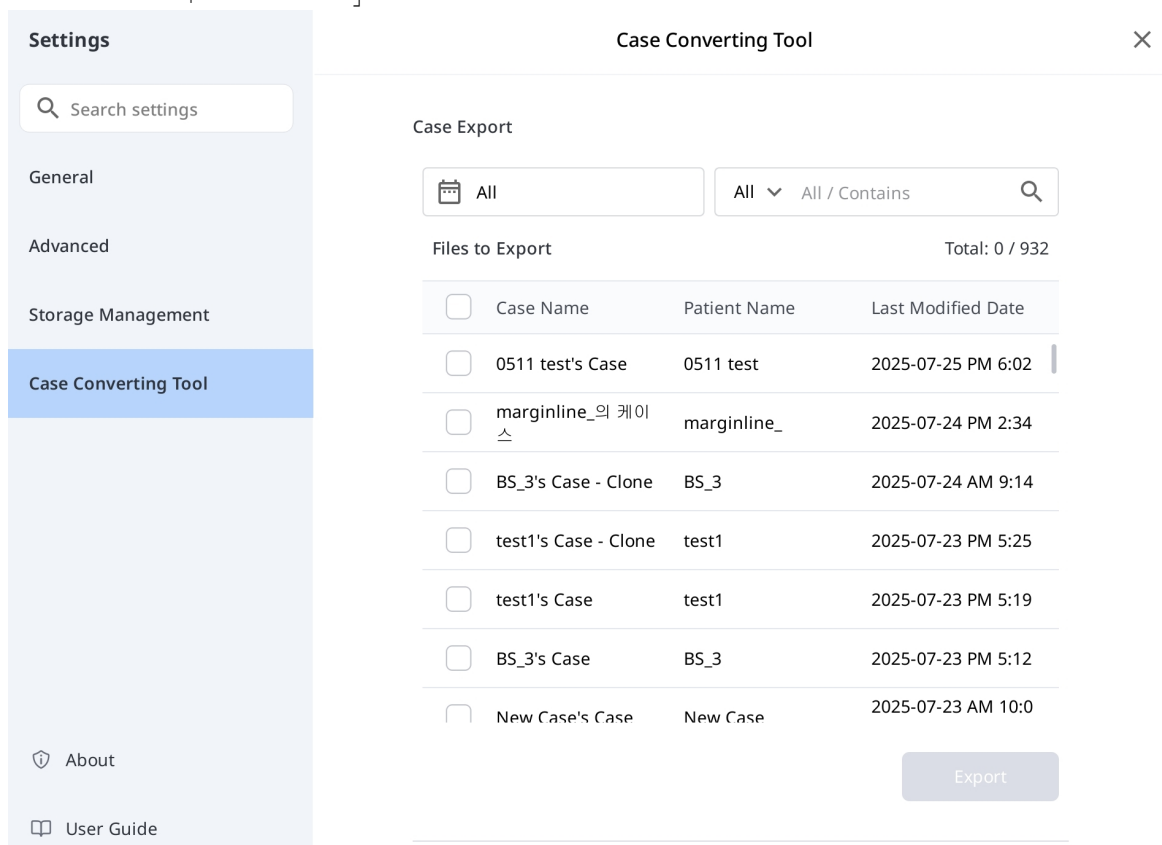
Question

Exporting successfully completed.
Do you want to open the folder containing the exported file?

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Settings

- General
- Advanced
- Storage Management
- Case Converting Tool**
- About
- User Guide

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Case Import

A new case is created for every imported case. It contains case data, patient information, a nd scan file. However, it does not contain order information.

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 Meditrial Srl
Via Po 9 00198, Rome Italy

 Meditrial Europe Ltd
Bahnhofstrasse 23 6300 Zug, Switzerland

 Meditrial UK Ltd
19 The Circle, Queen Elizabeth Street, London, SE1 2JE, United Kingdom

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
Merge all patient information to
Please select the patient you'd like to merge all data to.
After merging all data together, all information in other patient files will be deleted.

	Patient Name	Patient ID	Gender	Birth Date	Number of Cases	Date Registered	Recent Visit Date
☒	5-9_jScan_CAMOrd...				1	2020-09-18	2020-09-18
☐	test				2	2020-09-17	2020-09-17
☐	test				22	2019-07-26	2020-09-17

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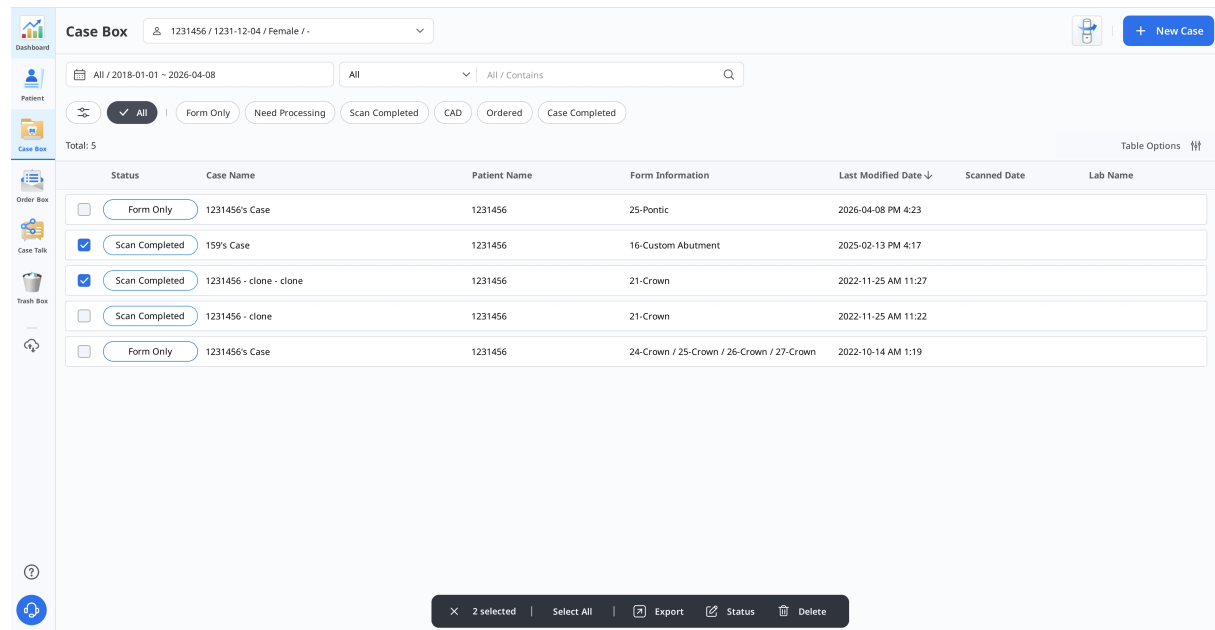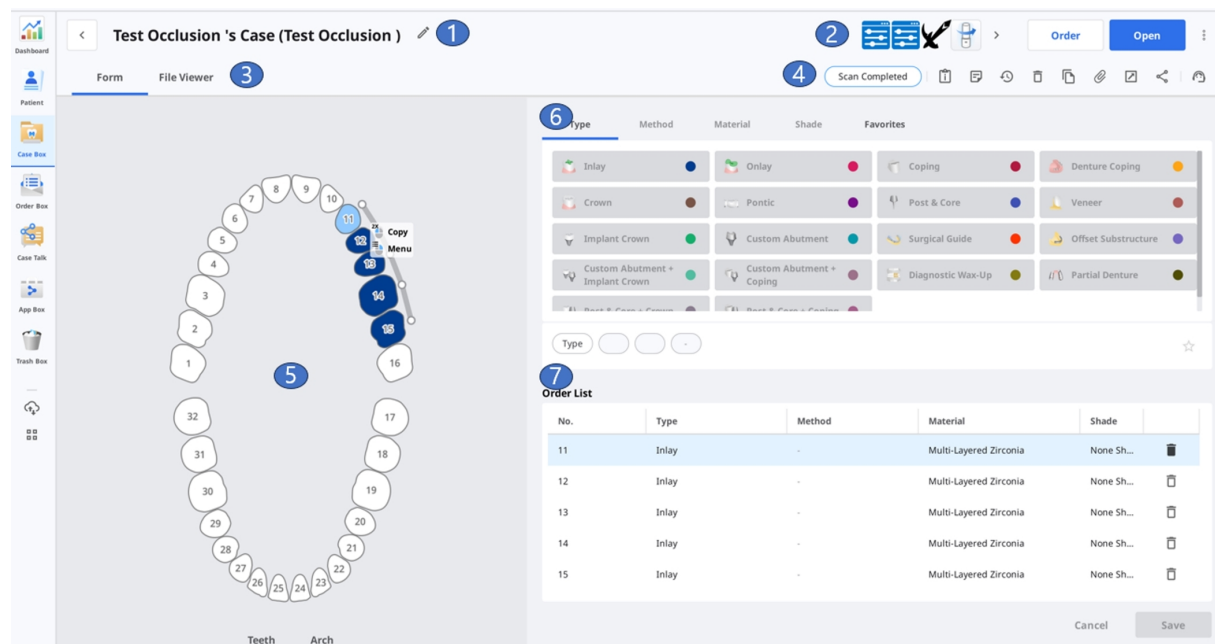
All / Contains

🔍

🔄 All | Pending | Canceled | Rejected | Accepted | Shipped | Case Completed

Total: 254

Table Options

	Status	Case Name	Patient Name	Order Date	Requested Delivery Date	Lab Name	Order ID
☐	Accepted	texture on's Case - Clone	texture on	2026-04-08 PM 5:40		MEDIT_LAB_STAGE	6012449
☒	Pending	Express Scan	Seol register	2026-04-02 PM 9:38	2026-04-09 AM 9:00	MEDIT_LAB_STAGE	6012242
☒	Case Completed	전설이의 케이스	전설이	2025-07-16 PM 5:46	2025-07-17 PM 3:00	MEDIT_LAB_STAGE	6004719
☒	Case Completed	27	1	2025-03-16 PM 3:11	2025-03-17 PM 3:00	MEDIT_LAB_STAGE	6002656
☐	Case Completed	3.3.0.136 TEST's Case	3.3.0.136 TEST	2024-04-12 PM 3:25	2024-04-13 PM 2:00	MEDIT_LAB_STAGE	5998987
☐	Case Completed	test의 케이스	Patient TEST	2024-04-09 PM 3:06	2024-04-10 PM 2:00	MEDIT_LAB_STAGE	5998955
☐	Case Completed	Mandibular Movement Case - Clone - Clone	MEDIT	2023-06-12 PM 12:12	2023-06-13 PM 2:00	MEDIT_LAB_STAGE	5996680
☐	Case Completed	Release New_2's Case	Release New_2	2023-04-25 PM 6:47	2023-04-25 PM 2:00	MEDIT_LAB_STAGE	5996180
☐	Case Completed	cX 문의's Case - Clone	cX 문의	2023-05-16 PM 5:26	2023-05-17 PM 2:00	MEDIT_LAB_STAGE	5996393
☐	Case Completed	cX 문의's Case	cX 문의	2023-05-15 PM 2:39	2023-05-16 PM 2:00	MEDIT_LAB_STAGE	5996362
☐	Case Completed	310.187's Case	310.187	2023-05-15 PM 1:33	2023-05-16 PM 2:00	MEDIT_LAB_STAGE	5996360
☐	Case Completed	310.187 test's Case - Clone		2023-05-03 PM 2:00		MEDIT_LAB_STAGE	5996306

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MEDIT Link

Add more information about your business. > seolahjeon_ English

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Memo

Tags

No.	Type	Method	Material	Shade
23	Crown	Anatomic	Multi-Layered Zi...	None

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Watch the Video >

The Doctors' Choice

Medit Intraoral Scanner Trusted by the Majority

Click to find out more >

Featured Apps (Top 5)

Diagnosis/Consultation

Order Placement

Medit Orthodontic Suite

Install

Diagnosis/Consultation

CAM

CAD

Medit ClinicCAD

CAD solution for creating a wide range of dental prostheses

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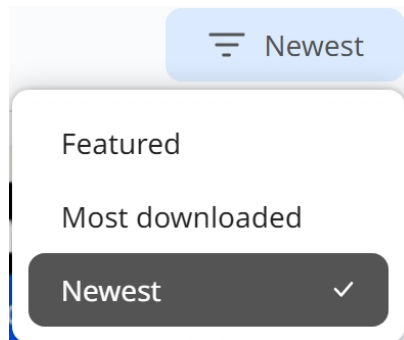
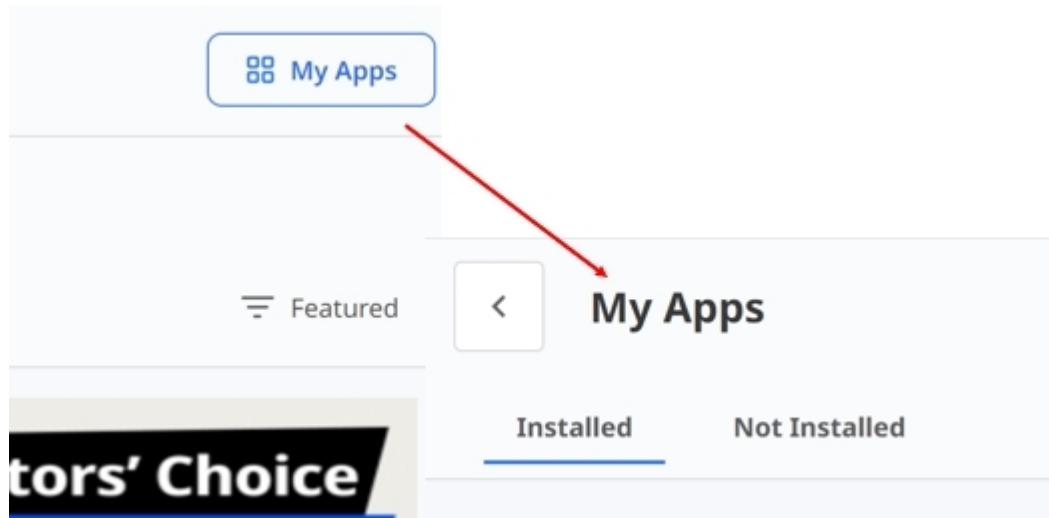
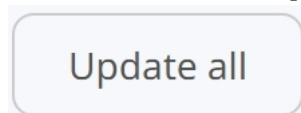
MEDIT Link 3.2.0 wonhoon choi

Trash Box

☐ Restore ☐ Selected Cases: 0 All / Contains

	Status	Case Name	Patient Name	Lab Name	Deletion Date ↓	Deleted By	Size
☐	Scan Completed	Jira, 임플란트 46, 47's Case - Clone - Clone	임플란트 46, 47		2023-04-26 오전 9:12	wonhoon choi	66,654 KB
☐	Form Only	G11598492304180810429	test test		2023-04-18 오후 5:15	wonhoon choi	0 KB
☐	Form Only	G11598492304180806313	test test		2023-04-18 오후 5:15	wonhoon choi	0 KB
☐	Form Only	G11598492304180811349	tes test		2023-04-18 오후 5:15	wonhoon choi	0 KB
☐	Form Only	G11598492304180754083	Glidewell_fastdesign		2023-04-18 오후 5:00	wonhoon choi	0 KB
☐	Form Only	G11598492304180733332	Glidewell_fastdesign		2023-04-18 오후 5:00	wonhoon choi	0 KB
☐	Form Only	G11598492304180648146	test test		2023-04-18 오후 5:00	wonhoon choi	0 KB
☐	Case Completed	Multi Occlusion Normal Case's Case - Clone	Multi Occlusion Normal Case		2023-04-18 오전 10:50	wonhoon choi	23,652 KB
☐	CAD	Multi Occlusion Normal Case's Case - Clone - Clone	Multi Occlusion Normal Case		2023-04-18 오전 10:50	wonhoon choi	23,652 KB
☐	Form Only	Occlusion Test's Case	Occlusion Test		2023-04-17 오후 6:31	wonhoon choi	0 KB
☐	Form Only	Test Occlusion 's Case	Test Occlusion		2023-04-17 오후 6:31	wonhoon choi	0 KB
☐	Form Only	Occlusion Test's Case	Occlusion Test		2023-04-17 오후 6:31	wonhoon choi	0 KB
☐	Form Only	Occlusion Test's Case	Occlusion Test		2023-04-17 오후 6:26	wonhoon choi	0 KB

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Total Cases: 210

☒ No Date
Data Type
All

☐ Case Name	Patient Name	Scanned Date↓	Size	Status
☐ BLATRIX Charline's Case	BLATRIX Charline	2025-01-06 PM 6:28	729.30 MB	Upload Required
☐ BLATRIX Charline's Case	BLATRIX Charline	2025-01-06 PM 6:28	729.30 MB	Upload Required
☐ BLATRIX Charline's Case - 복제	BLATRIX Charline	2025-01-06 PM 6:28	733.35 MB	Upload Required
☐ BLATRIX Charline's Case - Clone	BLATRIX Charline	2025-01-06 PM 6:28	734.01 MB	Upload Required
☐ test의 케이스	Patient TEST	2023-04-26 PM 7:58	1.25 GB	Upload Required
☐ i700 HD's Case	i700 HD	2023-04-24 PM 5:23	911.24 MB	Upload Required
☐ Mandibular Movement	MEDIT	2023-04-24 PM 1:59	1.88 GB	Upload Required
☐ Mandibular Movement Case	MEDIT	2023-04-24 PM 1:59	1.75 GB	Upload Required
☐ Mandibular Movement Case - Clo...	MEDIT	2023-04-24 PM 1:59	1.75 GB	Upload Required
☐ 서라_Mandibular Movement_TEST	Mandibular Movement_export_test	2023-04-10 PM 3:39	34.73 MB	Upload Required
☐ 서라_Mandibular Movement_TEST...	Mandibular Movement_export_test	2023-04-10 PM 3:39	731.43 MB	Upload Required

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
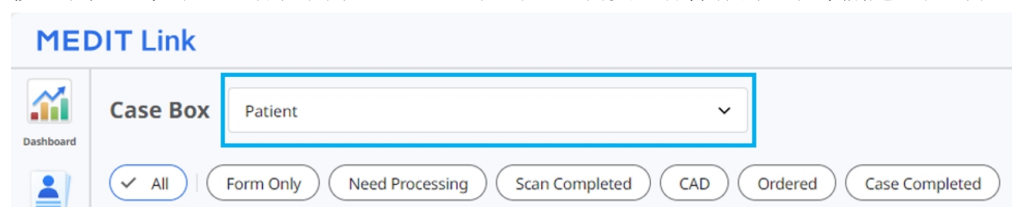
Case Box Patient

Dashboard

✓ All Form Only Need Processing Scan Completed CAD Ordered Case Completed

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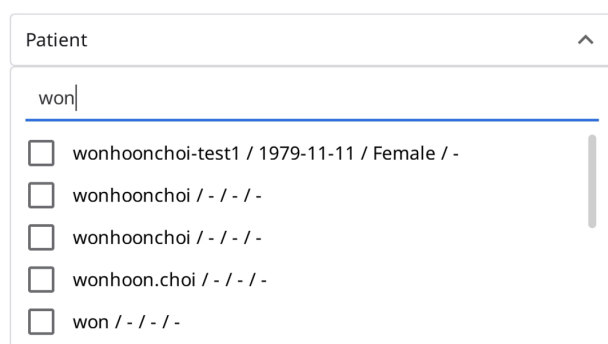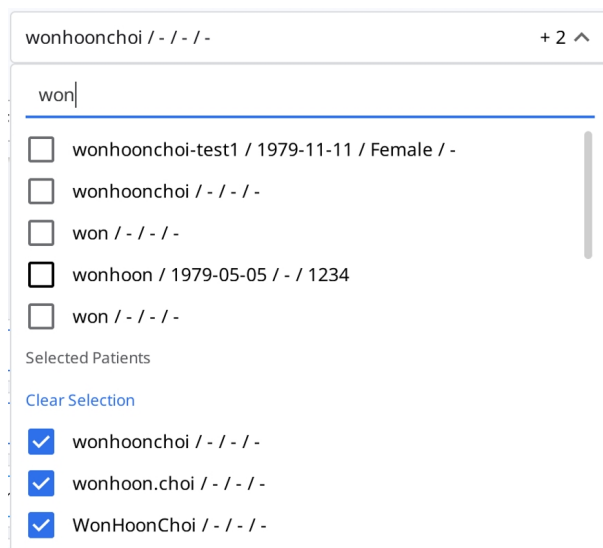
MEDIT Link 3.2.0 Wonhoon Choi_clinic

Dashboard
 Patient
 Case Box
 Order Box
 Case Talk
 App Box
 Trash Box
 More

Case Box

+ 2

All |
 Form Only |
 Need Processing |
 Scan Completed |
 CAD |
 Ordered |
 Case Completed

Total Cases: 20

2018-01-01~2023-10-24

All ▼ / Contains 🔍

Status	Patient Name	Case Name	Form Information	Last Modified Date ↴	Scan Date	Lab Name	Order Date
Form Only	wonhoonchoi	originalcase	-	2023-10-05 오전 8:10			
Form Only	wonhoonchoi	wonhoonchoi's Case	-	2023-10-05 오전 8:09			
Scan Completed	wonhoonchoi	FileVwerTest_1	8-Inlay / 10-Inlay / 12-Inlay / 13-Inlay / 14-Inl...	2023-06-15 오전 8:14	2023-02-09 오전 8:35		
Scan Completed	wonhoonchoi	FileVwerTest_1 - Clone	2-Inlay / 3-Inlay / 4-Inlay / 5-Inlay / 6-Inlay / 8...	2023-05-25 오후 2:05	2023-02-09 오전 8:35		
Form Only	wonhoonchoi	FileVwerTest_1 - Clone	8-Inlay / 10-Inlay / 12-Inlay / 13-Inlay / 14-Inl...	2023-03-15 오후 3:06			
Scan Completed	wonhoonchoi	FileVwerTest_1 - Clone - Clone	8-Inlay / 10-Inlay / 12-Inlay / 13-Inlay / 14-Inl...	2023-03-09 오후 4:41	2023-02-09 오전 8:35		
Need Processing	wonhoonchoi	FileVwerTest_1 - Clone	8-Inlay / 10-Inlay / 12-Inlay / 13-Inlay / 14-Inl...	2023-02-14 오전 8:13	2023-02-09 오전 8:35		
Scan Completed	wonhoonchoi	FileVwerTest_1 - Clone	8-Inlay / 10-Inlay / 12-Inlay / 13-Inlay / 14-Inl...	2023-02-10 오전 11:27	2023-02-09 오전 8:35		
Form Only	wonhoonchoi	testsamplecase	-	2023-02-10 오전 7:39			
Scan Completed	wonhoonchoi	Multi Occlusion	-	2023-02-09 오전 9:26	2023-02-09 오전 9:01		
Scan Completed	wonhoonchoi	FileVwerTest_1 - Clone - Clone	8-Inlay / 10-Inlay / 12-Inlay / 13-Inlay / 14-Inl...	2023-02-09 오전 8:50	2023-02-09 오전 8:35		
Case Completed	wonhoonchoi	wonhoonchoi's Case	12-Inlay / 13-Inlay / 14-Inlay	2022-10-13 오후 10:30	2022-10-11 오후 5:16		
Case Completed	wonhoonchoi	wonhoonchoi's Case	-	2022-10-13 오후 10:28			

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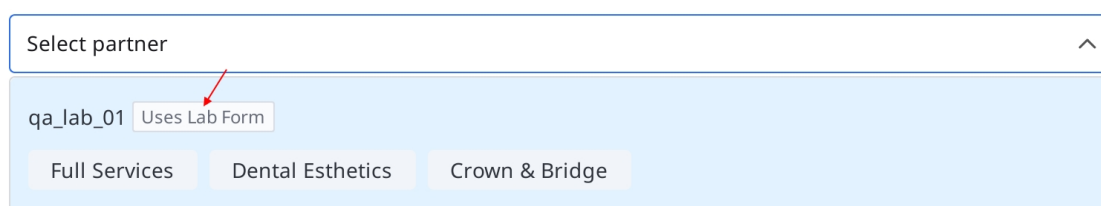
The screenshot shows a 'Select partner' dropdown menu. The selected partner is 'qa_lab_01', which has a tag that says 'Uses Lab Form'. Below the partner name, there are three buttons representing different service types: 'Full Services', 'Dental Esthetics', and 'Crown & Bridge'. A red arrow points to the 'Uses Lab Form' tag.

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The screenshot shows a dental arch diagram with teeth numbered 1-32. Tooth 26 is highlighted. The form on the right has tabs for Type, Method, Material, Shade, and Favorites. Under the Type tab, 'Crown' is selected. The Order List table shows one entry for tooth 26, type Crown, method Full-Contour, material Zirconia, and shade A1. A red arrow points to the 'Save' button at the bottom right.

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MEDIT Link 3.2.0

Case Box wonhoonchoi-test1 / 1979-11-11 / Female / - + 2 ▾

☒ All |
 ☐ Form Only |
 ☐ Need Processing |
 ☐ Scan Completed |
 ☐ CAD |
 ☐ Ordered |
 ☐ Case Completed

☒ **Export** ▾ |
 Status ▾ |
 ☐ |
 Selected Cases: 4

	Status	Patient Name	Case Name	Form Information
☒	Scan Completed	wonhoonchoi-test1	wonhoonchoi-test1's Case	-
☒	Form Only	wonhoonchoi-test1	wonhoonchoi-test1's Case	-
☒	Form Only	wonhoonchoi	originalcase	-
☒	Scan Completed	wonhoonchoi-test1	wonhoonchoi-test1's Case	-

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Export ▾

- Result Files**
- Case List (as shown)
- Case List (with all items)

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Share through Case Talk

Share with Anyone

Create a shared link and share it with anyone.

Expiration Date

2022-04-27

You can extend the period up to 90 days from today.

PIN Code

☐ PIN Code

Enter PIN Code

Set the PIN code using six digits for privacy protection.

☒ Share the patient's name

☒ Grant Download Permission Premium Feature

If you check this option, all users who have access to the shared link can download the files.

The premium feature is provided to all users for a limited time. Feature availability is subject to change depending on the membership plan.

☒ Allow comments

If you check this option, you will be able to communicate through the shared link.

Comment

Please enter your comment.

0/1024


By pressing the "Get Link" button, you agree that you understand personal information can be collected and provided to a third party (or shared in any public domain) ONLY with the consent of the person whose images on which you have worked. You also agree the consent of that person was truthfully obtained, and Medit is not responsible for any problems which may arise there after.

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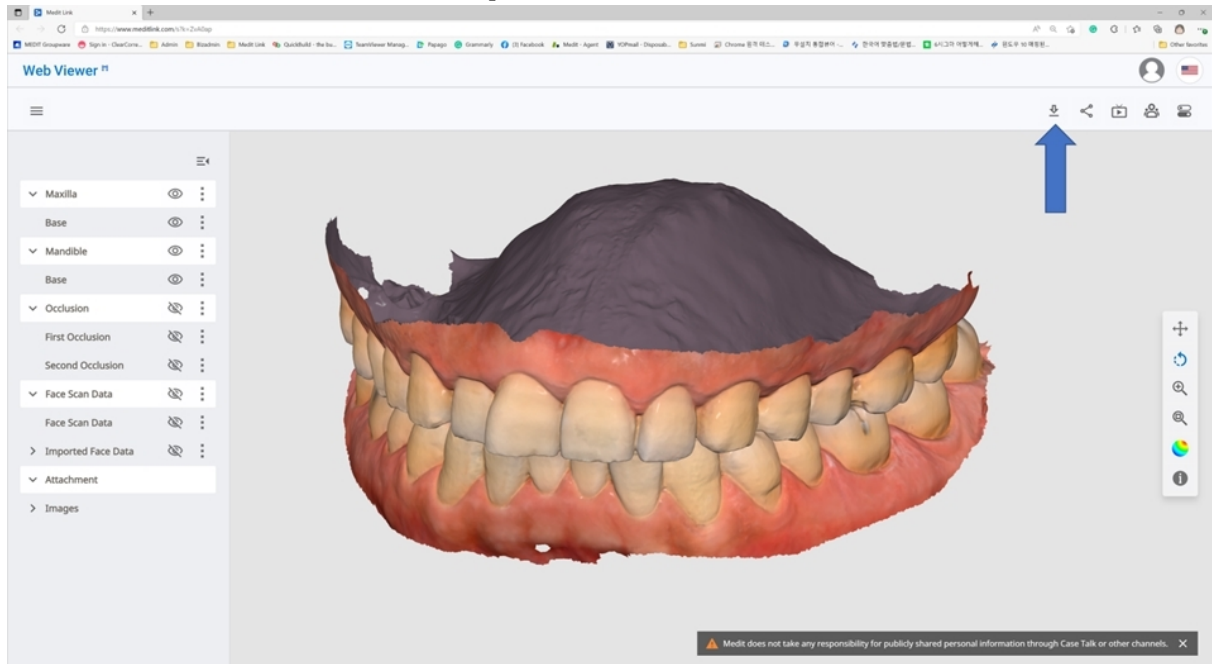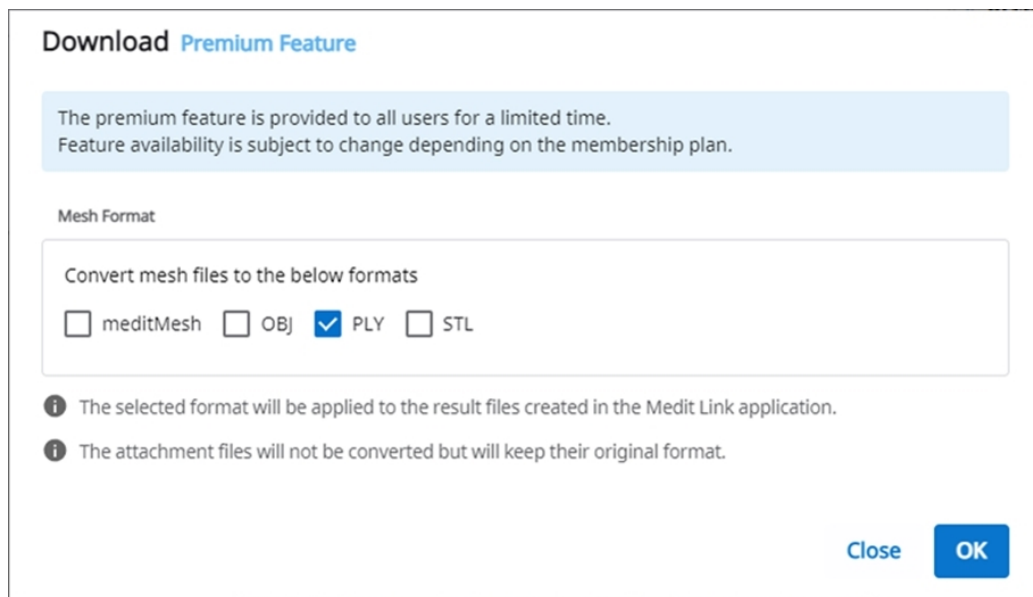
Partner

qa_lab_01

qa_lab_01 Uses Lab Form
San Diego, US
[92131] Residence Inn San Diego Rancho Bernardo
101111111

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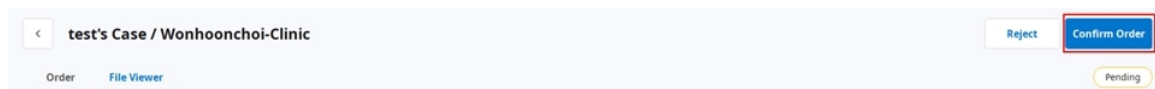
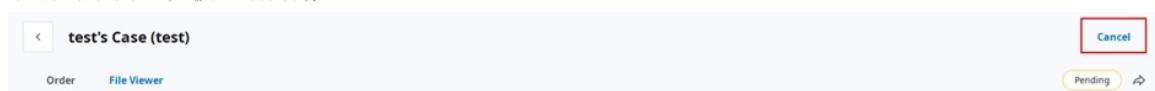
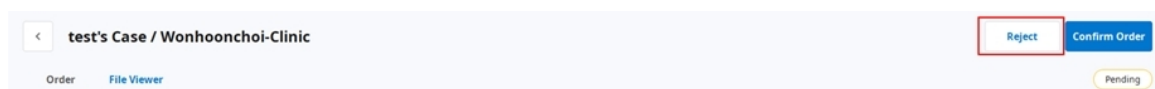
The screenshot shows a web interface for managing a case. The left sidebar contains icons for Dashboard, Patient, Case Box, Order Box, Case Talk, and Trash Box. The main area is titled 'test1's Case (test1)' and has tabs for 'Order' and 'File Viewer'. The 'Order' tab is active, showing a form with the following fields:

- Order ID: 6004858
- Assignee: (dropdown menu)
- Patient Name: test1
- Partner Name: Elysee Dental Test
- Scanned Date: (empty)
- Order Date: 2025-07-23 PM 5:17
- Requested Delivery Date: 2025-07-30 PM 05:00
- Expected Delivery Date: 2025-07-30 PM 05:00
- Memo: (empty text area)
- Tags: test

On the right, there is a 'Shipment' section with a table:

No.	Type	Method	Material	Shad
15	Crown	Anatomic	Zirconia	A4

At the top right of the main area, there is a 'Cancel Order' button and a 'Pending' status indicator.

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Shipment Registration

Select a delivery method and enter the required information to register.
The information you enter is saved in the Shipment tab and can be checked at any time.

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Case Name
mandibular movement - Clone - Clone

Patient Name
Test Occlusion

☒ Share patient information

Requested Delivery Date
2023-04-13

Delivery Time
오전 8:15 AM 8:00 PM 8:00

Memo

Tags
You can add the entered tag with Enter or Spacebar.

☒ Include all occlusion relationships in the order
☒ Include Tags In The Order

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Tags
tag

Order List Shipment

No.	Type	Method	Material	Shad
25	Crown	Anatomic	Multi-Layered Zirconia	A3.5
26	Crown	Anatomic	Multi-Layered Zirconia	A3.5

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Status	Case Name	Patient Name
☐ Scan Completed	0218's Case	0205 test2mn(456)
2024-12-23		
☐ Ordered	ORG_[Clinics][3.3.6] 01 Normal - Complete	ORG_[Clinics][3.3.6] 01 Normal
☒ Scan Completed	ORG_[Clinics][3.3.6] 01 Normal - Complete	ORG_[Clinics][3.3.6] 01 Normal
☒ Scan Completed	ORG_[Clinics][3.3.6] 01 Normal - Complete	ORG_[Clinics][3.3.6] 01 Normal
☒ Scan Completed	ORG_[Clinics][3.3.6] 01 Normal - Complete	ORG_[Clinics][3.3.6] 01 Normal

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BS9 25-Crown / 26-C

3 selected | Select All | Export | Status | D

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Patient Name*

Case Name*

 Creating patient information involves consent from the patient.
Patient information can be shared with 3rd parties while placing an order.

Tag

Doctorname × Specialmaterial × SpecialComment ×

 You can add up to 10 tags of 100 characters or less.

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Memo

Tags

Doctorname × Specialmaterial × SpecialComment ×

● You can add up to 10 tags of 100 characters or less.

☒ Include all occlusion relationships in the order

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