

MEDIT Link



ME-UG-702H
Revision 13 (29.09.2026)
SW version 3.5.0

Table of contents

Medit Link

Introduction	12
Medit Link	12
Service Overview	12
Cloud Service	13
Workflow Management Platform with Ordering System	14

Installation

System Requirements	15
System Requirements for Windows	15
System Requirements for macOS	15
System Requirements for Web Service	15
Recommended for Browsers	16
Minimum for Mobile Devices	16
Installation on Windows	17
Installation on macOS	20

Getting Started

Sign Up & Sign In	24
Medit Link Website	24
Sign Up	24

Login and Welcome Page	27
Medit Link Web	
Dashboard	29
Dashboard Overview	29
Medit Insights	30
Subscribe Medit Insights	30
Unsubscribe Medit Insights	32
Patient	33
Search Patient	33
Case Box/Work Box	34
Case Filtering	34
Filter by patient	34
Filter by status	36
Filter by duration	36
Filter by Search	37
Case Details	38
Order Box/In Box	39
Order Details	39
Case Talk	40
Creating a Case Talk	40
Setting the sharing permission type and sharing options	41
Case Talk details	43

Web Viewer	45
Annotation	46
Comments	48
Real-time screen sharing	48
Data Tree	49
Data Tree Controls	49
3D Model Controls	50
Partners	51
Search for Partners	51
My Partner List	51
Partner Invitation	52
Membership	53
Manage Plan	53
Managing an Active Plan	53
Subscribing to the 10 TB Plan	54
Canceling a Subscription	57
Membership	58
Medit Cloud Storage Plan	59
Cloud Storage	60
Payment Details	62
View Bill Details	64
Vouchers	64

Registering a Voucher	65
Activating a Voucher	65
My Info	67
My Account	67
My Clinic/My Lab (Admin Account Only)	67
My Scanner	68
Order Forms (Lab Admin Account Only)	70
Viewing Order Forms	70
Editing Order Forms	71
Managing Options	76
Viewing Change History	84
Other Menu	85
Notice & Update	85
Help & Resources	85
Download Medit Link App	86
Medit Link App	
Overview	87
The table below summarises what menus are available in each type of clinic and lab account.	87
Title Bar	89
My Account (My Information)	89
Notifications	90

Notice & Update	90
Settings	90
Network Status	91
Settings	93
Search Settings	93
General	93
Advanced	95
Export – Export Settings	96
Network – Proxy Server	99
Default App	100
User Experience Program	100
Storage Management – Overview & Cleanup	100
Set Period	101
Upgrade Storage	101
Storage Management – Advanced Settings	101
Case Converting Tool	102
Export Case	102
Import Case	104
About	105
User Guide	106
Menu	107
Dashboard	107

Patient	107
Add New Patient	108
Delete Patient	108
Merge Patients	109
Case Box/Work Box	110
Tabs in Case Box/Work Box	111
Functions in Case Box/Work Box	111
Order Box/In Box	114
Case Talk	115
App Box	116
Medit Apps Only	117
Sorting Options	117
My Apps	117
Trash Box	118
Case Synchronization Manager	118
Upload	119
Download	119
Help & Resources	120
File Viewer	121
Data Tree	121
Model Control Tools	122
3D Model Controls	122

2D Image Control	125
Function Tools	125
Image Capture and Memo	125
Ortho Capture	126
Fill Holes	127
Trimming	128
Optimize Occlusion Alignment	132
Mandibular Movement	134
Manage Multi Occlusion Data in File Viewer	135
Table Options	136
Status Filters	136
View Type	137
Grouping	138
Sorting	139
Filtering and Search	141
Filter by Patient	141
Filter by Date	142
Advanced Search Capabilities	142
Workflow	
Create/Edit/Clone/Delete Case	144
Creating Cases	144
Create Case in Case Box (Clinic)	144

Create Case in Work Box (Lab)	147
Editing Cases	148
Edit Case in Case Box (Clinic)	148
Edit Case in Work Box (Lab)	151
Cloning Cases	153
Clone Case in Case Box (Clinic)	153
Clone Case in Work Box (Lab)	153
Deleting Case	153
Delete Case (Clinic and Lab)	153
Form Information	155
Registering Form Information in Case Box	155
Enter on the Form Tab	155
Using a Lab-Defined Order Form	155
Registering Form Information in Work Box	156
Enter on the Form Tab	156
Connecting Teeth into Bridge	157
Setting Scanning Options	157
Medit Scan for Labs	158
Adding to Favorite	159
Add to Favorites	159
Manage Favorite Templates	159
Work on Scanning/CAD	160

Scan Process	160
CAD Process	160
Utilize Add-in Apps	160
Add-in Quick Launcher	160
App Launcher Management	160
Get Result Files	162
Export Scan Result Files	162
Export Result Files from Case List	163
Export Result Files from File Viewer	166
Export Result Files from Data Tree	167
Export Case List as CSV	170
Export Margin Line	171
Export Margin Line as a Color	171
Export Margin Line as a File	172
Download Files from Medit Link Web	173
Case Box/Work Box	174
Web Viewer	174
Order/Manage Orders	177
Ordering Cases	177
Confirm Order (In Box)	178
Cancel Order (Order Box)	178
Reject Order (In Box)	178

Checking Ordered Cases	179
Entering Shipment Information (In Box)	179
Managing Delivery Date	180
Setting Requested Delivery Date (Clinic account)	180
Setting Estimated Delivery Date (Lab account)	181
Completing Case	182
Complete (Work Box/In Box)	182
Complete (Case Box/Order Box)	184
Communicate Using Tags	185
Add Tags When Ordering	186
Check Tags on Case Details	187

Introduction

Medit Link

Medit Link is a collaborative tool designed to facilitate communication between clinics and labs to provide an optimized workflow for each party.

It offers differentiated features based on the type and permissions of each account. The program provides a user-friendly user interface for managing the workflow, from scanning with oral and model scanners, ordering cases, and designing and producing prostheses.

Medit Link consists of a web service and a computer application.

Service Overview

		Web	App
		• Register and manage the account • Search and manage partners • Share cases via Web Viewer with a sharable link • Access cases made from the App • Manage orders	• Enter form information • Scan • Manage orders • Link with CAD • Manage patients
Clinic	Admin	Search and manage partners and your clinics, and share cases. • Dashboard • Patient • Case Box • Order Box • Case Talk • Partners • Membership • My Info	Enter form information, scan, place orders, and manage patients. • Dashboard • Patient • Case Box • Order Box • Case Talk • App Box • Trash Box

		Web	App
	Member	• Patient • Case Box • Order Box • Case Talk • My Info	• Patient • Case Box • Case Talk • Order Box • App Box • Trash Box
Lab	Admin	Search and manage partners and your lab, and share cases. • Dashboard • In Box • Work Box • Case Talk • Partners • Membership • My Info	Enter form information, scan, and manage orders. • Dashboard • In Box • Work Box • Case Talk • App Box • Trash Box
	Member	• In Box • Work Box • Case Talk • My Info	• In Box • Work Box • Case Talk • App Box • Trash Box

Cloud Service

Medit Link supports the automatic upload/download of databases and files through AWS (Amazon Web Services).

The data stored in the cloud is automatically synchronized to the computers of multiple users in the organization so that they can access the same data from any computer at any time. Cloud-based backup and recovery keep all your data safe. It keeps your data safe even in case of data loss due to sudden computer failure.


The cloud service allows partnered dentists and labs to quickly share data without the hassle of emailing or manually transferring the scans. All data, including patient information, is securely encrypted, stored, and transmitted. Medit complies with medical laws and regulations like HIPAA and GDPR.

Workflow Management Platform with Ordering System

Medit Link offers its ordering system.

- Medit Link allows users to establish partnerships between dentists and dental laboratories and share data quickly and easily.
- You can check and manage the ordering status.

System Requirements

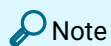

Please refer to the system requirements for Medit Scan for Clinics and Medit Scan for Labs if you want to use the scanning software with your scanner.

System Requirements for Windows

	Laptop & Desktop
CPU	Intel Core i5 2.6 GHz or faster
RAM	16 GB or more
Graphics	NVIDIA GeForce GTX 760 (2 GB) or higher, or equivalent AMD graphics card
OS	Windows 10 64 Bit (unavailable in 32 Bit) or later

System Requirements for macOS

	Mac
Processor	Apple M1 chip (8-core CPU, 7-core GPU, 16-core Neural Engine)
RAM	16 GB or more
OS	Sonoma 14 or later

System Requirements for Web Service

Recommended for Browsers

	Browsers
PC	Chrome, Firefox, Safari, Edge (Chromium)
Android	Chrome, Firefox
macOS	Safari, Chrome, Firefox


Minimum for Mobile Devices

	Android	iOS
Device	Galaxy S7 equivalent or higher	iPhone 7 equivalent or higher
OS	Android 7.1 Nougat	iOS 11.0
RAM	3 GB or higher	2 GB or higher
AP	Above Snapdragon 710	A10 Fusion


Installation on Windows

Medit Link is installed as a package with scan acquisition programs - Medit Scan for Clinics and Medit Scan for Labs.


1. Run Medit_LinkSetup.exe.
2. Select a language for installation.
3. Select the folder to install Medit Link. Please read and agree to the license terms and conditions and click "INSTALL."

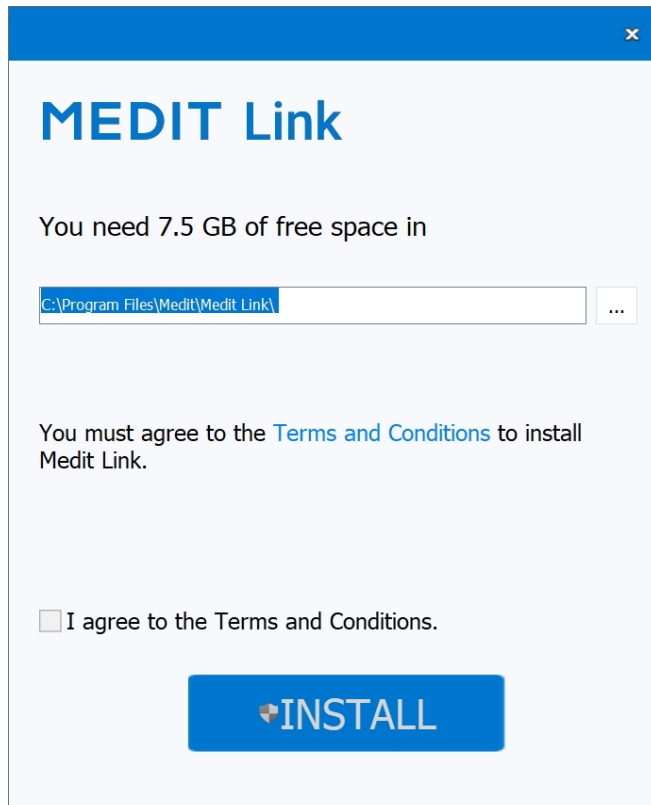


4. Make sure that the scanner cable is disconnected from the computer before proceeding to the next step. The installation may not proceed if the scanner is connected to your PC.
5. The program will be installed automatically and may take a few minutes. Do not turn off or unplug your PC until the installation is complete.

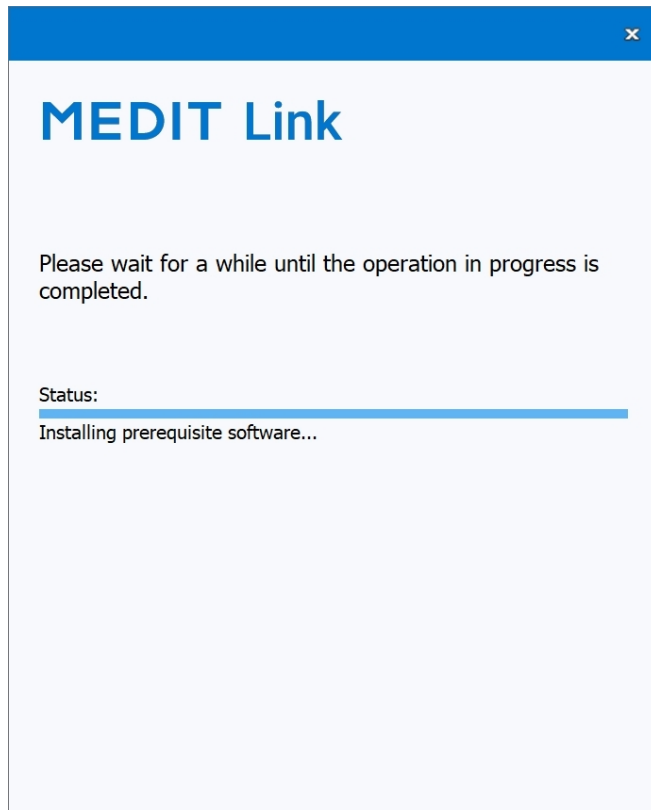

6. Please restart your PC after the installation is complete. You can proceed without restarting the PC if you do not have the scanner.

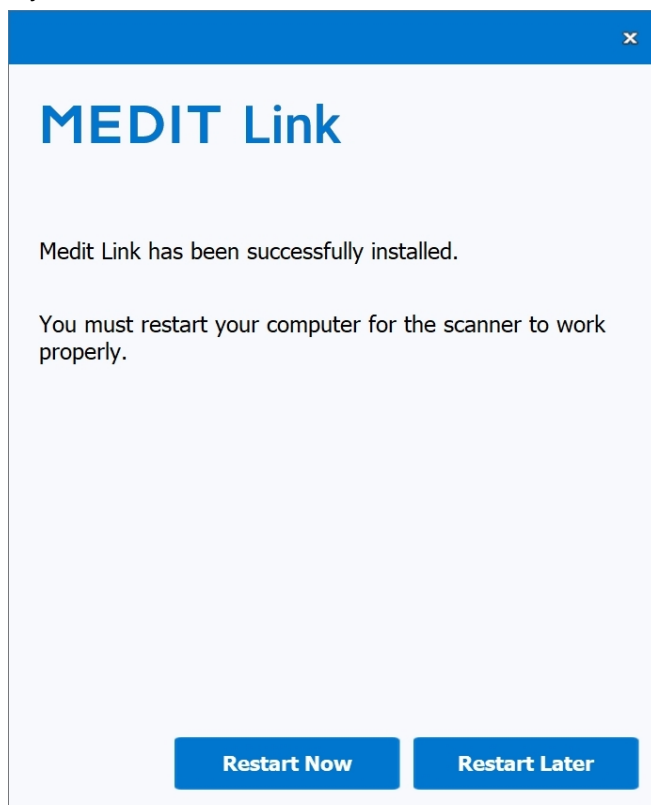

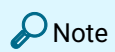


Note

If you do not restart your PC after the installation is complete, the scanner may not function properly.


Installation on macOS

Medit Link is installed as a package with Medit Scan for Clinics.


1. Run Medit_Link_x.x.x.pkg and click "Continue."

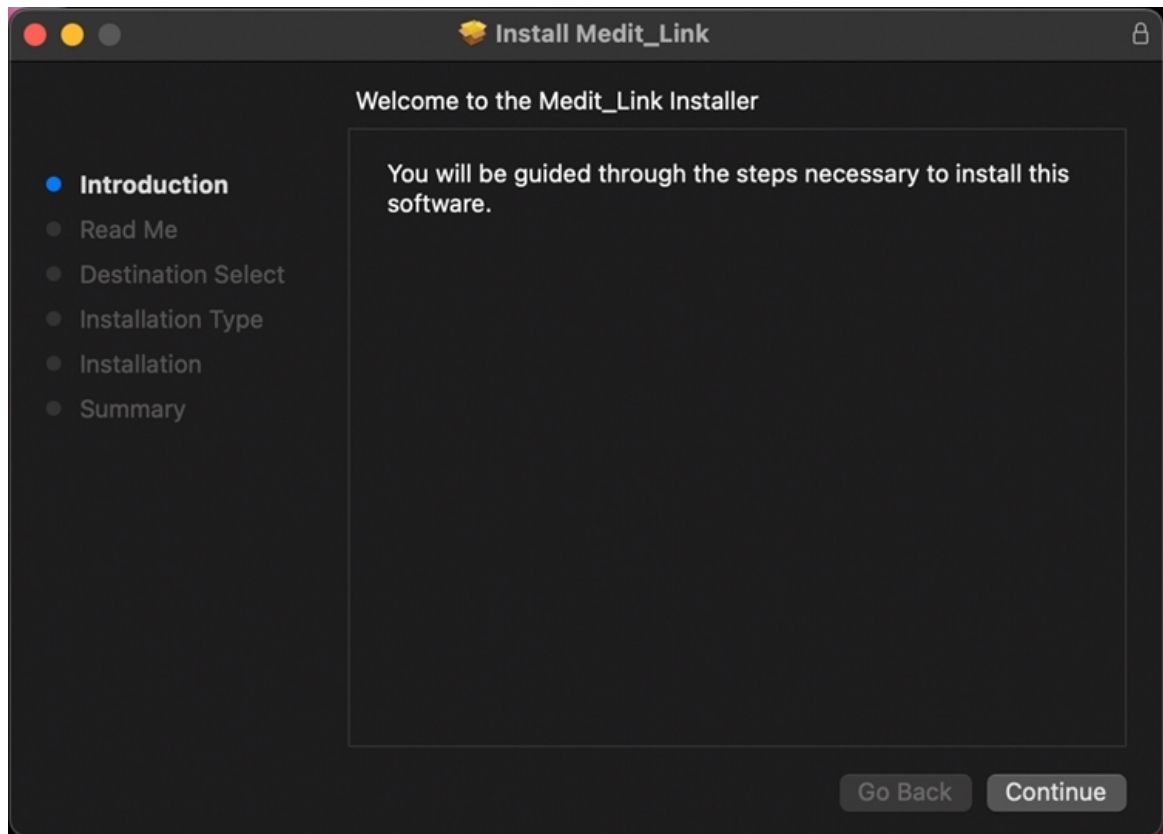


2. Read the End User License Agreement and click "Continue."

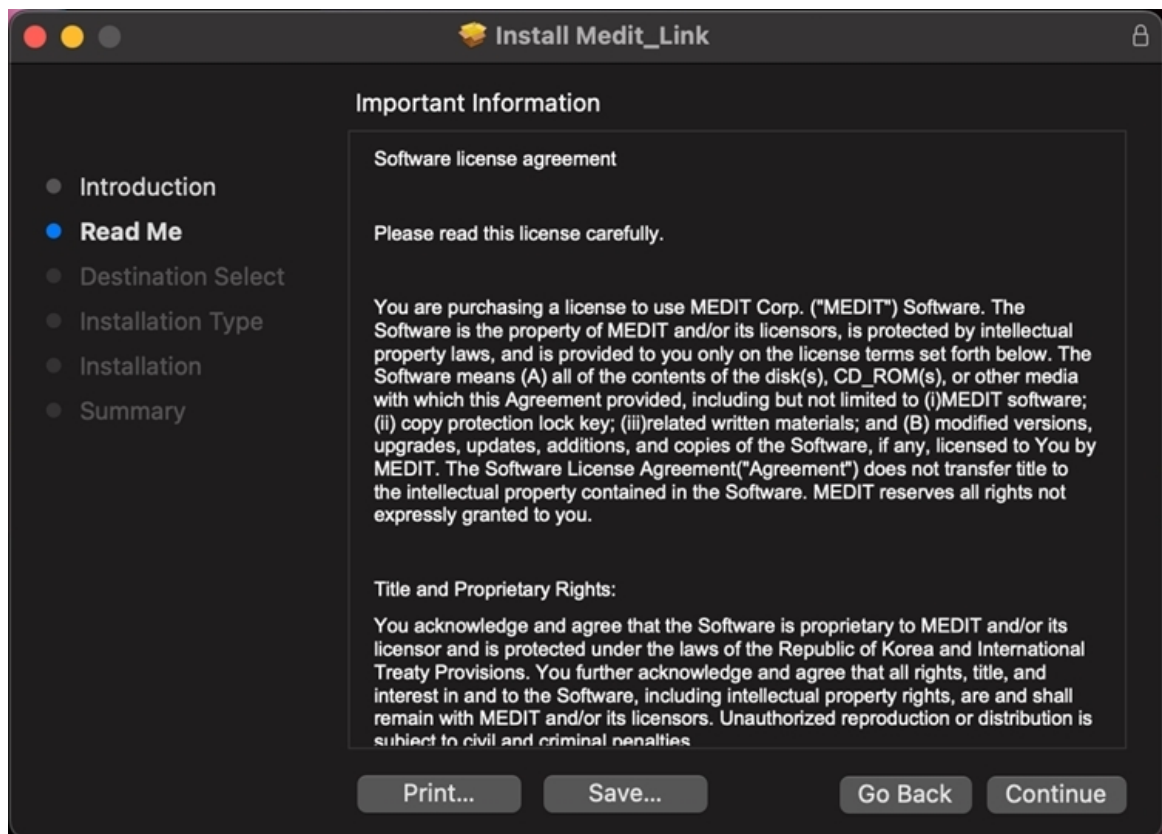

3. Select the disk to install and click "Continue."

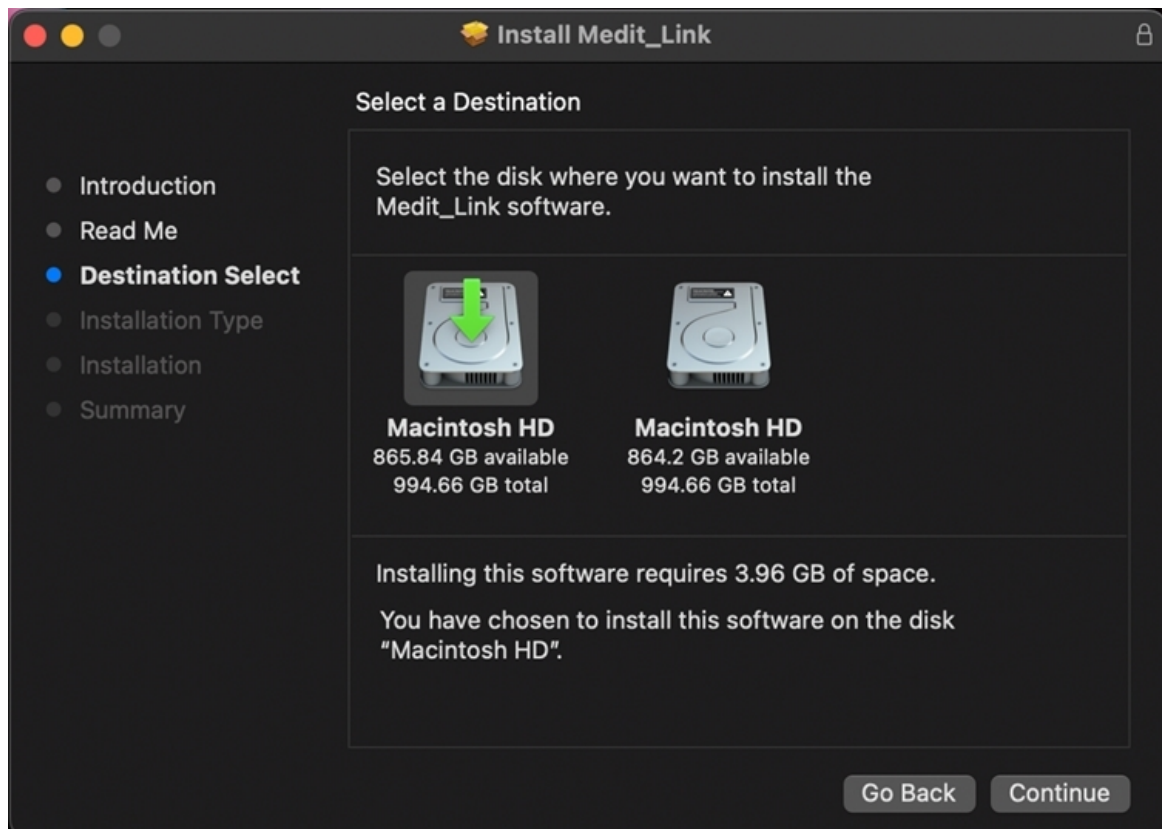


4. Click the "Change Install Location..." to change the installation destination.

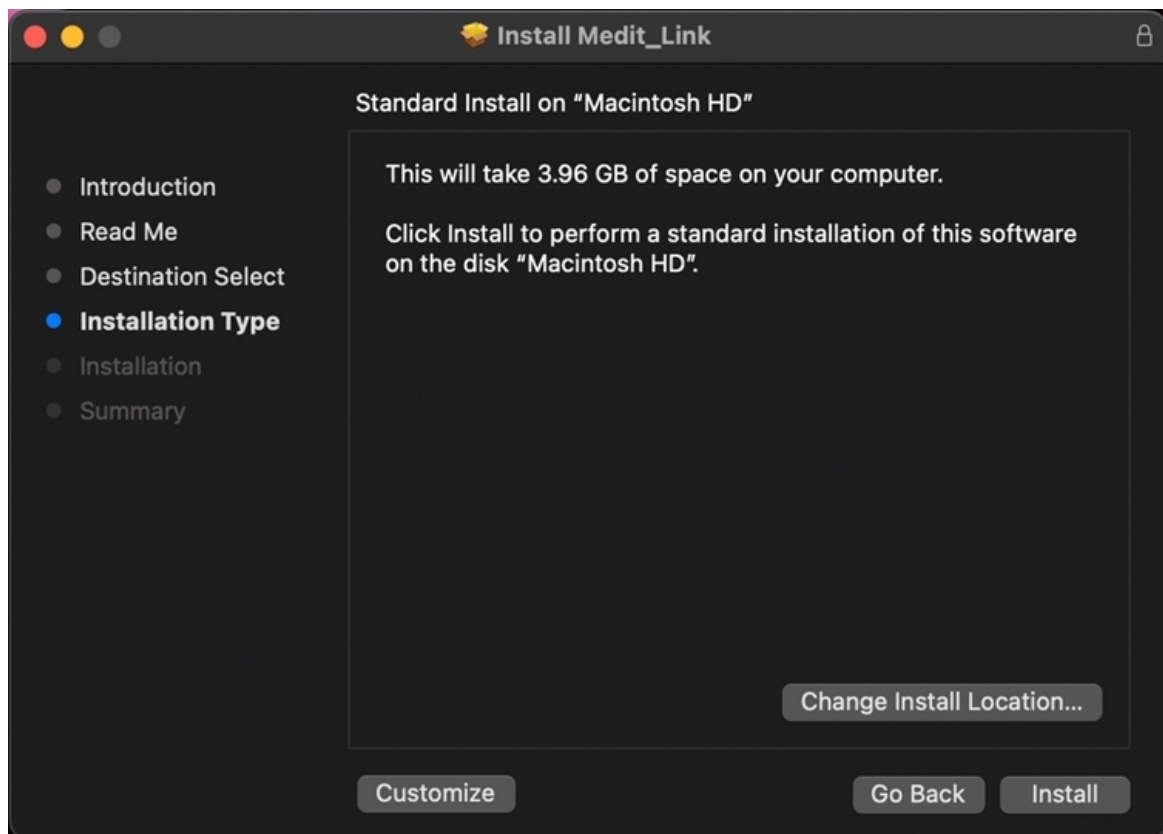


5. Click "Install" and type your Mac password to allow installation if requested.



6. Wait until the installation is completed.

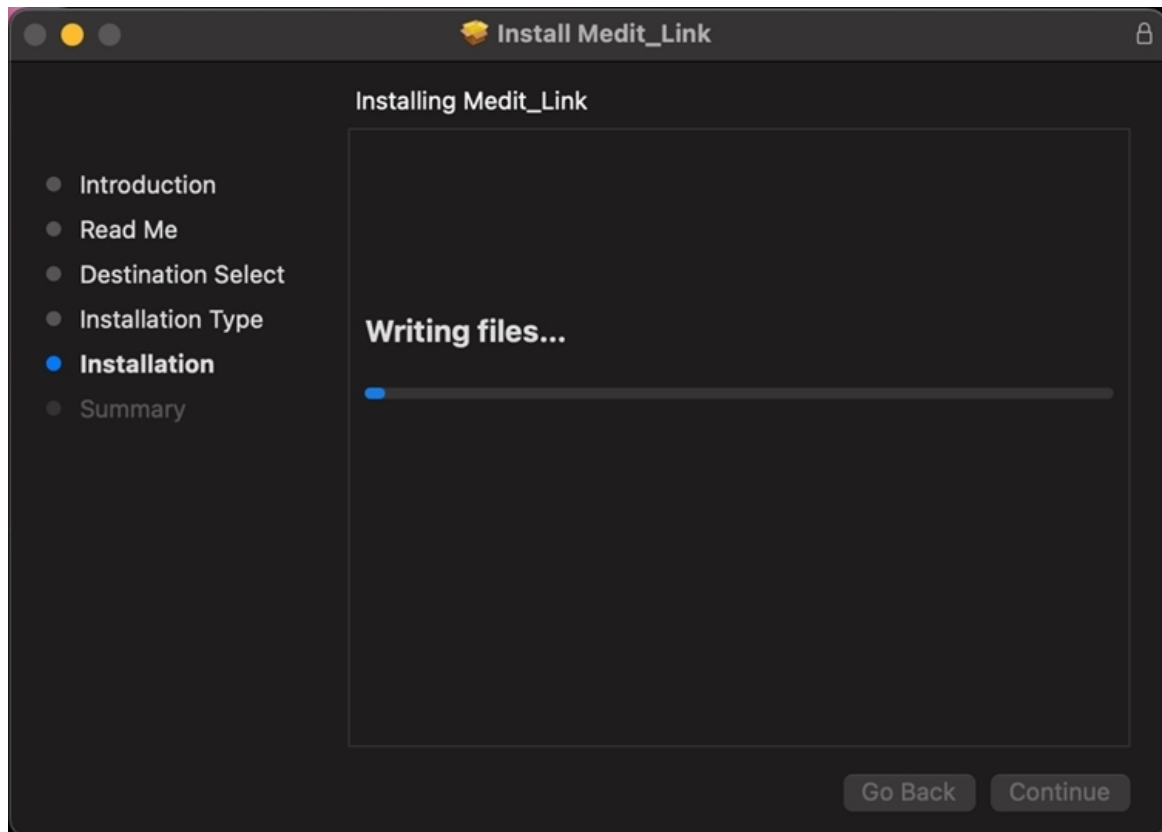



7. Click "Close" to exit the installation window.

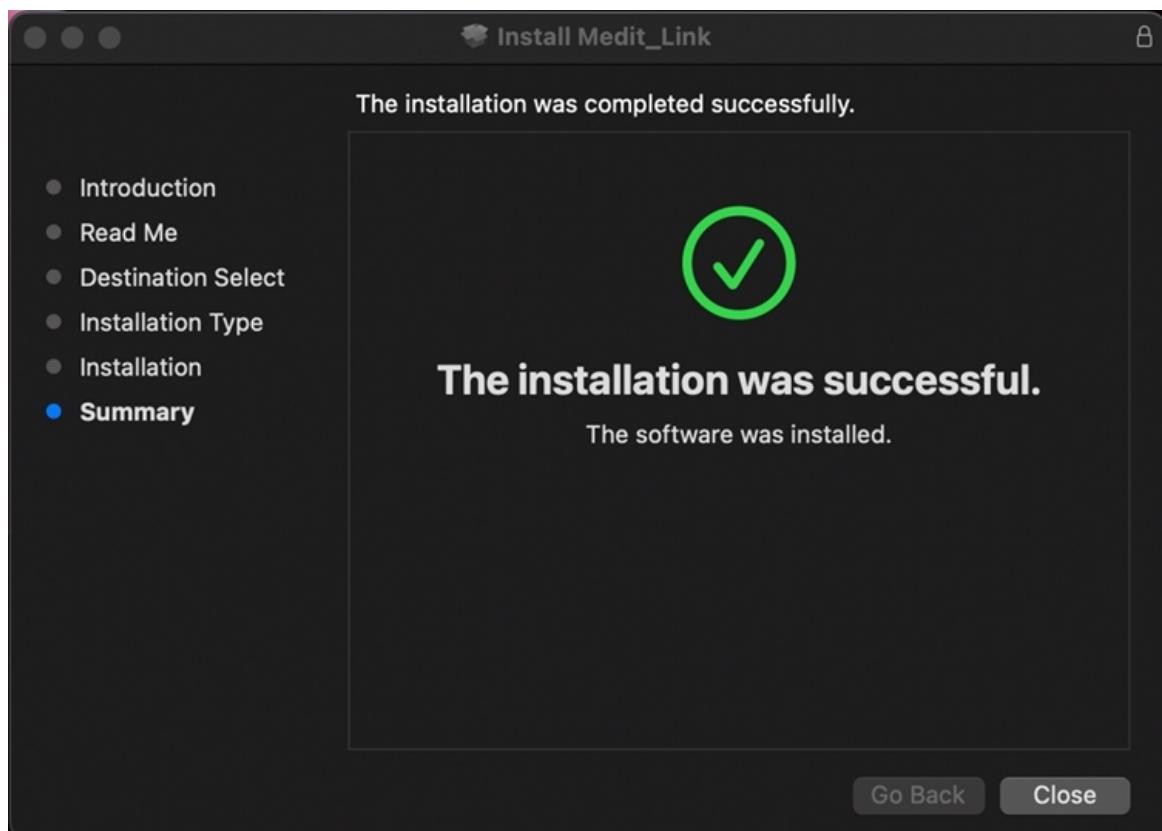

Sign Up & Sign In

Medit Link Website

Sign up for Medit Link service at www.meditlink.com.

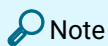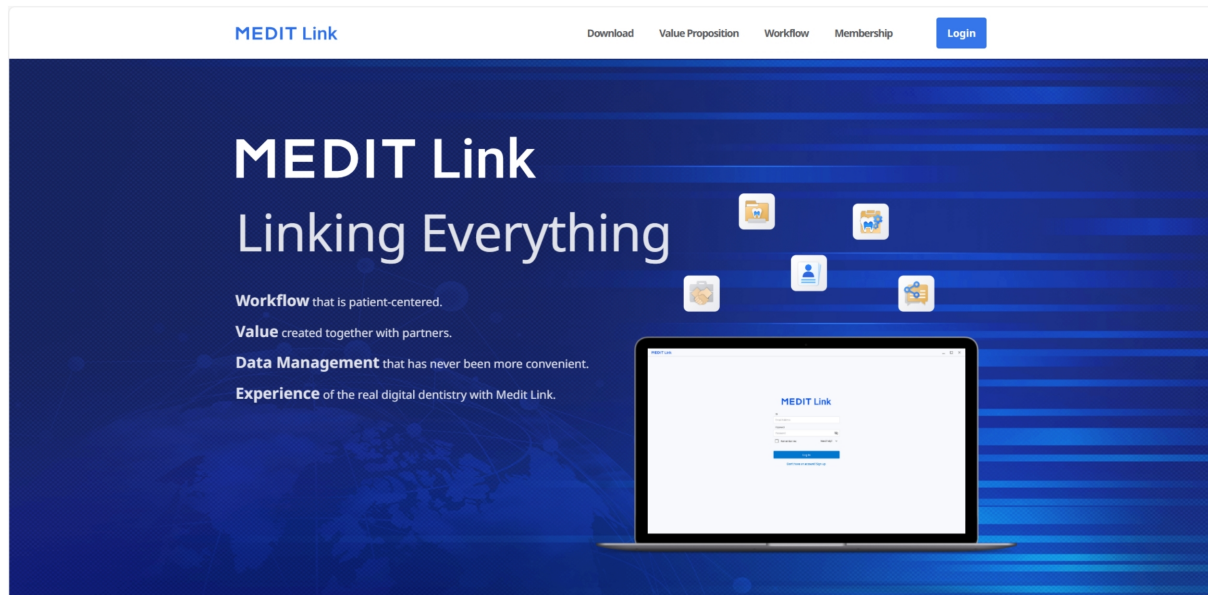

Note

Register for Medit Link China at www.meditlink.cn.

Sign Up

To use the Medit Link service, you must sign up for an account. You can utilize features provided in Medit Link after creating an account as a clinic or lab.

Users can select one of the following four account types to register for membership and use the service with full or limited features and privileges, depending on the account type.

- Clinic - Admin
- Clinic - Member
- Lab - Admin
- Lab - Member

The administrator is usually a representative of your organization and will have full access to all functions and services of Medit Link. In contrast, the member account has limited access.

To sign up for membership, follow these steps:

1. Click "Don't have an account yet? Sign up." on the login page.

MEDIT Link

ID

Password

☐ Keep me logged in. [Need help? -](#)

[Log In](#)

[Don't have an account? Sign up.](#)

 English (United States) 

2. After reading the provided descriptions, click the "Register" button for the desired account type.

Create Your Medit Link Account



Clinic

Administrator

Dentist, Manager

Have the authority to create a new dental clinic account and have full access to Medit Link.

- ✓ Scan/Order/Manage orders
- ✓ Manage members

[Register](#)


Member

Dental Hygienist, Dental Staff

Belong to the dental clinic their admin created and have limited access to Medit Link after the admin's approval.

- ✓ Scan/Order/Manage orders

[Register](#)

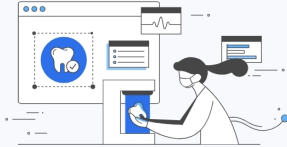

Lab

Administrator

Laboratory Director, Manager

Have the authority to create a new dental lab account and have full access to Medit Link.

- ✓ Scan/Manage ordered cases
- ✓ Manage members

[Register](#)

Member

Dental Technician, Laboratory Staff

Belong to the dental lab their admin created and have limited access to Medit Link after the admin's approval.

- ✓ Scan/Manage ordered cases

[Register](#)

3. The following registration form appears.

Create Your Medit Link Account

Clinic Name

Check Availability

Name

Email Address

Check

Please ensure your email address is active. Medit Link will send a verification email.

Registering with a disposable and/or fake e-mail address can make you vulnerable to security issues and make it difficult to receive information from Medit Link.

Password

Password

Password must be 8 to 16 characters long, including a combination of at least three letters, numbers, and special characters. (Passwords are accepted in English language only.)

Confirm Password

Confirm Password

Please type in the same password as above.

Select Country/Region

Select Country/Region

☐ I agree to all of the Terms and Conditions.

☐ Personal Information Collection Agreement (Mandatory)

☐ Terms & Conditions (Mandatory)

☐ Agree to the Marketing and Promotion Agreement (optional)

Check this option to receive marketing emails, which include product updates, software news, and promotional events.

View

View

View

Back

Register

4. Enter the clinic or lab name and then, depending on your desired account type, do the following:

- If you are registering as an administrator, click "Check Availability" to see if the entered name is available.
- If you are registering as a member, click "Search." The administrator must already register your organization before you sign up.

5. Fill out the rest of the registration form and click "Register."

 Note

If you agree to receive marketing information during registration, you will receive helpful information about the product and service from Medit.


6. A verification email will be sent to you.

26

7. To complete the registration process, verify your email address using the received link.


Login and Welcome Page

1. Enter your email address and password on the login page and click "Log In."

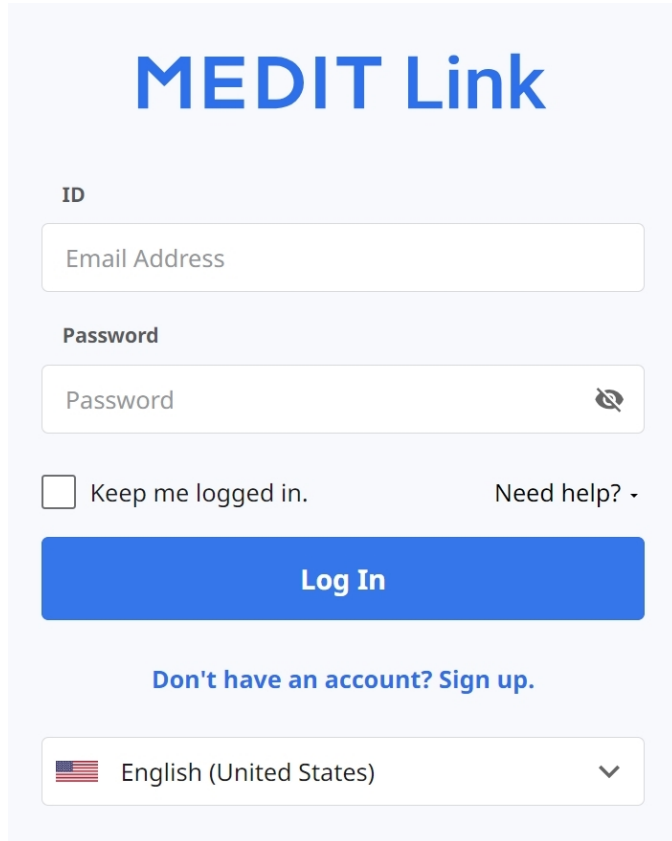The image shows a login form for 'MEDIT Link'. At the top, the text 'MEDIT Link' is displayed in a large, bold, blue font. Below this, there are two input fields: one labeled 'ID' with the placeholder text 'Email Address', and another labeled 'Password' with the placeholder text 'Password'. To the right of the password field is a small icon of an eye with a slash through it, indicating a toggle for password visibility. Below the input fields, there is a checkbox labeled 'Keep me logged in.' and a link 'Need help?' followed by a downward arrow. A large blue button with the text 'Log In' in white is positioned below these elements. Underneath the button, there is a link that says 'Don't have an account? Sign up.' At the bottom of the form, there is a dropdown menu showing a US flag icon and the text 'English (United States)' with a downward arrow.

2. The following pop-up will appear if you haven't filled out your business information.

Welcome, stage-test-wh!

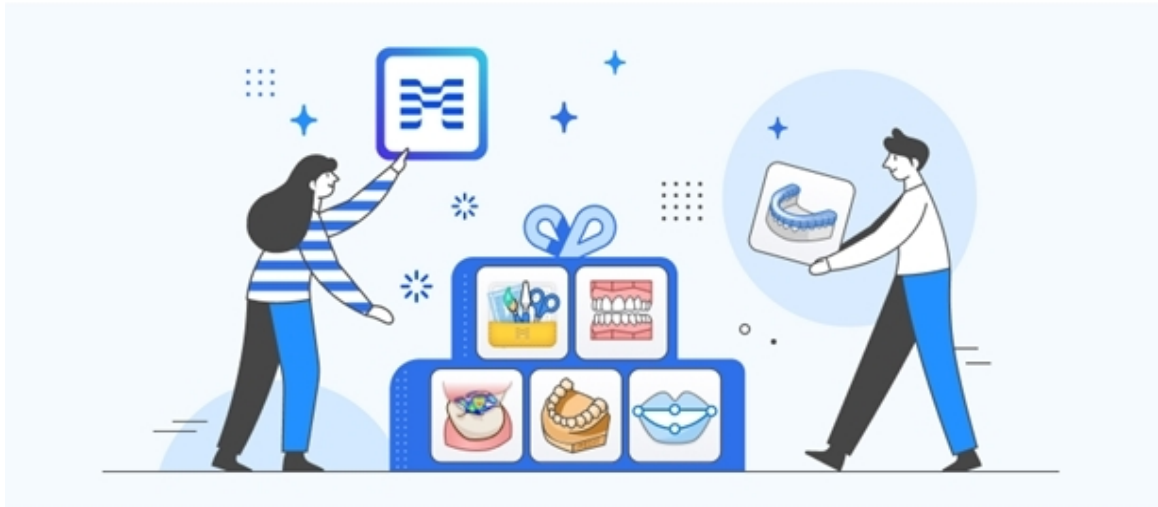

Utilize Various Apps

Download Medit Link to your desktop to utilize scanning programs and various apps.
All data generated is synced to the cloud and is also available on the web.

[Previous](#)

[Go to My Clinic](#)

Note


To use Medit Link services, your clinic or lab information should be entered using the administrator account.

3. Click the "Go to My Clinic" or "Go to My Lab" button and enter the clinic or lab information, including business name, cover/logo images, services you offer, office address, office number, etc.
4. The information entered on the welcome page can be edited anytime on My Info > My Clinic or My Lab tab in Medit Link Web.

Dashboard

Dashboard Overview

You can view the statistics for all the work done on Medit Link in visualized graphs and charts - both in the Medit Link App and Web.

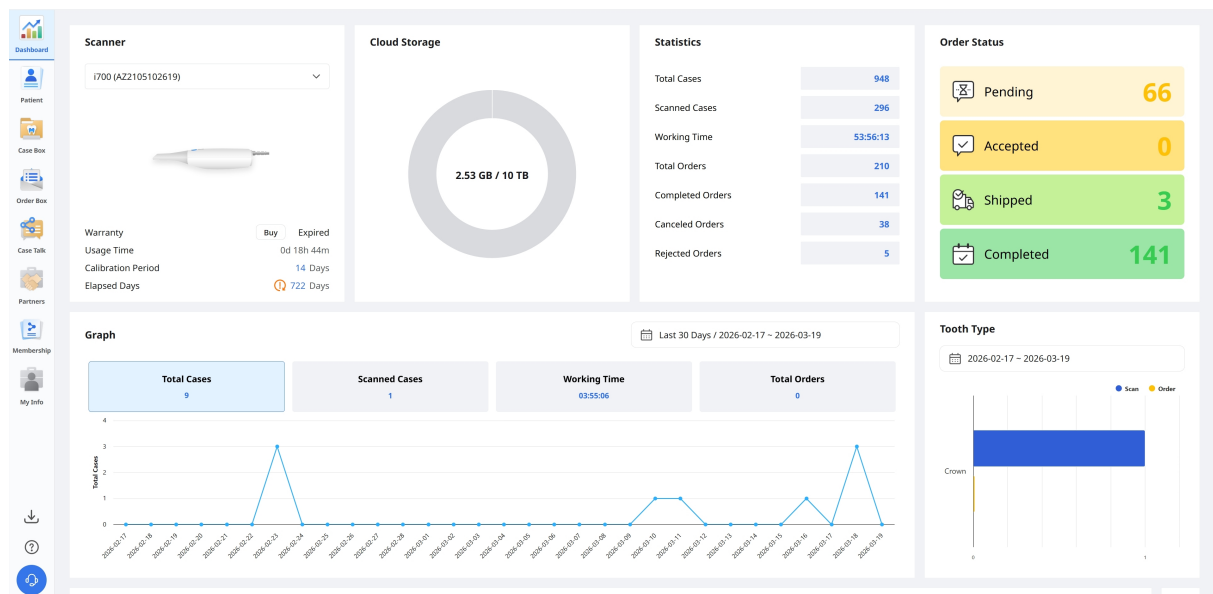


Graphs and charts in the dashboard show statistics for all the work done in Medit Link.

Scanner	Provide information about the active scanners, including warranty, scanner usage time, calibration interval, and calibration elapsed time.
Cloud Storage	Display the total capacity of the cloud and the amount in use.
Statistics	Provide statistics on cases based on job status.
Order Status	Provide statistics on cases based on order status.
Graph	Graph job information by period.
Tooth Type	Provide statistics on the number of scans and orders by treatment information.

Medit Insights

Medit Insights is a reporting service that helps you analyze your usage statistics.

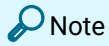

The following service is not available on Medit Link China.

- Medit Insights



Subscribe Medit Insights

You can receive a report email on the second day of every month. No additional sign-up process is required. Just enable "Subscribe Medit Insights."

If you don't use the email address you used for the Medit Link account, you can set another email address to receive Medit Insights.

Settings

Receive the report via

☒ Medit Link account


☐ Another email address

The report will be sent on the 2nd day of every month.

Cancel

OK

On the second of every month, you will receive an email with statistics, as shown below.

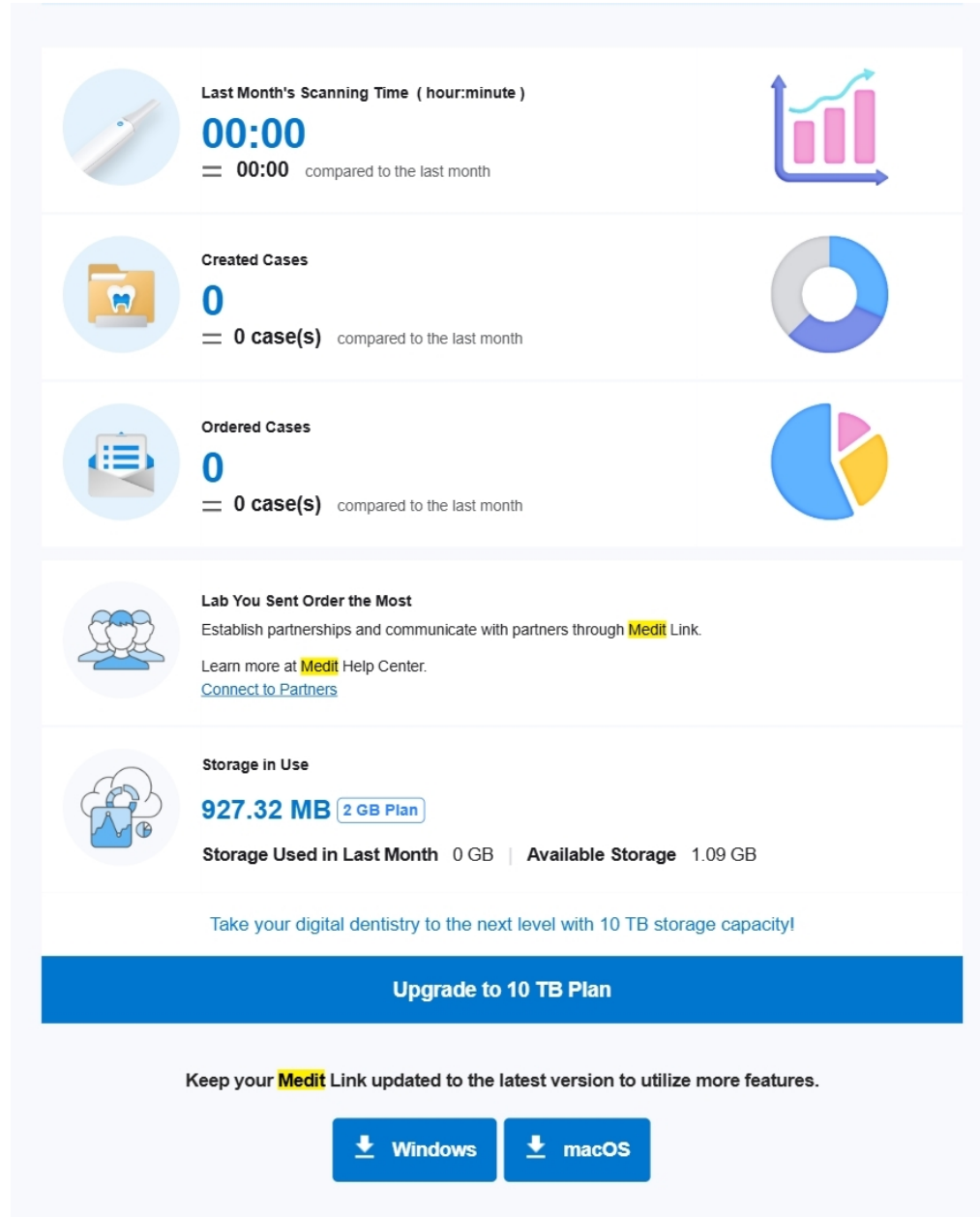

This is an automatically generated email. Please do not reply to this email.
Please contact us at insights@medit.com for more information.
[Unsubscribe](#). If you no longer wish to receive this email, click "Unsubscribe." You can change the email subscription settings in [Medit Link Dashboard](#).
Copyright © Medit Corp. All rights reserved.

Banner	Medit provides useful tips and information here. * If you don't want to see the banner, uncheck the Marketing and Promotion Agreement check box in My Info > My Account.
Statistical Information	The statistics information is displayed, such as scanning time, created cases, most ordered lab, etc.

Medit Link Download	The button provides the download link of the Medit Link.
Dashboard	All statistical information is based on your Dashboard information. You can see more details in the Dashboard.


Unsubscribe Medit Insights

Disable "Subscribe Medit Insights" in Dashboard or click "Unsubscribe" in the received email.

Patient

You can manage cases based on patient name and patient ID.

All patients registered through the Medit Link application are shown in a list, making it easier to manage patients and their cases.

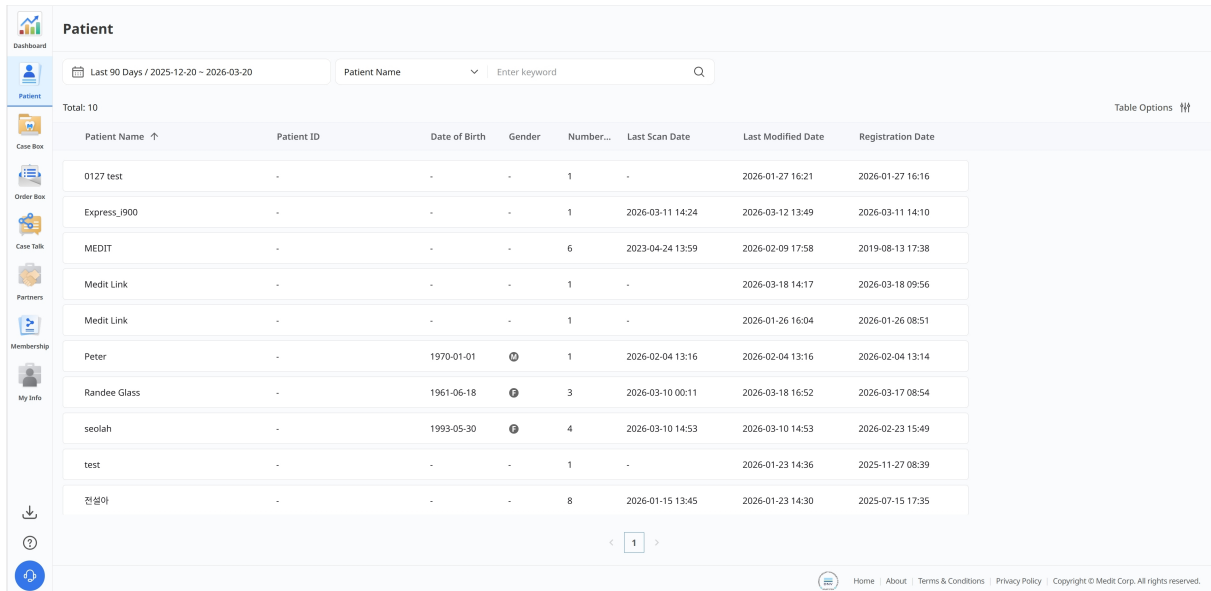


The screenshot shows the 'Patient' management interface. At the top, there's a search bar with a dropdown menu set to 'Patient Name' and a text input field 'Enter keyword'. Below the search bar, a table lists 10 patients. The table has columns: Patient Name, Patient ID, Date of Birth, Gender, Number..., Last Scan Date, Last Modified Date, and Registration Date. The patients listed are: 0127 test, Express_900, MEDIT, Medit Link, Medit Link, Peter, Randee Glass, seolah, test, and 전샘아. The interface also includes a sidebar with navigation icons and a footer with links to Home, About, Terms & Conditions, Privacy Policy, and Copyright.

Patient Name	Patient ID	Date of Birth	Gender	Number...	Last Scan Date	Last Modified Date	Registration Date
0127 test	-	-	-	1	-	2026-01-27 16:21	2026-01-27 16:16
Express_900	-	-	-	1	2026-03-11 14:24	2026-03-12 13:49	2026-03-11 14:10
MEDIT	-	-	-	6	2023-04-24 13:59	2026-02-09 17:58	2019-08-13 17:38
Medit Link	-	-	-	1	-	2026-03-18 14:17	2026-03-18 09:56
Medit Link	-	-	-	1	-	2026-01-26 16:04	2026-01-26 08:51
Peter	-	1970-01-01	♂	1	2026-02-04 13:16	2026-02-04 13:16	2026-02-04 13:14
Randee Glass	-	1961-06-18	♀	3	2026-03-10 00:11	2026-03-18 16:52	2026-03-17 08:54
seolah	-	1993-05-30	♀	4	2026-03-10 14:53	2026-03-10 14:53	2026-02-23 15:49
test	-	-	-	1	-	2026-01-23 14:36	2025-11-27 08:39
전샘아	-	-	-	8	2026-01-15 13:45	2026-01-23 14:30	2025-07-15 17:35

Search Patient

1. Choose Patient Name or Patient ID from the drop-down menu at the top center of the list.

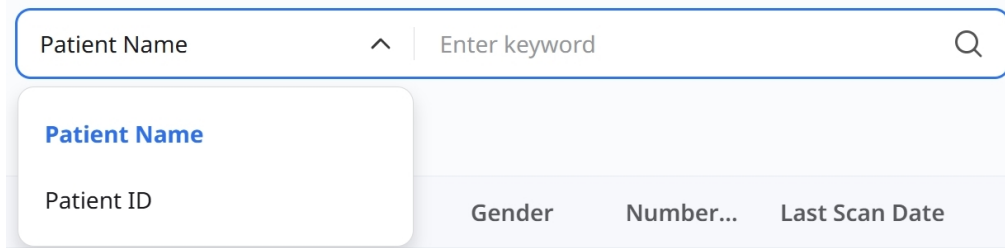

The screenshot shows the search patient interface. At the top, there's a search bar with a dropdown menu set to 'Patient Name' and a text input field 'Enter keyword'. Below the search bar, a dropdown menu is open, showing two options: 'Patient Name' and 'Patient ID'. The 'Patient Name' option is highlighted. The interface also includes a sidebar with navigation icons and a footer with links to Home, About, Terms & Conditions, Privacy Policy, and Copyright.

2. Enter the patient name or patient ID to search according to the selected search criteria.
3. Only patients with the entered patient name or ID are displayed in the list.

Case Box/Work Box

The Case Box (Clinic) or Work Box (Lab) shows all cases the user has created on the Medit Link application. You can see files via a web browser without installing a separate program.

The screenshot shows the Medit Link Case Box interface. At the top, there's a header with the Medit Link logo, a user profile (seolahjeon), and a language dropdown (English). Below the header, there's a sidebar with navigation icons for Dashboard, Patient, Case Box, Order Box, Case Talk, Partners, Membership, and My Info. The main area is titled 'Case Box' and features a dropdown menu set to 'Patient'. Below this, there's a date range filter ('Last 90 Days / 2025-12-20 ~ 2026-03-20') and a search bar ('Case Name' with 'Enter keyword'). A row of filter buttons includes 'All' (selected), 'Form Only', 'Need Processing', 'Scan Completed', 'CAD', 'Ordered', and 'Case Completed'. Below the filters, it says 'Total: 18' and 'Table Options'. A table lists cases with columns: Status, Case Name, Patient Name, Form Information, Last Modified..., Scanned Date, Order Date, Lab Name, and Assignee. The table contains 10 rows of case data. At the bottom, there's a pagination bar showing '1' and '2'.

Status	Case Name	Patient Name	Form Information	Last Modified...	Scanned Date	Order Date	Lab Name	Assignee
Scan Completed	Ortho Scan - Clone - Clo...	Randee Glass		2026-03-18 16:52	2026-03-10 00:11	-	-	Share
Scan Completed	Ortho Scan - Clone	Randee Glass		2026-03-18 14:19	2026-03-10 00:11	-	-	Share
Form Only	Medit Link의 케이스	Medit Link		2026-03-18 14:17	-	-	-	Share
Scan Completed	Ortho Scan	Randee Glass		2026-03-17 08:57	2026-03-10 00:11	-	-	Share
Need Processing	Express_i900의 케이스	Express_i900	23	2026-03-12 13:49	2026-03-11 14:24	-	-	Share
Form Only	Express Scan	seolah		2026-03-10 14:53	2026-03-10 14:53	-	-	Share
Form Only	Express Scan	seolah		2026-02-24 08:33	2026-02-24 08:33	-	-	Share
Need Processing	Express Scan	seolah		2026-02-23 16:38	2026-02-23 16:37	-	-	Share
Scan Completed	Express Scan	seolah	Maxilla	2026-02-23 15:50	2026-02-23 15:50	-	-	Share

Case Filtering

Filter by patient

You can filter cases by patients. It is useful when you find cases for specific patients.

The screenshot shows the Medit Link Case Box interface. The 'Case Box' dropdown menu is open, showing 'Patient' as the selected filter. Below the dropdown, there's a row of filter buttons: 'All' (selected), 'Form Only', 'Need Processing', 'Scan Completed', 'CAD', 'Ordered', and 'Case Completed'.

1. Enter a patient name in the search field.

Patient ^

test

☐ Occlusion Test / 0001-01-01 / - / -
☐ test_Initial / 0001-01-01 / - / -
☐ Test Occlusion / 0001-01-01 / - / -
☐ test / - / - / -
☐ Occlusion Test / 0001-01-01 / - / -

View More

2. Select a patient from the list of search results.

Occlusion Test / 0001-01-01 / - / - + 2 ^

test

☐ test / - / - / -
☐ Occlusion Test / 0001-01-01 / - / -
☐ Occlusion Test / 0001-01-01 / - / -
☐ test_Initial / 0001-01-01 / - / -
☐ Occlusion Test / 0001-01-01 / - / -


View More

Selected Patients

Clear Selection


☒ Occlusion Test / 0001-01-01 / - / -
☒ Test Occlusion / 0001-01-01 / - / -
☒ test / - / - / -

3. Only cases for the selected patients are shown in the list, making it easier to find the cases you want to work with.

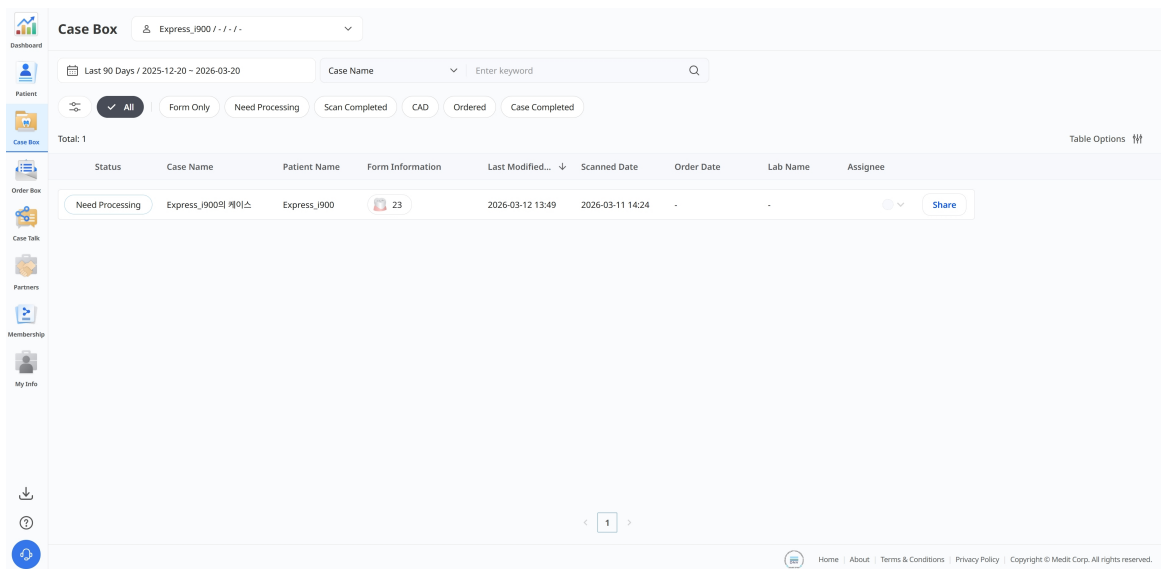

Filter by status

You can filter cases based on their status.

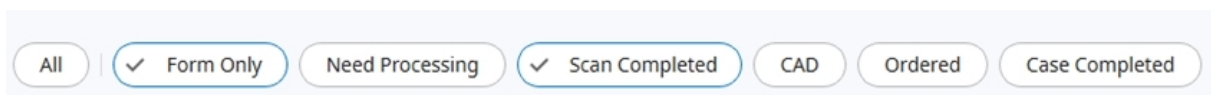

- The user can click the status they want to view. The browser will display only cases with the chosen status.
- The user can select multiple statuses at once.

The following status types are provided based on account type:

Case Box (Clinic)	All / Form Only / Need Processing / Scan Completed / CAD / Ordered / Case Completed
Work Box (Lab)	All / Form Only / Need Processing / Scan Completed / CAD / Done

Filter by duration

You can filter and display cases that were last modified within your set period.

1. Click the duration filter.



2. Select the range.

Date Range

All

Today

Yesterday

Last 7 Days

Last 30 Days

Last 90 Days

Current Month

Previous Month

Custom Range

- All: Show all files created since January 1, 2018, regardless of the selected date range.
- Custom Range: Select a start date (no earlier than Jan 1, 2018) and an end date to filter files.
- Today: Files created or modified on the current day only.
- Yesterday: Files created or modified one day before today.
- Last 7 Days: Files created or modified in the past 7 days, including today.
- Last 30 Days: Files created or modified in the past 30 days, including today.
- Last 90 Days: Files created or modified in the past 90 days, including today.
- Current Month: Files created or modified within the current calendar month.
- Previous Month: Files created or modified within the previous calendar month.

The case list will be updated in real time with the cases created or last modified within the selected duration.

Filter by Search

You can leverage search capabilities to search for the required information easily and quickly.

Case Name ▾

|


Case Name / Contains

Q

The following search options are provided based on account type:

Case Box (Clinic)	Case Name, Patient Name, Patient ID, Lab Name
Work Box (Lab)	Case Name, Patient Name, Clinic Name, Order ID

Case Details

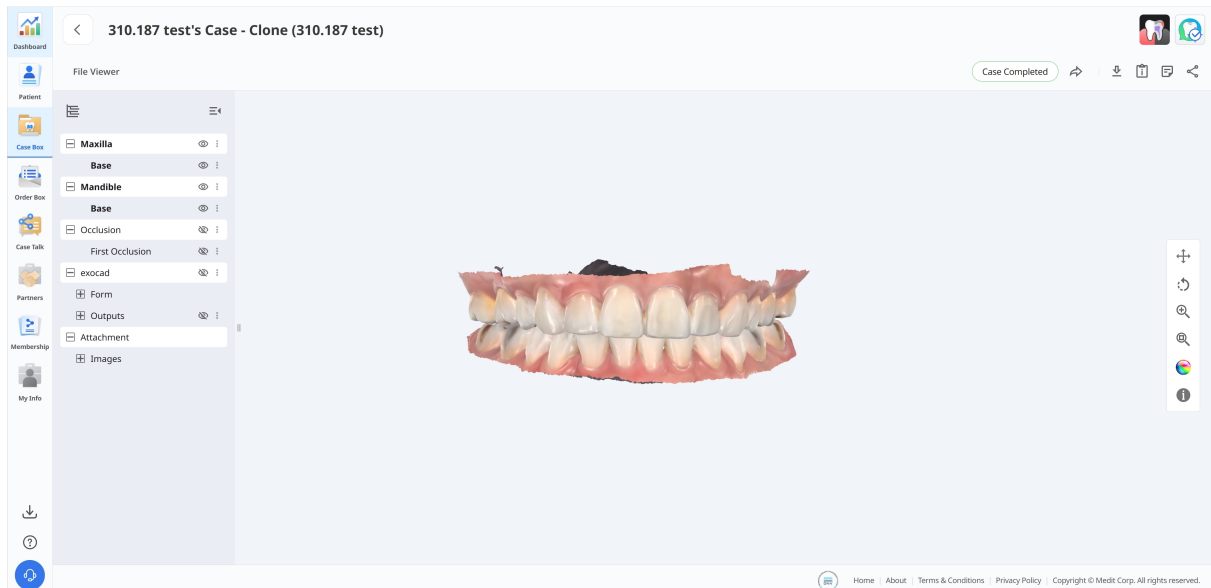


The following case management tools are provided in Case Box/Work Box.

	Case Status	Show the status of the case.
	Go to Order Box/In Box	Provide a shortcut for the same case in another box (Case Box ↔ Order Box; In Box ↔ Work Box).
	Download	Download files in the case to the local PC.
	Case Information	Review case details and save them as a PDF or print them out.
	Memo & Tags	Add a memo or tags when registering or ordering a case. Use them to communicate any additional or supplementary information about the case to your partners.
	Share	Create a link for sharing. Users with the shared link can access the 3D data in a browser such as Chrome, Firefox, or Safari without installing additional programs.

Order Box/In Box

The Order Box (Clinic) or In Box (Lab) shows ordered cases on the Medit Link application.

Order Details

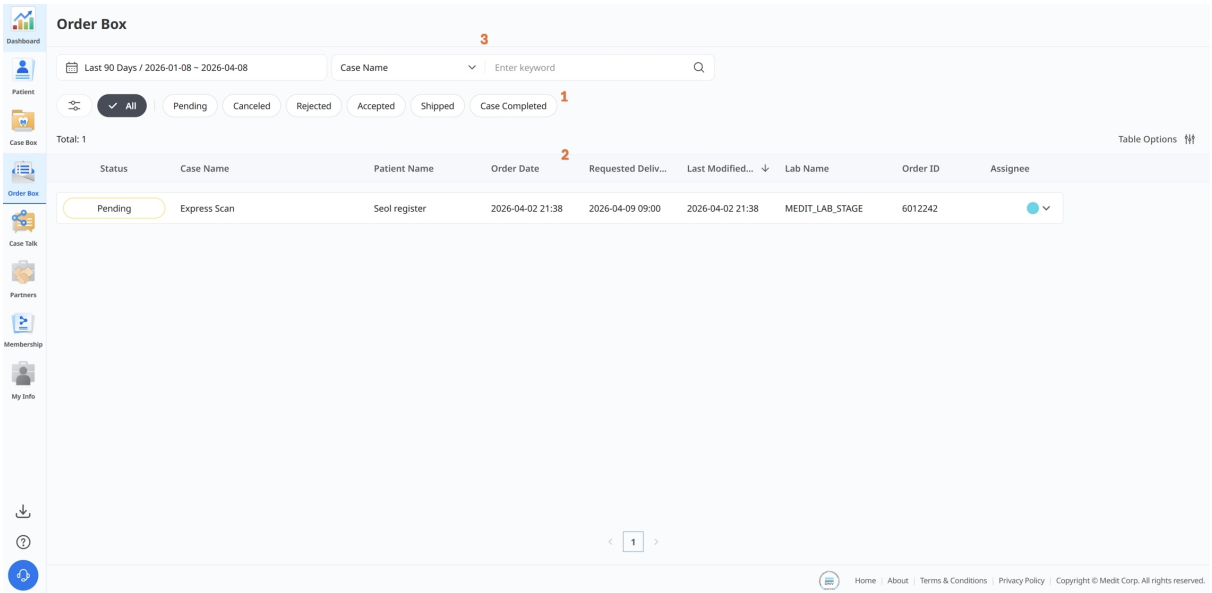


1	Filter cases by status such as Pending, Accepted, Shipped, Case Completed, Rejected, or Canceled.
2	See key details at a glance, including patient name, order date, delivery request date, lab name, and assigned technician.
3	Narrow down the list by case name or a specific date range.

Case Talk

Cases created in the Medit Link can be shared with specific people you specify, anyone with the shared link or registered members of Medit Link.

The user can see a list of all shared cases through "Case Talk" in Medit Link App and Medit Link Web.

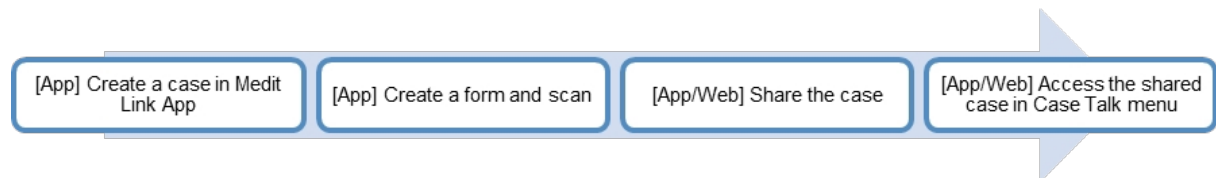


Shared cases can be viewed via a web browser without installing a separate program, enabling easy communication over the 3D data on various devices.


Creating a Case Talk

You can create a shared link to a case that was created in Case Box/Work Box.

1. Find a case you want to share and click the "Share" button on the right.

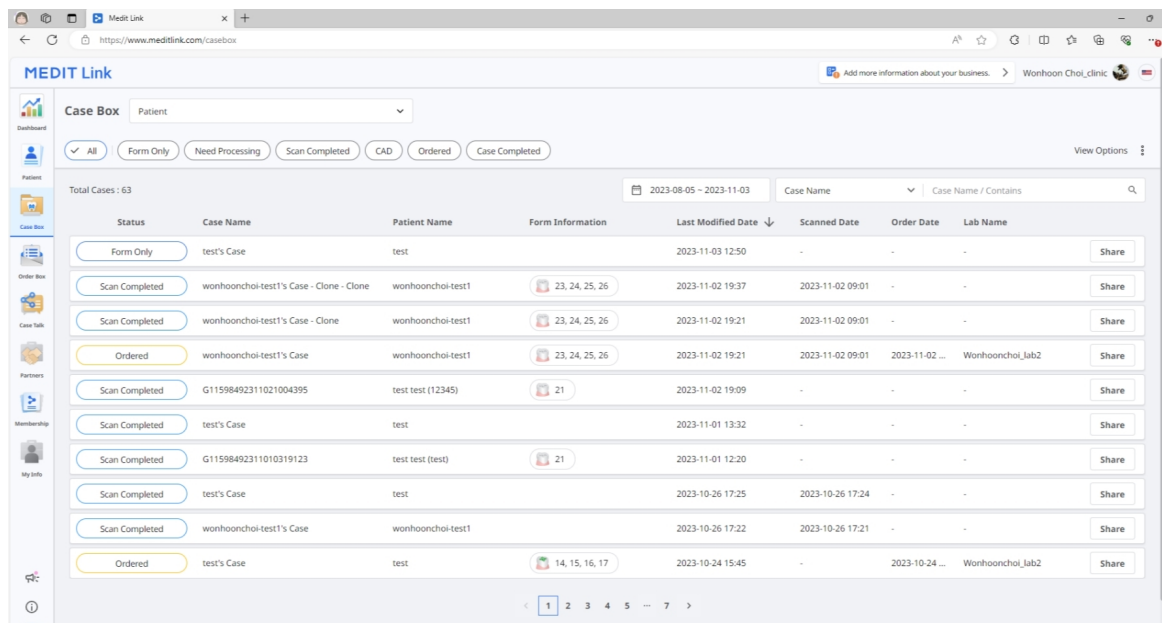

2. Select the files you want to share.

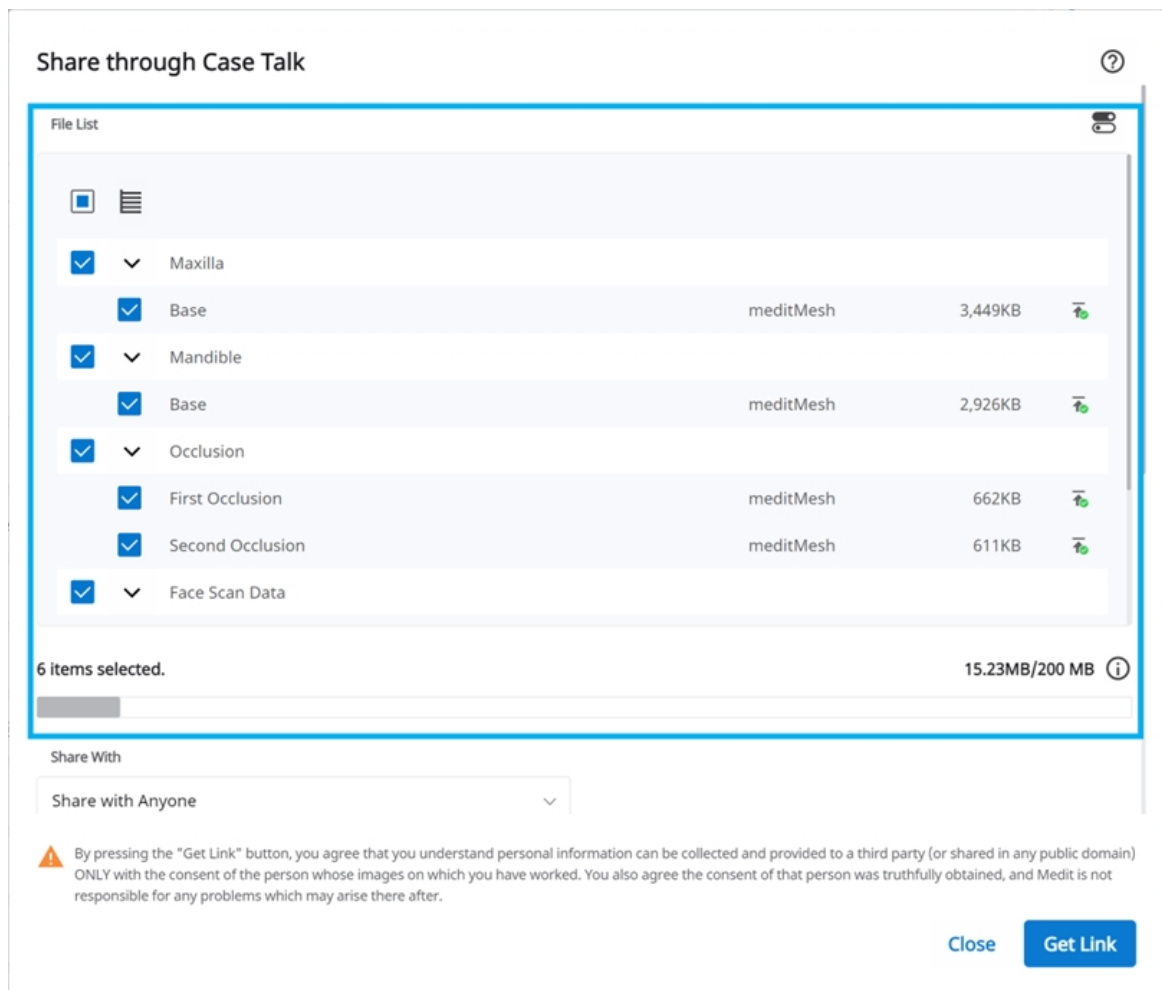

3. Set the sharing permission type and sharing options.
4. Click "Get Link" to create a shared link.

Setting the sharing permission type and sharing options

Share through Case Talk

?

Share With

Share with specific people

Create a shared link and share it with specific people.

Expiration Date

2022-06-28

You can extend the period up to 90 days from today.

Share via

PIN Code

PIN Code

Enter PIN Code

Set the PIN code using six digits for privacy protection.

✓

Share the patient's name

✓

Grant Download Permission Premium Feature

If you check this option, all users who have access to the shared link can download the files.

The premium feature is provided to all users for a limited time. Feature availability is subject to change depending on the membership plan.

✓

Allow comments

If you check this option, you will be able to communicate through the shared link.

Comment

Please enter your comment.

0/1024

By pressing the "Get Link" button, you agree that you understand personal information can be collected and provided to a third party (or shared in any public domain) ONLY with the consent of the person whose images on which you have worked. You also agree the consent of that person was truthfully obtained, and Medit is not responsible for any problems which may arise there after.

Close

Get Link


Share With	- Share with anyone Creates a shared link. You can share the files with anyone. - Share with specific people Creates and sends a shared link to specific people.
Share via	If you select "Share with specific people," you can share the link with specific people in various ways. - Email Address - SMS - Partners


42

Expiration Date	Set the expiration date for the link. After the link expires, the people you shared it with will no longer have access to files. You can extend the expiration period up to 90 days from today.
PIN Code	Set a 6-digit access code to ensure data security.
Share the patient's name	If you check this option, all people with the link will see the patient's name.
Grant Download Permission	If you check this option, all users with the link can download the files.
Allow comments	• If you check this option, you will be able to communicate through the shared link. • If you uncheck the option, the recipient can only view your comments but will not be able to add their own.
Comment	You can leave a message in the comment section.

Case Talk details

You can see files in the shared case. The following functions are available on the Case Talk details page.

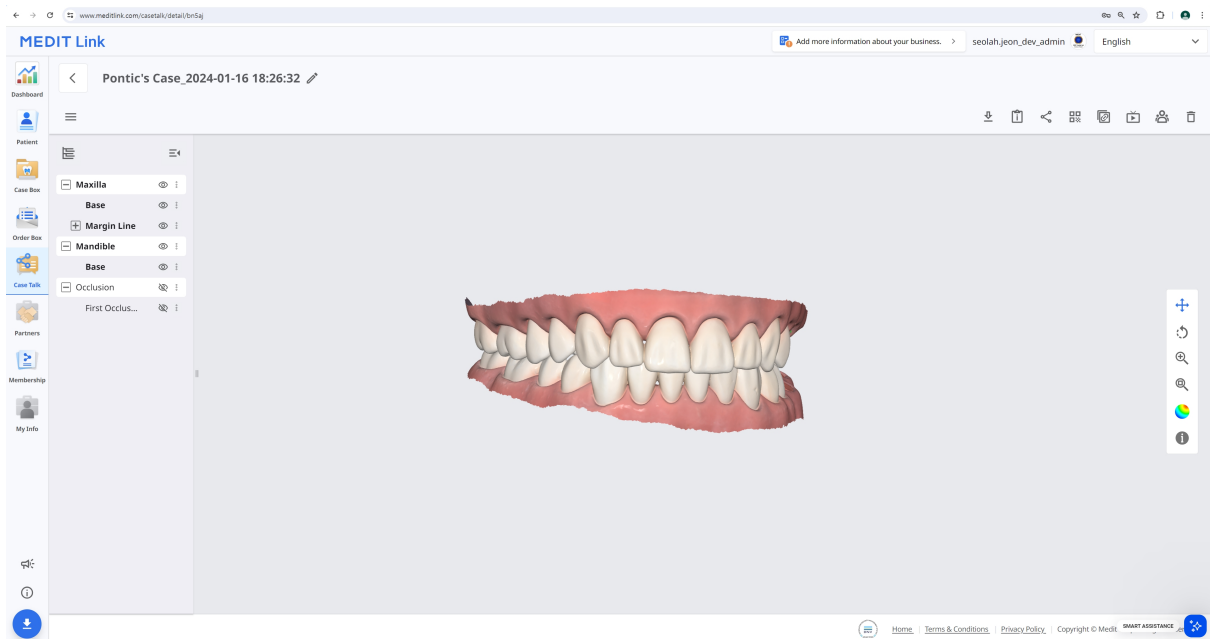

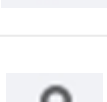	Download	Download all the files in the case to the local PC.
	Case Information	Show detailed information about the case.
	Share	Modify the share setting.
	QR Code	Share the case with a QR code.
	Copy Link	Copy the shared link to the clipboard. The copied link shows the shared files on your web browser.
	Live	Share your screen with the connected users. The screen sharing is synchronized with the connected users.
	Attendees	Show all users who are accessing the shared link.
	Delete	Delete the shared link. The deleted shared link is no longer accessible.

Note


The following services are not available on Medit Link China.

- Live
- Attendees

Web Viewer

You can view shared data through Case Talk in a web browser.

It is possible to communicate with the people who have access to the link using annotations and comments.

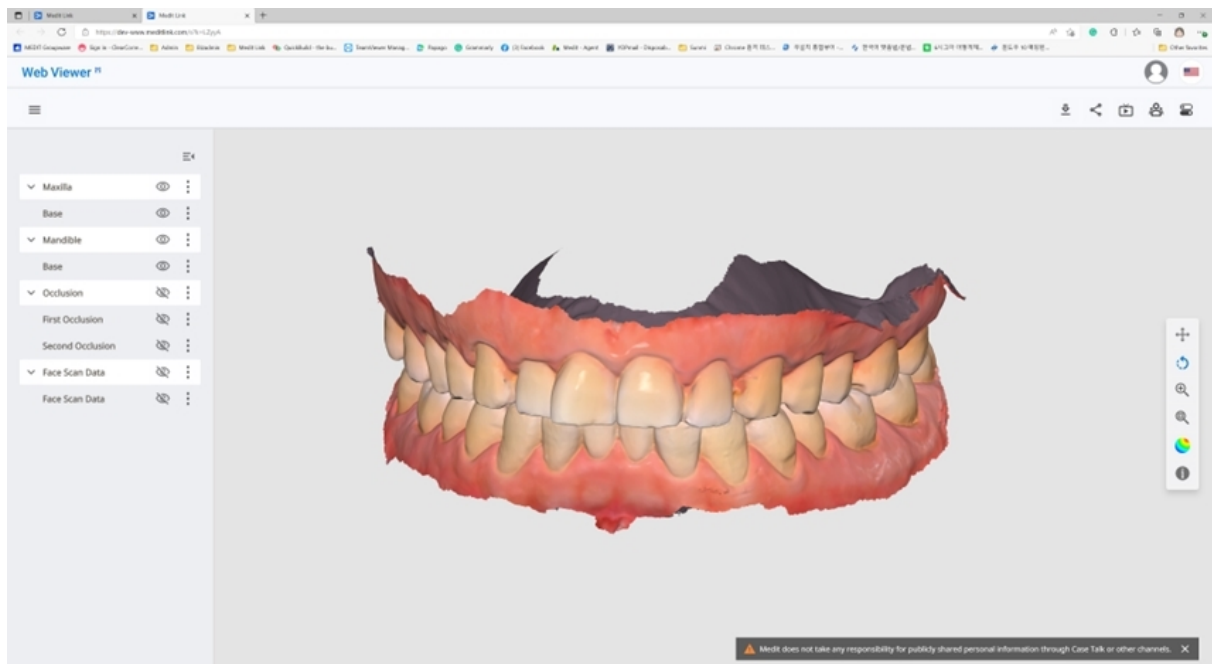

	Download	Download all the files in this case to the PC.
	Case Information	Show detailed information about the case.
	Share	Modify the share settings.

	QR Code	Share the case with a QR code.
	Live	Share your screen with the connected users. The screen sharing is synchronized with the connected users.
	Attendees	Show all users who are connected to the shared link.
	Tooth Labeling	Change the teeth labeling system.


Note

The following services are not available on Medit Link China.

- Live
- Attendees

Annotation

You can add annotations to a specified location on the 3D model. Anyone with the link can view the annotations you add.

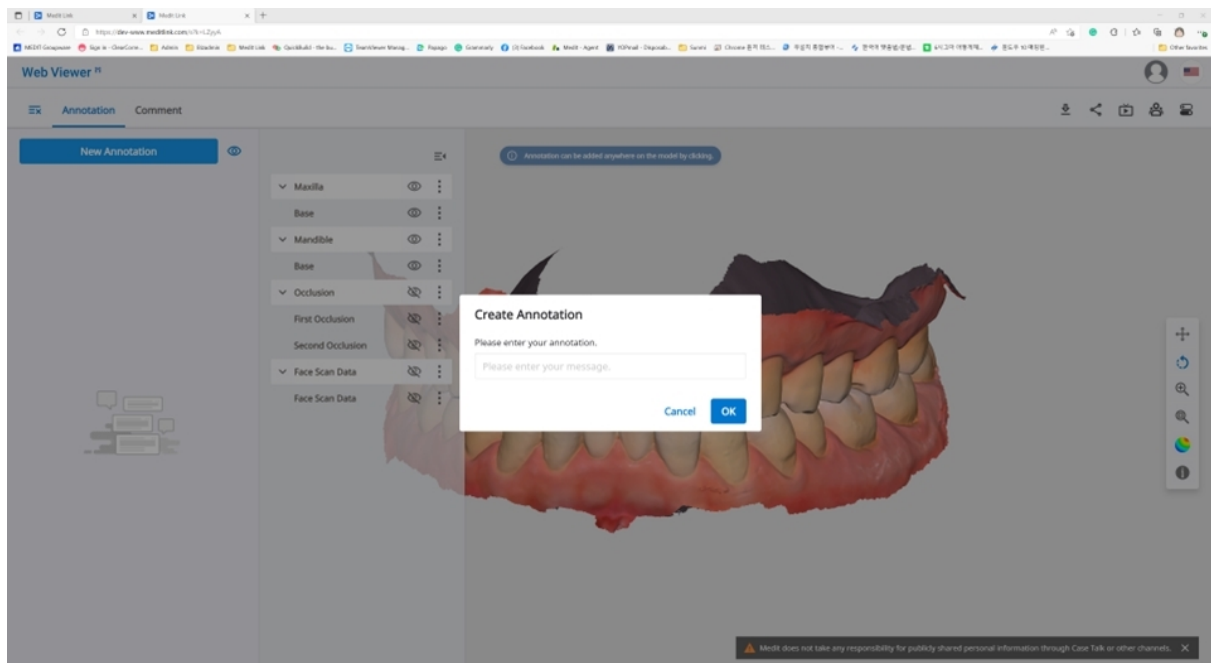


Note

The following service is not available on Medit Link China.

- Annotation

1. Find the menu icon in the top left corner and select the "Annotation" tab. Use the "New Annotation" button to add an annotation anywhere on the data with a click.
2. Type the message, then click "OK" to save it.

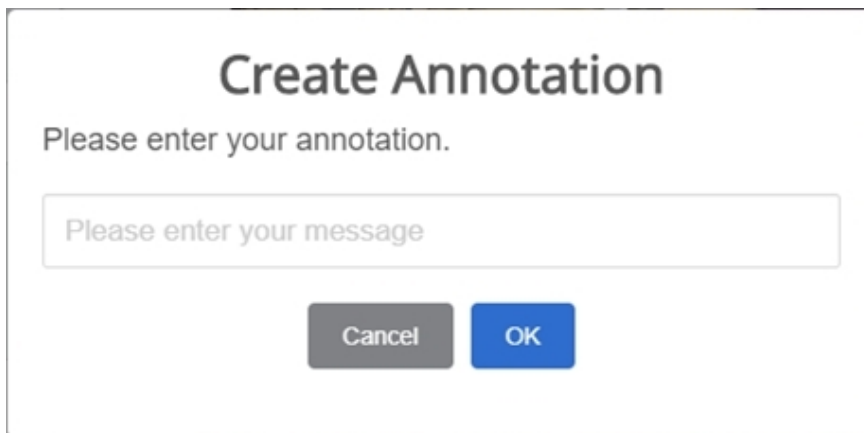A screenshot of the 'Create Annotation' dialog box. It has a title 'Create Annotation' and a prompt 'Please enter your annotation.' Below this is a text input field with the placeholder text 'Please enter your message'. At the bottom, there are two buttons: 'Cancel' and 'OK'.

3. You can use annotations as a simple communication channel.
Click the created annotation to enter a reply and hit "Send."

wonhoonchoi Beta wonhoonchoi

Annotation 1


1 09-09 20:31

Please enter your reply

Send

Comments

You can communicate with users who have access to the link through the "Comment" tab in the Web Viewer menu.

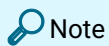


The following service is not available on Medit Link China.

- Comment

Real-time screen sharing


You can share your screen in real-time.

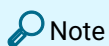

The following service is not available on Medit Link China.

- Live

1. Click "Live" to share your Web Viewer screen with the connected users.


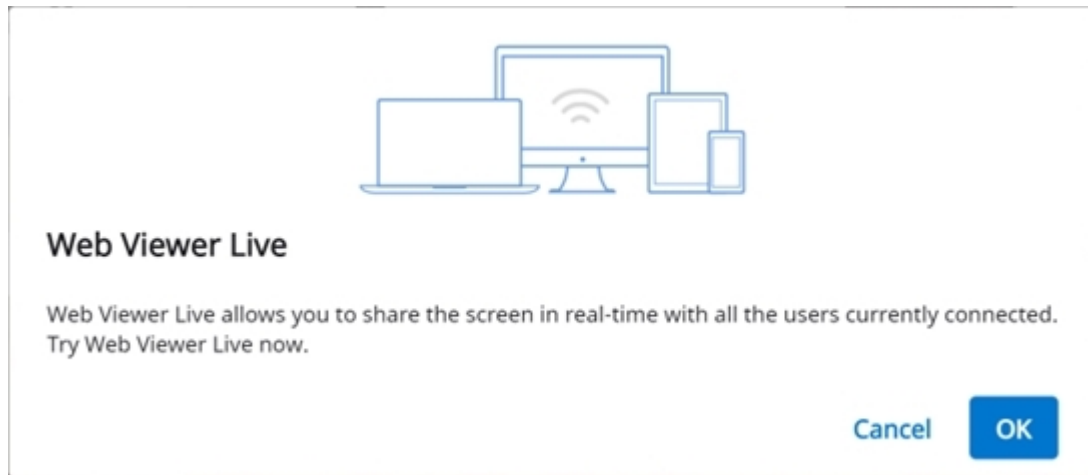


2. Click "OK" to start sharing your screen.
3. Other people who have accessed the link must accept the screen sharing.
4. The screen of the user who initiated screen sharing is displayed in sync with the screen of the person who received the request.


Data Tree

All data is grouped and displayed according to attributes in the data tree. You can show or hide data, adjust opacity, and so on.

Data Tree Controls

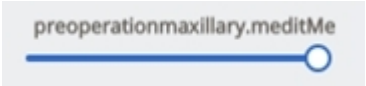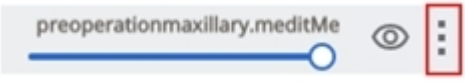Show/Hide	Click the 'eye' icon to show or hide data from the view. Control data visibility as a group or individual file.
Opacity	Hover the mouse over a visible element to see the opacity slide bar. Adjust the slider to change the transparency of your data. 
Extended Menu	The extended menu is provided on the right side of each Data Tree element. It includes the following options:  • Show This Only: Hide all other data except for the current one. • Hide: Hide the current data.

3D Model Controls

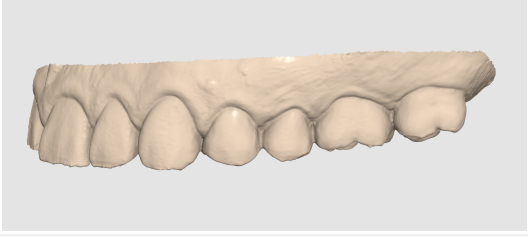	Pan		Move the model.
	Rotate		Rotate the model.
	Zoom In/Out		Zoom in and out on the model.
	Zoom Fit		Position the model in the center of the screen.
	Data Display Mode	Matte	See the data in color with matte finish (Phong rendering). 
		Monochrome	See the data in one color. 

Partners

You can search for potential partners to request a partnership and view the list of current partners.




Services using Google Maps are not available on Medit Link China.

Search for Partners


You can search for Medit Link users based on their business name and/or address.

In the search results, you can send a partnership request to a business you want to work with.

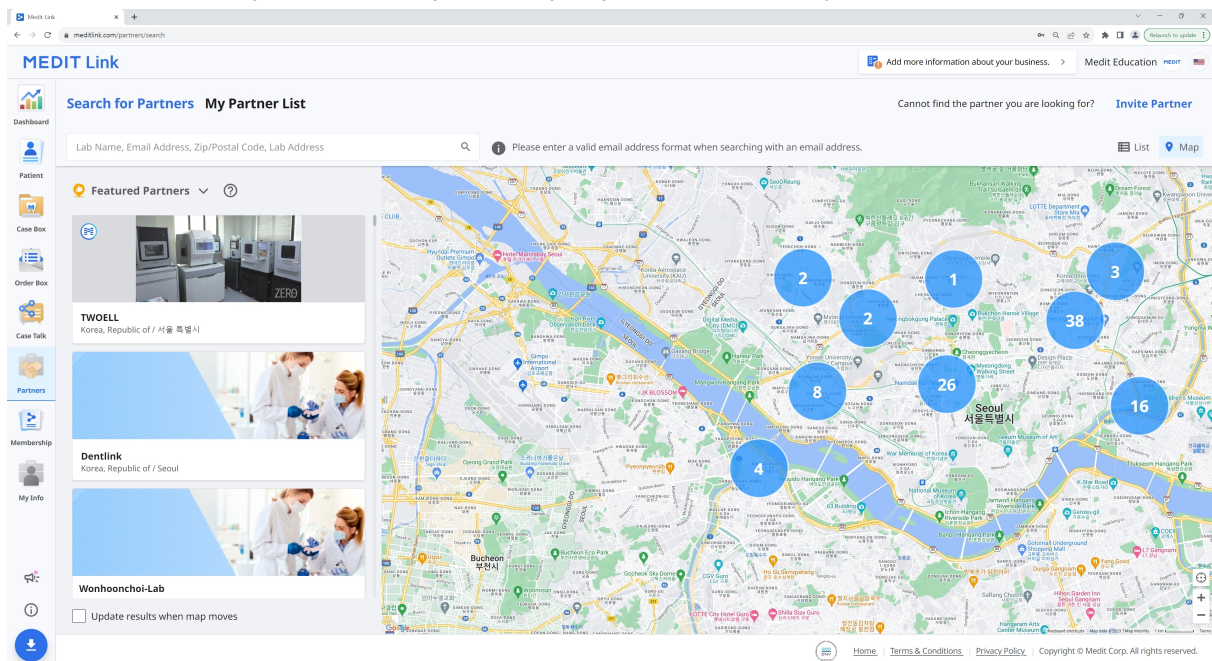

My Partner List

You can review your current partnerships with all your existing partner businesses in My Partner List.

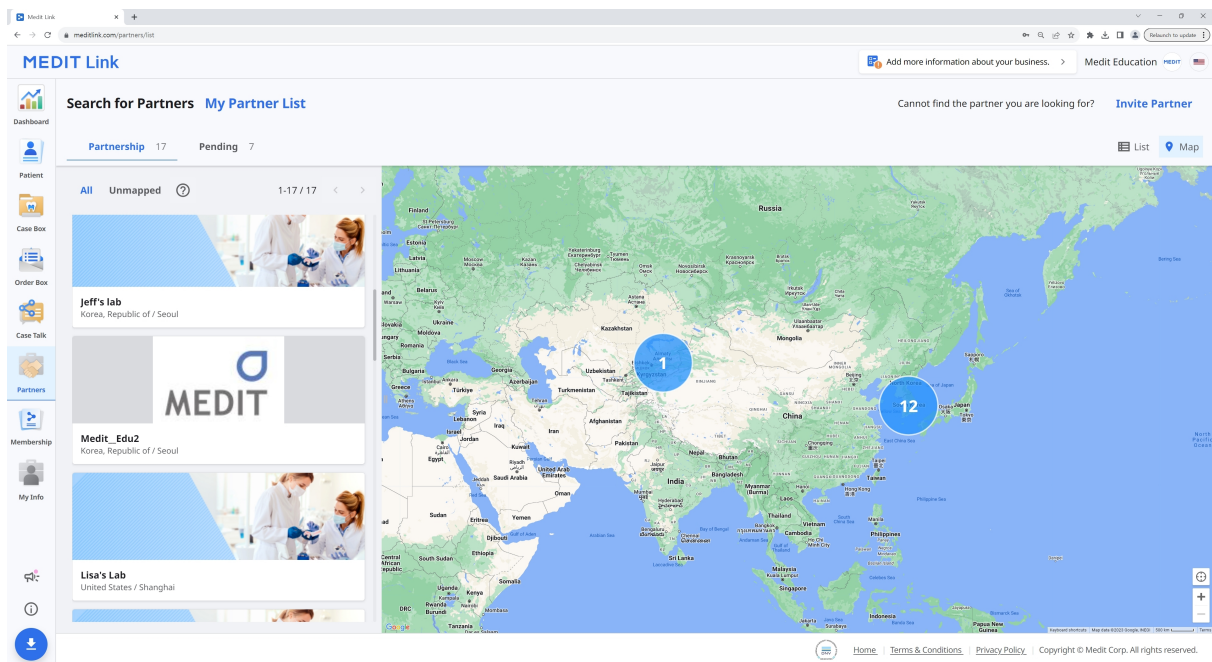

Partner Invitation

If you have a clinic or lab you want to work with, you can invite them to join Medit Link. The partnership will be automatically formed once they accept the invitation.

1. Click "Invite Partner" at the top right of the screen.
2. Enter the email address of a business you want to work with and send them an invitation with a partnership request.



Invite Partner

Invite your desired partner to join Medit Link. When the lab accepts your invitation, a partnership is automatically formed.


Cancel

Send Invitation

Membership

Manage Plan

On the **Manage Plan** tab, you can view and manage your active plans and browse products available for purchase or subscription.

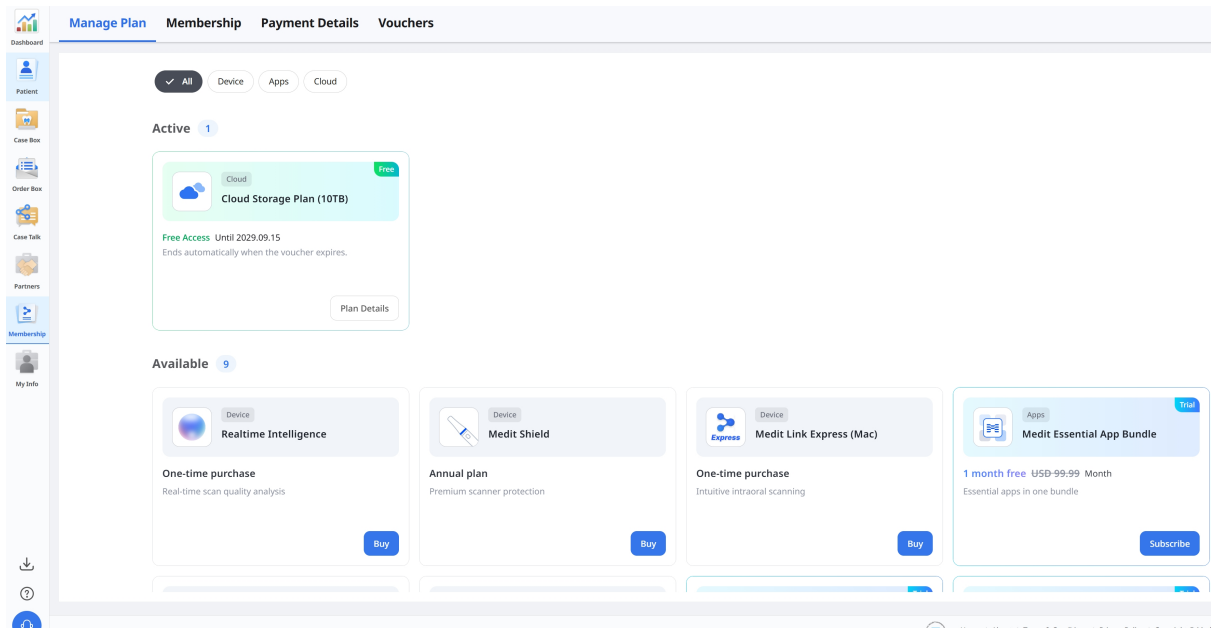

Managing an Active Plan

To view the details of an active plan, click **Plan Details** on the plan card.

[Manage Plan](#) [membership](#) [Payment Details](#) [vouchers](#)

✓ All

Device

Apps

Cloud

Active

1

Cloud

Free

 Cloud Storage Plan (10TB)

Free Access Until 2029.09.15
Ends automatically when the voucher expires.

Plan Details

Available

9

Device

Realtime Intelligence

Device

Me

Plan Details

×

Free Access

Your free access is active.

Product	Cloud Storage Plan (10TB)
Status	Free
Expiration Date	2029.09.15
Regular Price	USD 9.99/Month

⚠ Important

- You can't subscribe to this product while the voucher is active.
- Free access ends automatically when the voucher expires. A paid subscription is required to continue.

Confirm

Subscribing to the 10 TB Plan

If you are currently on the 2 GB Plan, you can subscribe to the 10 TB Plan. If you have free access through a voucher, you can subscribe after your free access ends.

1. Find the 10 TB Plan under **Available** and click **Subscribe**.

The screenshot displays a user's subscription management dashboard. At the top, there are four tabs: 'Manage Plan' (selected), 'Membership', 'Payment Details', and 'Vouchers'. Below the tabs, there are three filter buttons: 'All' (selected), 'Device', and 'Apps'. The 'Active' section shows 0 active subscriptions with a message: 'No active subscriptions. Browse our services below to get started.' The 'Available' section shows 8 available plans. A red arrow points from the 'Available' count to the 'Cloud Storage Plan (10TB)' card. The card features a 'Cloud' category tag, a 'Trial' badge, and the text '1 month free USD 9.99 Month' and 'Secure 10 TB cloud storage'. A blue 'Subscribe' button is located at the bottom right of the card.

Manage Plan Membership Payment Details Vouchers

✓ All Device Apps Cloud

Active 0

No active subscriptions.
Browse our services below to get started.

Available 8

Cloud Storage Plan (10TB)

1 month free ~~USD 9.99~~ Month

Secure 10 TB cloud storage

Subscribe

2. Review the plan details and click **Confirm**.

Subscribe



Cloud Storage Plan (10TB)

1 month free

~~USD 9.99~~ / Month

- ✓ Includes 10 TB of cloud storage
- ✓ Securely store large cases and other data
- ✓ Streamline your Medit Link data management in the cloud
- ✓ Store design files without separate local backups

- Actual charges may vary depending on your bank or card issuer's fees and policies.
- For more information about Cloud Storage Plan (10TB), see [here](#).

Trial Cloud Storage: First Month Free

1 month free >

☐ I agree to automatic renewal and recurring charges.

Cancel

Confirm

3. Enter the required payment information and complete your subscription.

The image shows a two-step subscription process. The top part is a payment summary screen with a blue header. It displays a coupon for \$0.00, valid until the coupon expires, after which the price will be \$9.99 per day. The summary includes the 'Medit Cloud Storage 10TB Plan for Clinics' at \$9.99, a subtotal of \$9.99, a 1-month free discount for new subscribers (-\$9.99), and a total due today of \$0.00. To the right is a 'Pay with card' form with fields for email (yhgam0829@gmail.com), card number (1234 1234 1234 1234), MM/YY, CVC, cardholder name, and country (South Korea). A 'Subscribe' button is at the bottom right. Below this is a 'Subscribe' modal window. It features a cloud storage icon and the message: 'Your Cloud Storage Plan (10TB) subscription is now active. Enjoy full access to all features.' It also lists payment details: Amount of Payment (Free), Next Payment Date (2026.10.23), and Last Payment Method (****). A 'Confirm' button is at the bottom right. A red arrow points from the 'Subscribe' button in the payment form to the 'Confirm' button in the modal.

Subscribe to Medit Cloud Storage 10TB Plan

\$0.00
Then \$9.99 per day after coupon expires

Medit Cloud Storage 10TB Plan for Clinics \$9.99
Medit Cloud storage plan
Billed daily

Subtotal \$9.99

1 Month Free for New Subscribers -\$9.99
100% off for a day

Tax \$0.00

Total due today \$0.00

Pay with card

Email yhgam0829@gmail.com

Card information

1234 1234 1234 1234 VISA

MM / YY CVC

Cardholder name

Full name on card

Country or region

South Korea

Subscribe

By confirming your subscription, you allow Medit USA Inc to charge you for future payments in accordance with their terms. You can always cancel your subscription.

Powered by stripe Terms Privacy

Subscribe

Your Cloud Storage Plan (10TB) subscription is now active.
Enjoy full access to all features.

Amount of Payment Free

Next Payment Date 2026.10.23

Last Payment Method ****

Confirm

Canceling a Subscription

1. Open **Plan Details** for the subscription you want to cancel and click **Unsubscribe**.

The image shows a 'Manage Plan' modal window. It has a title bar with a close button. Inside, there's a 'Promotion Applied' section stating 'Regular pricing will apply after the promotion ends.' Below this, the subscription details are listed: Product (Cloud Storage Plan (10TB)), Price (Free USD 9.99/Month), Expiration Date (2026.10.23), and Last Payment Method (****). An 'Edit' link is next to the payment method. At the bottom, there are two buttons: 'Close' and 'Unsubscribe'. A red arrow points to the 'Unsubscribe' button.

Manage Plan

Promotion Applied
Regular pricing will apply after the promotion ends.

Product Cloud Storage Plan (10TB)


Price Free USD 9.99/Month

Expiration Date 2026.10.23

Last Payment Method **** [Edit](#)

Close **Unsubscribe**

2. Review the cancellation details and click **Confirm Subscription Cancellation**.

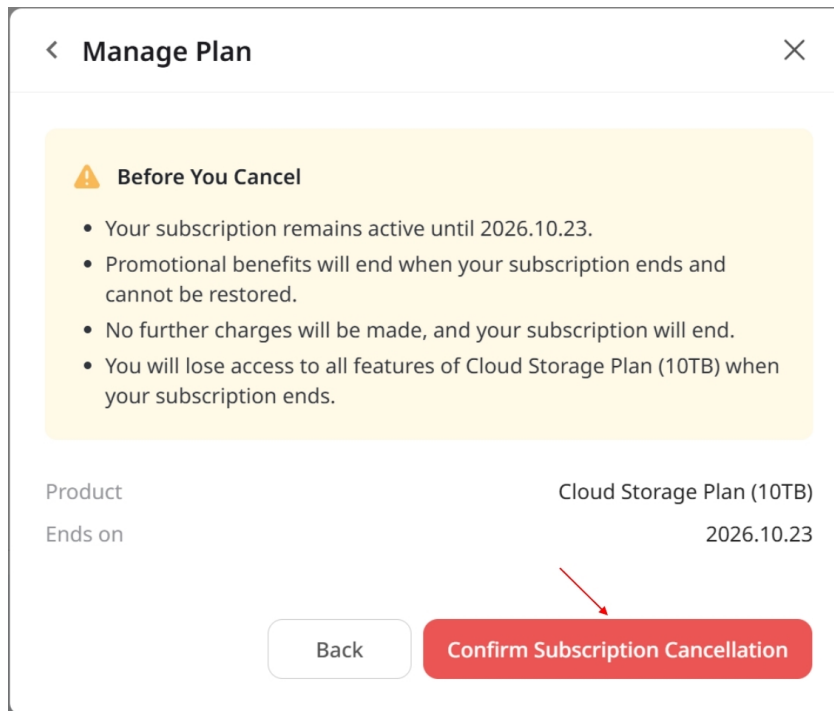

The screenshot shows a 'Manage Plan' dialog box with a close button (X) in the top right corner. Inside, there is a yellow warning box titled 'Before You Cancel' with a warning icon. It contains four bullet points: 'Your subscription remains active until 2026.10.23.', 'Promotional benefits will end when your subscription ends and cannot be restored.', 'No further charges will be made, and your subscription will end.', and 'You will lose access to all features of Cloud Storage Plan (10TB) when your subscription ends.' Below this, the 'Product' is listed as 'Cloud Storage Plan (10TB)' and 'Ends on' is '2026.10.23'. At the bottom, there are two buttons: a white 'Back' button and a red 'Confirm Subscription Cancellation' button. A red arrow points to the red button.

< **Manage Plan** ×


⚠ Before You Cancel

- Your subscription remains active until 2026.10.23.
- Promotional benefits will end when your subscription ends and cannot be restored.
- No further charges will be made, and your subscription will end.
- You will lose access to all features of Cloud Storage Plan (10TB) when your subscription ends.

Product Cloud Storage Plan (10TB)
Ends on 2026.10.23

Back Confirm Subscription Cancellation

3. When the cancellation is complete, click **Confirm**.

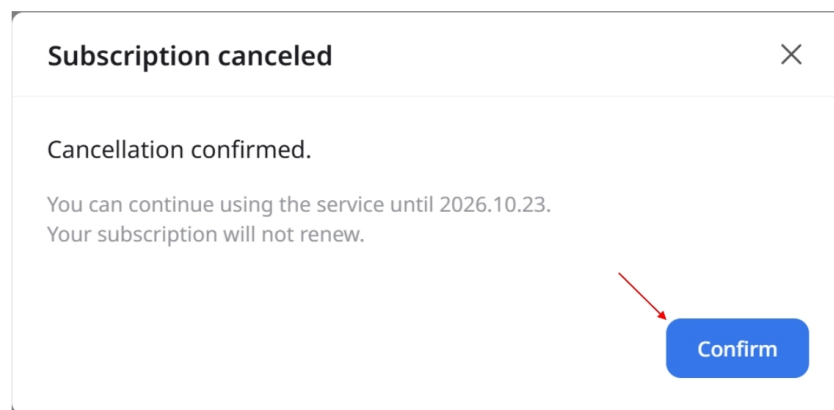

The screenshot shows a 'Subscription canceled' dialog box with a close button (X) in the top right corner. It contains the text 'Cancellation confirmed.' followed by 'You can continue using the service until 2026.10.23.' and 'Your subscription will not renew.' At the bottom right, there is a blue 'Confirm' button. A red arrow points to the blue button.


Subscription canceled ×

Cancellation confirmed.

You can continue using the service until 2026.10.23.
Your subscription will not renew.


Confirm

Membership

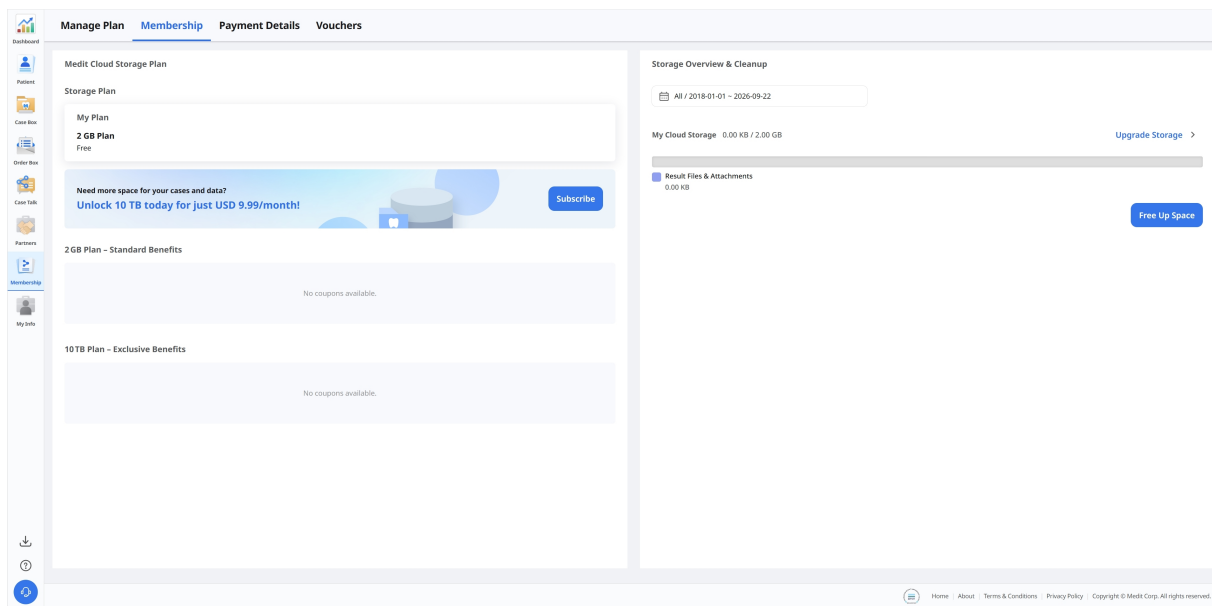

Medit Cloud Storage Plan

You can view your current cloud storage plan and the benefits available with each plan.

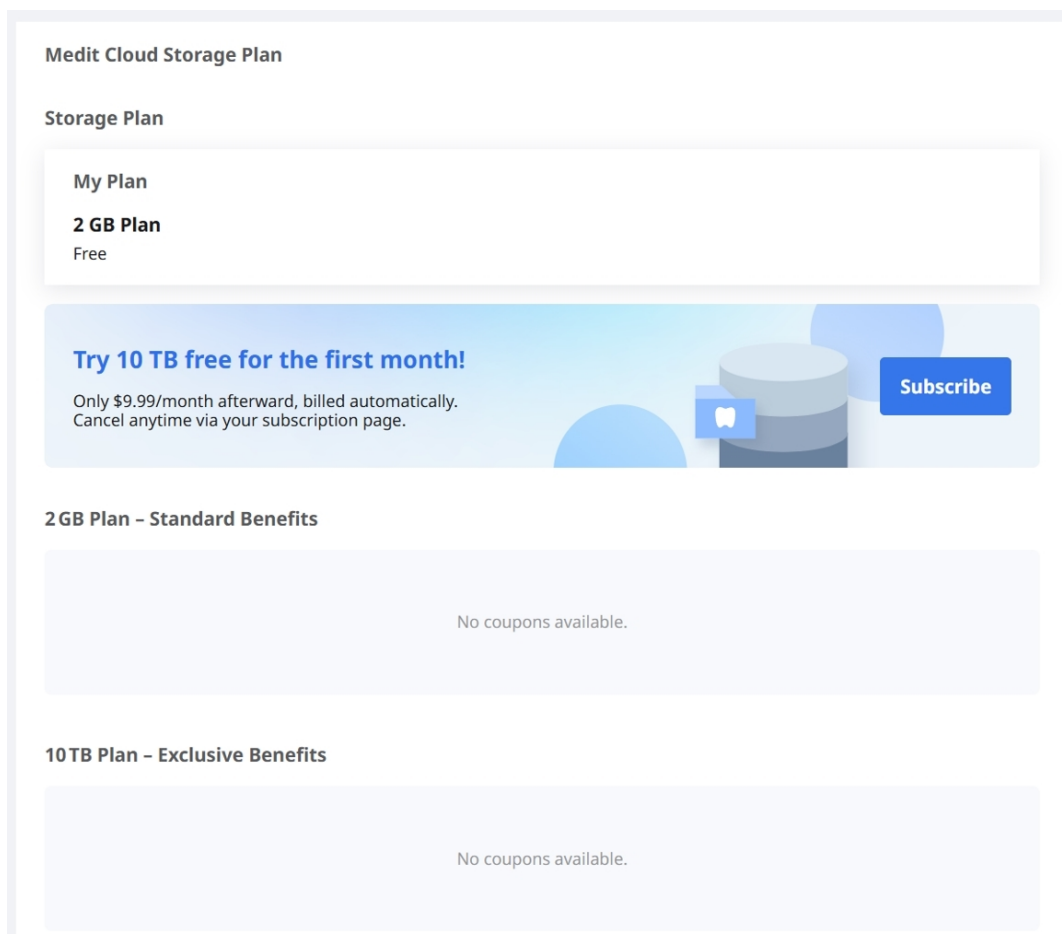

Medit Link provides the following cloud storage plans:

2 GB Plan

- Free standard plan for all signed-in users
- 2 GB cloud storage
- For new users

10 TB Plan

- Paid plan
- 10 TB cloud storage
- For users with many cases and various data
- One-month free trial

	2 GB Plan	10 TB Plan
Cloud Storage	2 GB	10 TB
Monthly Fee	Free of charge	\$9.99/month
* Cloud storage does not count raw data toward the used space.		

Note

The membership plans will be available on Medit Link China. You can enjoy the 10 TB Plan for free until the official release of the feature.

Cloud Storage

In the right-side section, you can manage your cloud storage.

Storage Overview & Cleanup

 2018-01-01 ~ 2025-11-05

My Cloud Storage 927.32 MB / 2.00 GB

[Upgrade Storage >](#)

Result Files & Attachments
927.32 MB

[Free Up Space](#)

Set Period 	Set the specific period to delete data created within that time frame. • All: Show all files created since January 1, 2018, regardless of the selected date range. • Custom Range: Select a start date (no earlier than Jan 1, 2018) and an end date to filter files. • Today: Files created or modified on the current day only. • Yesterday: Files created or modified one day before today. • Last 7 Days: Files created or modified in the past 7 days, including today. • Last 30 Days: Files created or modified in the past 30 days, including today. • Last 90 Days: Files created or modified in the past 90 days, including today. • Current Month: Files created or modified within the current calendar month. • Previous Month: Files created or modified within the previous calendar month.
Upgrade Storage	Click to upgrade to the 10 TB paid plan.
Result Files & Attachments	The result and attachment files take up space in the cloud. If the cloud is full, you can secure more storage space by removing unnecessary files.
Free Up Space	Remove selected files from the cloud storage to free up space.

Payment Details


You can view your subscription and payment history in the **Payment Details** tab.

Payment Date	Product	Payment Type	Duration	Details	Amount of Payment	Status
2026-09-21	Cloud Storage Plan (10TB)	Subscription	2026-09-21 - 2026-10-21	Monthly payment	USD 0.00	Payment Completed
2026-09-21	Cloud Storage Plan (10TB)	Subscription	2026-09-21 - 2026-10-21	Monthly payment	USD 0.00	Payment Completed
2026-09-21	Cloud Storage Plan (10TB)	Subscription	2026-09-21 - 2026-10-21	Monthly payment	USD 0.00	Payment Completed
2026-09-21	Cloud Storage Plan (10TB)	Subscription	2026-09-21 - 2026-10-21	Monthly payment	USD 0.00	Payment Completed
2026-09-21	Cloud Storage Plan (10TB)	Subscription	2026-09-21 - 2026-10-21	Monthly payment	USD 0.00	Payment Completed

Item	Description
Payment Date	Displays the date of payment.
Product	Displays the purchased plan or app.
Payment Type	Displays the type of payment: - Subscription - One-time payment
Duration	Displays the subscription period.
Details	Displays the detailed payment information.
Amount of Payment	Displays the amount for a subscription or one-time purchase.
Status	- Payment Completed: The payment has been processed. - Payment Failed: The payment was not processed. - Payment Incomplete: The payment has not been completed. Click "Pay Now" to proceed. - Waiting: The payment is pending for a subscription renewal.

View Bill Details

You can view your payment invoices and receipts by clicking the receipt icon on the right.

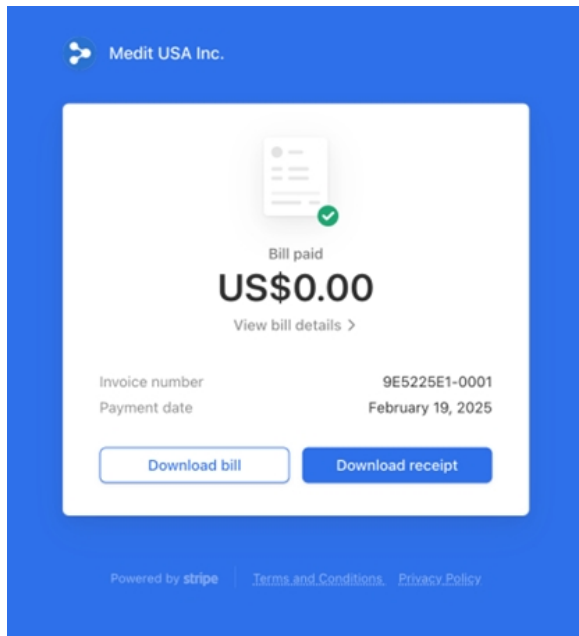

Click **View bill details** for payment details.

A screenshot of the same web interface as above, but with the "View bill details" link highlighted in blue. The rest of the content is identical.

Close Bill Details

Paid on February 19, 2025

SUMMATION

Target	yesoul.lee.medit.clinic@gmail.com
source	Medit USA Inc.
bill	#9E5225E1-0001

ITEM

FEBRUARY 19 - FEBRUARY 20, 2025	
Medit Cloud Storage 10TB Plan for Clinics	US\$0.60
Quantity 1 piece	
subtotal	US\$0.60
1 Month Free for New Subscribers(100% 할인)	-US\$0.60
Payment amount	US\$0.00
Amount paid	US\$0.00
Remaining amount	US\$0.00

Do you have any questions? [Contact Medit USA Inc.](#)

Vouchers

You can register promotion codes and manage your vouchers on the **Vouchers** tab.

Manage Plan Membership Payment Details **Vouchers**

Register Code

Enter voucher code Register

Vouchers 13

All

Product	Promotion	Code Number	Registration Date	Valid Until	From	To	Status
Cloud Storage Plan	Basic EJ_Basic_Cloud Storage Plan (10TB)	2QJWR0FP	2026-09-21	2026-12-10	2026-10-21	2026-11-21	Used
Cloud Storage Plan	Basic EJ_Basic_Cloud Storage Plan (10TB)	QQ3OXJ3T	2026-09-21	2026-12-10	-	-	Registered
Cloud Storage Plan	Basic EJ_Basic_Cloud Storage Plan (10TB)	RSP4GZOI	2026-09-21	2026-12-10	-	-	Registered
Cloud Storage Plan	Basic EJ_Basic_Cloud Storage Plan (10TB)	JYAYH1XH	2026-09-21	2026-12-10	-	-	Registered
Cloud Storage Plan	Basic EJ_Basic_Cloud Storage Plan (10TB)	RSMPS5A4	2026-09-21	2026-12-10	-	-	Registered
Cloud Storage Plan	Redeem EJ_Redeem_Cloud Storage Plan (10TB)	ZGYUDFRY	2026-09-21	2100-01-01	-	-	Registered
Cloud Storage Plan	Redeem EJ_Redeem_Cloud Storage Plan (10TB)	1MGYSVHF	2026-09-21	2100-01-01	2026-09-21	2026-10-21	Used
Cloud Storage Plan	Redeem EJ_Redeem_Expired1	3FWHQAA	2026-09-21	2026-09-21	-	-	Expired
Cloud Storage Plan	Redeem EJ_Redeem_Expired1	COUR8QKR	2026-09-21	2026-09-21	-	-	Expired
Cloud Storage Plan	Redeem EJ_Redeem_Terminated1	PZPJXH8	2026-09-21	2026-12-21	2026-09-21	2026-10-21	Used
Cloud Storage Plan	Redeem SQA_Redeem_Cloud Storage Plan_36M	GFZKBWJ5	2026-09-15	2100-01-01	2026-09-15	2029-09-15	Activated
Cloud Storage Plan	Redeem SQA_Redeem_Cloud Storage Plan_36M	VE00ZB46	2026-09-14	2100-01-01	2026-09-14	2029-09-14	Used

Registered vouchers are listed with their promotion details and current status.

Registering a Voucher

Enter the promotion code you received and click **Register**.

Register Code

Register

Once registered, the voucher appears in the list below.

Vouchers 2

All

Product	Promotion	Code Number
Cloud Storage Plan	Redeem SQA_Redeem_Cloud Storage Plan_36M	GFZ

Activating a Voucher

To use a **Redeem** voucher, click **Activate** next to the voucher.

Valid Until	From	To	Status	
2026-12-10	2026-10-21	2026-11-21	Used	<button>Remove</button>
2026-12-10	-	-	Registered	
2026-12-10	-	-	Registered	
2026-12-10	-	-	Registered	
2026-12-10	-	-	Registered	
2100-01-01	-	-	Registered	<button>Activate</button>
2100-01-01	2026-09-21	2026-10-21	Used	<button>Remove</button>
2026-09-21	-	-	Expired	<button>Remove</button>
2026-09-21	-	-	Expired	<button>Remove</button>
2026-12-21	2026-09-21	2026-10-21	Used	<button>Remove</button>

Review the voucher details and important information, and click **Activate**.

Activate Voucher

Voucher Code

ZG'

Product

Cloud Storage Plan

Free Access

30 Days

⚠ Important

- Once activated, the voucher cannot be canceled.
- You can't subscribe to this product while the voucher is active.


Cancel

Activate

My Info

My Account

In My Info, you can review and edit the information you provided during registration and check your membership status.

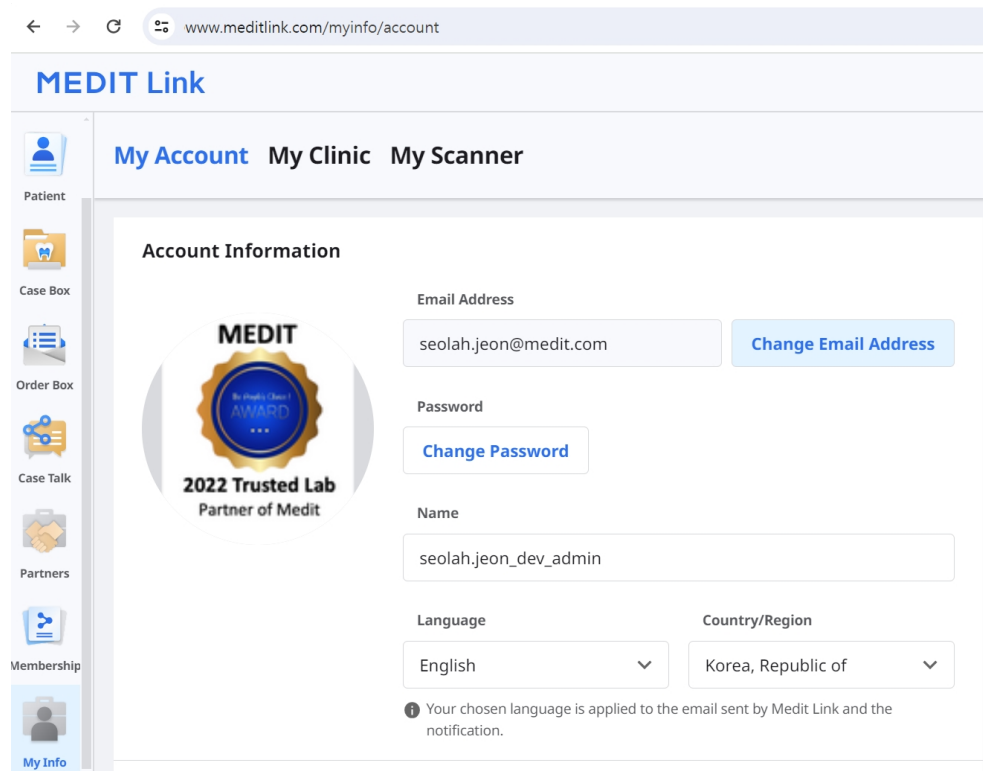

The screenshot shows a web browser window with the URL www.meditlink.com/myinfo/account. The page is titled "MEDIT Link" and has a navigation bar with "My Account", "My Clinic", and "My Scanner". A left sidebar contains icons for "Patient", "Case Box", "Order Box", "Case Talk", "Partners", "Membership", and "My Info" (which is highlighted). The main content area is titled "Account Information" and displays a circular logo for "2022 Trusted Lab Partner of Medit". Below the logo, there are input fields for "Email Address" (seolah.jeon@medit.com) with a "Change Email Address" button, "Password" with a "Change Password" button, and "Name" (seolah.jeon_dev_admin). There are also dropdown menus for "Language" (English) and "Country/Region" (Korea, Republic of). A note at the bottom states: "Your chosen language is applied to the email sent by Medit Link and the notification."

My Clinic/My Lab (Admin Account Only)

You can review and edit the information about your clinic or lab. This section is only available to admin accounts.

Dashboard

Patient

Case Box

Order Box

Case Talk

Partners

Membership

My Info

My Account

My Clinic

My Scanner

Clinic Name *

wonhoon-clinic

Edit

Cover and Logo

Looking for digitalized dental laboratories?

Find the most ordered Labs around your clinic with 'Featured Part

Services

Select

Specialties (choose up to 3)

Select

Select the services your business specializes in from the list.

About Us

e.g., Welcome to OOO Dental Clinic.
We are a licensed dental clinic situated in California specializing mainly in aesthetic restorations and implantology.

0/2000

Website URL

Add Website

Register up to 3 websites.

My Scanner

You can check the connected device information for the account.

No.	Type	Serial Number	Usage Time	Status	Last Connected	Last Calibrated	Warranty Start Date	Warranty Expiry
1	i700	AZ2105102619	0d 18h 44m	Unregistered	2024-04-12 15:19	2024-03-28 09:10	2024-02-27	2025-02-26
2	i700	AZ2105102622	0d 18h 39m	Unregistered	2023-04-26 20:30	2023-04-25 14:47	2021-07-09	2022-07-08
3	i700	AZ2305100410	0d 6h 38m	Unregistered	2025-04-03 13:35	2024-09-05 14:23	-	-
4	i600	BG2202100001	0d 0h 58m	Unregistered	2023-08-28 18:10	2023-07-13 15:00	2025-12-03	2026-12-02
5	i600	BG2202100026	0d 8h 32m	Unregistered	2022-10-28 17:03	2022-09-25 09:00	-	-
6	i600	BG2202100029	0d 2h 38m	Unregistered	-	2022-06-16 11:58	-	-
7	i900	CA2405100106	0d 11h 46m	Unregistered	2025-02-26 08:37	2024-11-07 10:11	-	-
8	i900	CA2512100037	0d 0h 0m	Registered	2026-04-08 17:35	2026-03-03 17:27	-	-
9	i900 Mobility	CG2502100039	0d 0h 19m	Unregistered	2026-04-06 15:34	2026-01-23 17:48	-	-
10	i900 Mobility	CG2506170009	0d 0h 0m	Unregistered	2026-03-12 12:08	2025-12-19 20:31	-	-

The scanner information is automatically registered once your scanner is connected to a Medit scan program.

Item	Description
Type	Displays the model name.
Serial Number	Displays the device's serial number.
Usage Time	Indicates how long the device has been used.
Status	Displays the device's registration status: • 'Registered' indicates that the scanner has been successfully connected and registered. • 'Unregistered' indicates that the scanner is no longer registered. If the scanner is registered with a different group, it will lose its registration in the current group.  Note To unregister your scanner: 1. Click the Registered button next to the scanner you want to unregister. 2. When the following pop-up appears, click Disconnect to confirm:  Do you want to unregister scanner I 4 from the current Clinic? Doing so will remove it from subscription benefits. Cancel Disconnect
Last Connected	Displays the date of the last connection.
Last Calibrated	Displays the date of the last calibration.
Warranty Start Date	Displays the date when the warranty was activated.

Warranty Expiry Date	Displays the date when the warranty expires.
Days Left	Displays the number of days remaining before the warranty expires.
Buy	Allows the user to purchase and extend the warranty.
Delete icon	Allows the user to remove the scanner from the list.

Order Forms (Lab Admin Account Only)

You can create and manage lab-defined order forms with the types, methods, materials, and shade systems your lab uses. Clinics can use your order form when your lab is assigned as the partner for a case.

In Medit Link Web, go to **My Info > Order Forms**.

MEDIT Link

My Account My Lab My Scanner **Order Forms**

Order Forms

Teeth 18 Arch 14

1. Inlay Anterior, Premolar, Molar

No.	Method	
1	-	N
		Z
		C
		R
		L
		C

Viewing Order Forms

Select **Teeth** or **Arch** to view the current order form for each category. Each type lists its available methods, materials, and shade systems. For teeth-based types, the allowed tooth groups are also shown.

Teeth 18 Arch 14

1. Inlay Anterior, Premolar, Molar

No.	Method	Material	Shade System
1	-	Multilayer Zirconia	VITA Classic / VITA 3D-Master / Ivoclar
		Zirconia	VITA Classic / VITA 3D-Master / Ivoclar
		Gold	-
		Resin	VITA Classic / VITA 3D-Master / Ivoclar
		LiSi	VITA Classic / VITA 3D-Master / Ivoclar
		Glass Ceramic	VITA Classic / VITA 3D-Master / Ivoclar
		Press Ceramic	VITA Classic / VITA 3D-Master / Ivoclar
		Hybrid Ceramic	VITA Classic / VITA 3D-Master / Ivoclar
		Composite Resin	VITA Classic / VITA 3D-Master / Ivoclar
		e.max	VITA Classic / VITA 3D-Master / Ivoclar
		Wax	-
		3D Print	VITA Classic / VITA 3D-Master / Ivoclar

Teeth 18 Arch 14

1. Denture

No.	Method	Material	Shade System
1	Full Denture	PMMA	VITA Classic / VITA 3D-Master / Ivoclar
		Resin	VITA Classic / VITA 3D-Master / Ivoclar
		Wax	-
2	Temporary	-	-
3	Denture Replica	Resin	VITA Classic / VITA 3D-Master / Ivoclar
4	Implant-Supported Denture	Resin	VITA Classic / VITA 3D-Master / Ivoclar
5	Wax Rim	Wax	-

2. Partial Denture

No.	Method	Material	Shade System
1	Framework	PMMA	VITA Classic / VITA 3D-Master / Ivoclar
		Metal	VITA Classic / VITA 3D-Master / Ivoclar

Editing Order Forms

1. Click **Edit** in the upper-right corner of the screen.

⌚ Change History **Edit**

2. Select **Teeth** or **Arch**, then select a type from the list.

Teeth 18 Arch 14

Manage Types

- ☐ Inlay
- ☐ Onlay
- ☐ Coping
- ☐ Denture Coping
- ☒ Crown
- ☐ Partial

Allowed Teeth

☒ Anterior ☒ Premolar ☒ Molar

3. For teeth-based types, select the applicable groups under **Allowed Teeth**.

Teeth 18 Arch 14

Manage Types

- ☐ Inlay
- ☐ Onlay
- ☐ Coping
- ☐ Denture Coping
- ☒ Crown
- ☐ Partial

Allowed Teeth

☒ Anterior ☒ Premolar ☒ Molar

4. Select the methods, materials, and shade systems to include.

Crown

Allowed Teeth

☒ Anterior ☒ Premolar ☒ Molar

Methods ¹

Search

Screw-Retained **Full-Contour** Offset Cement-Retained

SCRIP Temporary Cut-back Pontic Press Preform Screw-Retained Layered

Screw-Retained PFZ Screw-Retained PFG Screw-Retained PFM Cement-Retained Layered

[More](#)

☐ Do not use

Full-Contour

Materials ²

Search

Multilayer Zirconia Zirconia PMMA Gold Resin

LiSi Glass Ceramic Press Ceramic Hybrid Ceramic Composite Resin e.max

Titanium Metal Stainless Steel Co-Cr (Cobalt Chromium) Ni-Cr (Nickel-Chromium) Wax

PEEK 3D Print [+ Add Material](#)

[Hide](#)

☐ Do not use

Multilayer Zirconia

Shade Systems ³

[+ Add Shade System](#)

☒ VITA Classic 20 Shades [^](#)

None A1 A2 A3 A3.5 A4 B1 B2 B3 B4 C1 C2 C3 C4 D2 D3

D4 OM1 OM2 OM3

☐ VITA 3D-Master 82 Shades [v](#)

☐ Ivoclar 25 Shades [v](#)

☐ Material 5 Shades [v](#)

If a method or material is not required, select **Do not use** in the corresponding section. You can set this separately for methods and materials.

Crown

Allowed Teeth

☒ Anterior ☒ Premolar ☒ Molar

Methods

Search

Screw-Retained

Full-Contour

Offset

Cement-Retained

SCRIP

Temporary

Cut-back

Pontic

Press

Preform

Screw-Retained Layered

Screw-Retained PFZ

Screw-Retained PFG

Screw-Retained PFM

Cement-Retained Layered

☐ Do not use

Crown

Full-Contour

Allowed Teeth

☒ Anterior ☒ Premolar ☒ Molar

Materials

Search

LiSi

Glass C

Titanium

☐ Do not use

Multilayer Z

Shade Syst

☒ VITA Clas

☐ VITA 3D-I

☐ Ivoclar

☐ Material

Methods

☒ Do not use

Do not use 0

Materials

Search

Multilayer Zirconia

Zirconia

PMMA

Gold

Resin

LiSi

Glass Ceramic

Press Ceramic

Hybrid Ceramic

Composite Resin

e.max

Titanium

Metal

Stainless Steel

Co-Cr (Cobalt Chromium)

+ Add Material

☐ Do not use

Crown

Allowed Teeth

☒ Anterior ☒ Premolar ☒ Molar

Methods

☒ Do not use

Do not use 1


Materials

☒ Do not use

> Do not use

5. Repeat these steps for other types as needed.

6. Click **Save** in the upper-right corner, then click **Save** again in the confirmation dialog.

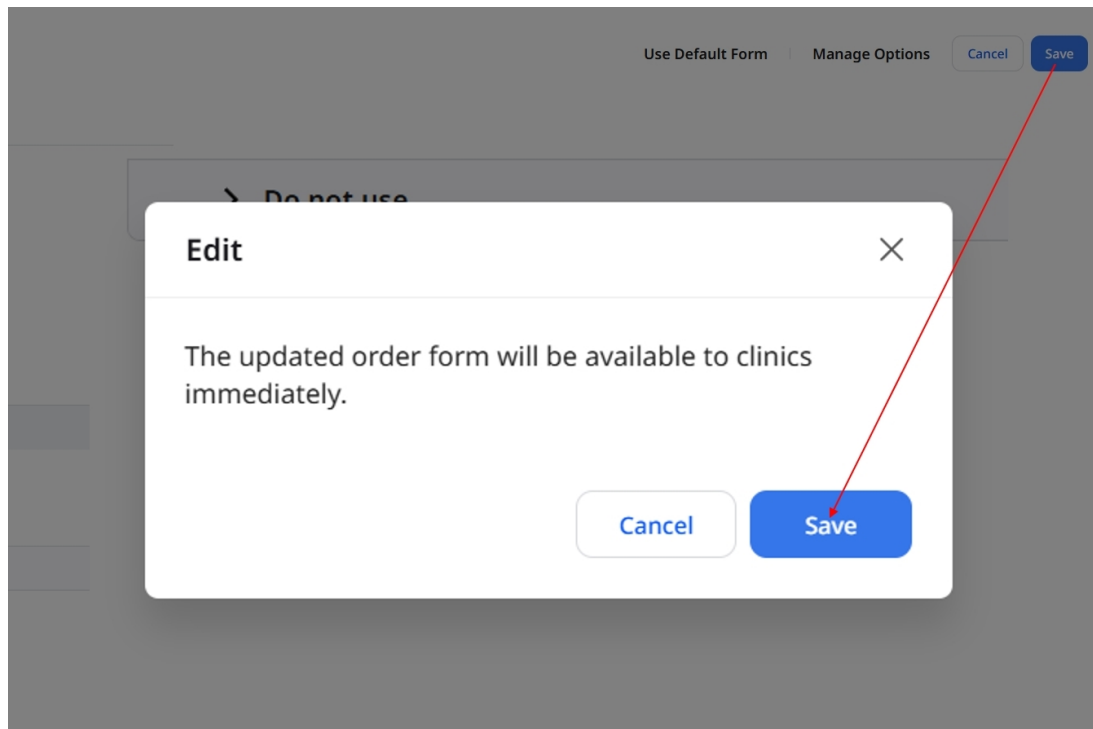


Note

For a combined type (a combination of two individual types), methods, materials, and shade systems follow the settings of its constituent types. Edit those types to change the options.

Allowed tooth groups can be set separately for the combined type.

Using the Default Form

1. On the order form editing screen, click **Use Default Form**.
2. Review the information in the dialog and click **Continue**.

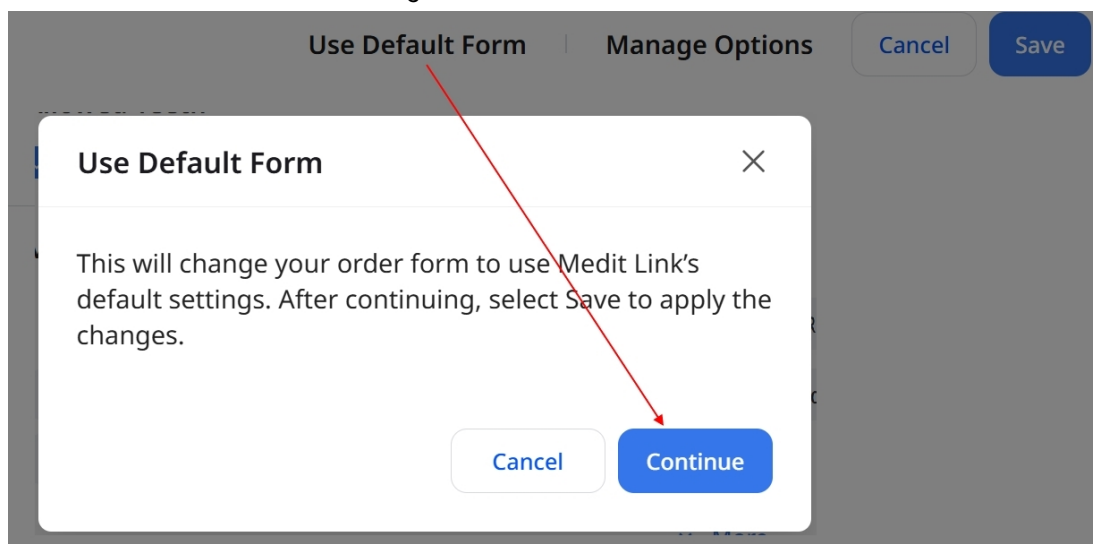

3. Click **Save** in the upper-right corner to apply the changes.

Use Default Form

Manage Options


Cancel

Save

Managing Options

Enabling or Disabling Options

1. On the order form editing screen, click **Manage Options**.

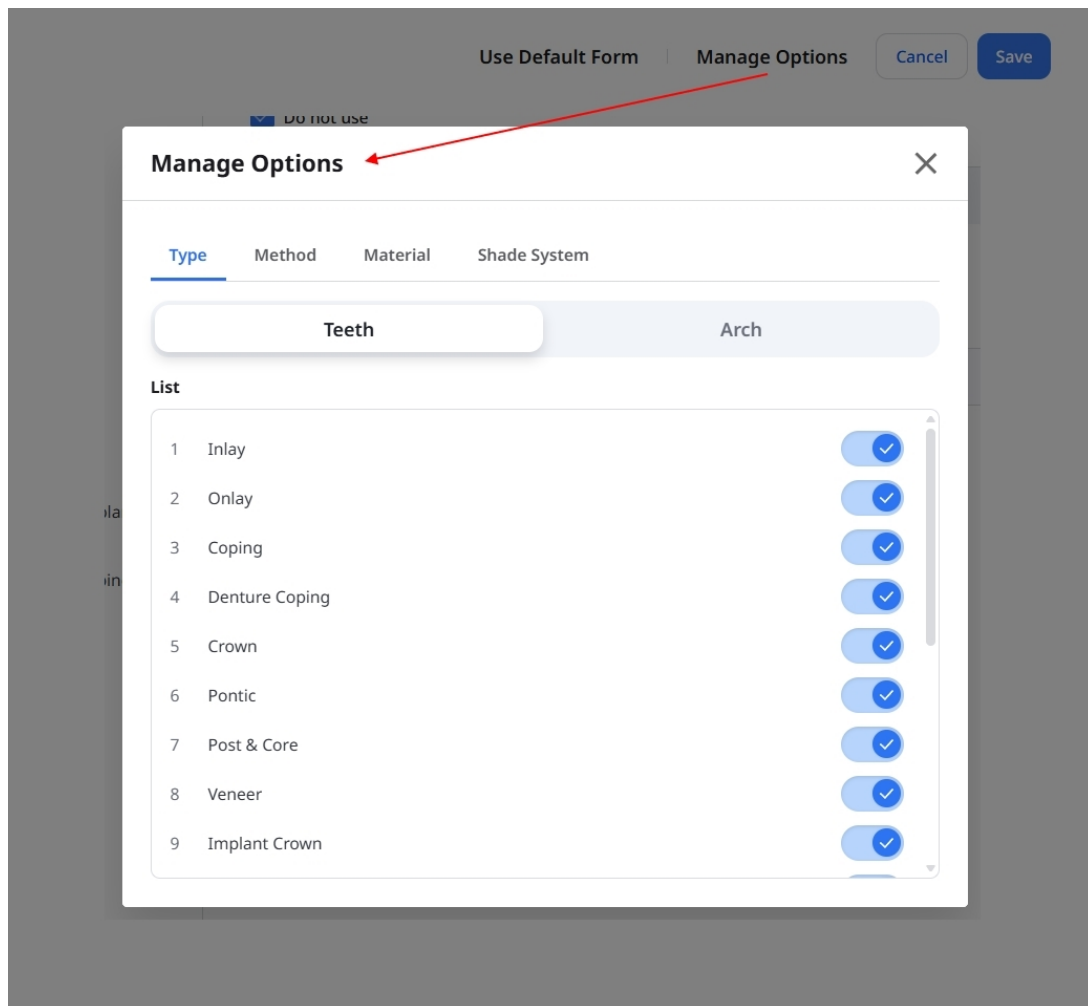

2. Select the **Type**, **Method**, **Material**, or **Shade System** tab.

- On the **Type** tab, select **Teeth** or **Arch**.

3. Turn the switch next to an option on or off.



4. Close **Manage Options** and review your form.
5. Click **Save** in the upper-right corner, then click **Save** again in the confirmation dialog.


Note

At least one type must remain enabled across the **Teeth** and **Arch** tabs.

Disabling a type also disables any combined types that include it. To re-enable a combined type, first enable the two types that make it up, then enable the combined type itself.

Adding Materials

1. On the order form editing screen, click **Manage Options** and select the **Material** tab.

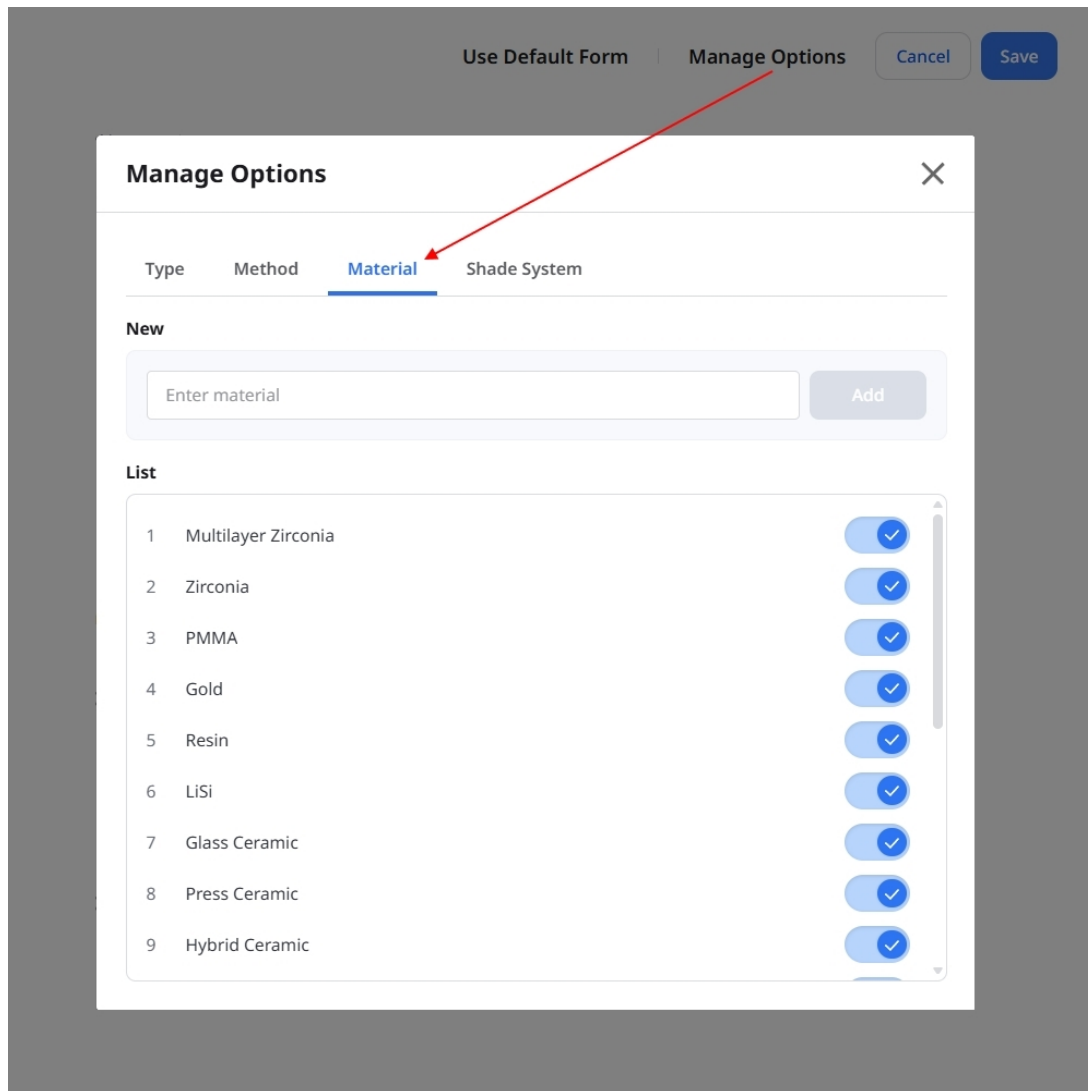

The screenshot shows a 'Manage Options' dialog box with a close button (X) in the top right corner. At the top of the dialog, there are four tabs: 'Type', 'Method', 'Material', and 'Shade System'. The 'Material' tab is currently selected and highlighted with a blue underline. A red arrow points from the 'Manage Options' button in the background to the 'Material' tab. Below the tabs, there is a 'New' section with a text input field labeled 'Enter material' and an 'Add' button. Below that is a 'List' section containing a scrollable list of materials, each with a checkbox and a checkmark icon. The materials listed are: 1 Multilayer Zirconia, 2 Zirconia, 3 PMMA, 4 Gold, 5 Resin, 6 LiSi, 7 Glass Ceramic, 8 Press Ceramic, and 9 Hybrid Ceramic. All checkboxes are currently checked.

Use Default Form | Manage Options | Cancel | Save

Manage Options X

Type | Method | **Material** | Shade System

New

Enter material Add

List

1	Multilayer Zirconia	☒
2	Zirconia	☒
3	PMMA	☒
4	Gold	☒
5	Resin	☒
6	LiSi	☒
7	Glass Ceramic	☒
8	Press Ceramic	☒
9	Hybrid Ceramic	☒

2. Under **New**, enter a material name and click **Add**. The new material appears in the list below.

Use Default Form | Manage Options | Cancel | Save

Manage Options

Type | Method | **Material** | Shade System

New

Zirconia HT

List

12	Titanium	☒
13	Metal	☒
14	Stainless Steel	☒
15	Co-Cr (Cobalt Chromium)	☒
16	Ni-Cr (Nickel-Chromium)	☒
17	Wax	☒
18	PEEK	☒
19	3D Print	☒
20	Zirconia HT edit delete	☒

3. Close **Manage Options** and select the new material for the types where you want to use it.

Crown

Allowed Teeth

☒ Anterior ☒ Premolar ☒ Molar

Methods

Search

Screw-Retained **Full-Contour** Offset Cement-Retained

SCRIP Temporary Cut-back Pontic Press Preform Screw-Retained Layered

Screw-Retained PFZ Screw-Retained PFG Screw-Retained PFM Cement-Retained Layered

More

☐ Do not use

Full-Contour 1

Materials

Search

Multilayer Zirconia Zirconia PMMA Gold Resin

LiSi Glass Ceramic Press Ceramic Hybrid Ceramic Composite Resin e.max

Titanium Metal Stainless Steel Co-Cr (Cobalt Chromium) Ni-Cr (Nickel-Chromium) Wax

PEEK 3D Print **Zirconia HT** + Add Material


Hide

☐ Do not use

4. Click **Save** in the upper-right corner, then click **Save** again in the confirmation dialog.

Adding Shade Systems


1. On the order form editing screen, click **Manage Options** and select the **Shade System** tab.

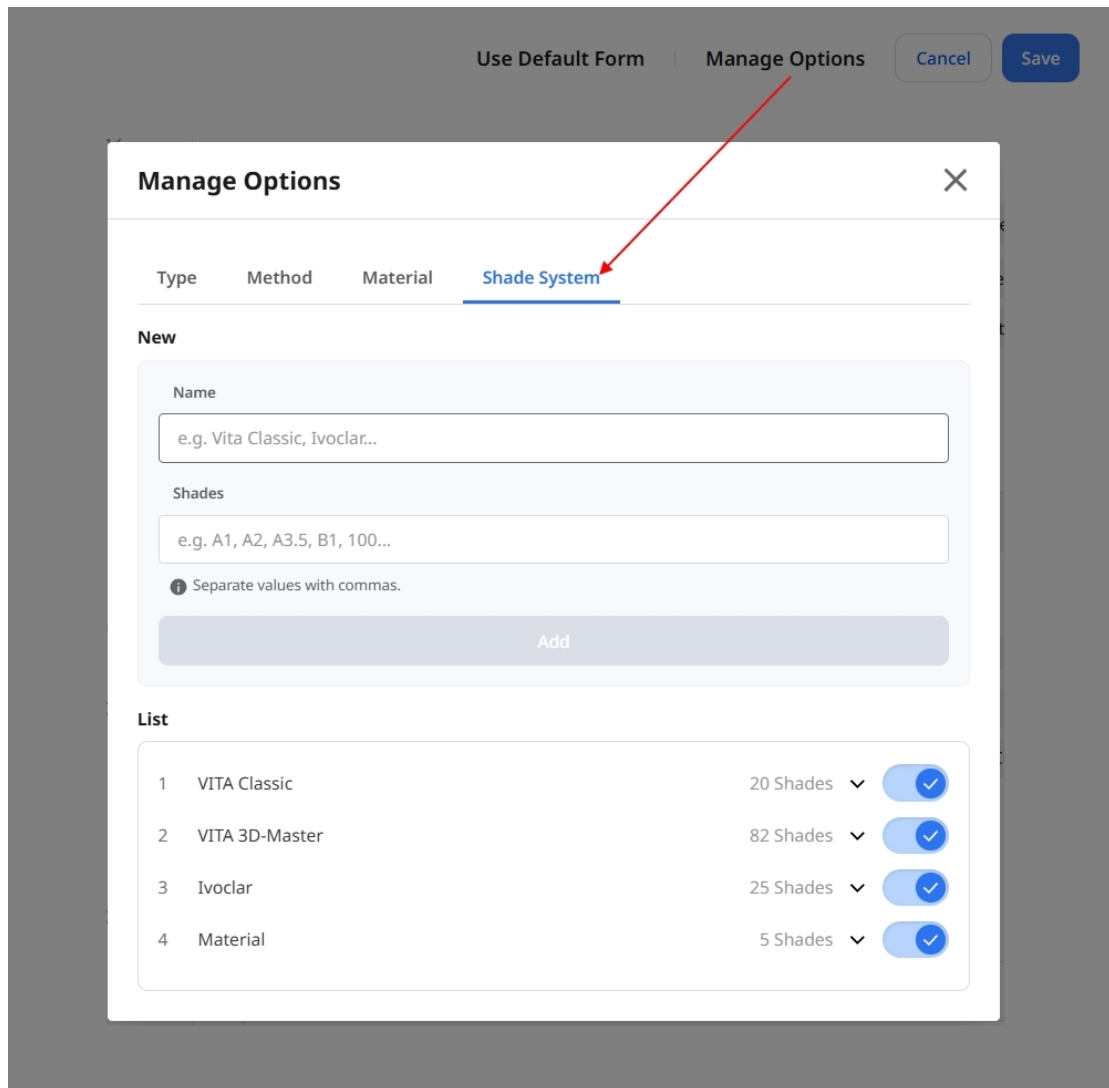

The screenshot shows the 'Manage Options' dialog box with the 'Shade System' tab selected. The dialog has a title bar with 'Use Default Form', 'Manage Options', 'Cancel', and 'Save' buttons. The 'Shade System' tab is highlighted with a red arrow. The 'New' section contains a 'Name' field with the placeholder 'e.g. Vita Classic, Ivoclar...' and a 'Shades' field with the placeholder 'e.g. A1, A2, A3.5, B1, 100...'. Below the 'Shades' field is a note 'Separate values with commas.' and an 'Add' button. The 'List' section shows a table of existing shade systems:


ID	Name	Shades	Active
1	VITA Classic	20 Shades	☒
2	VITA 3D-Master	82 Shades	☒
3	Ivoclar	25 Shades	☒
4	Material	5 Shades	☒

2. Enter the shade system name in the **Name** field. In the **Shades** field, enter the shade values separated by commas.

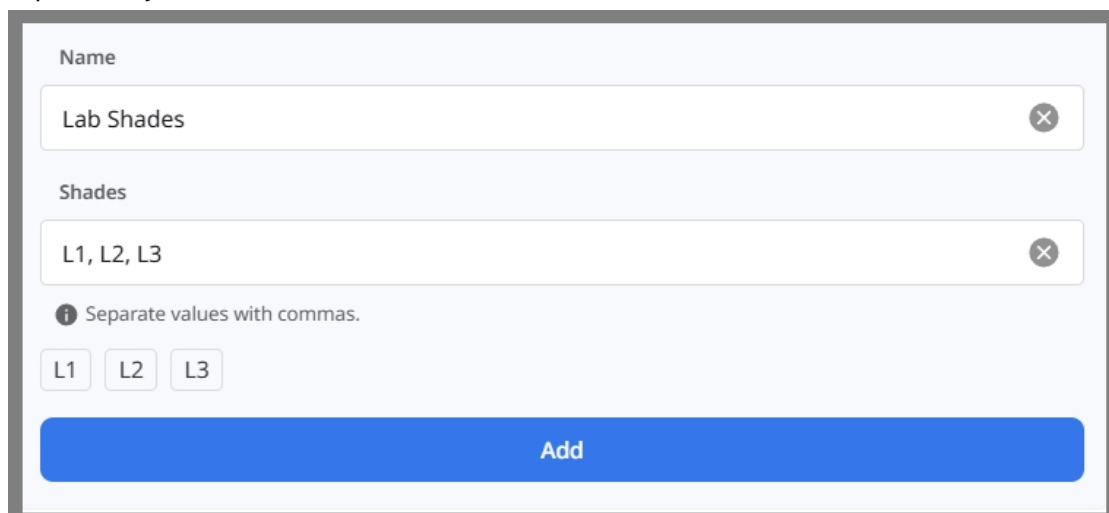

The screenshot shows the 'New' section of the 'Manage Options' dialog box. The 'Name' field contains 'Lab Shades' and the 'Shades' field contains 'L1, L2, L3'. Below the 'Shades' field is a note 'Separate values with commas.' and three buttons labeled 'L1', 'L2', and 'L3'. An 'Add' button is at the bottom.

3. Review the values and click **Add**.

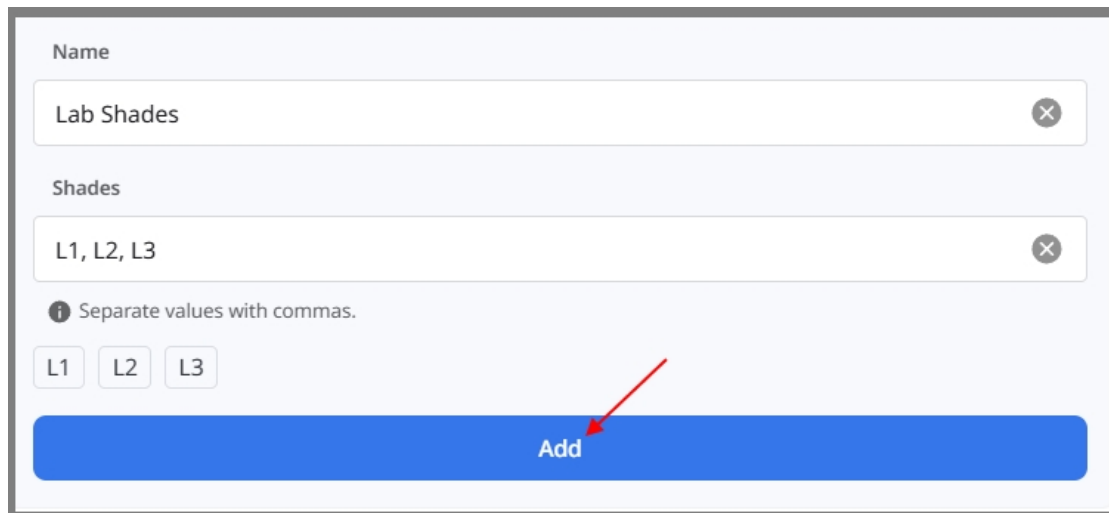

The form is titled "Name" and "Shades". The "Name" field contains "Lab Shades". The "Shades" field contains "L1, L2, L3". Below the "Shades" field, there is a hint: "i Separate values with commas." and three buttons labeled "L1", "L2", and "L3". At the bottom of the form is a large blue button labeled "Add". A red arrow points to the "Add" button.

Name

Lab Shades

Shades

L1, L2, L3

i Separate values with commas.

L1 L2 L3

Add

4. Close **Manage Options** and select the new shade system for the materials where you want to use it.

The screenshot shows the 'Crown' configuration window. Under 'Allowed Teeth', 'Anterior', 'Premolar', and 'Molar' are all checked. In the 'Methods' section, 'Full-Contour' is selected among various options like 'Screw-Retained', 'Offset', 'Cement-Retained', etc. Below this, a 'Do not use' checkbox is present. A section titled 'Full-Contour 1' contains a 'Materials' list where 'Zirconia HT' is highlighted. Another 'Do not use' checkbox is below the materials. The 'Zirconia HT' section is expanded, showing 'Shade Systems'. Here, 'Lab Shades' is checked, indicated by a red arrow, while others like 'VITA Classic' and 'Ivoclar' are unchecked. To the right of the shade systems, the number of shades (e.g., 20 Shades, 82 Shades) is listed with expand/collapse arrows. At the bottom of the shade systems list are buttons for 'L1', 'L2', and 'L3'.


5. Click **Save** in the upper-right corner, then click **Save** again in the confirmation dialog.

Note

New materials and shade systems are saved when you click **Add**. Canceling your order form edits does not remove these options.

Viewing Change History


On the **Order Forms** page, click **Change History** to review saved changes, with the most recent changes shown first.

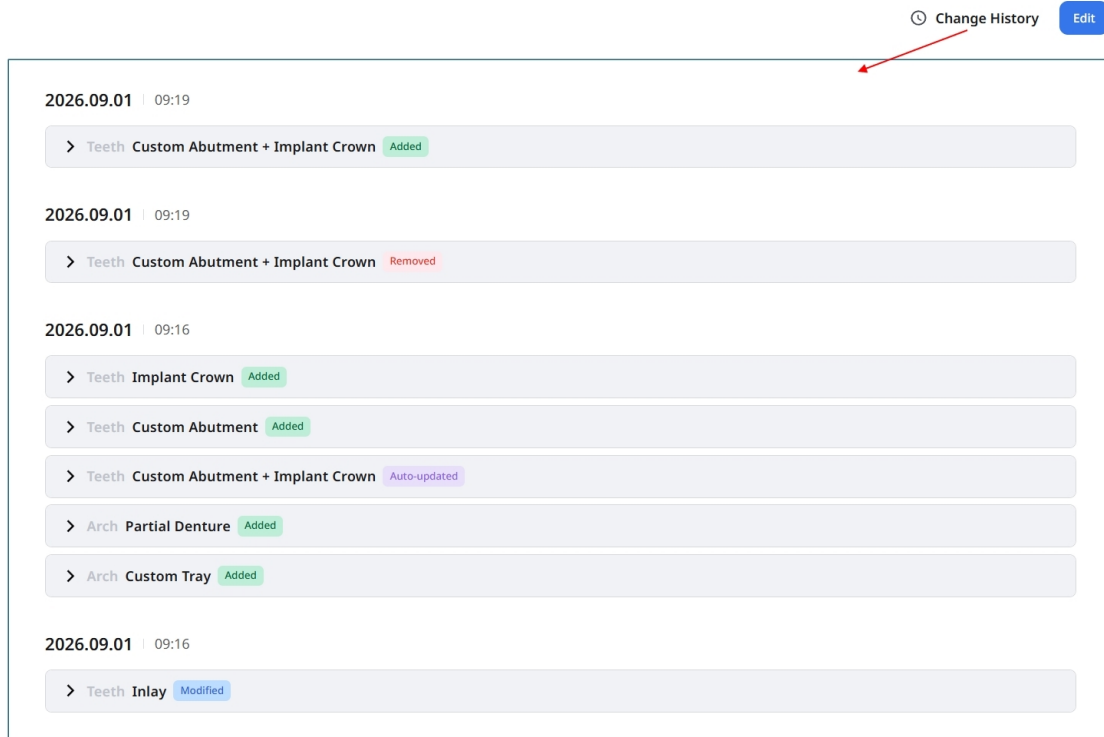

The screenshot shows the 'Change History' page with a list of changes. At the top right, there is a 'Change History' link with a clock icon and an 'Edit' button. A red arrow points to the 'Change History' link. The list of changes is as follows:

- 2026.09.01 | 09:19
 - > Teeth Custom Abutment + Implant Crown Added
- 2026.09.01 | 09:19
 - > Teeth Custom Abutment + Implant Crown Removed
- 2026.09.01 | 09:16
 - > Teeth Implant Crown Added
 - > Teeth Custom Abutment Added
 - > Teeth Custom Abutment + Implant Crown Auto-updated
 - > Arch Partial Denture Added
 - > Arch Custom Tray Added
- 2026.09.01 | 09:16
 - > Teeth Inlay Modified

Expand a type to view its details.

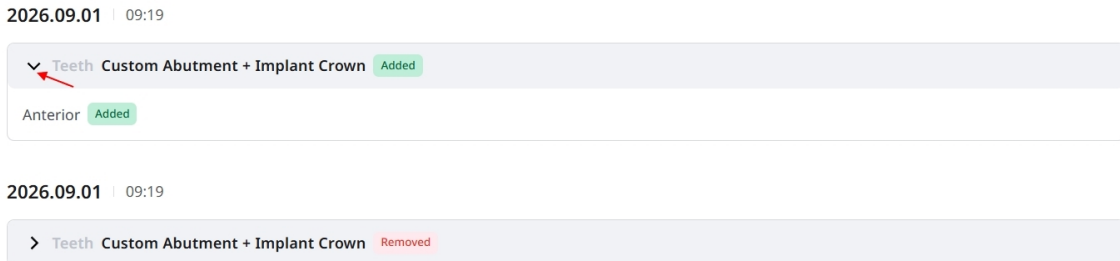

The screenshot shows the expanded details for the change 'Teeth Custom Abutment + Implant Crown' (Added) on 2026.09.01 at 09:19. A red arrow points to the expand/collapse icon (a downward arrow) next to the change type. The details show 'Anterior' (Added).

2026.09.01 | 09:19

> Teeth Custom Abutment + Implant Crown Added

Anterior Added

2026.09.01 | 09:19

> Teeth Custom Abutment + Implant Crown Removed

- **Added** and **Removed** identify options added to or removed from the form.
- **Modified** identifies changes that include both additions and removals.
- **Auto-updated** identifies changes to a combined type resulting from changes to its constituent types.

Note

Only changes applied to the order form appear in Change History. Previous versions cannot be restored from this page.

Other Menu

Notice & Update

You can check the important notices from Medit and program updates of Medit software.

Notice & Update

Notice
NFT Service Termina...
2023-04-19

Notice
End of Support Notic...
2023-04-14

Update
Release Note for Me...
2023-03-29

Notice
[Medit] We've updat...
2023-03-29

Update
Release Note for Me...
2023-03-17

Update
Medit Occlusion Ana...
2023-03-08

Notice
NFT Service Termination Notice
2023-04-19

To our users,

We want to thank all Medit Link users who have supported and used our NFT service.

Unfortunately, we must inform you that this service will be terminated following the upcoming program update.

Please see the details below for more information:

- **Feature in Question**
 - Share via NFT (in Medit Link Case Talk)
- **Effective Date**
 - Starting from the update in the first half of 2023 (scheduled for April 27)

※ The update schedule is subject to change.


We'll continuously work to provide better service.

Medit Link Team

Help & Resources

Click the Help & Resources icon to open the menu, then select a resource.



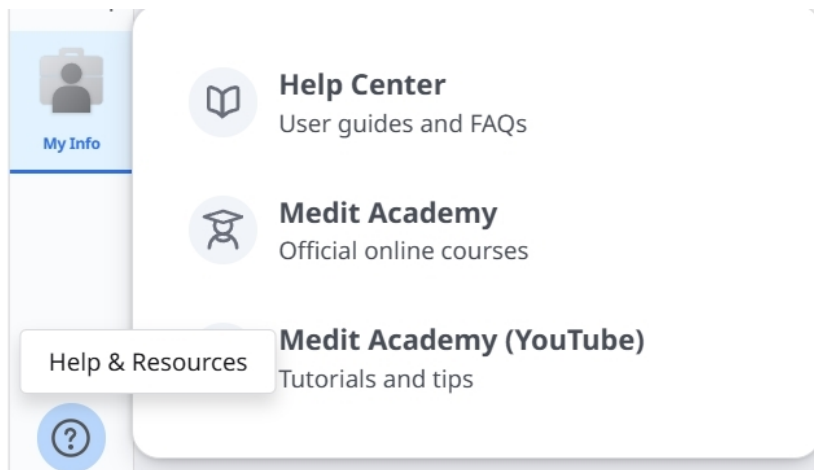

Download Medit Link App

You can download the latest version of the program installer from the web. This menu is helpful when you lose the installer provided in the product package, or you have an outdated installer.

Note

Medit Link updates automatically when you have an internet connection.



Overview

The screenshot displays the 'Patient' section of the MEDIT Link application. The top bar shows the title 'MEDIT Link' and version '3.4.10'. The sidebar on the left contains icons for Dashboard, Patient, Case Box, Order Box, Case Talk, and Trash Box. The main area shows a table of patient cases with columns: Patient Name, Patient ID, Gender, Date of Birth, Number of Cases, and Registration Date. The table is filtered by 'All' and shows a total of 640 cases. The table lists several cases, including '0511 test', '0127 test', '19.12.09 test', '150', '123454321', '1231456', '123', '111', '1', '1', '1', and '1'. Each case has a 'Scan' button and a vertical ellipsis menu icon.

1	Title Bar	Refer to Medit Link App > Title Bar for detailed information.
2	Menu	Refer to Medit Link App > Menu for detailed information.
3	Tools and Functions	Refer to Medit Link App > File Viewer and Medit Link App > View Options for detailed information.
4	Case List	Refer to Medit Link App > Filtering and Search for detailed information.

The table below summarises what menus are available in each type of clinic and lab account.

Account Type		Patient	Case Box	Order Box	In Box	Work Box	Case Talk
Clinic	Admin	0	0	0	X	X	

Account Type		Patient	Case Box	Order Box	In Box	Work Box	Case Talk
	Member	0	0	0	X	X	0
Lab Accounts	Admin	X	X	X	0	0	
	Member	X	X	X	0	0	0

Trash Box	App Box
0	0
0 * Restore feature only	0
0	0
0 * Restore feature only	0


Title Bar

The title bar provides quick access to features such as account information, notifications, notices and updates, settings, and network status.



My Account (My Information)

You can access your account, membership, and partner settings. Select My Info, Membership, or Partners to open the web app.

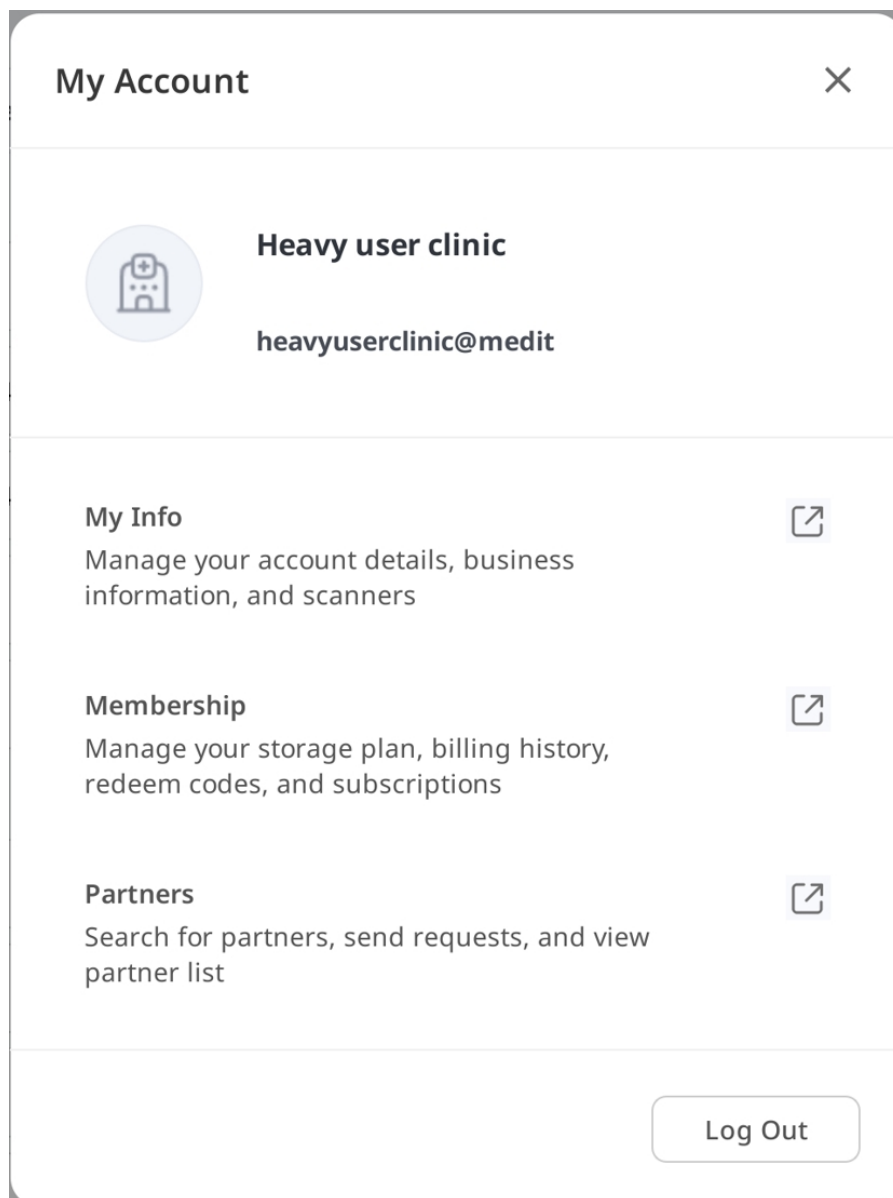

Notifications

You can view notifications about your account activity, such as partnerships, orders, and deliveries.

Notifications

All

Unread





Promotion

1일 클리닉 Premium 프로모션 프로모션이 종료되었습니다.

2026-03-14 AM 12:50

>

Delivery

"heavyuserlab@meditqa.com"에서 주문에 대한 배송을 시작하였습니다.

2026-03-13 PM 03:56

>

Order

"heavyuserlab@meditqa.com"에서 주문을 거절했습니다.

2026-03-13 PM 03:45




Notice & Update

You can check the important notices from Medit and program updates of Medit software.

Settings

The Settings dialog allows users to configure the program.

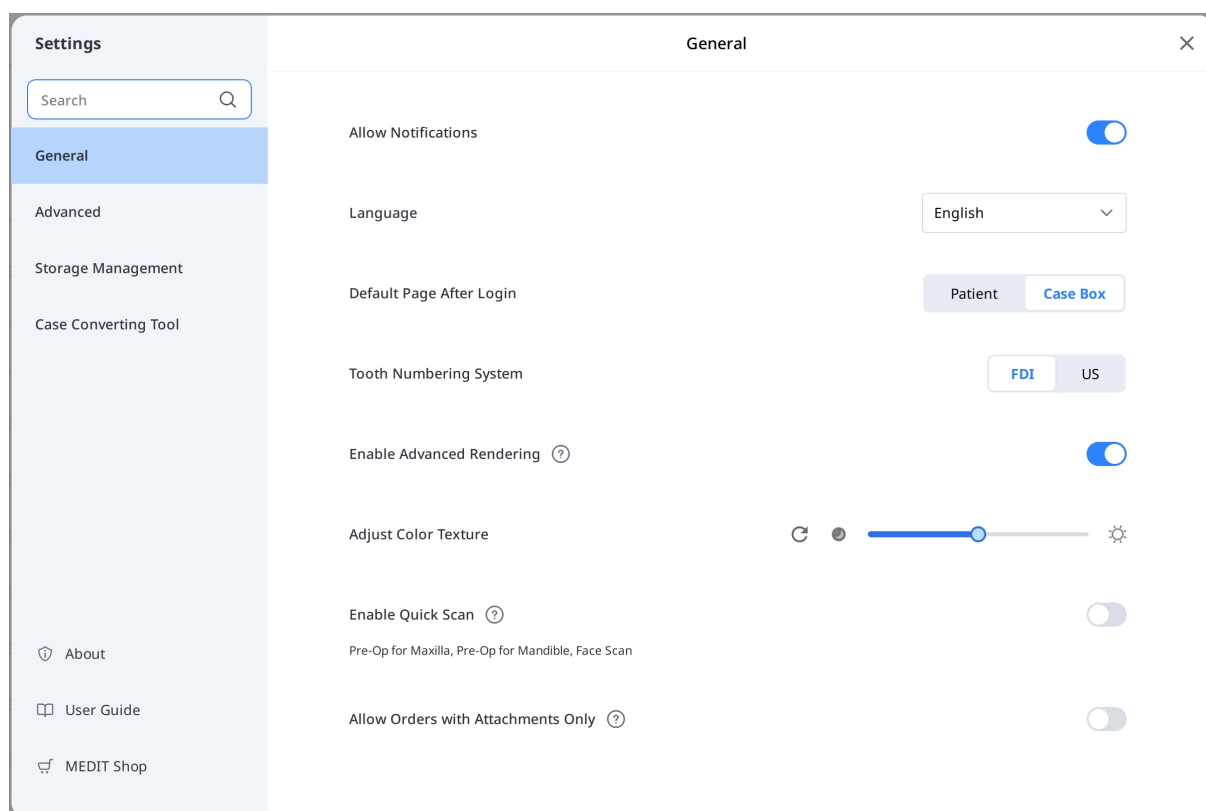

Please refer to the [Medit Link App > Settings](#) for detailed information about the setting options.

Network Status

Medit Link enables you to use all the features seamlessly when connected to a network.

To help you check the status of your network connection in real time, the network status icon is displayed in the title bar.

Network Status



Network speed is good.



Network speed is fast enough.



Network speed is poor. Can be switched to offline mode automatically.



Either network is not connected or Medit Link is in offline mode.



Synchronization is in progress. Some case lists may not match.



Synchronization failed. The case lists may not match.



Work Offline

Close

If the network connection is not good, you may experience the following restrictions:

- Some functions will be restricted, such as ordering a case, creating a case, sharing via Case Talk, and cloud synchronization.
- The overall performance of the service can be reduced due to the impact of the network environment.

In this case, you can enter the Offline Mode to perform basic tasks and then use the features that require internet once you are connected.



Note


The network symbol indicates the response speed between the current area and AWS (Amazon Web Service). It does not mean your network speed, so sometimes it may show the status as poor or bad, even though your network speed is fast enough.

If the other process uses the same network, the speed may show as poor or bad, and the overall performance of Medit Link could be bad. In this case, please do not share the network with other processes.

Settings


Search Settings

Enables you to find a feature name across all tabs in the **Settings** window. Selecting a result automatically opens the related tab, helping you locate options more quickly.

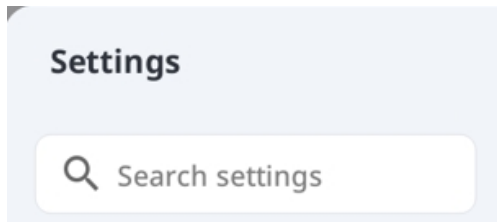

General

Provides options for overall system setup.

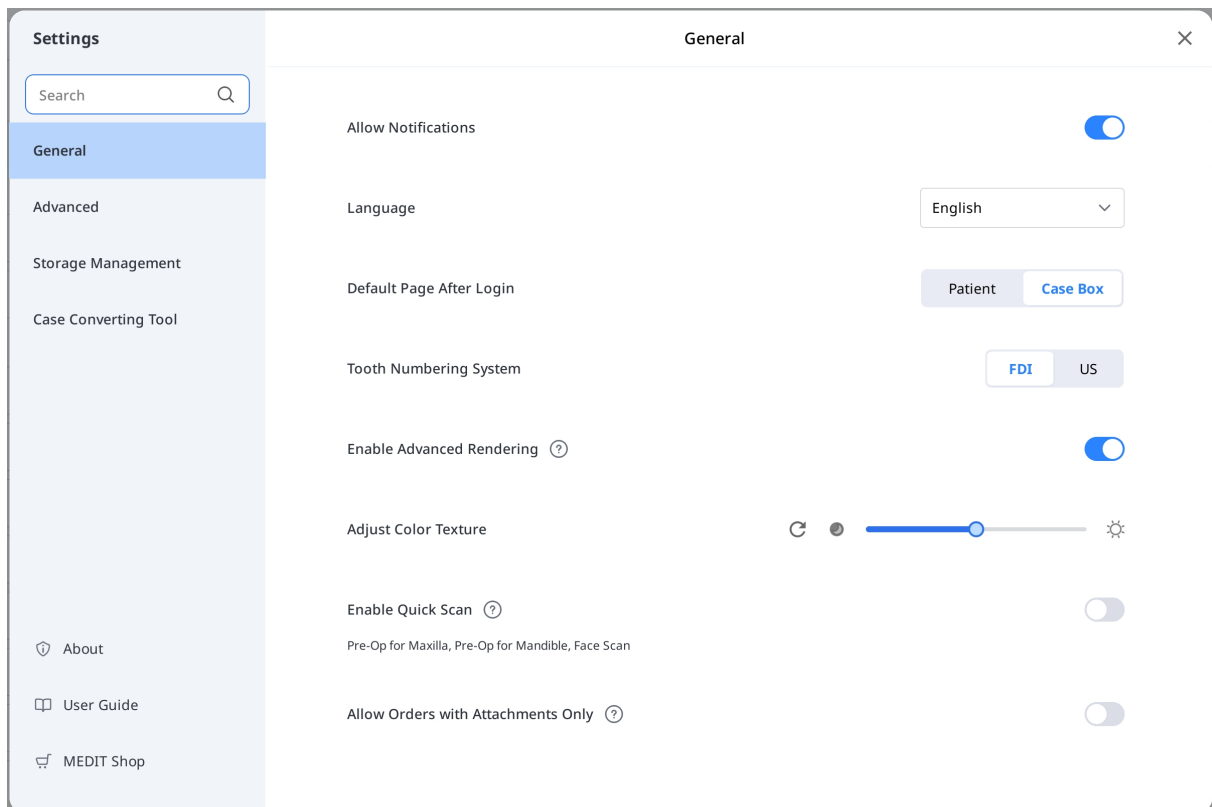

Language	Select the language for the program.
Tooth Numbering System	Select a teeth labeling system between FDI and US.

Allow Notifications	When checked, you will receive notifications for program updates.	
Enable Quick Scan	Choose if you want to have the scan stages below automatically appear in Medit Scan for Clinics without having to fill in the form information. • Pre-op for Maxilla • Pre-op for Mandible • Face Scan	
Enable Advanced Rendering	Advanced Rendering	When checked, the 3D model is displayed in a more vivid and realistic rendering.
	Adjust Color Texture	Adjust the brightness of the 3D model. * The resulting colors may slightly differ from Medit Link when you open the same data with other software.
Allow Orders with Attachments Only	When disabled, orders with only attachments cannot be placed. This helps prevent unintended submissions. Settings saved in the app are also applied to the web app.	

Allow Orders with Attachments Only	When disabled, orders with only attachments cannot be placed. This helps prevent unintended submissions. Settings saved in the app are also applied to the web app.	
Default Page After Login	Select the default menu to display when starting the program. • Clinic accounts can choose between Patient and Case Box • Lab accounts can choose between In Box and Work Box	

Advanced

Configure advanced options such as network proxy, default app, and experience sharing settings.

Settings

Search settings

General

Advanced

Storage Management

Case Converting Tool

About

User Guide

Advanced

Export

Export Settings

Network

Proxy Server

Proxy Server IP (or DNS)

Proxy Server Port

0

Proxy Server Scheme

HTTPHTTPS

Proxy Server Username (Option)

Proxy Server Password (Option)

Connect

Default App

Medit Scan for Clinics

C:/Program Files/Medit/Medit Link/Medit Link/Medit_iScan.exe

User Experience Program

Join Medit User Experience Improvement Program

Export – Export Settings

Set how to export files from Medit Link.

Settings

Q Search settings

General

Advanced

Storage Management

Case Converting Tool

About

User Guide

Q Search settings

General

Advanced

Storage Management

Case Converting Tool

About

User Guide

<

Export Settings

×

Type

Medit

Export Files To

C:/Users/medit/Downloads/

If duplicate files exist ?

Ask every time

Mesh Format ?

✓ MEDITMesh

OBJ

✓ PLY

STL

✓

 Do not convert the attached files

Image Format ?

No image

Compress Files

Combine Individual Mesh

Rule for Folder Name

✓ Case Creation Date

Scanned Date

Order Date

✓ Patient Name

Case Name

Order ID

Partner Name

Patient ID

Assignee

Preview Folder Name

Case Creation Date - Patient Name

Rule for File Name

✓ Case Creation Date

Scanned Date

Order Date

✓ Patient Name

Case Name

Order ID

Partner Name

Patient ID

Assignee

Preview File Name

Type	Select the CAD program you want to use for designing. The program will adjust the results based on your selection, allowing you to proceed with the design without any additional axis alignment.
Export Files To	Specify the path where you want to store the results.
If duplicate files exist	Choose the default action for cases when the same folder or file name exists in the target path. • Always ask when exporting. • Export a new file and overwrite the existing files. The file currently stored will be lost. • Create a folder and export the files
Mesh Format	Specify the format of the final mesh file.
Do not convert the attached files	When enabled, attached files will be exported in their original formats without conversion.
Image Format	Specify the format of the exported image files.
Compress & Combine	Compress Files: Compress all results into a single file. Combine Individual Mesh: Combine individual mesh files for the maxilla or mandible into a single mesh file. Also keep individual mesh files: Keep the individual mesh files in addition to the combined mesh file. This option is available when Combine Individual Mesh is enabled.

Compress & Combine	Compress Files: Compress all results into a single file. Combine Individual Mesh: Combine individual mesh files for the maxilla or mandible into a single mesh file. Also keep individual mesh files: Keep the individual mesh files in addition to the combined mesh file. This option is available when Combine Individual Mesh is enabled.
Rule for Folder Name	Define rules for creating a folder name when exporting files. The set rules will be applied to folder naming, and the exported file will be saved in the created folder. The below information can be used to set the rule. • Clinic account Case Creation Date, Scan Date, Order Date, Patient Name, Case Name, Partner Name, Patient ID, Order ID, Assignee • Lab account Case Creation Date, Scan Date, Order Date, Patient Name, Case Name, Partner Name, Order ID, Assignee
Rule for File Name	Define rules for the exported file name. The below information can be used to set the rule. • Clinic account Case Creation Date, Scan Date, Order Date, Patient Name, Case Name, Partner Name, Patient ID, Order ID, Assignee • Lab account Case Creation Date, Scan Date, Order Date, Patient Name, Case Name, Partner Name, Order ID, Assignee

Network – Proxy Server

Set a proxy server for Medit Link using a custom IP address, port number, and optional credentials.

Default App

Manage the programs that support integration in Medit Link.

You can connect the following programs:

	Clinic Account	Lab Account
Scan Program	Medit Scan for Clinics	Medit Scan for Clinics Medit Scan for Labs



Integration of the following services is not available for Medit Link for Mac.

- Medit Scan for Labs

User Experience Program

Enable this option to help improve Medit products by sharing anonymous usage data.

Storage Management – Overview & Cleanup

This section provides a visual breakdown of your cloud and local storage usage. It helps you understand how your storage is being used and which file types consume the most space. Click "Free Up Space" to delete unnecessary files from cloud or local storage within the selected date range.

Settings

Search

General

Advanced

Storage Management

Case Converting Tool

Storage Management

All / 2018-01-01 ~ 2026-04-09

My Cloud Storage 3.59 GB / 2.00 GB

Result Files & Attachments 3.59 GB

Local Storage 329.16 GB / 475.87 GB

Raw Data 0.00 KB

Result Files & Attachments 0.00 KB

Others 329.16 GB

Free Up Space

Raw data refers to scan data captured during scanning, while result files and attachments are processed scan results and additional files attached to the case. [Learn More](#)

Set Period

Set the specific period to delete data created within the time frame below:

- All: Show all files created since January 1, 2018, regardless of the selected date range.
- Custom Range: Select a start date (no earlier than Jan 1, 2018) and an end date to filter files.
- Today: Files created or modified on the current day only.
- Yesterday: Files created or modified one day before today.
- Last 7 Days: Files created or modified in the past 7 days, including today.
- Last 30 Days: Files created or modified in the past 30 days, including today.
- Current Month: Files created or modified within the current calendar month.
- Previous Month: Files created or modified within the previous calendar month.

Upgrade Storage

Select Upgrade Storage to subscribe to a paid storage plan. Depending on your available storage, an upgrade message may appear.



Your cloud storage is full.

[Upgrade Storage](#)

Storage Management – Advanced Settings

Advanced Settings

Folder Path

C:/Users/medit/Downloads/



Apps Installation Folder

C:/Program Files/Medit/



Auto Backup 



Register Offline Data

[Register](#)

Finds and registers offline data on this device that is not registered in Medit Link.

Folder Path	Specify a folder to store the data results generated in Medit Link and connected programs.
-------------	--

Apps Installation Folder	Specify an installation path for Medit Apps and partner apps available in the App Box.
Auto Backup	Enable this option to automatically back up your scanned data to the cloud without using your cloud storage quota.
Register Local Cases	Find and register locally stored cases that are not yet in Medit Link.


Case Converting Tool

You can move cases between different accounts by exporting existing cases of the account or importing cases exported from other accounts.

Export Case

Export cases of the currently logged-in account to the local PC.

1. Open Medit Link and click the Settings icon.

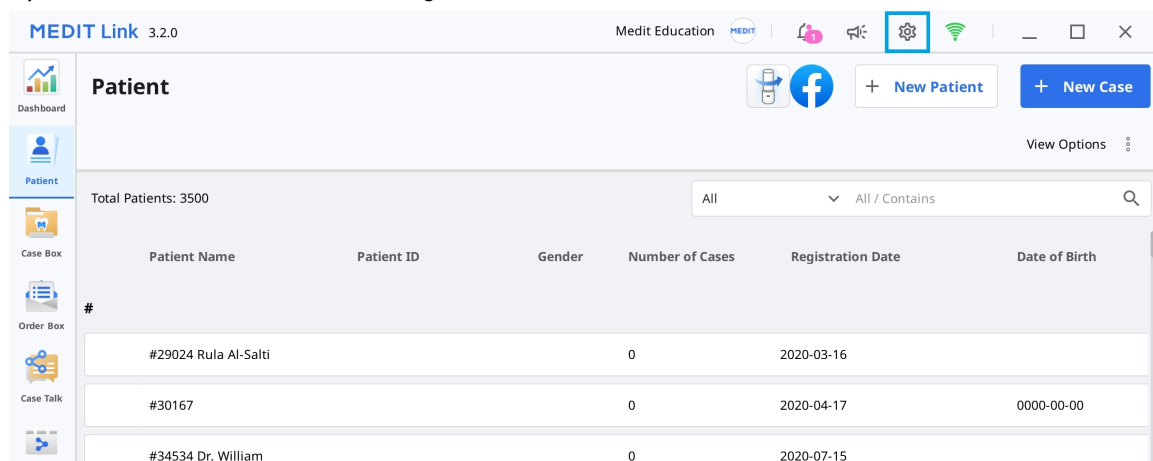

2. Choose the "Case Converting Tool" section in the left-side menu.

Settings

Search settings

General

Advanced

Storage Management

Case Converting Tool

About

User Guide

Case Converting Tool

Case Export

All

All ▾ All / Contains

Files to Export Total: 0 / 932

☐	Case Name	Patient Name	Last Modified Date
☐	0511 test's Case	0511 test	2025-07-25 PM 6:02
☐	marginline_의 케이스	marginline_	2025-07-24 PM 2:34
☐	BS_3's Case - Clone	BS_3	2025-07-24 AM 9:14
☐	test1's Case - Clone	test1	2025-07-23 PM 5:25
☐	test1's Case	test1	2025-07-23 PM 5:19
☐	BS_3's Case	BS_3	2025-07-23 PM 5:12
☐	New Case's Case	New Case	2025-07-23 AM 10:0

Export

3. In the Case Export section, set a period or search for cases to export if necessary.

All

All ▾ All / Contains

All ▾

From

2018-01-01

To

2025-02-21

0 / 767

Patient Name	Last Modified Date
	2025-02-07 PM 3:34
	2025-02-07 PM 3:34

4. Choose cases you want to export from the list in the "Files to Export" section.

Case Export

All

All ▾ All / Contains

Files to Export Total: 932 / 932

☒	Case Name	Patient Name	Last Modified Date
☒	0511 test's Case	0511 test	2025-07-25 PM 6:02
☒	marginline_의 케이스	marginline_	2025-07-24 PM 2:34

- Set the file name by editing the case name in the File Name field.
- Set the destination folder to export the selected cases.
- Click "Export."
- Click "Yes" to open the folder where the file is exported.

Question



Exporting successfully completed.
Do you want to open the folder containing the exported file?


No

Yes

9. The exported cases are saved in .meditex file format.

Import Case

1. Open Medit Link and click the Settings icon.

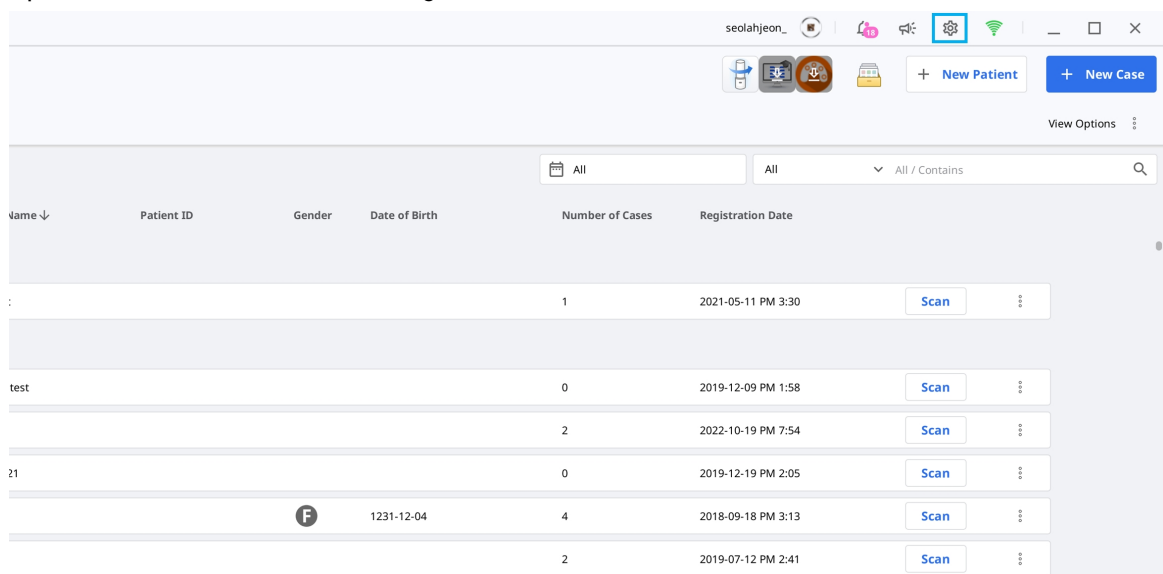

2. Choose the "Case Converting Tool" section in the left-side menu.

Settings

Search settings

General

Advanced

Storage Management

Case Converting Tool

About

User Guide

Case Converting Tool

Case Export

All

All / Contains

Files to Export

Total: 0 / 932

☐	Case Name	Patient Name	Last Modified Date
☐	0511 test's Case	0511 test	2025-07-25 PM 6:02
☐	marginline_의 케이 스	marginline_	2025-07-24 PM 2:34
☐	BS_3's Case - Clone	BS_3	2025-07-24 AM 9:14
☐	test1's Case - Clone	test1	2025-07-23 PM 5:25
☐	test1's Case	test1	2025-07-23 PM 5:19
☐	BS_3's Case	BS_3	2025-07-23 PM 5:12
☐	New Case's Case	New Case	2025-07-23 AM 10:0

Export

- In the Case Import section, click the folder icon to browse your PC for the .meditex file you want to import.

Case Import

A new case is created for every imported case. It contains case data, patient information, and scan file. However, it does not contain order information.

C:/Users/medit/Downloads

- Click the "Import" button and check the files to import in the list.
- The imported case will also appear in your Case Box (Clinic account) or Work Box (Lab account).

About

Show detailed information about Medit Link.

About

MEDIT Link

Product Name: Medit Link

Version: 3.5.0

Build Number: 101

Server: Stage Korea

UDI DI: (01)8800026700449

UDI PI: (10)3.5.0

 **Medit Corp.**

9F, 10F, 13F, 14F, 16F, 8, Yangpyeong-ro 25-gil, Yeongdeungpo-gu,
Seoul, 07207, Republic of Korea
Tel: +82-2-2193-9600

 2026.08

[Terms & Conditions](#)

[Open Source License](#)



EC REP Meditrial Srl
Via Po 9 00198, Rome Italy

CH REP Meditrial Europe Ltd
Bahnhofstrasse 23 6300 Zug, Switzerland

UK REP Meditrial UK Ltd
19 The Circle, Queen Elizabeth Street, London, SE1 2JE, United Kingdom

☒ Send Anonymous Usage Statistics

Close

User Guide


Get the user guide for Medit Link.

Opening the user guide requires an Internet connection.

Menu

Dashboard

You can view the statistics for the entire work done on Medit Link in visualized graphs and charts – both in the Medit Link App and Web.

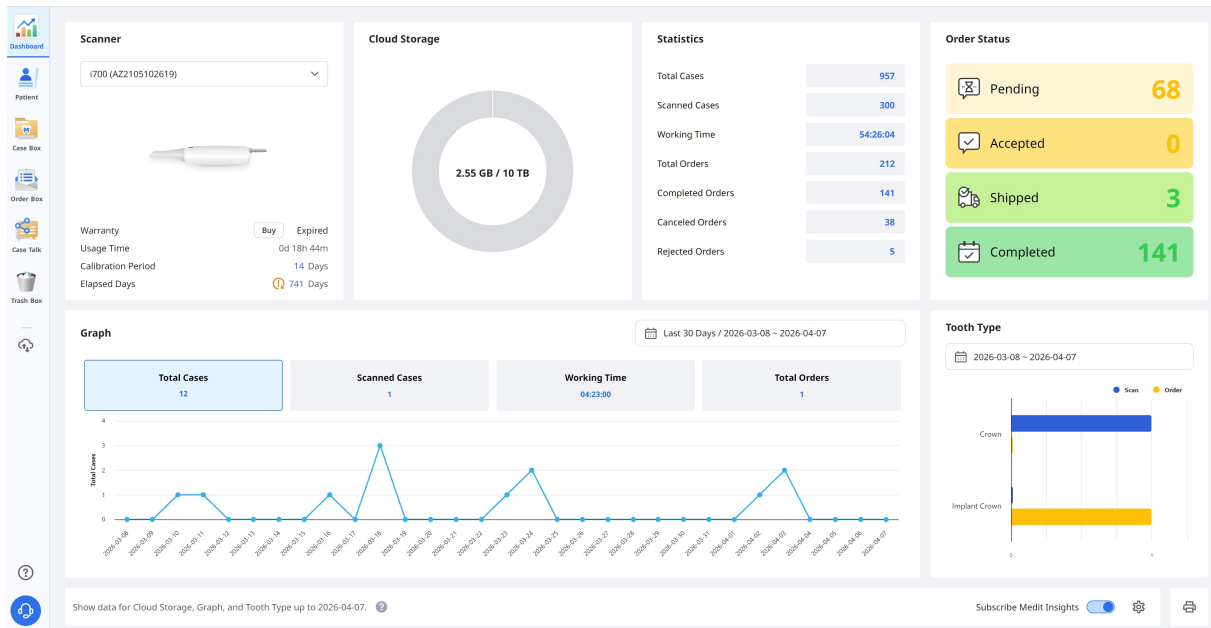


Note

Please refer to the [Medit Link Web > Dashboard](#) for more information on the dashboard and Medit Insights.

Patient

You can add, edit, delete, or merge patient records in the Patient menu.

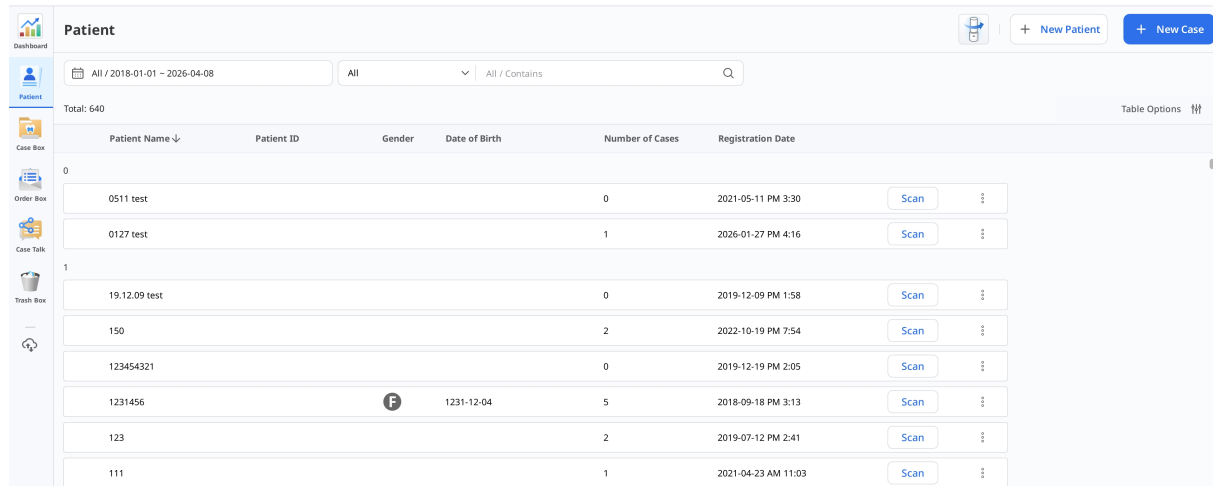


The screenshot shows the 'Patient' menu in a software application. At the top, there's a header with 'Patient' and buttons for '+ New Patient' and '+ New Case'. Below the header, there's a search bar with a date range 'All / 2018-01-01 - 2026-04-08' and a filter dropdown set to 'All'. A 'Total: 640' is displayed. The main area contains a table with columns: Patient Name, Patient ID, Gender, Date of Birth, Number of Cases, and Registration Date. The table lists several patient records, each with a 'Scan' button and a vertical ellipsis menu icon. On the left side, there's a sidebar with icons for Dashboard, Patient, Case Box, Order Box, Case Task, and Trash Box.

Patient Name	Patient ID	Gender	Date of Birth	Number of Cases	Registration Date
0511 test				0	2021-05-11 PM 3:30
0127 test				1	2026-01-27 PM 4:16
19.12.09 test				0	2019-12-09 PM 1:58
150				2	2022-10-19 PM 7:54
123454321				0	2019-12-19 PM 2:05
1231456			1231-12-04	5	2018-09-18 PM 3:13
123				2	2019-07-12 PM 2:41
111				1	2021-04-23 AM 11:03

Add New Patient

1. Click "New Patient."

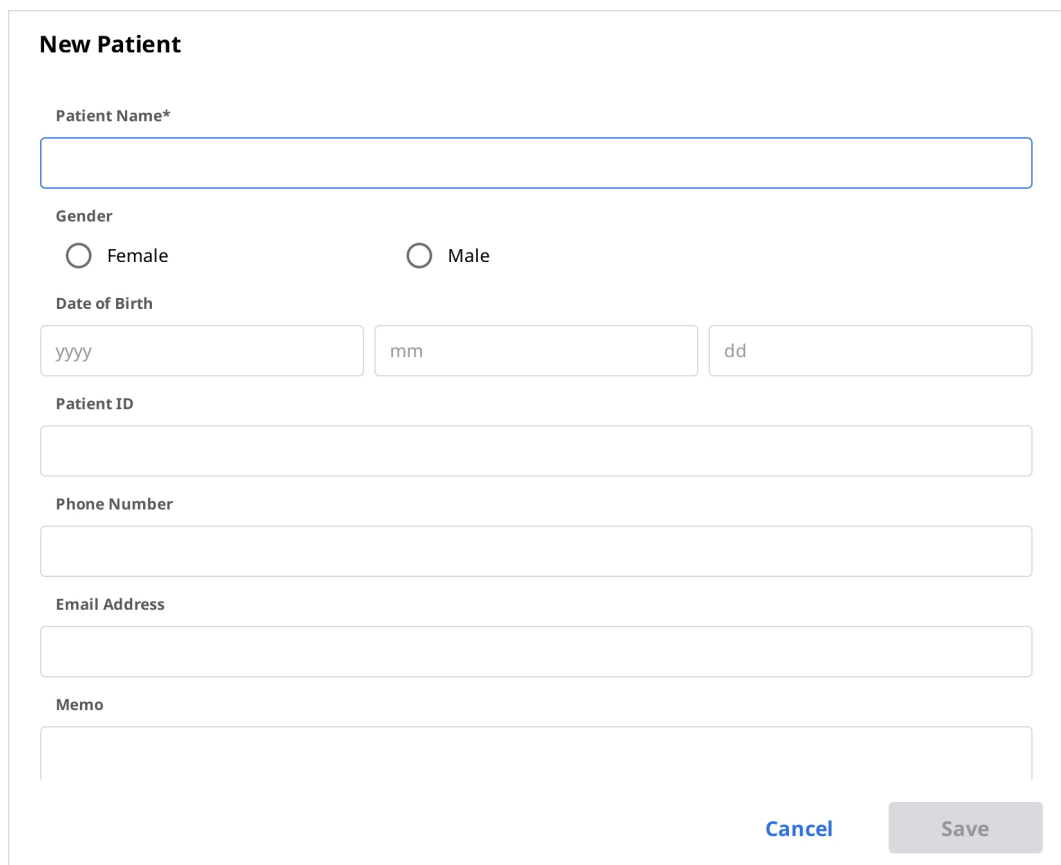


The 'New Patient' dialog form contains the following fields and options:

- Patient Name***: A text input field.
- Gender**: Radio buttons for ☐ Female and ☐ Male.
- Date of Birth**: Three input fields for year (yyyy), month (mm), and day (dd).
- Patient ID**: A text input field.
- Phone Number**: A text input field.
- Email Address**: A text input field.
- Memo**: A text area.
- Buttons**: 'Cancel' and 'Save' buttons at the bottom right.

2. Fill in the required information in the New Patient dialog, then click "Save."

Delete Patient

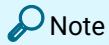

Note

You can delete a patient only when there are no cases for the patient.

1. Select one or more patients from the list.

	Patient Name ↓	Patient ID	Gender	Date of Birth	Number of Cases
☐	BS4				1
☒	BLATRIX Charline	d82cb01f1e1d4cc7b4bc		1986-01-05	1
☒	BLATRIX Charline	d82cb01f1e1d4cc7b4bc		1986-01-05	3

2. Click "Delete."



Merge Patients

You can merge all cases of multiple patients into one patient.

1. Select the patients you want to merge from the list.
2. Click "Merge Patients."



3. Select the target patient you would like to merge all cases into.
4. All cases are merged into the target patient and the patient information of the target patient maintains.

Merge all patient information to


Please select the patient you'd like to merge all data to.
After merging all data together, all information in other patient files will be deleted.

	Patient Name	Patient ID	Gender	Birth Date	Number of Cases	Date Registered	Recent Visit Date
☒	5-9_iScan_CAMOrd...				1	2020-09-18	2020-09-18
☐	test				2	2020-09-17	2020-09-17
☐	test				22	2019-07-26	2020-09-17


Cancel Apply

5. Click "Apply" to proceed with the merge.

Case Box/Work Box

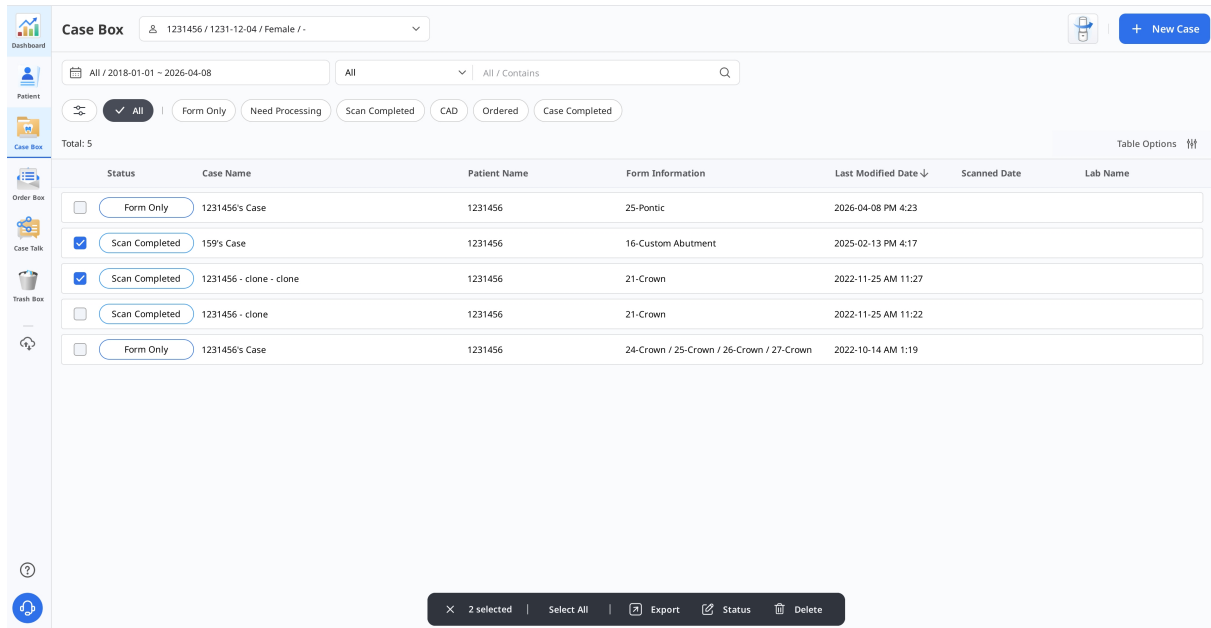

The action bar appears at the bottom of the screen when one or more cases are selected.

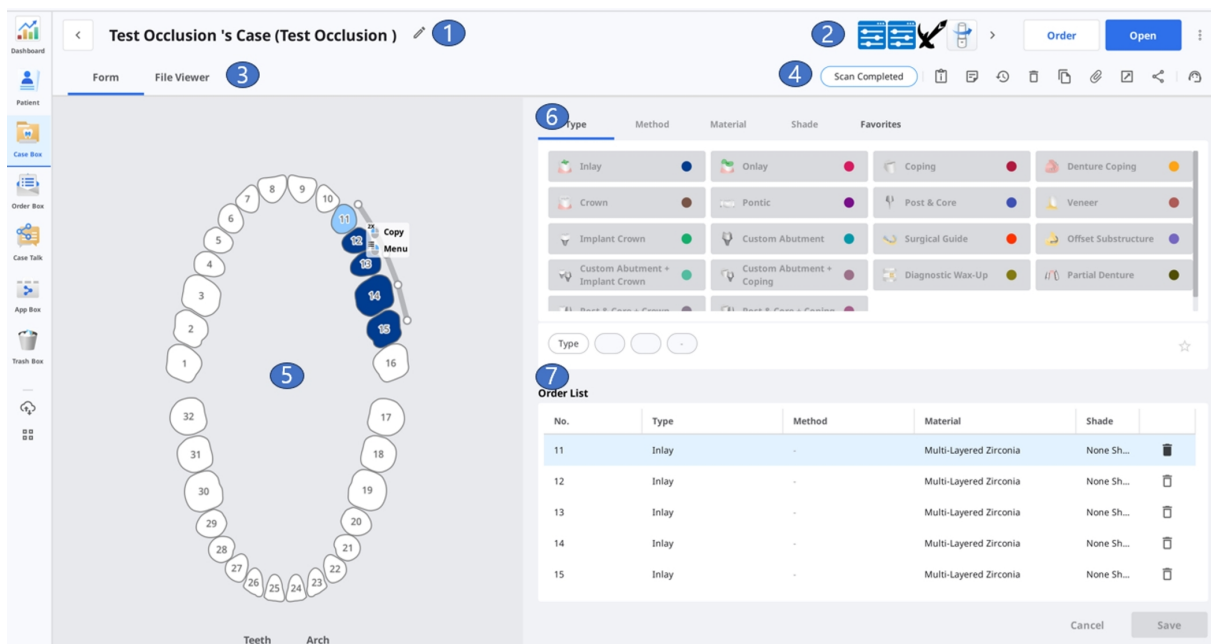

1	Show case name and patient name.
2	Provide key features available: Scan or Open button, Order button, app icons.
3	Show available tabs: Form, File Viewer.
4	Show functions for managing cases and files.

5	Show the form registration UI for the selected scanning program.
6	Provide the type/method/material/shade information for users to choose for the selected tooth number.
7	Show the type/method/material/shade information applied to the selected teeth.

Tabs in Case Box/Work Box

Form	File Viewer
The form tab provides a UI for entering basic information for scanning work, such as type, method, material for teeth numbers, and scan options. The information you enter is passed to the connected scanning program and used as order information when placing an order with your partner.	View result files generated by scanning programs and CAD. The file viewer displays and manages the acquired or attached 3D models and 2D images.

Functions in Case Box/Work Box

The following functions are provided in Case Box/Work Box.

	Case Status
	Go to Order Box/In Box
	Case Information
	Memo & Tags
	Case History

	Delete
	Clone
	Attach
	Export
	Share
	Submit Support Request

Shows the status of the case.

Provides a shortcut for the same case in another box (Case Box ↔ Order Box; In Box ↔ Work Box).

Review case details and save them as a PDF or print them out.

Add a memo or tags when registering or ordering a case. Use them to communicate any additional or supplementary information about the case to your partners.

Provides detailed information about the case history.

History

Last Modified Date	3/29/2022 1:57 PM
Scanning Date	4/12/2021 2:06 PM
Creation Date	3/22/2022 1:17 PM
Case Size	435MB

[Close](#)

Delete the case.

You cannot delete ordered cases or cases that were created offline.

Create a duplicate of the current case in Case Box. All data will be copied to the cloned case. This feature helps replicate treatment information for the same patient or make an additional order with the same information or files.

Attach files in various formats.

Export scan or CAD result files to the local PC. You can choose your preferred file format when exporting files.

The exported files can be used in CAD and CAM software for further work.

Create a link for sharing.

Users with the shared link can access the 3D data in a browser such as Chrome, Firefox, or Safari without installing additional programs.

You can share the link freely with anybody or specific people only. A PIN code or a link expiration date can be set to secure private information.

Land on a Medit Help Center page to submit the support request.

Dashboard

Patient

Case Box

Order Box

Case Talk

Trash Box

Order Box

All / 2018-01-01 ~ 2026-04-08

All

All / Contains

All

Pending

Canceled

Rejected

Accepted

Shipped

Case Completed

Total: 254

Table Options

	Status	Case Name	Patient Name	Order Date	Requested Delivery Date	Lab Name	Order ID
☐	Accepted	texture on's Case - Clone	texture on	2026-04-08 PM 5:40		MEDIT_LAB_STAGE	6012449
☒	Pending	Express Scan	Seol register	2026-04-02 PM 9:38	2026-04-09 AM 9:00	MEDIT_LAB_STAGE	6012242
☒	Case Completed	전설이의 케이스	전설이	2025-07-16 PM 5:46	2025-07-17 PM 3:00	MEDIT_LAB_STAGE	6004719
☒	Case Completed	27	1	2025-03-16 PM 3:11	2025-03-17 PM 3:00	MEDIT_LAB_STAGE	6002656
☐	Case Completed	3.3.0.136 TEST's Case	3.3.0.136 TEST	2024-04-12 PM 3:25	2024-04-13 PM 2:00	MEDIT_LAB_STAGE	5998987
☐	Case Completed	test의 케이스	Patient TEST	2024-04-09 PM 3:06	2024-04-10 PM 2:00	MEDIT_LAB_STAGE	5998955
☐	Case Completed	Mandibular Movement Case - Clone - Clone	MEDIT	2023-06-12 PM 12:12	2023-06-13 PM 2:00	MEDIT_LAB_STAGE	5996680
☐	Case Completed	Release New_2's Case	Release New_2	2023-04-25 PM 6:47	2023-04-25 PM 2:00	MEDIT_LAB_STAGE	5996180
☐	Case Completed	cX 문익's Case - Clone	cX 문익	2023-05-16 PM 5:26	2023-05-17 PM 2:00	MEDIT_LAB_STAGE	5996393
☐	Case Completed	cX 문익's Case	cX 문익	2023-05-15 PM 2:39	2023-05-16 PM 2:00	MEDIT_LAB_STAGE	5996362
☐	Case Completed	310.187's Case	310.187	2023-05-15 PM 1:33	2023-05-16 PM 2:00	MEDIT_LAB_STAGE	5996360
☐	Case Completed	310.187 test's Case - Clone			2023-05-03 PM 2:00	MEDIT_LAB_STAGE	5996306

3 selected

Select All

Export

Delete

1	Show case name and patient name.
2	Provide functions for order management (accepting/canceling; delivering/receiving).
3	Show available tabs: Order, File Viewer
4	Show functions for managing cases and files: case status, shortcuts for detailed information, case information, memo & tags, etc.

5	Show the detailed information for the order including the order number, patient name, partner name, scanning date, order date, delivery date, and memo & tags.
6	Show the ordered product information and delivery status.
7	Show the assignee for the case.


Tabs in Order Box/In Box

Order	File Viewer
Check the detailed information related to the order, such as order number, assignee, patient name, account name, scanning date, order date, delivery date, memo, and tags.	View files generated by scan/CAD operations. The file viewer displays and manages the acquired or attached 3D models and 2D images.

Case Talk

Cases created in the Medit Link can be shared with specific people, anyone with the shared link or registered members of Medit Link.

The user can see a list of all shared cases through Case Talk in the Medit Link App and Medit Link Web.

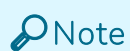

Case Talk service requires an internet connection.

Dashboard

Patient

Case Box

Order Box

Case Talk

App Box

Trash Box

Case Talk

All / 2018-01-01 ~ 2026-04-10

Case Talk Name

Enter keyword

Total: 9

Case Talk Name	Patient Name	Form Information	Scanned Date	Shared Date	Expiration Date	
20260313_new's Case - Clone_2026-03-19 14:36:08	20260313_new222	11	2026-03-13 15:27	2026-03-19 14:36	2026-06-17 23:59	
20260313_new's Case - Clone_2026-03-19 14:34:48	20260313_new222	11	2026-03-13 15:27	2026-03-19 14:34	2026-06-17 23:59	
크라온 보릿지 44-46's Case_2026-03-13 11:05:36	크라온 보릿지 44-46	44, 46 45	2020-10-20 11:46	2026-03-13 11:05	2026-06-11 23:59	
크라온 보릿지 44-46's Case_2026-03-13 11:05:12	크라온 보릿지 44-46	44, 46 45	2020-10-20 11:46	2026-03-13 11:05	2026-06-11 23:59	
크라온 보릿지 44-46's Case_2026-03-13 10:57:30	크라온 보릿지 44-46	44, 46 45	2020-10-20 11:46	2026-03-13 10:57	2026-06-11 23:59	
크라온 보릿지 44-46's Case_2026-03-13 10:53:08	크라온 보릿지 44-46	44, 46 45	2020-10-20 11:46	2026-03-13 10:53	2026-06-11 23:59	
크라온 보릿지 44-46's Case_2026-03-13 10:49:01	크라온 보릿지 44-46	44, 46 45	2020-10-20 11:46	2026-03-13 10:49	2026-06-11 23:59	
ORG_[Clinics][3.3.6] 01 Normal - Complete_2026-03-12 02:00:00	ORG_[Clinics][3.3.6] 01 Normal	16 11 26 31, 32	2024-12-23 15:11	2026-03-12 18:00	2026-06-11 15:59	
크라온 보릿지 44-46's Case_2026-03-12 17:19:49	크라온 보릿지 44-46	44, 46 45	2020-10-20 11:46	2026-03-12 17:19	2026-06-10 23:59	
20260313_new's Case - Clone_2026-03-19 14:36:08	20260313_new222	11	2026-03-13 15:27	2026-03-19 14:36	2026-06-17 23:59	

1

Shared cases can be viewed via a web browser without installing a separate program, enabling easy communication over the 3D data on various devices.

App Box

You can browse, install, and manage apps by category, including Imaging, CAD/CAM, Diagnosis/Consultation, and Management.

Dashboard

Patient

Case Box

Order Box

Case Talk

App Box

Trash Box

Help

App Box

My Apps

Search by App Information

Medit Apps Only

All

Management

Imaging

Diagnosis/Consultation

CAD

CAM

Utilities

Order Placement

Featured

Medit : Transforming Digital Dentistry

Featured on CBS News

Watch the Video >

The Doctors' Choice

Medit Intraoral Scanner Trusted by the Majority

Click to find out more >

Featured Apps (Top 5)

Diagnosis/Consultation

Order Placement

Medit Orthodontic Suite

Install

Diagnosis/Consultation

CAM

CAD

Medit ClinicCAD

CAD solution for creating a wide range of dental prostheses

Install

Time to Unwrap!

Diagnosis/Consultation

Medit AuraVue

Management

Andcom Integration

Management

Open Dental Integration

Utilities

1114_Standard Apps

Management

국가 제함 App TEST


Medit Apps Only



Use this toggle to display only Medit Apps. When enabled, partner apps are hidden.

Sorting Options


Use the sorting options to change how apps are displayed. The app list and section titles update accordingly.

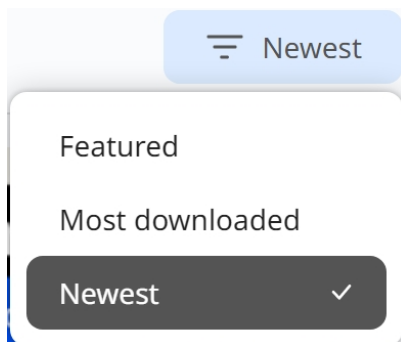


- **Featured:** Displays selected apps recommended by Medit.
- **Most downloaded:** Displays apps based on download popularity.
- **Newest:** Displays recently added apps.

My Apps

Click "My Apps" to access a dedicated page where you can manage your apps.

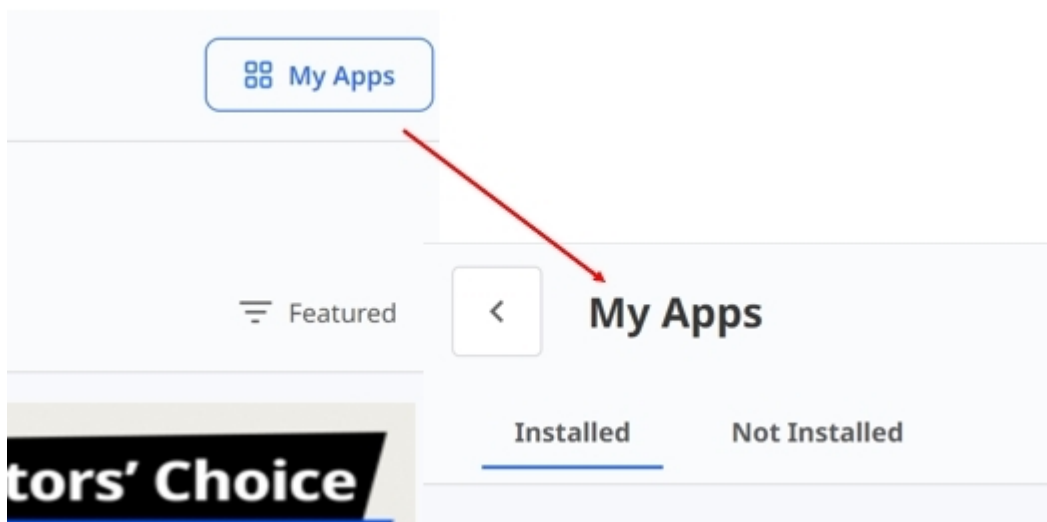

- **Installed:** Displays and updates apps currently installed on your system.
- **Update All:** Click "Update all" to update all installed apps with available updates at once.

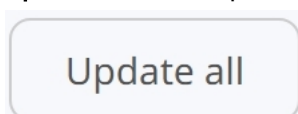

- **Not Installed:** Displays apps that are available but not yet installed.

Trash Box

The Trash Box temporarily stores the deleted cases. Any case in the Trash Box can be either restored or permanently deleted.

⚠️ Caution

The data deleted from the Trash Box will be permanently deleted from your PC and the server and cannot be recovered.

Status	Case Name	Patient Name	Lab Name	Deletion Date ↓	Deleted By	Size
☐ Scan Completed	Jira_임플란트 46, 47's Case - Clone - Clone	임플란트 46, 47		2023-04-26 오전 9:12	wonhoon choi	66,654 KB
☐ Form Only	G11598492304180810429	test test		2023-04-18 오후 5:15	wonhoon choi	0 KB
☐ Form Only	G11598492304180806313	test test		2023-04-18 오후 5:15	wonhoon choi	0 KB
☐ Form Only	G11598492304180811349	tes test		2023-04-18 오후 5:15	wonhoon choi	0 KB
☐ Form Only	G11598492304180754083	Gildewell_fastdesign		2023-04-18 오후 5:00	wonhoon choi	0 KB
☐ Form Only	G11598492304180733332	Gildewell_fastdesign		2023-04-18 오후 5:00	wonhoon choi	0 KB
☐ Form Only	G11598492304180648146	test test		2023-04-18 오후 5:00	wonhoon choi	0 KB
☐ Case Completed	Multi Occlusion Normal Case's Case - Clone	Multi Occlusion Normal Case		2023-04-18 오전 10:50	wonhoon choi	23,652 KB
☐ CAD	Multi Occlusion Normal Case's Case - Clone - Clone	Multi Occlusion Normal Case		2023-04-18 오전 10:50	wonhoon choi	23,652 KB
☐ Form Only	Occlusion Test's Case	Occlusion Test		2023-04-17 오후 6:31	wonhoon choi	0 KB
☐ Form Only	Test Occlusion 's Case	Test Occlusion		2023-04-17 오후 6:31	wonhoon choi	0 KB
☐ Form Only	Occlusion Test's Case	Occlusion Test		2023-04-17 오후 6:31	wonhoon choi	0 KB
☐ Form Only	Occlusion Test's Case	Occlusion Test		2023-04-17 오후 6:26	wonhoon choi	0 KB

Case Synchronization Manager

The Case Synchronization Manager helps you to upload and download cases efficiently. You can upload multiple cases to the cloud or download them to your local PC. To identify cases requiring synchronization, simply search for them or apply date filters.

Case Synchronization Manager

Upload
Download

All

Case Name
 Case Name / Contains

Total Cases: 210

☒ No Date
 Data Type
 All

☐	Case Name	Patient Name	Scanned Date↓	Size	Status
☐	BLATRIX Charline's Case	BLATRIX Charline	2025-01-06 PM 6:28	729.30 MB	Upload Required
☐	BLATRIX Charline's Case	BLATRIX Charline	2025-01-06 PM 6:28	729.30 MB	Upload Required
☐	BLATRIX Charline's Case - 복제	BLATRIX Charline	2025-01-06 PM 6:28	733.35 MB	Upload Required
☐	BLATRIX Charline's Case - Clone	BLATRIX Charline	2025-01-06 PM 6:28	734.01 MB	Upload Required
☐	test의 케이스	Patient TEST	2023-04-26 PM 7:58	1.25 GB	Upload Required
☐	i700 HD's Case	i700 HD	2023-04-24 PM 5:23	911.24 MB	Upload Required
☐	Mandibular Movement	MEDIT	2023-04-24 PM 1:59	1.88 GB	Upload Required
☐	Mandibular Movement Case	MEDIT	2023-04-24 PM 1:59	1.75 GB	Upload Required
☐	Mandibular Movement Case - Clo...	MEDIT	2023-04-24 PM 1:59	1.75 GB	Upload Required
☐	서라_Mandibular Movement_TEST	Mandibular Movement_export_test	2023-04-10 PM 3:39	34.73 MB	Upload Required
☐	서라_Mandibular Movement_TEST...	Mandibular Movement_export_test	2023-04-10 PM 3:39	731.43 MB	Upload Required

Close

Upload


Upload	Upload the selected cases.
Cancel	Cancel the upload for selected cases.
Retry	Retry uploading cases if the upload fails.

Download

 Download	Download the selected cases.
 Cancel	Cancel the download for selected cases.
 Retry	Retry downloading cases if the download fails.

Help & Resources

Click the Help & Resources icon to open the menu, then select a resource.

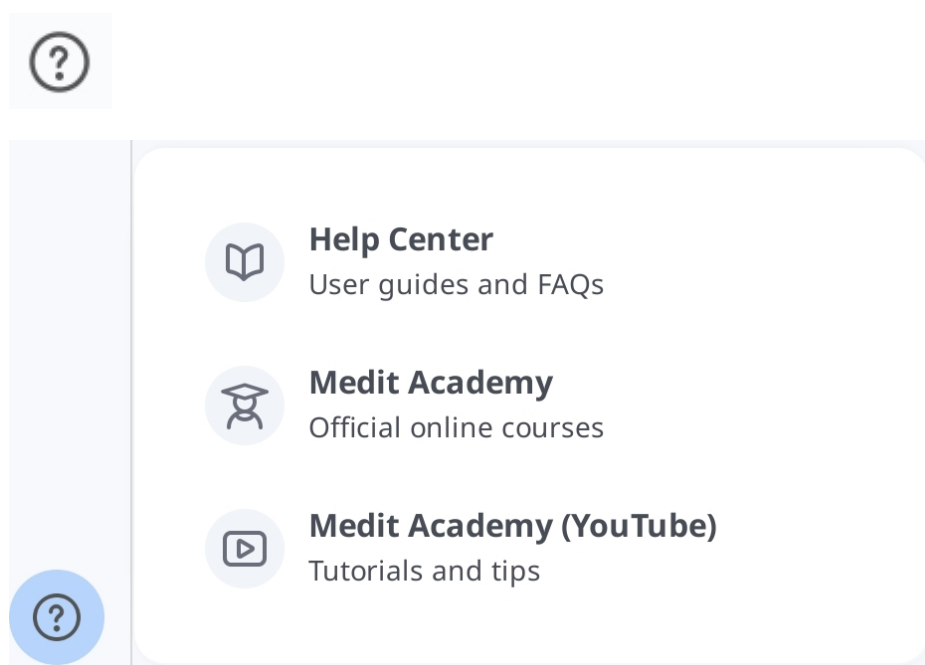

File Viewer

File Viewer groups and shows various data, including the data acquired from programs linked to Medit Link and data attached by users.

It provides basic management functions such as showing/hiding data, adjusting opacity, and deleting data.


It supports previews of files in the following formats:


- 3D data: STL, OBJ, PLY
- 2D images: PNG, JPG, JPEG, BMP, GIF, TIF, TIFF
- Video files: MKV, MP4, AVI, WMV

Data Tree

You can group and display all data for the case in the data tree. You can also show/hide data, adjust transparency, and delete data from the data tree.

	Show/Hide	You can show or hide the data.
	Opacity	The slide bar appears when you click one of the data on the data tree. You can adjust the opacity of the shown data by moving the slider bar. 

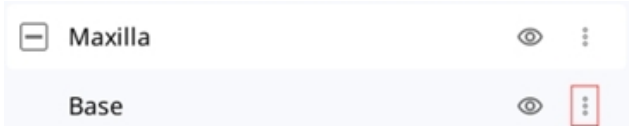	Opacity	The slide bar appears when you click one of the data on the data tree. You can adjust the opacity of the shown data by moving the slider bar. 
	More	Clicking the More icon shows additional features available.  • Show This Only: Hide all other data but display only the corresponding data. • Hide: Hide the data. • Duplicate: Copy the data and add it to the attachment group. • Export: Export the selected files. • Rename: Change the name of the data. • Delete: Delete the data.

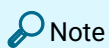

Note

You can rename data when online.

Model Control Tools

3D Model Controls

	Pan
---	-----

	Rotate	
	Zoom In/Out	
	Zoom Fit	
	Data Display Mode	
		Matte
		Monochrome
	Grid Settings (mm)	
		Grid Off



Overlay On

Move the model.

Rotate the model.

Zoom in and out on the model.


Position the model in the center of the screen.

Glossy

See the data in color with matte finish (Phong rendering).




See the data in one color.

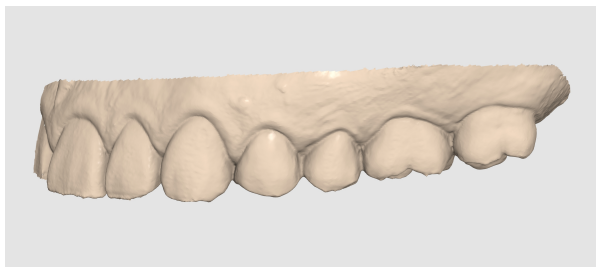


Grid On

Hide the grid in the background.

Overlay the grid over the model.

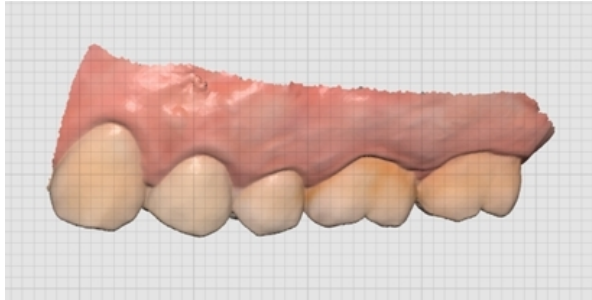

2D Image Control

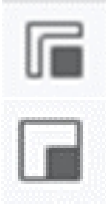
	Next/Previous Image	Move among images.
	Edit	Edit the image.
	Pinch/Paste Window	You can detach the image window and adjust it to the desired position. You can resize the window. When the window is restored to its original location, it returns to its default size.

Function Tools

Image Capture and Memo


Capture the current screen as a 2D image, draw on it, or add notes.

You can save it as a file with your notes and send it as a reference.

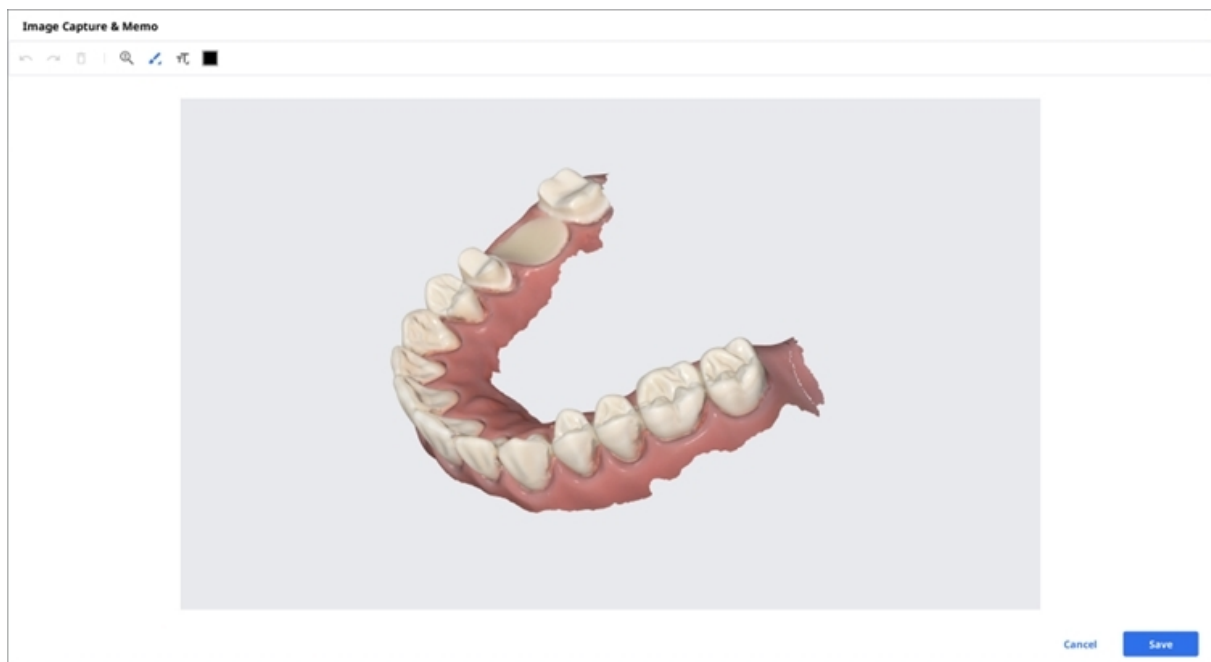


Ortho Capture

Capture the images in the current viewing direction and top/bottom/left/right from the current viewing direction.

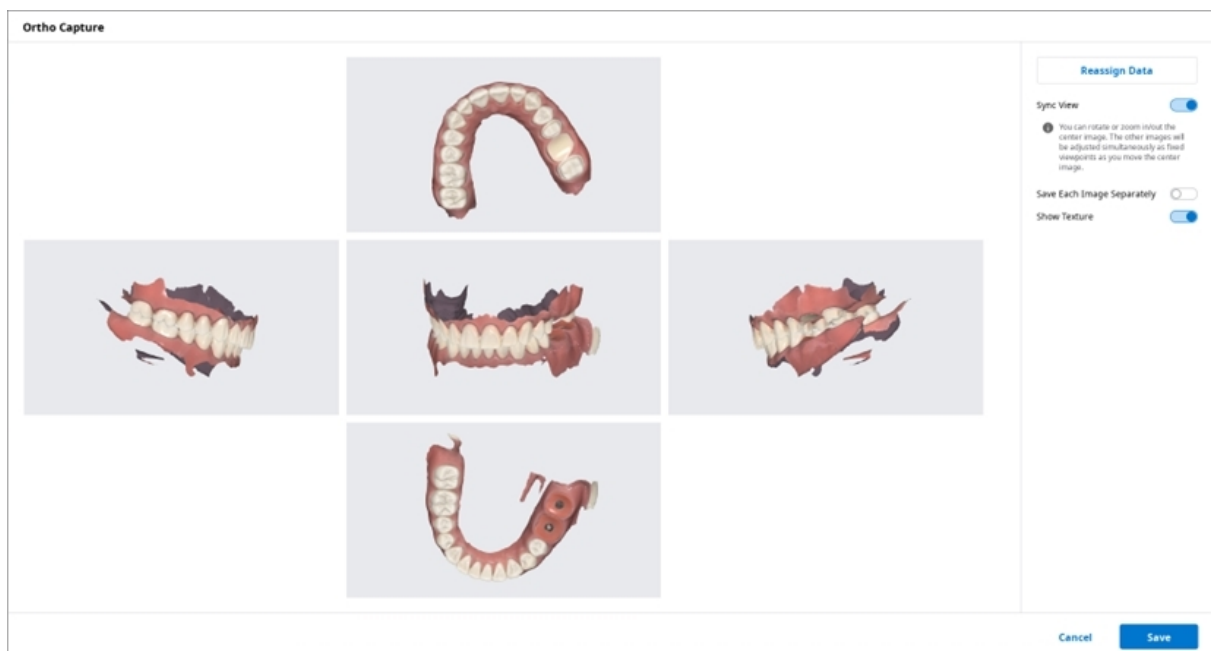


Sync View	When on, you can rotate or zoom in/out the center image. The other images will be adjusted simultaneously as fixed viewpoints as you move the center image. When off, you can pan, rotate, or zoom in/out all images individually.
-----------	--

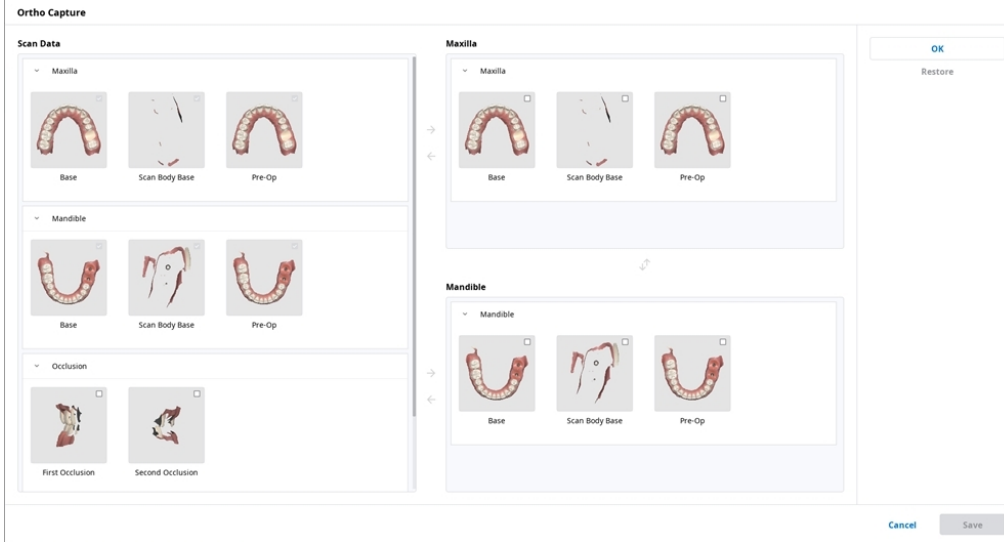Save Each Image Separately	Save each image in high resolution.
Show Texture	Turn on/off the texture for the model.
Reassign Data	Change the data for maxilla/mandible with the Reassign Data feature. You can use the changed data to create and utilize the result images of Ortho Capture.  • Scan Data: Scan data for the case are all shown. • Maxilla: The data you selected for the maxilla are displayed. • Mandible: The data you selected for the mandible are displayed.


Fill Holes

Fill holes in 3D mesh data.

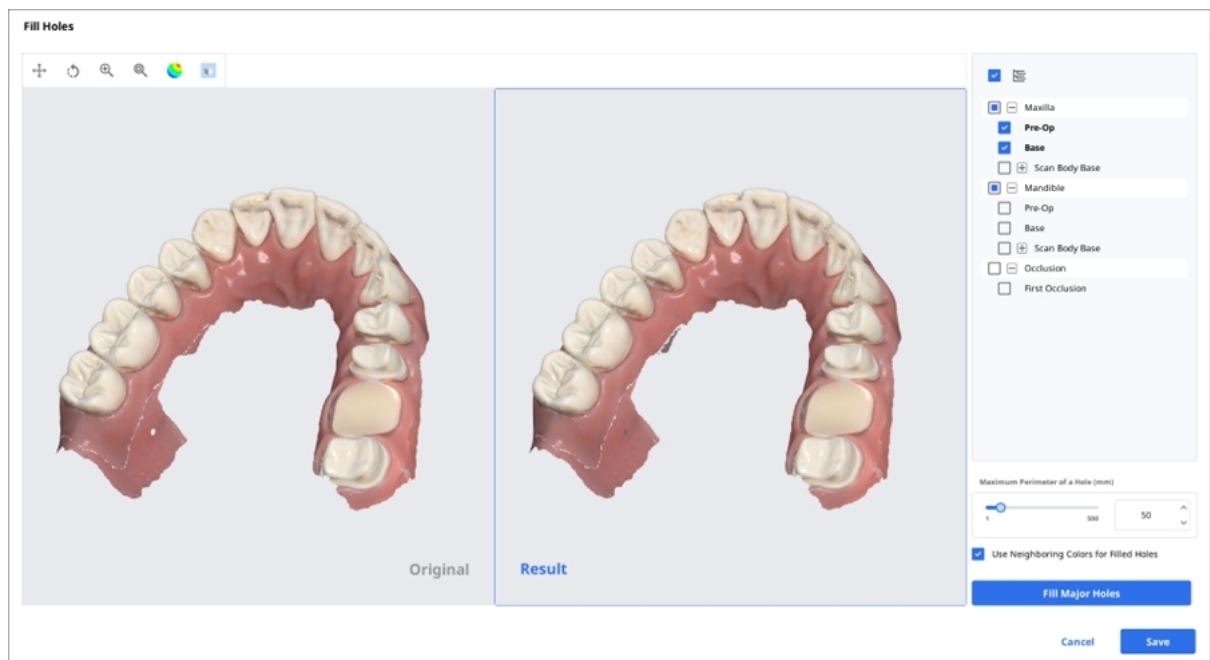



Trimming

Trimming provides various tools to edit data.

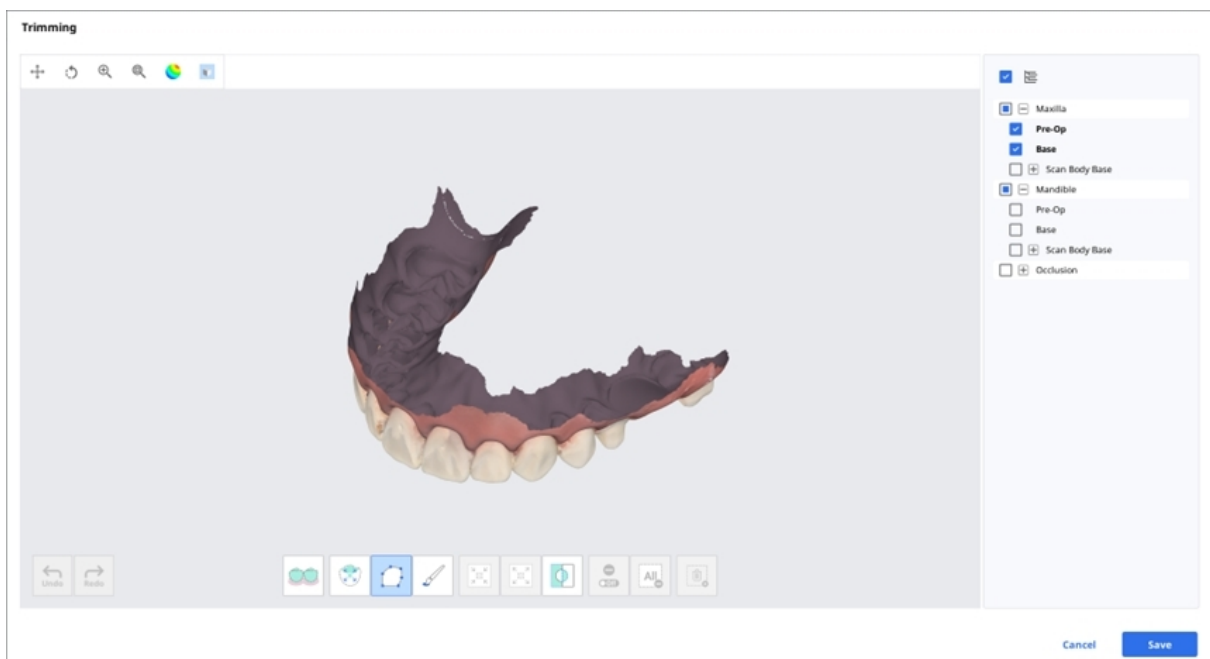

The selected 3D mesh data can be edited using different tools located at the bottom of the window.

Selection		Smart Teeth Selection	
-----------	---	-----------------------	--

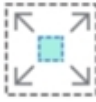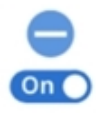		Smart Single Tooth Selection	Automatically select the area of a single tooth with a click. You can click or drag on the tooth.
		Polyline Selection	Select all entities within a polyline shape drawn on the screen.
		Brush Selection	Select all entities on a freehand-drawn path on the screen. Only the front face will be selected. The brush comes in three different sizes.
Adjustment of Selected Area		Shrink Selected Area	
		Expand Selected Area	Expand the selected area each time you press the button.
		Invert Selected Area	Invert the selection.
Deselection		Selection/Deselection	
		Clear All Selection	Clear all selected areas.

		Delete Selected Area	Delete the data from the selected area.
--	---	----------------------	---

Automatically select all teeth of the arch, leaving out gingiva parts.


Note

This function is only available for the scan data acquired using Medit Scan for Clinics with the "Use GPU" option on.

Reduce the selected area each time you press the button.

When on, deselect the area using various tools.

1. Select the area you want to trim with the area selection tools.

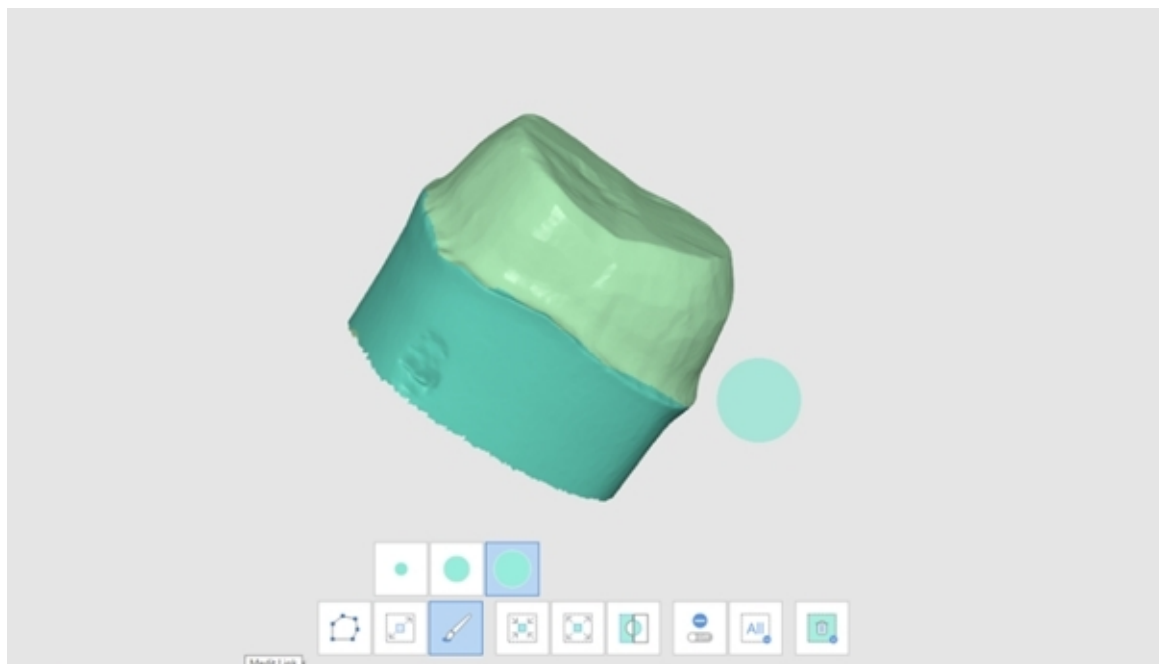


Note


The polyline trimming tool selects all entities within a polyline shape drawn on the screen, while the brush trimming tool selects only the front face data.

2. You can invert, shrink, or expand the selection with the following three tools.


- Invert Selected Area

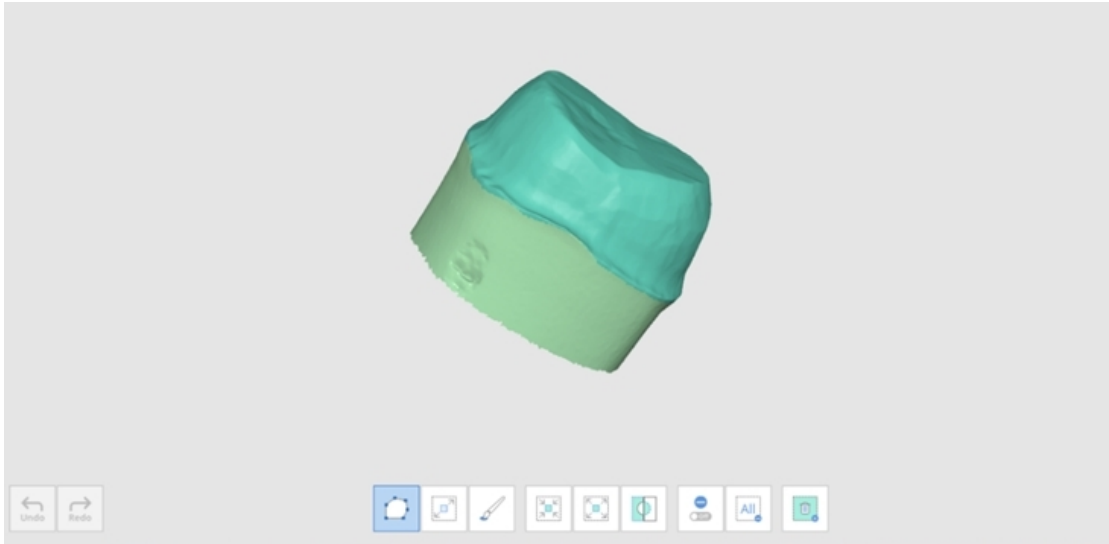

- Shrink Selected Area

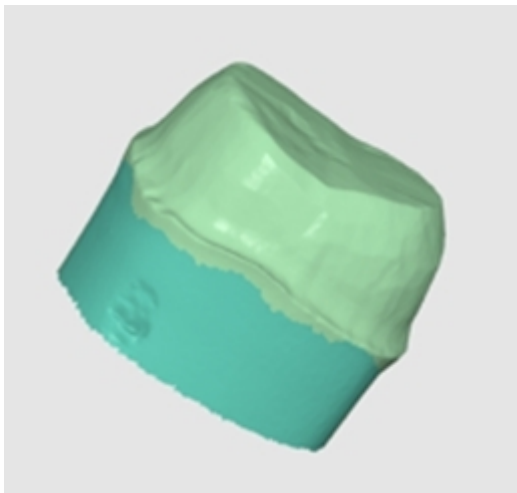


- Expand Selected Area

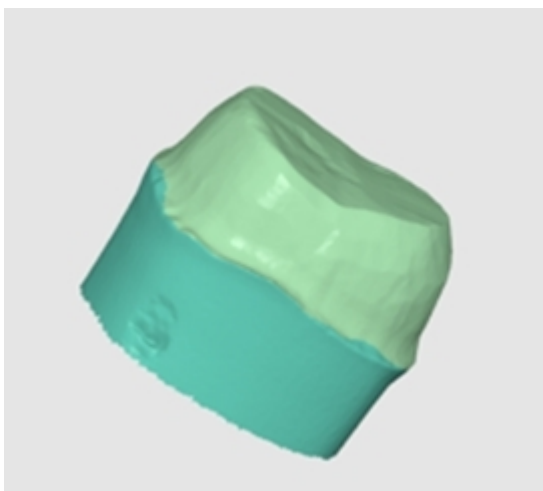

3. You can deselect the selected areas with the Delete Selected Area tool.




Optimize Occlusion Alignment

Optimize the occlusion alignment. You can choose the target data and adjust the occlusion alignment.

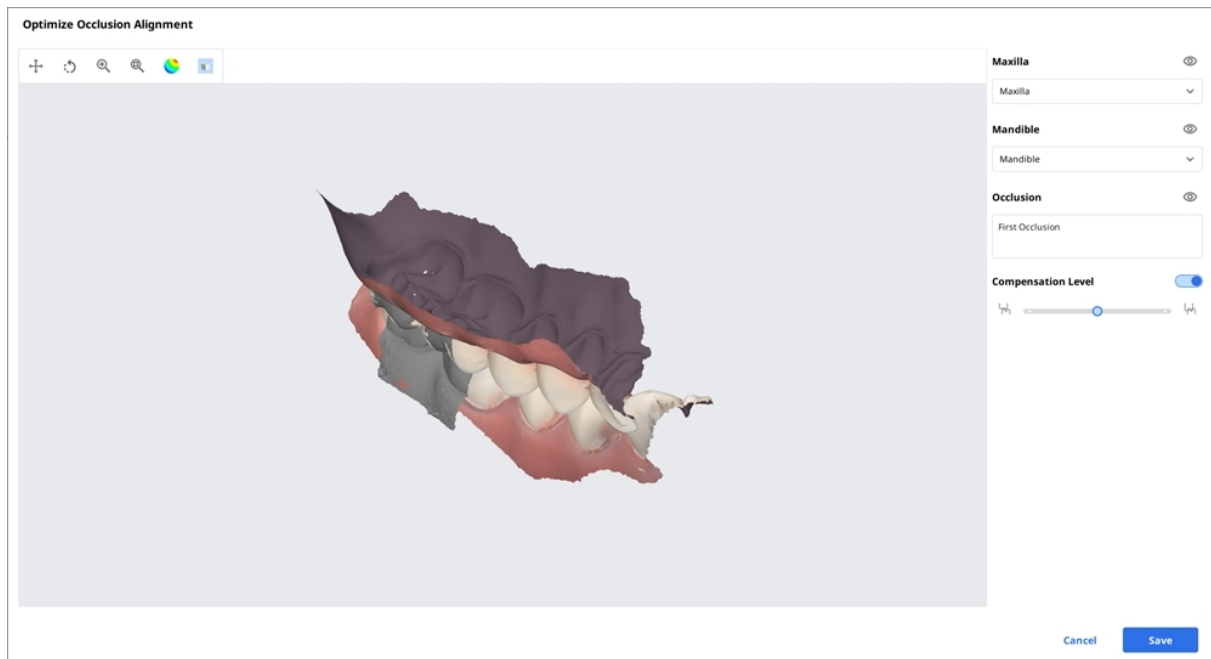

Compensation Level


Low	Compensation Level  
Mid	Compensation Level  
High	Compensation Level  
Off	Compensation Level  

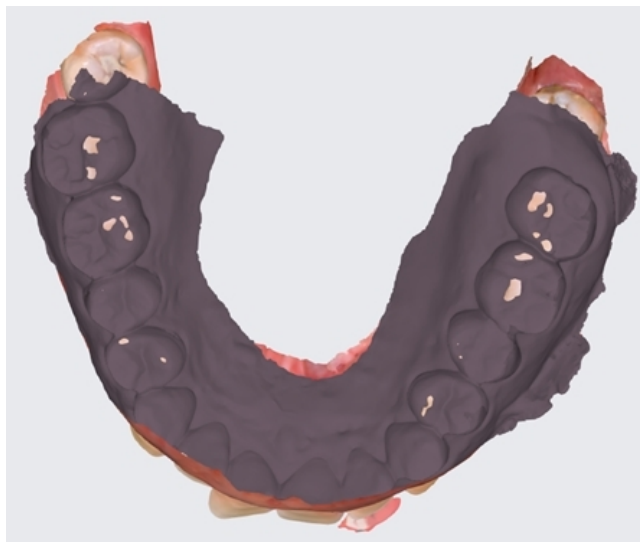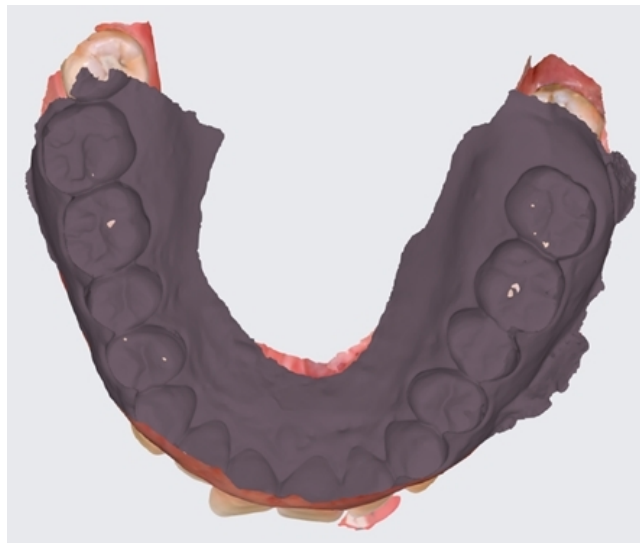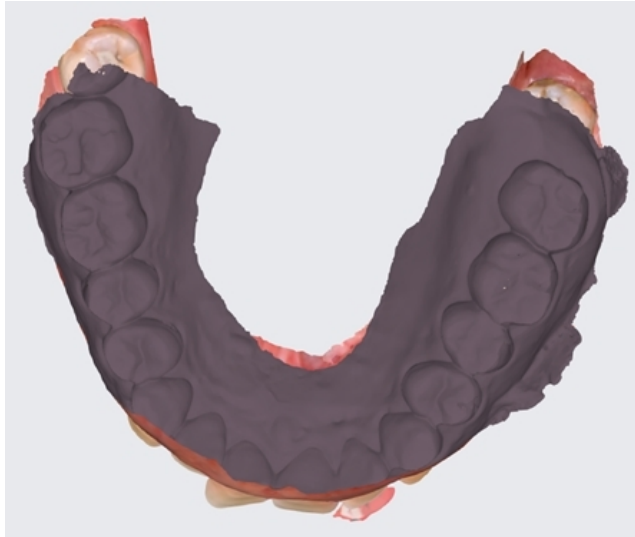



Mandibular Movement

You can play the mandibular movements recorded in Medit Scan for Clinics to see the dynamic occlusal contacts between the upper and lower dentitions.

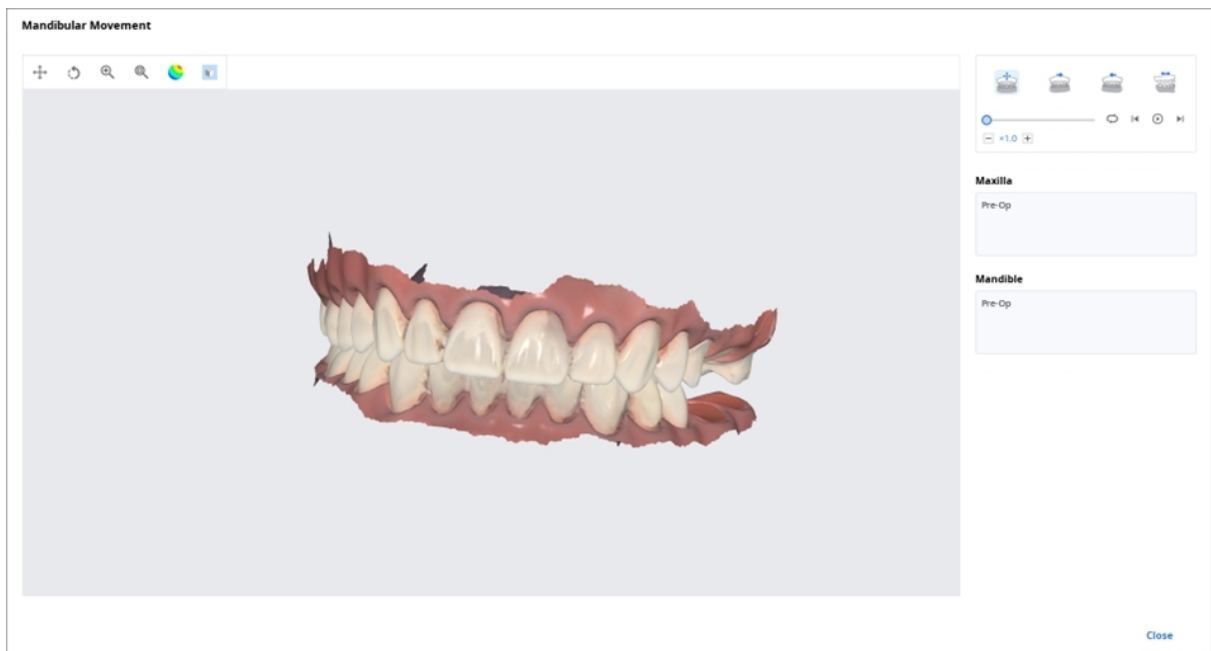

	Free Movement	Play the recording of free mandibular movement.
	Left Lateral Movement	Play the recording of the mandible moving to the left.
	Right Lateral Movement	Play the recording of the mandible moving to the right.

	Protrusive Movement	Play the recording of the mandible moving forward.
---	----------------------------	--

Manage Multi Occlusion Data in File Viewer

You can get multiple results of occlusion data from Medit Scan Application.

Medit Link provides the function to manage the results through the "File Viewer."

How to manage multiple results

When the user gets multiple occlusion results from Medit Scan for Clinics, Medit Link provides the option to select a result.

The list of results is listed by name, and you can select one item.

The selected model will be shown in the "Data Tree" and "Model Display Area."



You can collapse the list and change the item using the arrow icons.

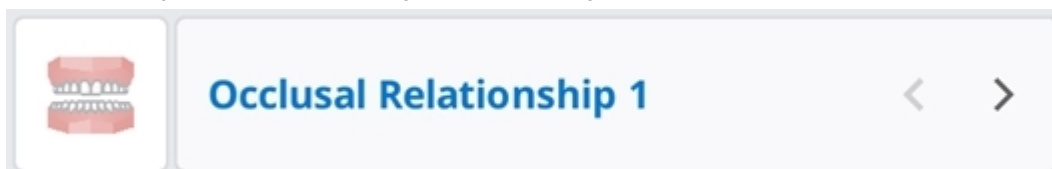


Table Options

You can change the Table Options to customize the layouts of the case list in the Patient/Case Box/Work Box/In Box/Order Box.

Table Options	×
Table Items	8 in use >
View Type	List View >
Grouping	None >

Status Filters

Show or hide each filter representing the status of the case. The hidden status filters are not displayed on the screen in the Case Box/Work Box/In Box/Order Box.

You can easily find the case you need with status filters.

- You can click the status you want to view to display only cases in that state.
- You can select multiple statuses at a time.

For each box, the following filtering options are available:

Box	Status Filters
Case Box	All / Form Only / Need Processing / Scan Completed / CAD / Ordered / Case Completed
Order Box	All / Pending / Canceled / Rejected / Accepted / Shipped / Case Completed
In Box	All / Pending / Canceled / Rejected / Accepted / Ready / Shipped / Case Completed
Work Box	All / Form Only / Need Processing / Scan Completed / CAD / Done

Table Items

Select the items to display in the case list in the Patient/Case Box/Work Box/In Box/Order Box. You can change the order of table items or adjust the table column width in the list view to best optimize for the user environment.

For each box, the following table items are available:

Box	Table Items
Patient	Patient Name / Patient ID / Gender / Number of Cases / Registration Date / Date of Birth / Last Scanned Date / Last Modified Date
Case Box	Status / Case Name / Date of Birth / Form Information / Last Modified Date / Scanned Date / Order Date / Lab Name / Modified By / Case Creation Date / Created By / Patient Name / Assignee
Order Box	Status / Case Name / Patient Name / Scanned Date / Order Date / Lab Name / Requested Delivery Date / Order ID / Last Modified Date / Assignee
In Box	Status / Case Name / Order Date / Requested Delivery Date / Last Modified Date / Clinic Name
Work Box	Status / Case Name / Patient Name / Form Information / Last Modified Date / Requested Delivery Date / Scanned Date / Clinic Name / Order Date

View Type

Change the view type between the Gallery View and List View on the Case Box.

For each box, the following view type options are available:


Box	View Type Options
Patient	List View
Case Box	List View / Gallery View
Order Box	List View
In Box	List View
Work Box	List View

Grouping

Select a grouping option to classify and rearrange cases in the list by the selected option.

For each box, the following grouping options are available:

Box	Grouping Options
Case Box	None / Patient Name / Scanned Date / Order Date
Order Box	None / Scanned Date / Order Date / Requested Delivery Date
In Box	None / Order Date / Requested Delivery Date
Work Box	None / Requested Delivery Date

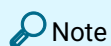

Note

The grouping options are unavailable in the Patient, Case Talk, and Trash Box.

Sorting

Select a sorting option to change the order of cases in the list according to the selected option.

For each box, the following sorting options are available:

Box	Grouping	Sorting Options
Patient	-	Patient Name / Patient ID / Date of Birth / Registration Date / Number of Cases / Last Scanned Date / Last Modified Date
Case Box	None	Patient Name / Last Modified Date / Scanned Date/ Order Date / Lab Name / Created By / Case Creation Date/ Modified By
	Patient Name	Last Modified Date / Scanned Date / Order Date / Lab Name / Created By / Case Creation Date/ Modified By
	Scan Date	Patient Name / Last Modified Date / Order Date / Scanned Date / Lab Name / Created By / Case Creation Date/ Modified By
	Order Date	Patient Name / Last Modified Date / Order Date / Scanned Date / Lab Name / Created By / Case Creation Date/ Modified By
Order Box	None	Patient Name / Scanned Date / Order Date/ Requested Delivery Date / Lab Name / Last Modified Date
	Scan Date	Patient Name / Order Date / Requested Delivery Date / Lab Name / Scanned Date / Last Modified Date
	Order Date	Patient Name / Order Date / Scanned Date / Requested Delivery Date / Lab Name / Last Modified Date
	Requested Delivery Date	Patient Name / Scanned Date / Order Date / Lab Name / Requested Delivery Date / Last Modified Date
In Box	None	Patient Name / Order Date / Requested Delivery Date / Clinic Name

Box	Grouping	Sorting Options
	Order Date	Patient Name / Scanned Date / Requested Delivery Date / Clinic Name
	Requested Delivery Date	Patient Name / Order Date / Clinic Name
Work Box	None	Patient Name / Last Modified Date / Requested Delivery Date / Clinic Name
	Requested Delivery Date	Patient Name / Last Modified Date / Clinic Name


Note

The sorting options are unavailable in the Case Talk and Trash Box.

Filtering and Search

Filter by Patient

Filter and display cases by patient name. You can create a list by filtering only the patients you need to work on.

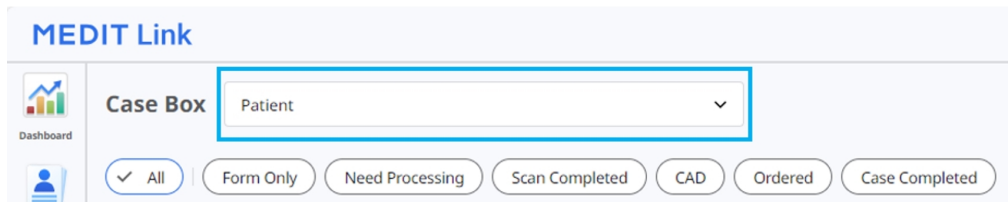


MEDIT Link

Case Box Patient

Dashboard

✓ All Form Only Need Processing Scan Completed CAD Ordered Case Completed

1. Enter a patient name.

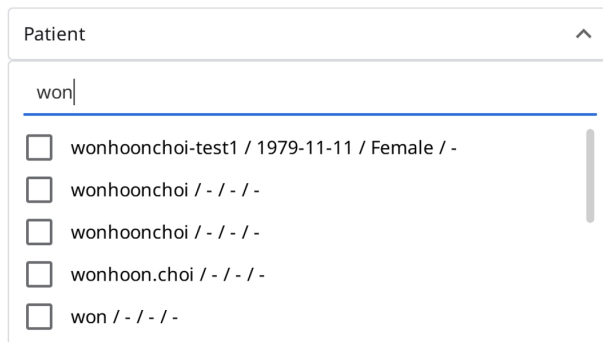


Patient

won

- ☐ wonhoonchoi-test1 / 1979-11-11 / Female / -
- ☐ wonhoonchoi / - / - / -
- ☐ wonhoonchoi / - / - / -
- ☐ wonhoon.choi / - / - / -
- ☐ won / - / - / -

2. Select the patients among the search results below.

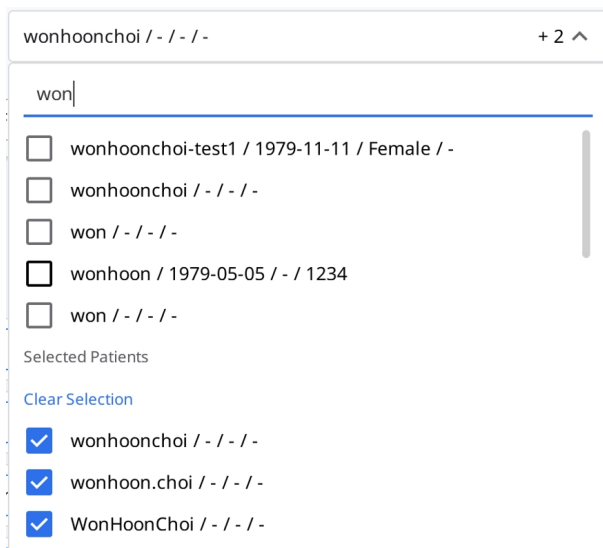

wonhoonchoi / - / - / - + 2

won

- ☐ wonhoonchoi-test1 / 1979-11-11 / Female / -
- ☐ wonhoonchoi / - / - / -
- ☐ won / - / - / -
- ☐ wonhoon / 1979-05-05 / - / 1234
- ☐ won / - / - / -

Selected Patients

[Clear Selection](#)

- ☒ wonhoonchoi / - / - / -
- ☒ wonhoon.choi / - / - / -
- ☒ WonHoonChoi / - / - / -

3. Only cases of the selected patients are displayed on the list.

Status	Patient Name	Case Name	Form Information	Last Modified Date	Scan Date	Lab Name	Order Date
Form Only	wonhoonchoi	originalcase	-	2023-10-05 오전 8:10			
Form Only	wonhoonchoi	wonhoonchoi's Case	-	2023-10-05 오전 8:09			
Scan Completed	wonhoonchoi	FileViverTest_1	8-Inlay / 10-Inlay / 12-Inlay / 13-Inlay / 14-Inl...	2023-06-15 오전 8:14	2023-02-09 오전 8:35		
Scan Completed	wonhoonchoi	FileViverTest_1 - Clone	2-Inlay / 3-Inlay / 4-Inlay / 5-Inlay / 6-Inlay / 8...	2023-05-25 오후 2:05	2023-02-09 오전 8:35		
Form Only	wonhoonchoi	FileViverTest_1 - Clone	8-Inlay / 10-Inlay / 12-Inlay / 13-Inlay / 14-Inl...	2023-03-15 오후 3:06			
Scan Completed	wonhoonchoi	FileViverTest_1 - Clone - Clone	8-Inlay / 10-Inlay / 12-Inlay / 13-Inlay / 14-Inl...	2023-03-09 오후 4:41	2023-02-09 오전 8:35		
Need Processing	wonhoonchoi	FileViverTest_1 - Clone	8-Inlay / 10-Inlay / 12-Inlay / 13-Inlay / 14-Inl...	2023-02-14 오전 8:13	2023-02-09 오전 8:35		
Scan Completed	wonhoonchoi	FileViverTest_1 - Clone	8-Inlay / 10-Inlay / 12-Inlay / 13-Inlay / 14-Inl...	2023-02-10 오전 11:27	2023-02-09 오전 8:35		
Form Only	wonhoonchoi	testsamplecase	-	2023-02-10 오전 7:39			
Scan Completed	wonhoonchoi	Multi Occlusion	-	2023-02-09 오전 9:26	2023-02-09 오전 9:01		
Scan Completed	wonhoonchoi	FileViverTest_1 - Clone - Clone	8-Inlay / 10-Inlay / 12-Inlay / 13-Inlay / 14-Inl...	2023-02-09 오전 8:50	2023-02-09 오전 8:35		
Case Completed	wonhoonchoi	wonhoonchoi's Case	12-Inlay / 13-Inlay / 14-Inlay	2022-10-13 오후 10:30	2022-10-11 오후 5:16		
Case Completed	wonhoonchoi	wonhoonchoi's Case	-	2022-10-13 오후 10:28			

Filter by Date

Filter and display cases last modified within the set duration.

1. Click the duration filter.

 All / 2018-01-01 ~

2. Select the duration from the list.

- All: Shows all cases.
- Custom Range: Specify the start and end dates to display cases last modified within the period.
- Today: Show the cases last modified today.
- Yesterday: Show the cases last modified yesterday.
- Last 7 days: Show the cases last modified during the previous 7 days excluding today.
- Last 30 days: Show the cases last modified during the previous 30 days excluding today.
- Last 90 days: Show the cases last modified during the previous 90 days excluding today.
- Current Month: Show the cases last modified during this month.
- Previous Month: Show the cases modified during the last month.

The case list will be updated in real time with the cases created or last modified within the selected duration.

Advanced Search Capabilities

Leverage search capabilities to quickly find the information you want.

All

▼

All / Contains

Q

The following Searching Criteria are available for each box:

Patient	All, Patient Name, Patient ID
Case Box	All, Case Name, Patient Name, Patient ID, Form Information, Lab Name
Order Box	All, Case Name, Patient Name, Patient ID, Lab Name, Order ID
In Box	All, Case Name, Patient Name, Clinic Name, Order ID
Work Box	All, Case Name, Patient Name, Form Information, Clinic Name, Order ID
Case Talk	Case Talk Name (Contains), Patient Name (Equals)
Trash Box (Clinic)	All, Case Name, Patient Name, Lab Name, Deleted By
Trash Box (Lab)	All, Case Name, Patient Name, Clinic Name, Deleted By

You can select the Searching Criteria and choose between the Contains and Equals logic.

All

▼ All / Contains

🔍

Searching Criteria

☒ All
 ☐ Case Name
 ☐ Patient Name
 ☐ Patient ID
 ☐ Form Information
 ☐ Lab Name

☒ Contains
 ☐ Equals

Create/Edit/Clone/Delete Case

Creating Cases

Create Case in Case Box (Clinic)

1. Click the "New Case" button in the upper right corner.

New Case Information

Patient Name*



Case Name*

 Creating patient information involves consent from the patient.
Patient information can be shared with 3rd parties while placing an order.

Partner

Select partner

▼

Assignee

Select assignees

▼



 You can select up to 5 assignees.

Tags

You can add the entered tag with Enter or Spacebar.

 You can add up to 10 tags of 100 characters or less.

[Terms & Conditions](#)

Cancel

Register

Register & Scan

2. Enter basic information, such as the patient's name and case name.
 - If you enter a name not on the patient list, a new patient will be added.
 - If you enter a name that is in the patient list, you can select it from the list.

- If registering a case for an existing patient, click the Patient List icon to select an existing patient.



- In the Patient List dialog, you can register a new patient or edit the information of an existing patient.

Patient List

Select a patient in the list to set the patient of the case. You can add or manage the patient information in Medit Link.

ALL Search / Contains + New Patient

	Patient Name	Patient ID	Gender	Birth Date	Number of Cases	Registered Date	Recent Visit Date	
☐	test				0	2020-08-28	2020-08-28	
☐	test				0	2020-08-28	2020-08-28	
☐	test				13	2020-08-23	2020-08-28	
☒	test				85	2020-08-23	2020-05-08	
☐	Wonhoon Choi				1	2020-04-01	2020-04-01	
☐	test				0	2020-08-18		
☐	13				0	2020-08-20		
☐	14				0	2020-08-20		
☐	15				0	2020-08-20		

Close Select

3. The case name will be filled in automatically, but you can still change it.

4. Select a partner for the case.

- Click the Partner field and select a partner from the list.

Partner

Select partner ^

heavyuserlab@meditqa.com


Full Services

- To use a lab-defined order form, select a partner marked Uses Lab Form. For details, see [Using a Lab-Defined Order Form](#).

qa_lab_01 Uses Lab Form

Full Services Dental Esthetics Crown & Bridge

- Partner selection is available to admin accounts only.
- Partner selection is not available in offline mode.
- If no partners are available, click "Find Partners" to find and connect with a partner.

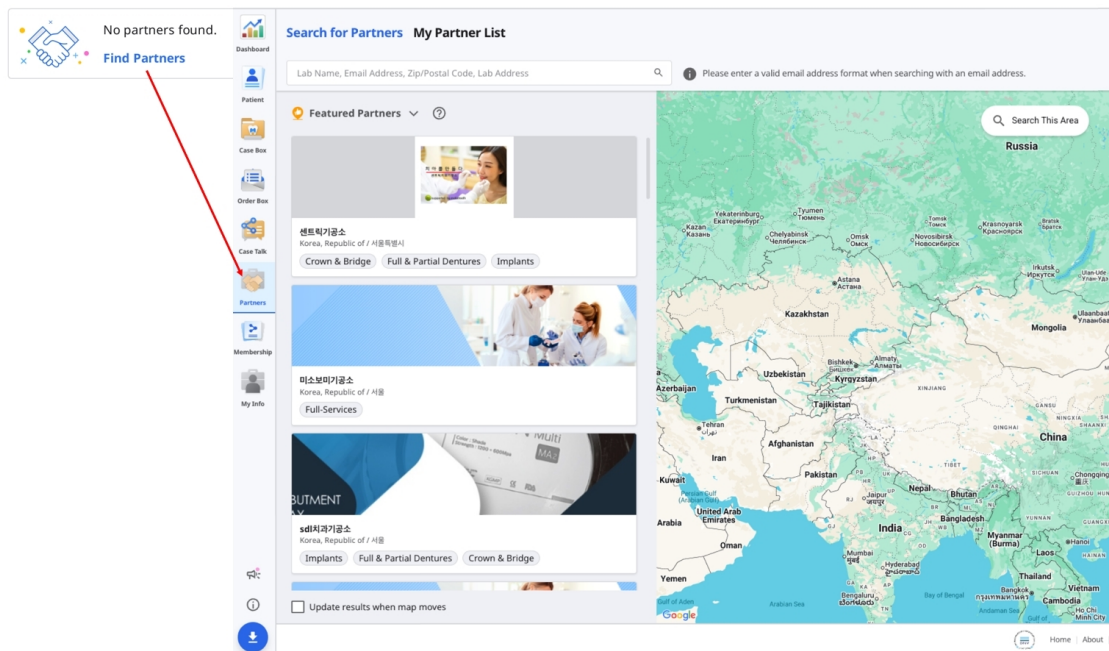

5. Select or add assignees for the case in the Assignee field.

- You can search for and select up to 5 assignees per case using the dropdown list.

Assignee

Select assignees

Search assignees

Dentist

☐ seolah.jeon_2

Dental Technician

☐ kisoo.jeon

Dental Hygienist

☐ rana.jeon

- To add a new assignee:

- Click the Assignee Management icon in the top-right corner of the dropdown.

Assignee

Select assignees

Search assignees

Dentist

☐ seolah.jeon_2

- In the window that opens, click Add Assignee.

Assignee Management

Add and manage assignees per case.

Total 3 ■ Dentist 1 ■ Dental Technician 1 ■ Dental Hygienist 1

+ Add Assignee

Name	Role
☐ rana.jeon	■ Dental Hygienist
☐ kisoo.jeon	■ Dental Technician
☐ seolah.jeon_2	■ Dentist

- Enter a name (up to 50 characters) and select a role for the new assignee. Each assignee can have only one role, chosen from: Dentist, Dental Technician, or Dental Hygienist.

6. Click "Register & Scan" to complete the case registration and proceed to scan.

- **Register:** You will be redirected to the main page. This is useful when registering multiple cases.
- **Register & Scan:** Move to the Form tab to register form information for the case.

Create Case in Work Box (Lab)

1. Click the "New Case" button in the upper right corner.

New Case Information

Clinic Name*

▼

Patient Name*

Case Name*

Assignee*

Add assignees

▼

+

ℹ You can select up to 5 assignees.

Requested Delivery Date

📅

2025-07-30

Tags

You can add the entered tag with Enter or Spacebar.

ℹ You can add up to 10 tags of 100 characters or less.

Cancel

Register

2. Enter or select the clinic name.
3. Enter the patient's name and case name. The case name will be filled in automatically, but you can still change it.
4. Choose an assignee for this case.
5. Enter the estimated delivery date considering the completion date of the work.
6. After registering the case, you will be redirected to the main page.

Editing Cases

Edit Case in Case Box (Clinic)

1. Click the "Edit" icon next to the case name.

Edit Case Information

Patient Name*

Test



 "Test Male 19700114" is selected.

Case Name*

Test's Case

 Creating patient information involves consent from the patient.
Patient information can be shared with 3rd parties while placing an order.

Partner

Select partner

 No partners found.
[Find Partners](#)

Assignee

Add assignees



 You can select up to 5 assignees.

Tags

You can add the entered tag with Enter or Spacebar.

 You can add up to 10 tags of 100 characters or less.

[Terms & Conditions](#)

Cancel

Apply

- Click the "Search" icon to change the patient's name. Note that it cannot be edited directly.; instead, you can find a patient name you would like to replace the current one with and can edit the patient information from the list.

Patient List

Select a patient in the list to set the patient of the case. You can add or manage the patient information in Medit Link.

All
Search / Contains
+ New Patient

	Patient Name	Patient ID	Gender	Birth Date	Number of Cases	Registered Date	Recent Visit Date	
☐	test				0	2020-08-28	2020-08-28	
☐	test				0	2020-08-28	2020-08-28	
☒	test				13	2020-08-23	2020-08-28	
☐	test				85	2020-08-23	2020-05-08	
☐	Wonhoon Choi				1	2020-04-01	2020-04-01	
☐	ttttt				0	2020-08-18		
☒	13				0	2020-08-20		
☐	14				0	2020-08-20		
☐	15				0	2020-08-20		

Close
Select

- Check the patient's name you want to replace the current one with.
- To create a new patient, press the "New Patient" button.
- Use the "Edit" button to change the registered information about the patient.
- To change the partner for a case, click the Partner field, select a different partner, then click Apply.

Partner

Elysee Dental Test
^

Elysee Dental Test
Lab_Hari_Standard
MEDIT_LAB_STAGE
Dental Esthetics
Crown & Bridge
Implant Prosthetics

You can select up to 5 assignees.

Tags

You can add the entered tag with Enter or Spacebar.

You can add up to 10 tags of 100 characters or less.

[Terms & Conditions](#)

Cancel

Apply

- To modify an existing assignee:

- In the Assignee Management modal, click the three-dot icon (⋮) next to the assignee you want to edit.

Assignee Management
Add and manage assignees per case.

Total 3 ■ Dentist 1 ■ Dental Technician 1 ■ Dental Hygienist 1 + Add Assignee

Name	Role	
☐ ranaj.jeon	■ Dental Hygienist	⋮
☐ kisoo.jeon	■ Dental Technician	⋮
☒ seolah.jeon_2	■ Dentist	⋮

- You can update both the name and the role.

Modify Assignee
Update the name and role of the assignee.


Name*

13/50

Role*

■ Dentist

- Dentist
- Dental Technician
- Dental Hygienist

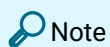

Note

Updates to assignee information are automatically applied to all existing cases. However, these changes won't be reflected in cases that have already been ordered.

If you remove an assignee, they will also be removed from any cases to which they were assigned.

8. Click "Apply" to finish the process.

Edit Case in Work Box (Lab)

1. Click the "Edit" icon next to the case name.

Edit Case Information

Clinic Name*

seolah

Patient Name*

3.4.2_ScanBody_Labs

Case Name*

Scan for Labs 3.4.2_ScanBody_Labs's Case

Assignee*

Add assignees

+

i You can select up to 5 assignees.

Requested Delivery Date

2025-04-10

Tags

You can add the entered tag with Enter or Spacebar.

i You can add up to 10 tags of 100 characters or less.

Cancel

Apply

2. Click "Apply" after the editing is done.

Note

You cannot edit the information of the cases received through the order in In Box.

You can add or edit assignees for ordered cases.

Cloning Cases

Clone Case in Case Box (Clinic)

1. Click the "Clone" icon at the top.



2. You will be asked what kind of information you would like to clone to the new case.

Information
 Select a clone type.

Cancel Form Data Only All Data

3. Choose whether to replicate only form information or all data.
4. When cloning is complete, a case named "[Case name] – Clone" will be created.

Clone Case in Work Box (Lab)

1. Click the "Clone" icon at the top.



2. Duplicate all data to create a new case.
3. When cloning is complete, a case named "[Case name] – Clone" will be created.


Deleting Case

Delete Case (Clinic and Lab)

1. Click the "Delete" icon.



2. Click "OK" to delete the case.

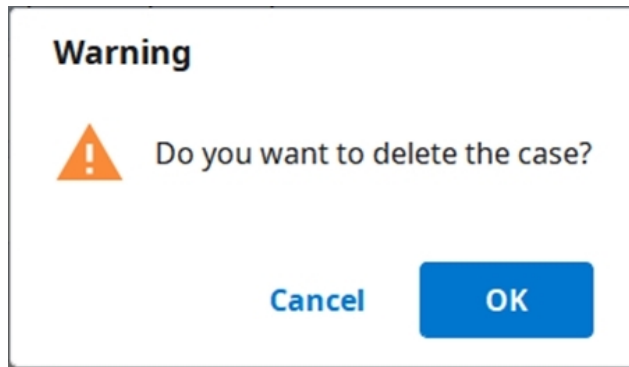

3. Deleted cases are moved to the Trash Box and can be recovered if necessary.

Form Information

You can enter the form information based on the connected program and the treatment type.

The information you enter will be used as basic information for working on scanning or CAD programs and ordering.

Registering Form Information in Case Box

Enter on the Form Tab


1. Select whether you are going to input treatment information by teeth or by arch.
2. Select Type – Method – Material – Shade.
 - If you entered information incorrectly, you go back to select it again.
 - The selected information is automatically registered to the order list.
3. Click "Save" to save your registration.

Using a Lab-Defined Order Form

Some labs provide their own order forms. These labs are marked **Uses Lab Form**.

When you select one of these labs as the partner for a case, only the options available in the lab's order form can be selected.

1. Select a partner marked **Uses Lab Form** when creating or editing a case.

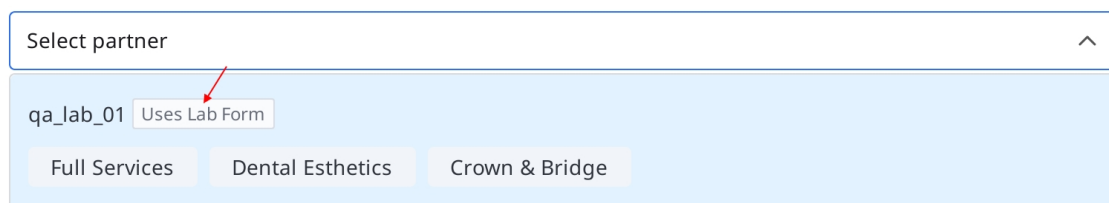

The screenshot shows a web interface for selecting a partner. At the top is a search bar labeled 'Select partner' with an upward arrow icon on the right. Below the search bar, a dropdown menu is open, displaying a list of partners. The first partner is 'qa_lab_01', which has a small tag next to it that says 'Uses Lab Form'. A red arrow points to this tag. Below the partner list, there are three buttons: 'Full Services', 'Dental Esthetics', and 'Crown & Bridge'.

- After creating or updating the case, go to the **Form** tab, select the teeth or arch, complete the required form information, and click **Save**.

The screenshot shows the 'Form' tab in a dental software interface. The main workspace displays a dental arch diagram with teeth numbered 1 through 48. A red arrow points to the 'Form' tab at the top. Another red arrow points to the 'Teeth' label at the bottom of the arch diagram. On the right side, there is a panel titled 'Order' with a grid of options categorized by Type, Method, Material, Shade, and Favorites. Below this is an 'Order List' table.

Type	Method	Material	Shade	Favorites
Inlay	Onlay	Coping	Denture Coping	
Crown	Pontic	Post & Core	Veneer	
Implant Crown	Custom Abutment	Surgical Guide	Offset Substructure	
Custom Abutment + Impla...	Custom Abutment + Coping	Diagnostic Wax-Up	Partial Denture	

No.	Type	Method	Material	Shade	
26	Crown	Full-Contour	Zirconia	A1	

At the bottom right, there are 'Cancel' and 'Save' buttons, with a red arrow pointing to the 'Save' button.

To place the order, see [Ordering Cases](#).

Registering Form Information in Work Box

Enter on the Form Tab

- Select the scan program you want to use.
- Select one of the following options that will define your scanning strategy: Teeth/Flexible Multi-Die/Multi-Die.
- Methods may differ depending on the selected scan program.

	Medit Scan for Clinics	Medit Scan for Labs
Teeth	0	0

Arch	O	O
Flexible Multi-Die	X	O
Multi-Die	X	O

4. Select teeth, arch, or die numbers for the scanning process.
5. Select Type – Method – Material – Shade.
 - If you entered information incorrectly, you go back to select it again.
 - The selected information is automatically registered to the list.
6. Click "Save" to save your registration.

Connecting Teeth into Bridge

You can connect multiple teeth into a bridge when the same material is applied to multiple adjacent teeth. Click the bridge icon beside the selected teeth to connect or disconnect them.



Setting Scanning Options

Depending on the program you are working with, various scan options are provided.

The selected scan options are passed to each scanning program.

Medit Scan for Labs

General options

	Adjacent Teeth	Mark the prepared and adjacent teeth and pass this information to the scanning program.
	Antagonist	Automatically mark the antagonist and passes this information to the scanning program.
	Articulator	Display options for selecting the articulator in the scanning program.
	Movement Marker	Allow using the "Jaw Motion" scan strategy in the scanning program.

Options by teeth

	Scan Body	Allow setting a strategy for scanning scan bodies in the scanning program.
	Pre-Op Scan	Allow setting a strategy for scanning a pre-op model in the scanning program.
	Wax-Up Scan	Allow setting a strategy for scanning wax-ups in the scanning program.

Adding to Favorite

You can quickly select a product by registering it as a favorite.

The list of registered favorites is managed as a template and can be shared through the export and import features.

Add to Favorites

1. Select the tooth number.
2. Select Type – Method – Material – Shade.
 - Type – Method – Material – Shades can be registered as favorites even if not all are selected.
3. Click the "Star" icon.



4. A combination of the selected parameters, excluding tooth number, is on your favorite search list.

Manage Favorite Templates

You can register, manage, and share different types of templates.

Add	Add a new template.
Rename	Rename the current template.
Delete	Delete the current template.
Export	Export the current template to a file.
Import	Import a template file.

Work on Scanning/CAD

Scan Process

Based on the entered information, you proceed with scanning.

1. Enter/save treatment information and select the necessary scan options.
2. Click the "Scan" button in the upper-right corner of the window.

CAD Process


Medit Link provides the DentalCAD integration app to connect Medit Link with DentalCAD.

1. Install the DentalCAD Integration app in the App Box.
2. Go to the case details and click on the DentalCAD Integration app icon on Quick Launcher.
3. Specify the path to the DentalCAD app installed on your PC in the pop-up dialog.
4. Enter the case details of the scanned case and click the DentalCAD Integration app icon, which will pop up a dialog to create a dentalProject file.
5. Register form information for DentalCAD and click "Open" to open the CAD program.
6. If you install the DentalDB Integration app, the CAD results will be automatically delivered to the specified save path in exocad. (Lab accounts only)


Utilize Add-in Apps

Add-in Quick Launcher

You can directly run various add-in apps by clicking on the installed app icons in the Quick Launcher. You can change the order of app icons by dragging them.

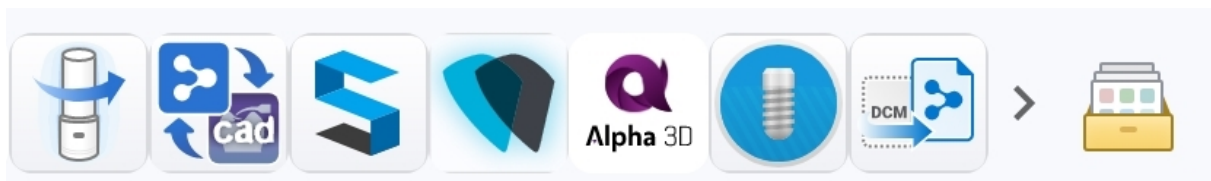

App Launcher Management

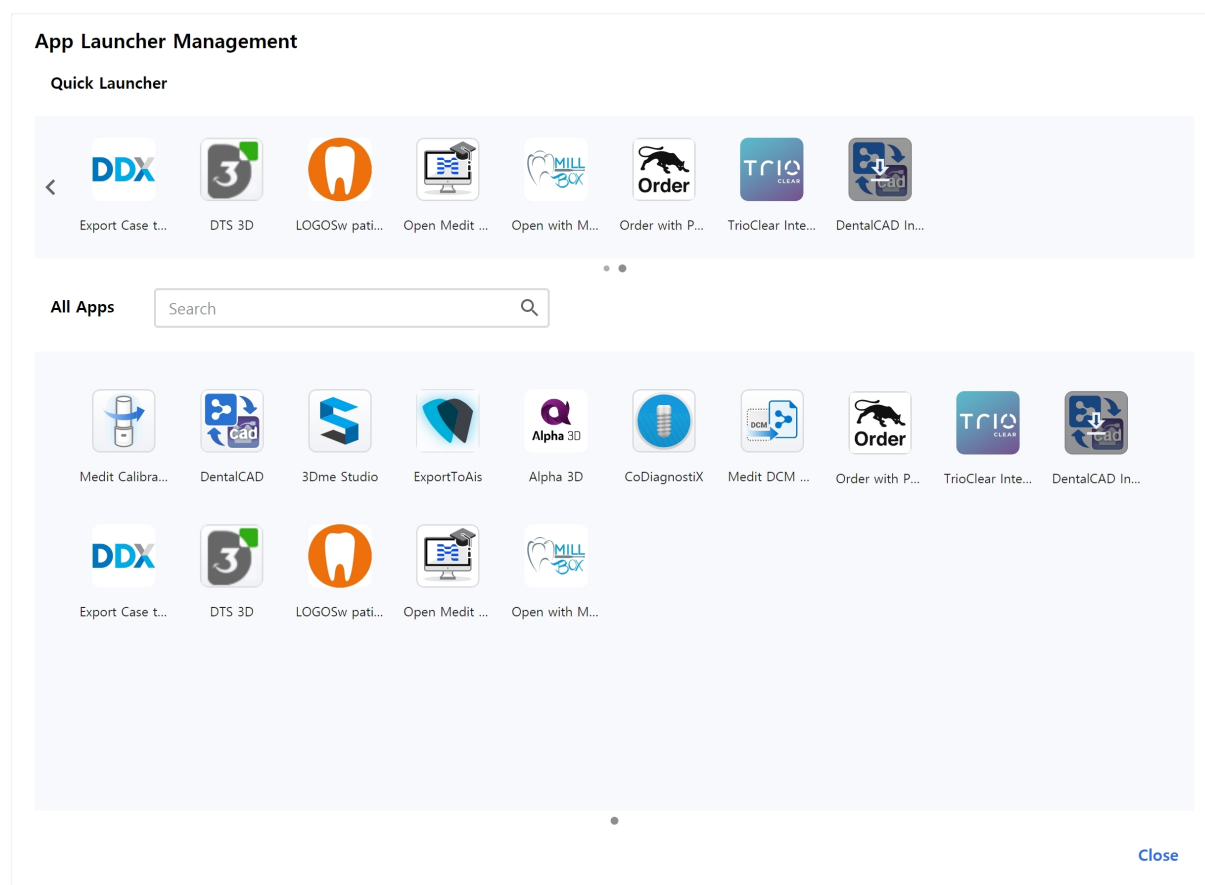

Quick Launcher

You can define which apps appear in the Quick Launcher. You can change the order of app icons by dragging them, or you can create folders to group similar apps together to organize them.

All Apps

All apps you have installed are shown. You can right-click on an app icon to add the app to the Quick Launcher. From the right-click the context menu, you can also select to run, install, or uninstall the app.

Get Result Files

Export Scan Result Files

In addition to sharing results with partners through orders, Medit Link allows you to store your result files in the desired file format using the export function.

Export Files in This Case

Options

File List

Export Files To

C:/ReleaseBuild/

Folder Name

5162022-Crown Fit Sample data

File Name

5162022-Crown Fit Sample data

Type

exocad

☐ Enabled "automatic .dentalProject import" in DentalDB

Mesh Format

Convert mesh files to the below formats

☒ meditMesh

☐ OBJ

☐ PLY

☐ STL

☐ Do not convert the attached files

Image Format

Convert images to the below format

☒ Convert images to the below format

☐ JPG

☒ JPEG

☐ PNG

☐ BMP

☒ Compress Files

☒ Combine Individual Mesh

☐ Save Case Information as PDF

[Preview](#)

Close


Export

Export Files To	Specify the path where you want to store the result files.
-----------------	--


Folder Name	The system creates a new folder and stores all result files in it. The default folder name is set by the naming rule. But you can change the folder name directly. You can define the naming rule on Settings > Export Options > Rules for Folder Name.
File Name	Set the file name for the result. The naming rule sets the default file name. But you can change the file name directly. You can define the naming rule through Settings > Export Options > Rules for File Name
Type	Choose the design program you want to use. Medit Link will automatically adjust the files according to your selection, enabling you to proceed with the design without the need for additional axis alignment.
Mesh Format	Specify the format of the mesh file.
Image Format	Specify the format of the exported image files.
Compress Files	Compress all results into a single file.
Combine Individual Mesh	Combine the result of the maxilla or mandible with the individual tooth. It is useful when importing several files into a program is not possible.
Export All Occlusal Relationship	When multiple occlusal results are generated from the scan app, you can export all occlusal results at once.
Save Case Information as PDF	Save the case information as a PDF with the result.

Export Result Files from Case List


You can export result files of a case or multiple cases simultaneously with the "Export" option above the case list in Case Box and Work Box.

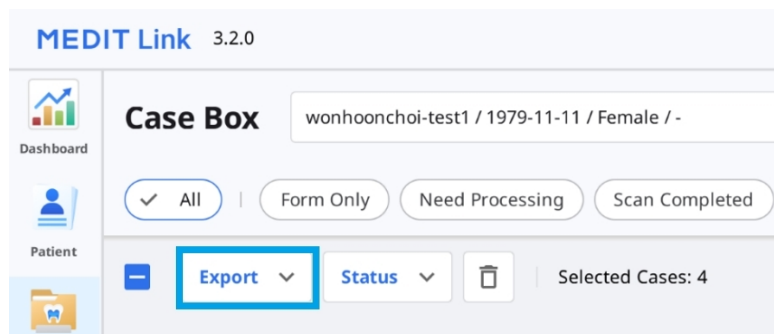

1. Go to Case Box or Work Box.
2. Select all cases you want to export result files from the list.

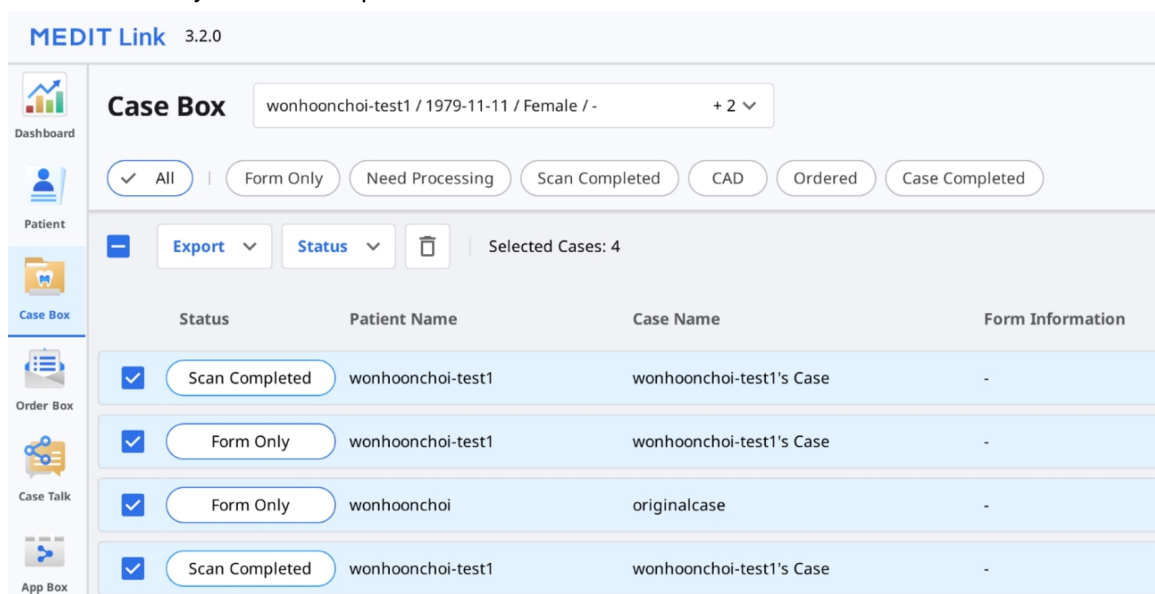

3. Click "Export" and select "Result Files."

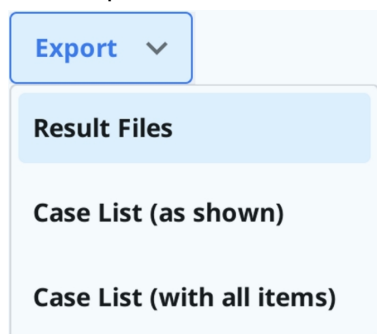

4. Set the destination folder, file type, mesh format, and image format, and select export options on the "Options" tab.

Export Files in Selected Cases

Options

Export Files To

C:/Users/Wonhoon Choi Laptop/Desktop/ExportTest/3/

Folder Name

Case Creation Date - Case Name

File Name

Case Creation Date - Case Name

Type

Medit

Mesh Format

Convert mesh files to the below formats

☐ meditMesh
 ☐ OBJ
 ☒ PLY
 ☐ STL

☐ Do not convert the attached files

Image Format

☐ Convert images to the below format

☐ JPG
 ☐ JPEG
 ☒ PNG
 ☐ BMP

☐ Compress Files

☐ Combine Individual Mesh

☐ Export All Occlusal Relationship

☐ Save Case Information as PDF

Close Export

Note

The "Folder Name" and "File Name" fields are not editable when exporting files for multiple cases.

The folder name and file name are automatically created according to the rules you set for the "Rule for Folder Name" and "Rule for File Name" options in Settings > Export Options.

- Click "Export."
- Click "Yes" to open the folder where the files are saved.

Question




Exporting successfully completed.
Do you want to open the folder containing the exported file?

No


Yes

Export Result Files from File Viewer

You can export scan result files for a case with the "Export" icon on the top right corner of the File Viewer tab.

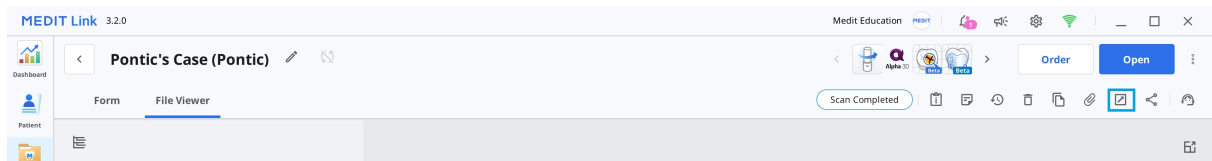

1. Go to Case Box or Work Box.
2. Click to open a case from the list.
3. Click the "Export" icon on the top right corner of the File Viewer tab.

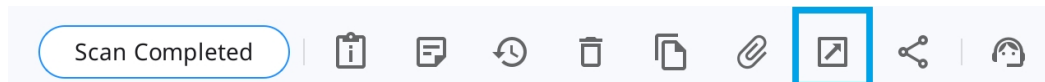

4. Set the destination folder, folder name, file name, file type, mesh format, and image format, and select export options on the "Options" tab.

Export files in case

Options

File List

Export Files To

C:/Users/Medit-Wonhoonchoi/Desktop/test/

Folder Name

Multi Occlusion - Scan - Before Scanning

File Name

Multi Occlusion - Scan - Before Scanning

CAD Type

Maestro

Mesh Format

Convert mesh files to the below formats

meditMesh

OBJ

☒ PLY

STL

☒ Do not convert the attached files

Image Format

☐ Convert images to the below format

JPG

JPEG

☒ PNG

BMP

☐ Compress Files

☐ Combine Individual Mesh


☐ Export All Occlusal Relationship

☐ Save Case Information as PDF [Preview](#)

Close


Export

5. In the "File List" tab, you can select the required files for the case by deselecting unnecessary ones.

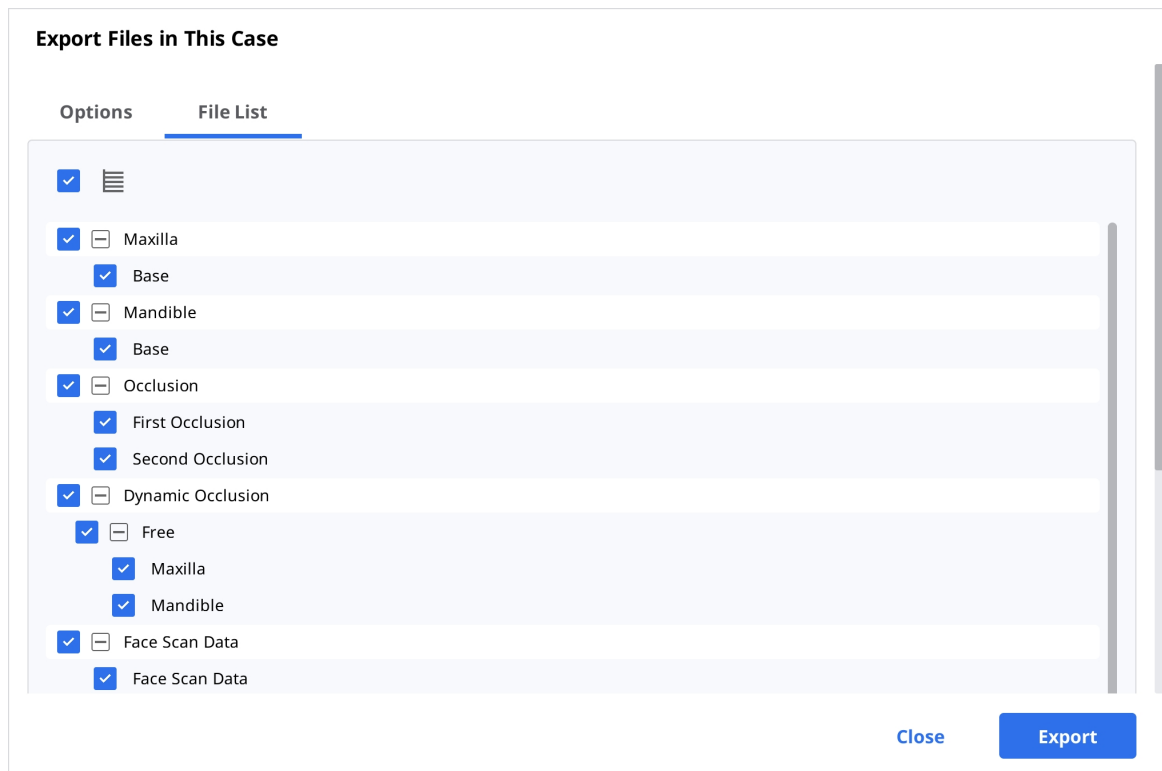


6. Click "Export."

7. Click "Yes" to open the folder where the files are saved.

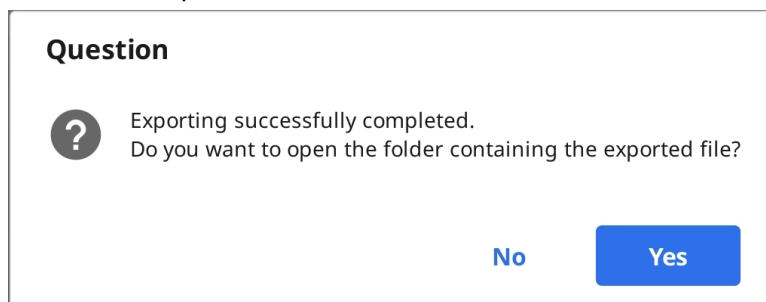

Export Result Files from Data Tree


You can export a file or files in a specific data group from the expanded menu in the file viewer.

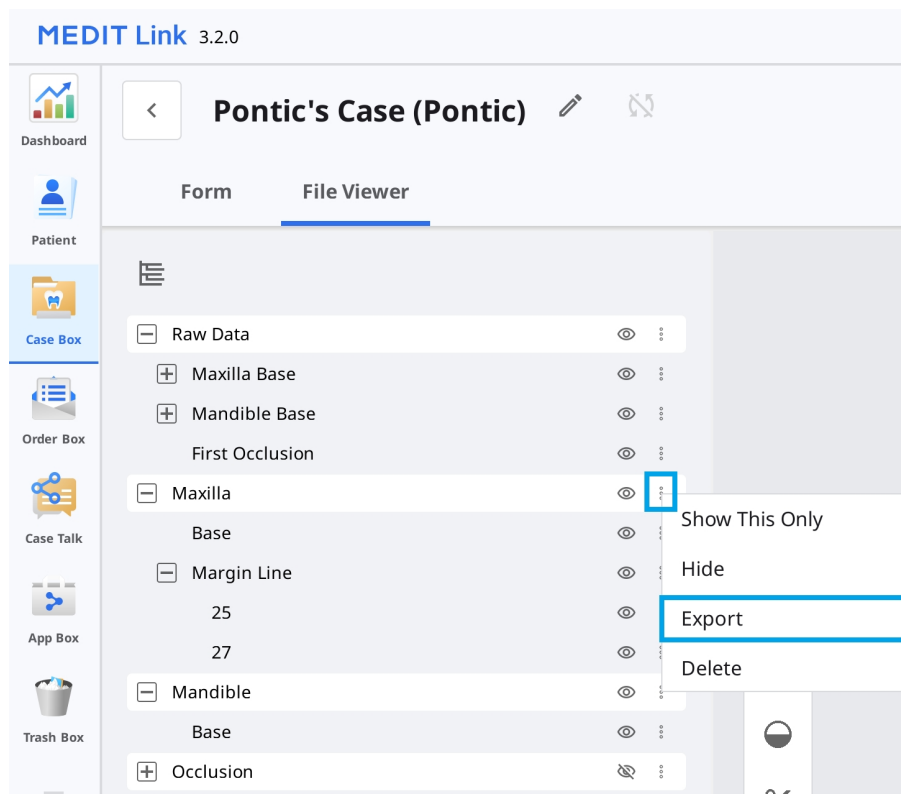

1. Go to Case Box or Work Box.
2. Click to open a case from the list.
3. Click the three dots icon next to the data group or an entity in the data tree in the File Viewer tab.



4. Then click "Export."

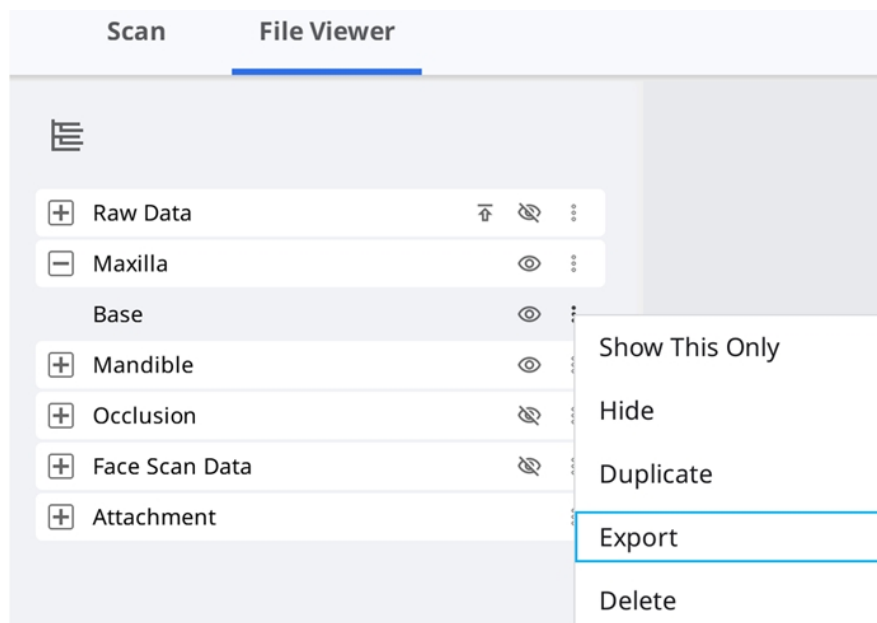

5. Set the destination folder, folder name, file name, file type, mesh format, and image format, and select export options on the "Options" tab.

Export files in case

Options

File List

Export Files To

C:/Users/Medit-Wonhoonchoi/Desktop/test/

Folder Name

Multi Occlusion - Scan - Before Scanning

File Name

Multi Occlusion - Scan - Before Scanning

CAD Type

Maestro

Mesh Format

Convert mesh files to the below formats

☐ meditMesh

☐ OBJ

☒ PLY

☐ STL

☒ Do not convert the attached files

Image Format

☐ Convert images to the below format

☐ JPG

☐ JPEG

☒ PNG

☐ BMP

☐ Compress Files

☐ Combine Individual Mesh

☐ Export All Occlusal Relationship

☐ Save Case Information as PDF [Preview](#)

Close

Export

6. You can see only the files you chose in the data tree are selected in the File List.

Export Files in This Case

Options **File List**

☐ Maxilla
 ☐ Base
 ☒ Mandible
 ☒ Base
 ☐ Occlusion
 ☐ First Occlusion
 ☐ Second Occlusion
 ☐ Dynamic Occlusion
 ☐ Free
 ☐ Maxilla
 ☐ Mandible
 ☐ Face Scan Data
 ☐ Face Scan Data

Close **Export**

- Click "Export."
- Click "Yes" to open the folder where the files are saved.

Question

? Exporting successfully completed.
Do you want to open the folder containing the exported file?

No **Yes**

Export Case List as CSV

Medit Link provides the function to export the case list in CSV format.

You can export the case list with basic information and manage cases easily.

- Go to Case Box, Order Box, In Box, or Work Box.
- Select the cases that you want to export.

Case Box Patient



All / 2018-01-01 ~ 2026-09-23 All / Contains

☒ All
 ☐ Form Only
 ☐ Need Processing
 ☐ Scan Completed
 ☐ CAD
 ☐ Ordered
 ☐ Case Completed

Total: 1093

Status	Case Name	Patient Name	Form Information	Last Modified Date	Scanned Date	Lab Name	Date of Birth
☒ Ordered	B59_3의 케이스	B59	25-Crown / 26-Crown / 27-Crown	2023-04-17 PM 2:49	2023-03-29 PM 3:19	heavyuserlab@m...	
☒ Case Completed	B59_2의 케이스	B59	25-Crown / 26-Crown / 27-Crown	2023-03-29 PM 3:55	2023-03-29 PM 3:19	heavyuserlab@m...	
☒ Case Completed	B59_1의 케이스	B59	25-Crown / 26-Crown / 27-Crown	2023-03-29 PM 3:55	2023-03-29 PM 3:19	heavyuserlab@m...	

- Click "Export" in the action bar at the bottom of the screen and select "Case List (as shown)" or "Case List (with all items)."

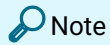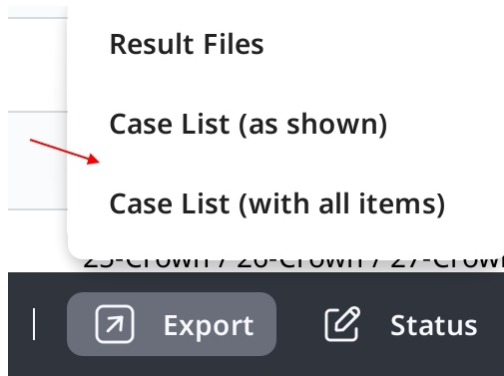

Note

Case List (as shown): Export the case list with the currently displayed items on the screen. Go to View Options > Table Items to show or hide the items as you desire.

Case List (with all items): Export the case list with all items, including Status, Case Name, Date of Birth, Form Information, Last Modified date, Scanned Date, Order Date, and Lab Name.

4. Define the path and file name, then save the list.
5. The case list is saved as a CSV file.

Export Margin Line


Medit Scan for Clinics provides a function to create the margin line, and Medit Link can transfer the margin line to the CAD program.

Export Margin Line as a Color

In general, Medit Link converts the margin line as color and exports the result with maxilla or mandible together. You can create the margin line in their cad referencing the color information.

Example:

Maxilla

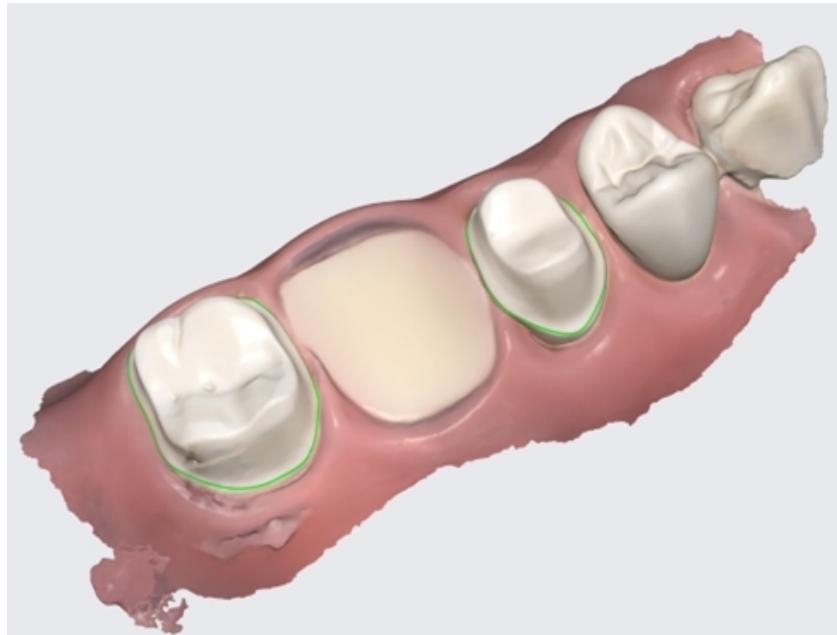

Result with the margin lines as color information

Note

This feature is available for formats that include color information, such as meditMesh, OBJ, and PLY. The margin line does not convert into STL.

Note


Type > 3Shape does not support this feature.

Export Margin Line as a File

In the case of Exocad, Medit Link export the margin line to "~margin.xyz" which can be imported into exocad directly.

You can use this file in exocad without additional work.

Example:

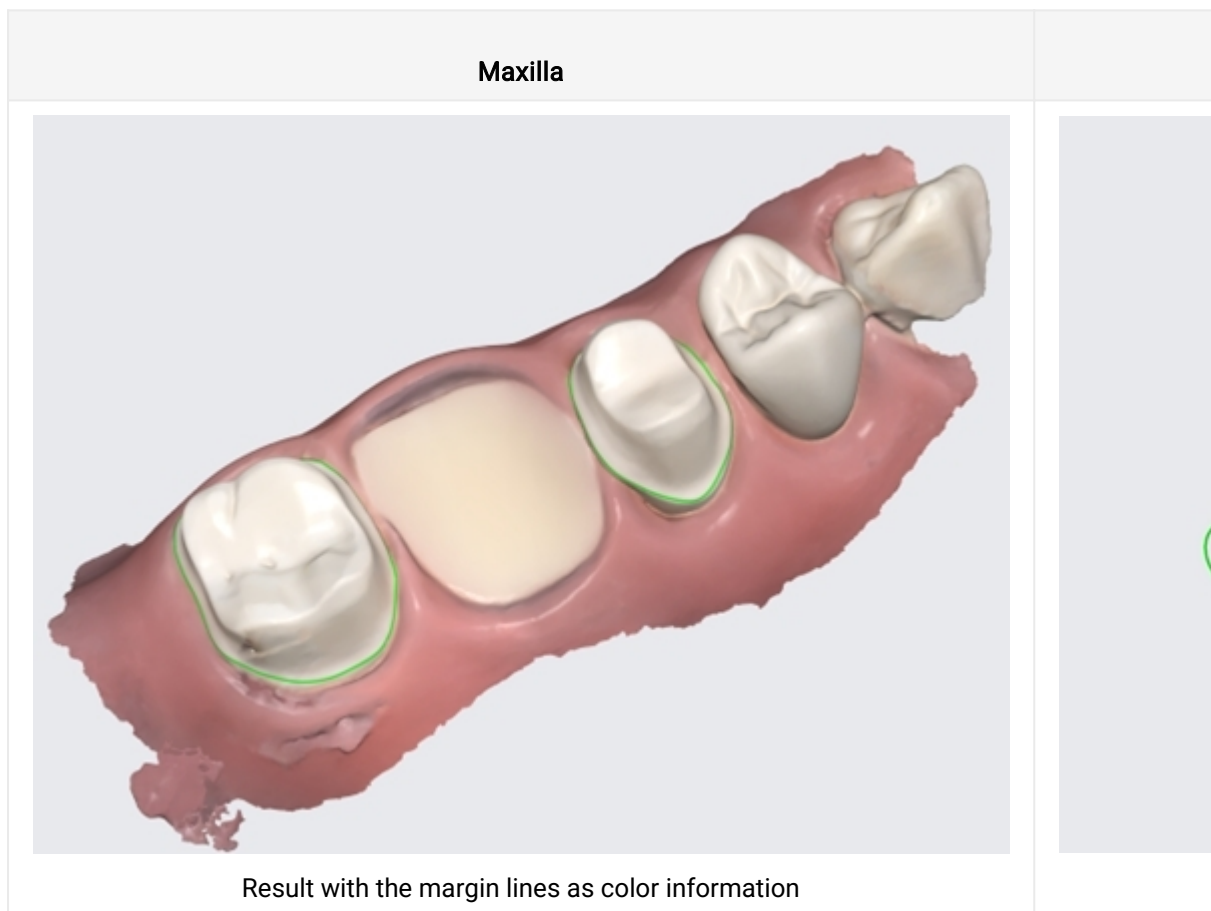

Download Files from Medit Link Web

Medit Link provides a download function.

If you want to get the result files without installing the Medit Link application, you can download the files from the Medit Link web service.

Note


The data should be uploaded on the Medit Cloud to download the files from the web.

The download function is available in the following menu of the Medit Link web service.

- Case Box
- Work Box
- Web Viewer

Case Box/Work Box

You can find the "Download" button at the top right.

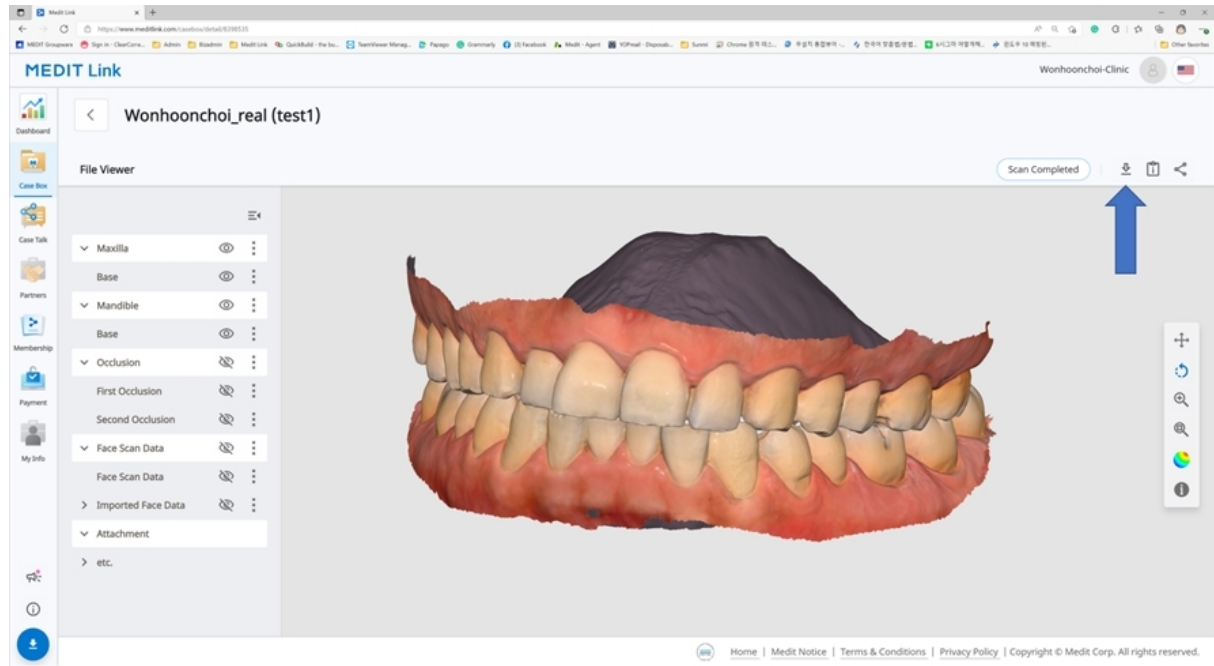

Web Viewer

When you share the case using Case Talk, you can set whether to allow downloading or not.

If you check “Grant Download Permission” and share the case, all users who connect to the shared link can download the files.

Medit Link

Share through Case Talk

Share with Anyone

Create a shared link and share it with anyone.

Expiration Date

2022-04-27

You can extend the period up to 90 days from today.

PIN Code

☐ PIN Code

Enter PIN Code

Set the PIN code using six digits for privacy protection.

☒ Share the patient's name

☒ Grant Download Permission [Premium Feature](#)

If you check this option, all users who have access to the shared link can download the files.

The premium feature is provided to all users for a limited time. Feature availability is subject to change depending on the membership plan.

☒ Allow comments

If you check this option, you will be able to communicate through the shared link.

Comment

Please enter your comment.

0/1024


By pressing the "Get Link" button, you agree that you understand personal information can be collected and provided to a third party (or shared in any public domain) ONLY with the consent of the person whose images on which you have worked. You also agree the consent of that person was truthfully obtained, and Medit is not responsible for any problems which may arise there after.

Close

Get Link

175


In the Web Viewer, you can find the "Download" icon at the top right.

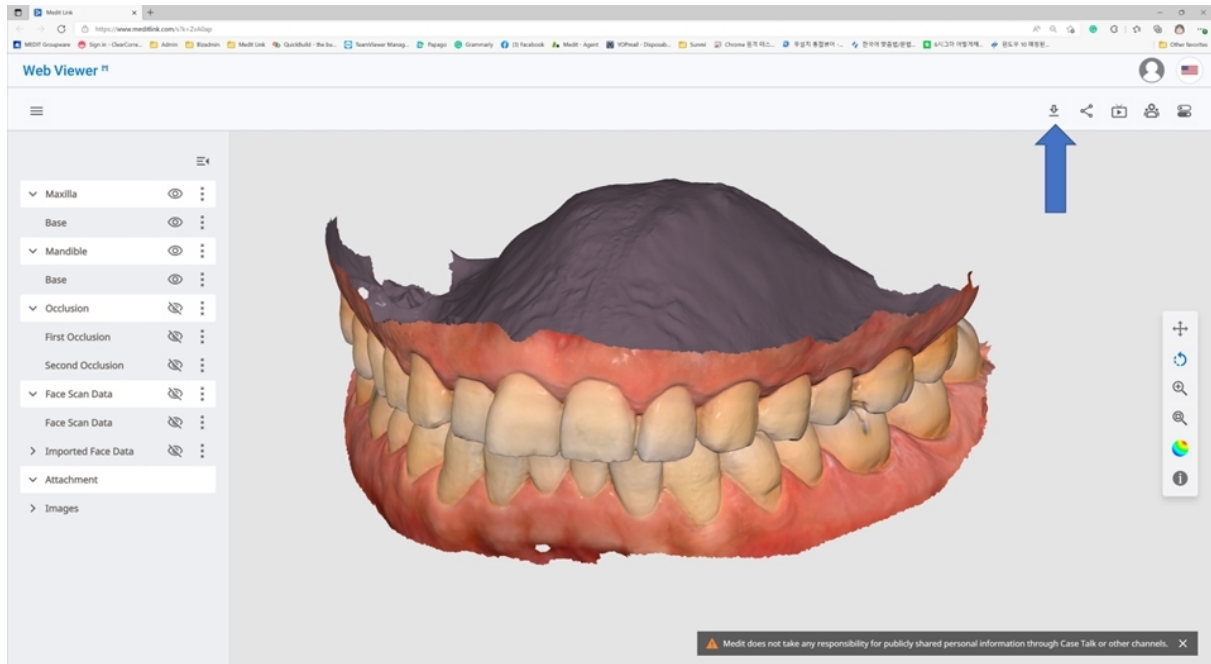

1. Click "Download."



2. Select the mesh format from the download dialog.

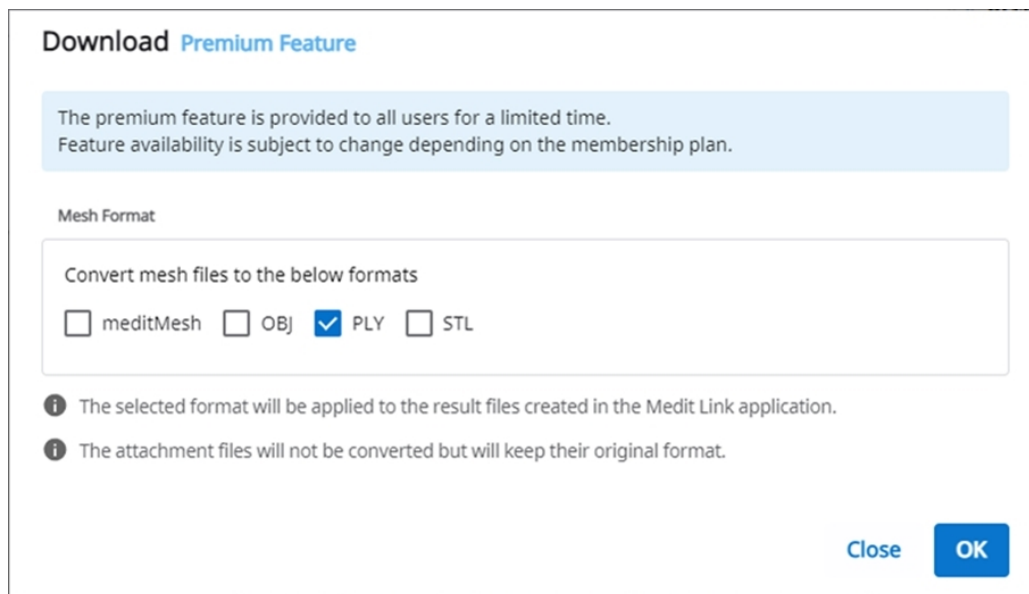

3. Then click "OK."

Order/Manage Orders

You can easily order a case with data acquired from connected programs and attached by the user.

Ordering Cases

You can place an order for your partners in Case Box.

1. Click the "Order" button.
2. Check the information required for your order.

Order

General

Order List

File List

Partner

qa_lab_01



qa_lab_01

Uses Lab Form

San Diego, US

[92131] Residence Inn San Diego Rancho Bernardo

101111111

Case Name

0923 Medit's Case

Patient Name

0923 Medit

☒ Share patient information

Requested Delivery Date

2026-09-30

Delivery Time

오전 10:30

AM 8:00

PM 8:00

Assignee

가나다라마바사아... ✕

☒ Share assignee information

Memo

0/2048

Tags

☒ Include tags in the order

Close

OK

3. Review the selected partner. To change it, select another partner from the dropdown list.


Note

If you change the partner, any form information not supported by the new partner will be cleared.

4. Check the product information on the Order List tab.
5. Check the files to order on the File List tab.
6. Choose whether to share patient information.
7. Enter the requested delivery date and time.
 - The default delivery date is set to 7 days after the order date.
8. Select assignees for the case. If you don't want to share assignee information with the lab, clear the "Share assignee information" checkbox.
9. Add a memo or tags as needed to provide additional information to the partner.
10. When multiple occlusal results are generated from the scan app, you can send all occlusal results at once by selecting the "Order All Occlusal Relationship" checkbox.
11. Proceed with your order.


Confirm Order (In Box)

- Accept the order upon receiving it.
- Accepted orders will be forwarded to the Work Box.

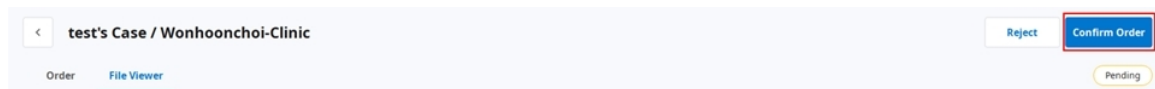


Cancel Order (Order Box)

- You can cancel an order before your partner accepts it.
- Canceled cases can be re-ordered.

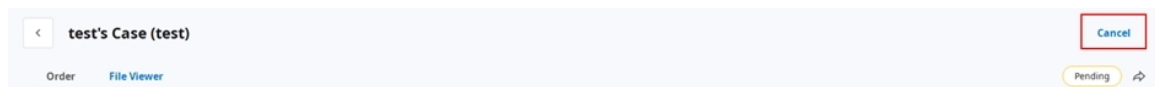

Reject Order (In Box)

- You can reject a received order.
- Rejected cases can be re-ordered.

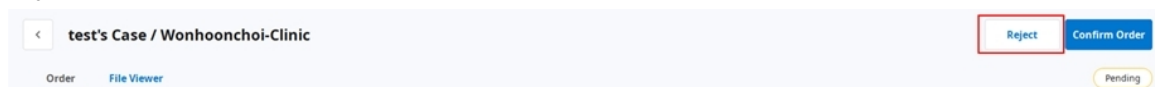

Checking Ordered Cases

You can check the details of the case you ordered from the Order Box in the clinic account. You can also check the case you received from the In Box in the lab account.

The screenshot shows a web application interface for managing dental cases. The main header is 'test1's Case (test1)' with a 'Cancel Order' link. Below the header, there are two tabs: 'Order' and 'File Viewer'. The 'Order' tab is active, displaying a form with the following fields:

- Order ID:** 6004858
- Assignee:** (dropdown menu)
- Patient Name:** test1
- Partner Name:** Elysee Dental Test
- Scanned Date:** (calendar icon)
- Order Date:** 2025-07-23 PM 5:17
- Requested Delivery Date:** 2025-07-30 PM 05:00
- Expected Delivery Date:** 2025-07-30 PM 05:00
- Memo:** (text area)
- Tags:** test

On the right side, there is a table titled 'Order List' and 'Shipment'.

No.	Type	Method	Material	Shad
15	Crown	Anatomic	Zirconia	A4

Order Info	The information that can be included in the order is as follows: order ID/assignee/patient name/account name/scan date/order date/delivery date/delivery due date/memo. You can check the details of the ordered product and delivery information.
Order List	The order list displays the details of the ordered product.
Shipment	Enter/check the shipping information.

Entering Shipment Information (In Box)

Enter the shipping information of the completed case and start shipping.

1. Click the Shipment tab.
2. Enter the shipment information and click "Register."

Shipment Registration

Select a delivery method and enter the required information to register.
The information you enter is saved in the Shipment tab and can be checked at any time.

☒ Shipment ☐ Self-Delivery

Company

Tracking Number

Memo

Cancel

Register

3. The entered information is sent to your partner.

Managing Delivery Date

You can set the estimated delivery date when placing an order, and this information will be forwarded to your partner.

Setting Requested Delivery Date (Clinic account)

- When ordering from Case Box, you can enter the requested delivery date.

Order

General Order List File List

Partner
WonhoonChoi_Lab_Stage

 **WonhoonChoi_Lab_Stage**
Seoul, KR
[02835] 23 Incheon-ro 22-gil, Seongbuk-gu, Seoul, Korea
01026088546

Case Name
mandibular movement - Clone - Clone

Patient Name
Test Occlusion

☒ Share patient information

Requested Delivery Date
2023-04-13

Delivery Time
오전 8:15 AM 8:00 PM 8:00

Memo

Tags
You can add the entered tag with Enter or Spacebar.

☒ Include all occlusion relationships in the order
☒ Include Tags In The Order

Cancel OK

- The default value is 7 days after the date of the order.
- Enter the date/time you wish to receive the products and proceed with the order.
- You can set the time easily using the slider.

Setting Estimated Delivery Date (Lab account)

- You can enter the estimated delivery date for the case ordered from In Box.

Dashboard < 전철아의 케이스 / Meditcompany

Reject Confirm Order

Pending

Order File Viewer

Order ID
6004719

Assignee
No one has been assigned.

Patient Name

Partner Name
Meditcompany

Scanned Date

Order Date
2025-07-16 PM 5:46

Requested Delivery Date
2025-07-17 PM 03:00

Expected Delivery Date
2025-07-17 PM 03:00

Memo

Tags
tag

Order List Shipment

No.	Type	Method	Material	Shad
25	Crown	Anatomic	Multi-Layered Zirconia	A3.5
26	Crown	Anatomic	Multi-Layered Zirconia	A3.5

1. Check the requested delivery date with the partner and reschedule using the "Edit" button.



2. Enter the estimated delivery date considering the actual working day. If required, you can leave a comment also.
3. Click "Apply."

The screenshot shows a date selection interface for March 2022. The calendar grid has the 17th highlighted with a blue border. Below the calendar, there is a 'Delivery Time' section with a clock icon and a time input field showing '11:45 AM'. A horizontal slider below the input field ranges from 'AM 8:00' to 'PM 8:00', with a blue dot indicating the current time. Below the time section is a 'Comment' text area. At the bottom of the form are two buttons: 'Cancel' and 'Apply'.

4. The estimated delivery date and the memo will be forwarded together with the order.
 - You can check the date and memos in the Order Box.
 - You can change the estimated delivery date, and the changes will be forwarded to your partner.
5. You can check the case history using the "History" button.



Completing Case

Complete (Work Box/In Box)


Once all tasks are done, click "Done" to complete the case.

Complete a case from the case detail

You can find the three dots icon on the right from the "Open" or "Rescan" button.



Complete the work by clicking the icon and selecting "Done."

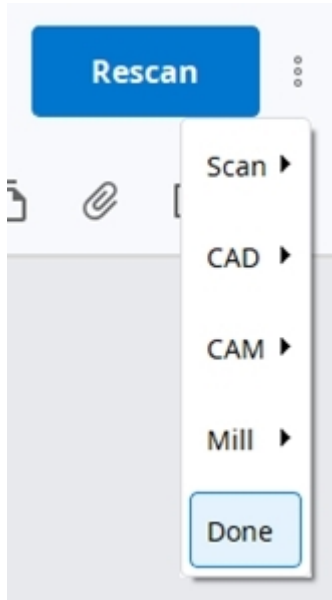

Complete the multiple cases from the case list

1. Go to Work Box.
2. Select all cases to complete.

MEDIT Link 3.2.0

Work Box

Dashboard | In Box | **Work Box** | Case Talk | App Box | Trash Box

✓ All | Form Only | Need Processing | Scan Completed | CAD | CAM

Export ▾ Status ▾ Selected Cases: 4

Done

Status	Case Name	Patient Name	Form Information
✓ Form Only	tsstt's Case	tsstt	-
✓ Form Only	test's Case	test	-
✓ CAD	Maxilla Crown	Scan for Clinics	14-Crown
✓ CAD	Maxilla Crown	Scan for Clinics	14-Crown


3. Click "Status" and select "Done."

Complete (Case Box/Order Box)

After the case is completed on the lab side, the clinic needs to confirm the case to complete the whole order process.

The clinic can find the "Confirm Delivery" button in the Order Box if the lab completes its work.

Click this, the whole order process completes, and the case status will be changed to "Case Completed."

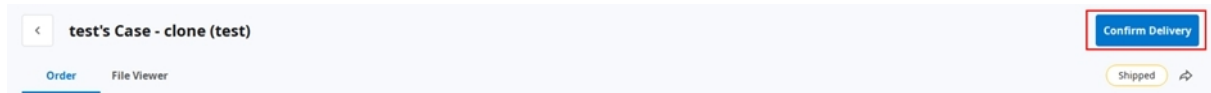

Complete a non-ordered case in the Case Box

For cases that are not ordered to the Lab, you can complete these cases from the Case Box.

1. Go to the Case Box.
2. Select the case that you want to complete.

A screenshot of the "Case Box" interface. On the left is a sidebar with icons for "Case Box", "Order Box", "Case Talk", "App Box", and "Trash Box". The "Case Box" is selected. At the top, there are filter buttons: "All" (selected), "Form Only", "Need Processing", "Scan Completed", "CAD", "Ordered", and "Case Completed". Below the filters, it says "Total: 1093". A table lists cases with columns for "Status", "Case Name", and "Patient Name". The first row shows a case with status "Scan Completed" and name "0218's Case". Below this is a date separator "2024-12-23". The following rows show cases with status "Ordered" or "Scan Completed" and names starting with "ORG_[Clinics][3.3.6] 01 Normal - Complete".

Status	Case Name	Patient Name
Scan Completed	0218's Case	0205 test2mn(456)
2024-12-23		
Ordered	ORG_[Clinics][3.3.6] 01 Normal - Complete	ORG_[Clinics][3.3.6] 01 Normal
Scan Completed	ORG_[Clinics][3.3.6] 01 Normal - Complete	ORG_[Clinics][3.3.6] 01 Normal
Scan Completed	ORG_[Clinics][3.3.6] 01 Normal - Complete	ORG_[Clinics][3.3.6] 01 Normal
Scan Completed	ORG_[Clinics][3.3.6] 01 Normal - Complete	ORG_[Clinics][3.3.6] 01 Normal

3. Click "Status" in the action bar at the bottom of the screen and select "Case Completed."

A screenshot of the bottom action bar of the application. It contains several buttons: "3 selected", "Select All", "Export", "Status", and a trash icon. The "Status" button is highlighted, and a dropdown menu is open, showing the option "Case Completed" with a red arrow pointing to it. Above the action bar, the text "BS9" and "25-Crown / 26-C" is visible.

BS9 25-Crown / 26-C


3 selected | Select All | Export | Status | D

Case Completed

Communicate Using Tags

You can add additional information for your case using tags and send them to your partners.

Add Tags When Creating a New Case You can enter tags when you create a new case or edit case information.

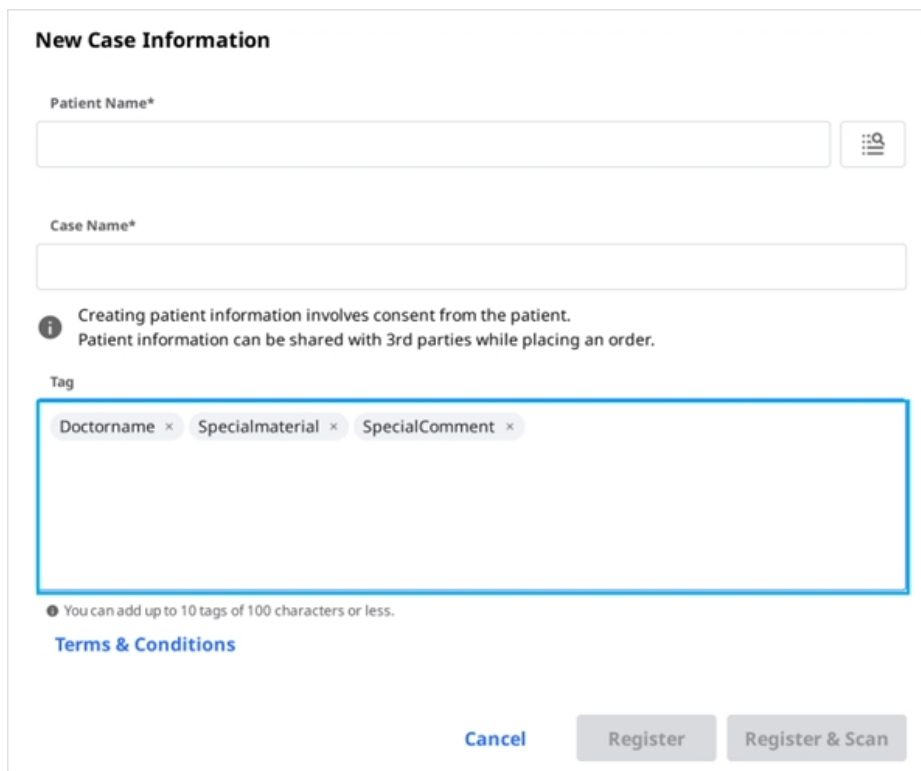

New Case Information

Patient Name*

Case Name*

*Creating patient information involves consent from the patient.
Patient information can be shared with 3rd parties while placing an order.*

Tag

Doctorname × Specialmaterial × SpecialComment ×

You can add up to 10 tags of 100 characters or less.

[Terms & Conditions](#)

[Cancel](#) [Register](#) [Register & Scan](#)

You can easily search and select tags you've used before.

New Case Information

Patient Name*

⌵

Case Name*

i

Creating patient information involves consent from the patient.
Patient information can be shared with 3rd parties while placing an order.

Tag

t

t (new)

Suggestions

tet

tagtest

tag

tat2

Add Tags When Ordering

186

You can check and edit tags in the Order dialog while ordering.

Order

General Order List File List

Partner
WonhoonChoi_Lab_Stage

 **WonhoonChoi_Lab_Stage**
Seoul, KR
[02835] 23 Incheon-ro 22-gil, Seongbuk-gu, Seoul, Korea
01026088546

Case Name
mandibular movement - Clone - Clone

Patient Name
Test Occlusion

☒ Share patient information

Requested Delivery Date
2023-04-13

Delivery Time
오전 8:15 AM 8:00 PM 8:00

Memo

Tags
Doctorname × Specialmaterial × SpecialComment ×

● You can add up to 10 tags of 100 characters or less.

☒ Include all occlusion relationships in the order

☒ Include Tags In The Order

Cancel OK

Check Tags on Case Details

MEDIT Link 3.2.0

wonhoon choi

Order Open

Scan Completed

Occlusal Relationship 1

Memo Edit

Tags Cancel Save

Dr_Name × Worker_Name ×

● You can add up to 10 tags of 100 characters or less.

+Z

eIFU download link:

<https://support.medit.com/hc/en-us/articles/53571022051737-Medit-Apps-PDF>

Medit webpage:

<https://www.medit.com>



Meditrial Srl

Via Po 9 00198, Rome Italy

ecrep@meditrial.eu



Meditrial Europe Ltd

Bahnhofstrasse 23 6300 Zug, Switzerland



Medit Corp.

9F, 10F, 13F, 14F, 16F, 8, Yangpyeong-ro 25-gil, Yeongdeungpo-gu, Seoul, 07207, Republic of Korea

Tel: +82-2-2193-9600

Contact for Product Support

Email: support@medit.com

Tel: +82-2-2193-9600