MEDIT Link

ME-UG-702H
Revision 13 (29.09.2026)
SW version 3.5.0

Table of contents

Medit Link

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Administrator

Dentist, Manager

Have the authority to create a new dental clinic account and have full access to Medit Link.

- ✓ Scan/Order/Manage orders
- ✓ Manage members

Register

Lab

Administrator

Laboratory Director, Manager

Have the authority to create a new dental lab account and have full access to Medit Link.

- ✓ Scan/Manage ordered cases
- ✓ Manage members

Register

Member

Dental Hygienist, Dental Staff

Belong to the dental clinic their admin created and have limited access to Medit Link after the admin's approval.

- ✓ Scan/Order/Manage orders

Register

Member

Dental Technician, Laboratory Staff

Belong to the dental lab their admin created and have limited access to Medit Link after the admin's approval.

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Email Address

Check

 Please ensure your email address is active. Medit Link will send a verification email.

 Registering with a disposable and/or fake e-mail address can make you vulnerable to security issues and make it difficult to receive information from Medit Link.

Password

 Password must be 8 to 16 characters long, including a combination of at least three letters, numbers, and special characters. (Passwords are accepted in English language only.)

Confirm Password

 Please type in the same password as above.

Select Country/Region

☐ I agree to all of the Terms and Conditions.

☐ Personal Information Collection Agreement (Mandatory)

View

☐ Terms & Conditions (Mandatory)

View

☐ Agree to the Marketing and Promotion Agreement (optional)

View

 Check this option to receive marketing emails, which include product updates, software news, and promotional events.

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Settings

Receive the report via

- ☒ Medit Link account
- ☐ Another email address

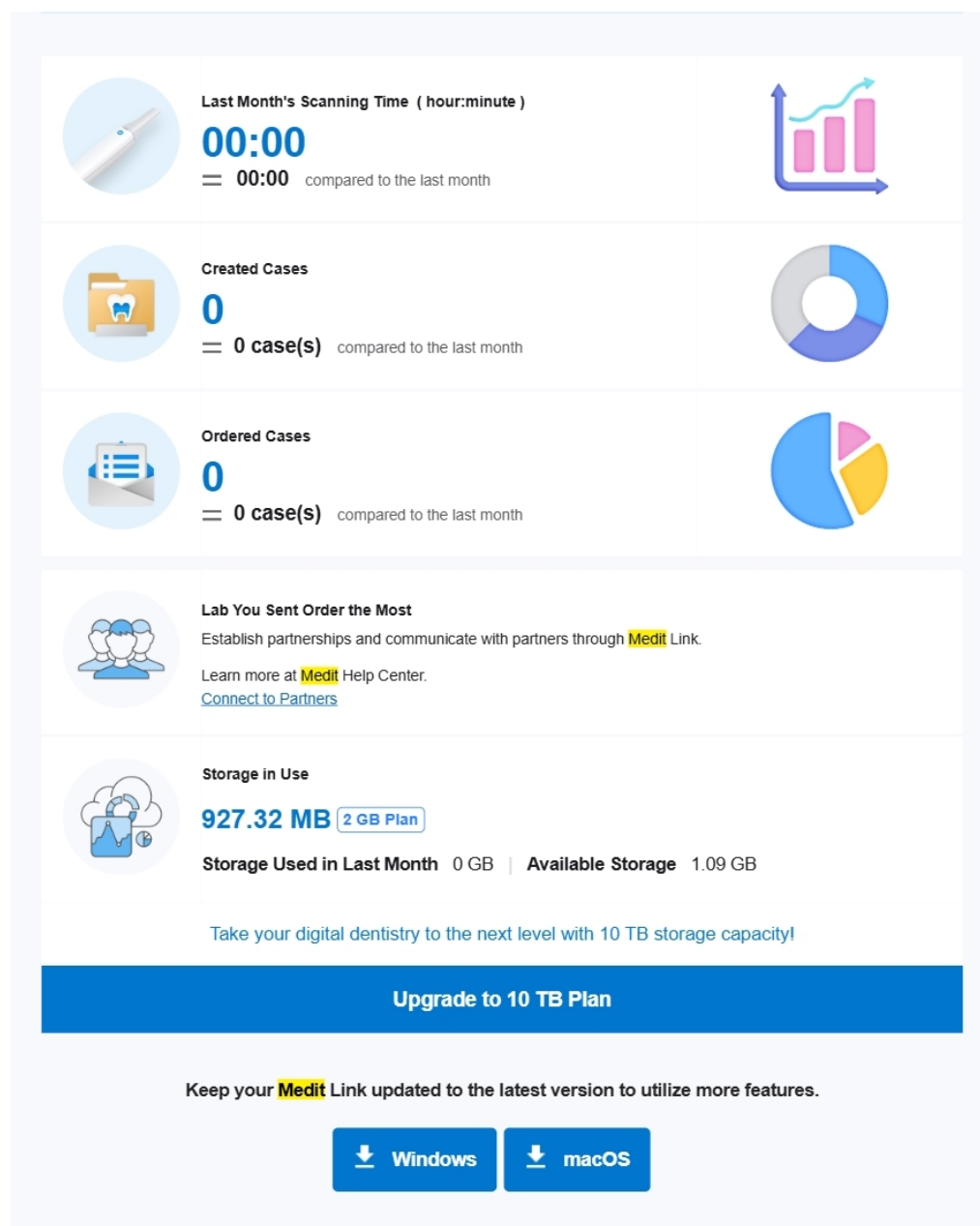
Enter another email address to receive the report email.

The report will be sent on the 2nd day of every month.

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This is an automatically generated email. Please do not reply to this email.

Please contact us at insights@medit.com for more information.

[Unsubscribe](#) If you no longer wish to receive this email, click "Unsubscribe." You can change the email subscription settings in [Medit Link Dashboard](#).

Copyright © Medit Corp. All rights reserved.

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The screenshot shows the Medit Link Case Box interface. At the top, there's a header with the Medit Link logo, a user profile (seolahjeon), and a language dropdown (English). Below the header, the 'Case Box' section is active, showing a filter for 'Patient'. A date range 'Last 90 Days / 2025-12-20 ~ 2026-03-20' and a search bar 'Case Name' with a keyword input field are visible. A row of filter buttons includes 'All' (selected), 'Form Only', 'Need Processing', 'Scan Completed', 'CAD', 'Ordered', and 'Case Completed'. Below these, it says 'Total: 18'. A table with 10 columns (Status, Case Name, Patient Name, Form Information, Last Modified..., Scanned Date, Order Date, Lab Name, Assignee) displays a list of cases. Each row includes a status button (e.g., 'Scan Completed', 'Form Only', 'Need Processing') and a 'Share' button. The footer contains a navigation menu with 'Home', 'About', 'Terms & Conditions', 'Privacy Policy', and 'Copyright © Medit Corp. All rights reserved.'

Status	Case Name	Patient Name	Form Information	Last Modified...	Scanned Date	Order Date	Lab Name	Assignee
Scan Completed	Ortho Scan - Clone - Clo...	Randee Glass		2026-03-18 16:52	2026-03-10 00:11	-	-	
Scan Completed	Ortho Scan - Clone	Randee Glass		2026-03-18 14:19	2026-03-10 00:11	-	-	
Form Only	Medit Link의 케이스	Medit Link		2026-03-18 14:17	-	-	-	
Scan Completed	Ortho Scan	Randee Glass		2026-03-17 08:57	2026-03-10 00:11	-	-	
Need Processing	Express_i900의 케이스	Express_i900	23	2026-03-12 13:49	2026-03-11 14:24	-	-	
Form Only	Express Scan	seolah		2026-03-10 14:53	2026-03-10 14:53	-	-	
Form Only	Express Scan	seolah		2026-02-24 08:33	2026-02-24 08:33	-	-	
Need Processing	Express Scan	seolah		2026-02-23 16:38	2026-02-23 16:37	-	-	
Scan Completed	Express Scan	seolah	Maxilla	2026-02-23 15:50	2026-02-23 15:50	-	-	

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The screenshot displays the 'Order Box' section of the Medit Link application. It features a sidebar with navigation icons for Dashboard, Patient, Case Box, Order Box, Case Talk, Partners, Membership, and My Info. The main area shows a table of orders with columns: Status, Case Name, Patient Name, Order Date, Requested Deliv..., Last Modified..., Lab Name, Order ID, and Assignee. A filter bar at the top allows searching by Case Name and filtering by status (All, Pending, Canceled, Rejected, Accepted, Shipped, Case Completed). A red '1' points to the status filter buttons, a red '2' points to the table headers, and a red '3' points to the search input field.

Status	Case Name	Patient Name	Order Date	Requested Deliv...	Last Modified...	Lab Name	Order ID	Assignee
Pending	Express Scan	Seol register	2026-04-02 21:38	2026-04-09 09:00	2026-04-02 21:38	MEDIT_LAB_STAGE	6012242	

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Invite Partner

Invite your desired partner to join Medit Link. When the lab accepts your invitation, a partnership is automatically formed.

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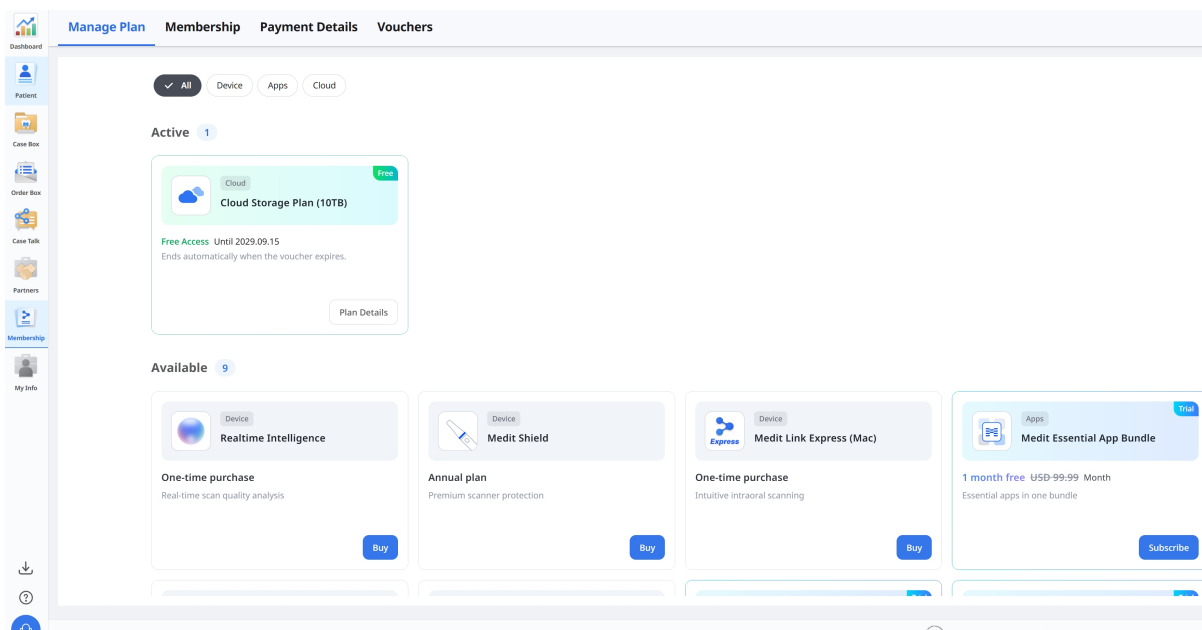
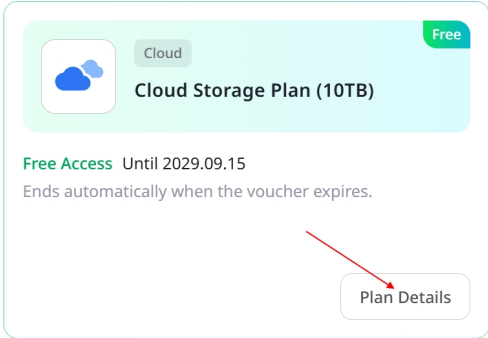
manage plan **membership** **Payment Details** **vouchers**

✓ All Device Apps Cloud

Active 1

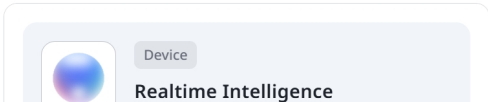
Cloud Storage Plan (10TB)

Free Access Until 2029.09.15
Ends automatically when the voucher expires.

Plan Details

Available 9

Realtime Intelligence

Plan Details

Free Access
Your free access is active.

Product	Cloud Storage Plan (10TB)
Status	Free
Expiration Date	2029.09.15
Regular Price	USD 9.99/Month

⚠ Important

- You can't subscribe to this product while the voucher is active.
- Free access ends automatically when the voucher expires. A paid subscription is required to continue.

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The screenshot displays a user's subscription management dashboard. At the top, there are four tabs: 'Manage Plan', 'Membership', 'Payment Details', and 'Vouchers'. Below these tabs are three filter buttons: 'All' (selected), 'Device', and 'Apps'. The 'Active' section shows 0 subscriptions with a message: 'No active subscriptions. Browse our services below to get started.' The 'Available' section shows 8 subscriptions. A red arrow points from the 'Available' count to a specific plan card. The card is for the 'Cloud Storage Plan (10TB)' under the 'Cloud' category. It features a 'Trial' badge, a price of '1 month free' (with 'USD 9.99 Month' crossed out), and the text 'Secure 10 TB cloud storage'. A blue 'Subscribe' button is located at the bottom right of the card.

Manage Plan Membership Payment Details Vouchers

✓ All Device Apps Cloud

Active 0

No active subscriptions.
Browse our services below to get started.

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~~USD 9.99~~ Month

✓ Includes 10 TB of cloud storage

✓ Securely store large cases and other data

✓ Streamline your Medit Link data management in the cloud

✓ Store design files without separate local backups

• Actual charges may vary depending on your bank or card issuer's fees and policies.

• For more information about Cloud Storage Plan (10TB), see [here](#).

Trial

Cloud Storage: First Month Free

1 month free >

☐ I agree to automatic renewal and recurring charges.

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The image shows two parts of a subscription process. The top part is a blue payment summary screen titled 'Subscribe to Medit Cloud Storage 10TB Plan'. It shows a price of \$0.00 with a note 'Then \$9.99 per day after coupon expires'. Below this, it lists 'Medit Cloud Storage 10TB Plan for Clinics' at \$9.99, billed daily. A 'Subtotal' of \$9.99 is shown, followed by a '1 Month Free for New Subscribers' coupon for -\$9.99. The 'Tax' is \$0.00, and the 'Total due today' is \$0.00. The bottom part is a 'Pay with card' form with fields for Email (yhgam0829@gmail.com), Card information (1234 1234 1234 1234), MM / YY, CVC, Cardholder name (Full name on card), and Country or region (South Korea). A 'Subscribe' button is at the bottom. A red arrow points from the 'Subscribe' button to the 'Confirm' button in the modal below.

Subscribe

Your Cloud Storage Plan (10TB) subscription is now active.
Enjoy full access to all features.

Amount of Payment	Free
Next Payment Date	2026.10.23
Last Payment Method	****

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The image shows a 'Manage Plan' modal. It has a title bar with 'Manage Plan' and a close button. Below the title bar, there is a 'Promotion Applied' section with the text 'Regular pricing will apply after the promotion ends.' Below this, there is a table with subscription details: Product (Cloud Storage Plan (10TB)), Price (Free ~~USD 9.99~~/Month), Expiration Date (2026.10.23), and Last Payment Method (**** [Edit](#)). At the bottom, there are two buttons: 'Close' and 'Unsubscribe'. A red arrow points to the 'Unsubscribe' button.

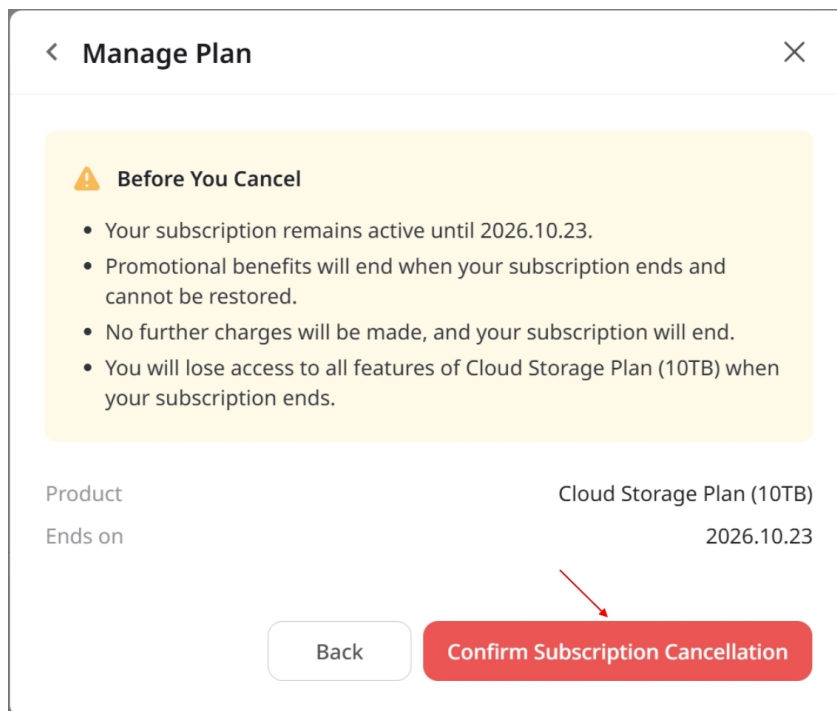
Manage Plan

Promotion Applied
Regular pricing will apply after the promotion ends.

Product	Cloud Storage Plan (10TB)
Price	Free USD 9.99 /Month
Expiration Date	2026.10.23
Last Payment Method	**** Edit

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The screenshot shows a 'Manage Plan' dialog box with a close button (X) in the top right corner. Inside, there is a yellow warning box titled 'Before You Cancel' with a warning icon. It contains four bullet points: 'Your subscription remains active until 2026.10.23.', 'Promotional benefits will end when your subscription ends and cannot be restored.', 'No further charges will be made, and your subscription will end.', and 'You will lose access to all features of Cloud Storage Plan (10TB) when your subscription ends.' Below this, the 'Product' is listed as 'Cloud Storage Plan (10TB)' and 'Ends on' is '2026.10.23'. At the bottom, there are two buttons: 'Back' and 'Confirm Subscription Cancellation'. A red arrow points to the 'Confirm Subscription Cancellation' button.

< Manage Plan X

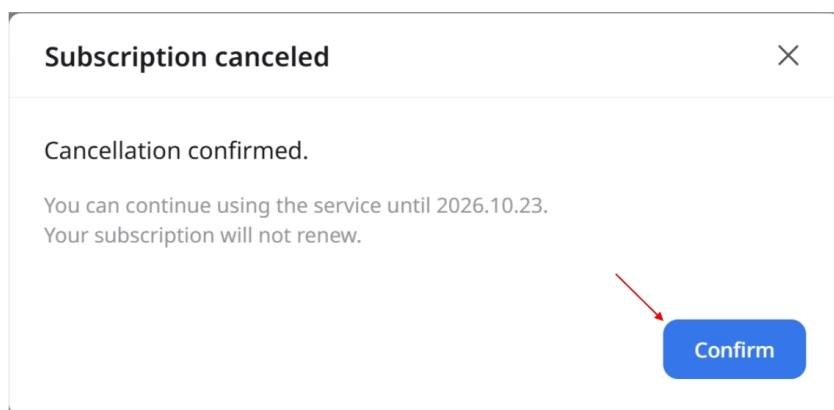
⚠ Before You Cancel

- Your subscription remains active until 2026.10.23.
- Promotional benefits will end when your subscription ends and cannot be restored.
- No further charges will be made, and your subscription will end.
- You will lose access to all features of Cloud Storage Plan (10TB) when your subscription ends.

Product Cloud Storage Plan (10TB)
Ends on 2026.10.23

Back Confirm Subscription Cancellation

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The screenshot shows a 'Subscription canceled' dialog box with a close button (X) in the top right corner. The text inside says 'Cancellation confirmed.' followed by 'You can continue using the service until 2026.10.23.' and 'Your subscription will not renew.' At the bottom right, there is a blue 'Confirm' button. A red arrow points to the 'Confirm' button.

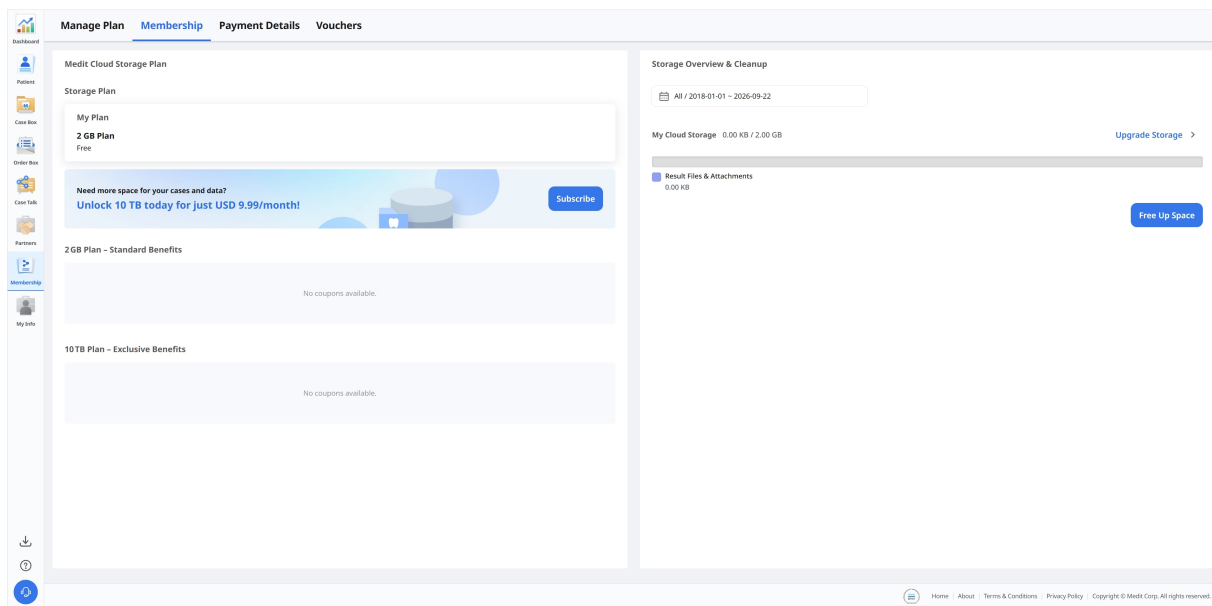
Subscription canceled X

Cancellation confirmed.

You can continue using the service until 2026.10.23.
Your subscription will not renew.

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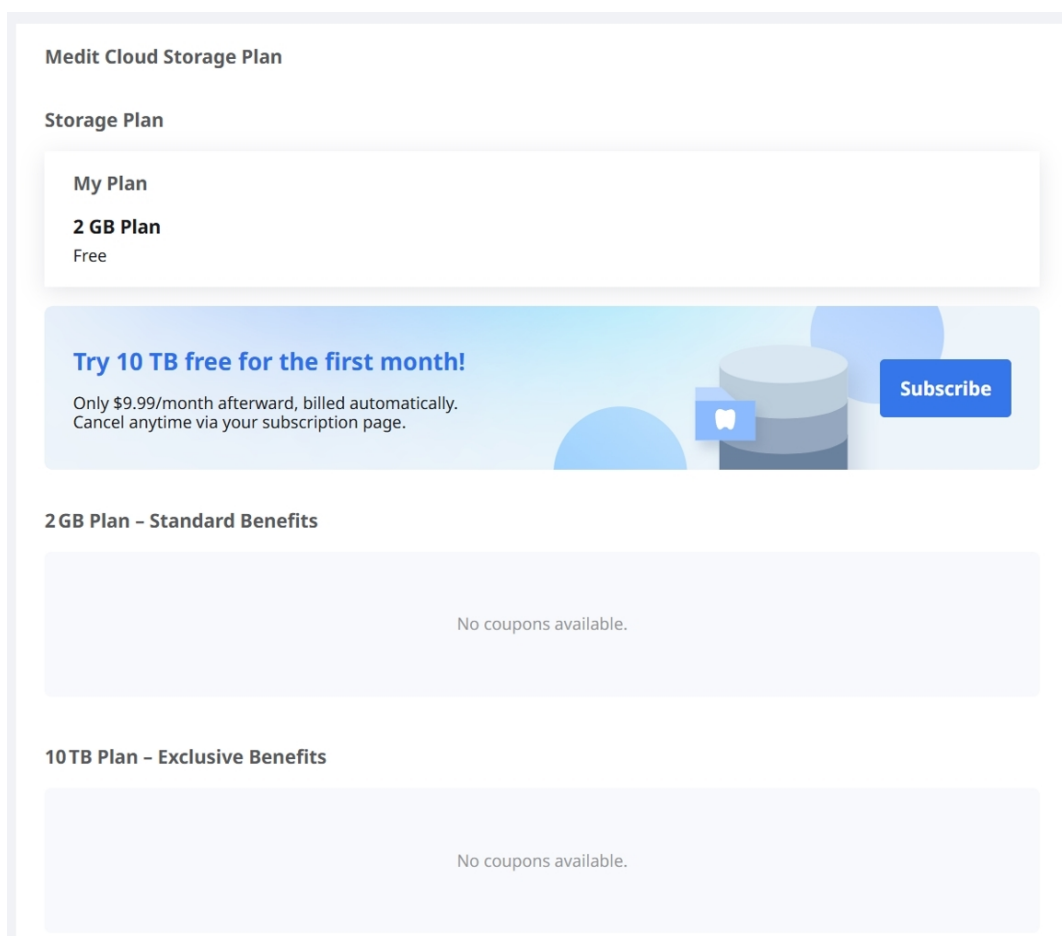
Storage Overview & Cleanup

 2018-01-01 ~ 2025-11-05

My Cloud Storage 927.32 MB / 2.00 GB [Upgrade Storage >](#)

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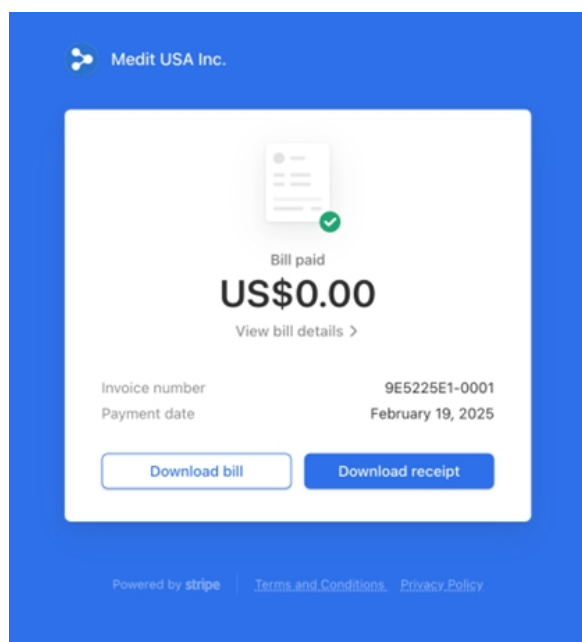
Payment Date	Product	Payment Type	Duration	Details	Amount of Payment	Status
2026-09-21	Cloud Storage Plan (10TB)	Subscription	2026-09-21 - 2026-10-21	Monthly payment	USD 0.00	Payment Completed
2026-09-21	Cloud Storage Plan (10TB)	Subscription	2026-09-21 - 2026-10-21	Monthly payment	USD 0.00	Payment Completed
2026-09-21	Cloud Storage Plan (10TB)	Subscription	2026-09-21 - 2026-10-21	Monthly payment	USD 0.00	Payment Completed
2026-09-21	Cloud Storage Plan (10TB)	Subscription	2026-09-21 - 2026-10-21	Monthly payment	USD 0.00	Payment Completed
2026-09-21	Cloud Storage Plan (10TB)	Subscription	2026-09-21 - 2026-10-21	Monthly payment	USD 0.00	Payment Completed

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ITEM

FEBRUARY 19 - FEBRUARY 20, 2025

Medit Cloud Storage 10TB Plan for Clinics

Quantity 1 piece

US\$0.60

subtotal

US\$0.60

1 Month Free for New Subscribers(100% 할인)

-US\$0.60

Payment amount

US\$0.00

Remaining amount

US\$0.00

Do you have any questions?[Contact Medit USA Inc.](#)

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All

▼

Product	Promotion	Code Number
Cloud Storage Plan	Redeem SQA_Redeem_Cloud Storage Plan_36M	GFZl

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30 Days

⚠ Important

- Once activated, the voucher cannot be canceled.
- You can't subscribe to this product while the voucher is active.

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The screenshot shows a web browser window with the URL www.meditlink.com/myinfo/account. The page is titled "MEDIT Link" and has a navigation bar with "My Account", "My Clinic", and "My Scanner". On the left is a sidebar with icons for Patient, Case Box, Order Box, Case Talk, Partners, Membership, and My Info (which is highlighted). The main content area is titled "Account Information" and features a circular MEDIT logo with the text "2022 Trusted Lab Partner of MEDIT". To the right of the logo are several form fields: "Email Address" (seolah.jeon@medit.com) with a "Change Email Address" button, "Password" with a "Change Password" button, "Name" (seolah.jeon_dev_admin), "Language" (English), and "Country/Region" (Korea, Republic of). A note at the bottom states: "Your chosen language is applied to the email sent by MEDIT Link and the notification."

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Edit

Cover and Logo

Looking for digitalized dental laboratories?

Find the most ordered Labs around your clinic with 'Featured Partners'

123

Services

Select

Specialties (choose up to 3)

Select

Select the services your business specializes in from the list.

About Us

e.g., Welcome to OOO Dental Clinic.
We are a licensed dental clinic situated in California specializing mainly in aesthetic restorations and implantology.

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The scanner information is automatically registered once your scanner is connected to a Medit scan program.

Home | About | Terms & Conditions | Privacy Policy | Copyright © Medit Corp. All rights reserved.

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MEDIT Link

My Account My Lab My Scanner **Order Forms**

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Methods ¹

Search

Screw-Retained Full-Contour Offset Cement-Retained

SCRP Temporary Cut-back Pontic Press Preform Screw-Retained Layered

Screw-Retained PFZ Screw-Retained PFG Screw-Retained PFM Cement-Retained Layered

[More](#)

☐ Do not use

Full-Contour

Materials ²

Search

Multilayer Zirconia Zirconia PMMA Gold Resin

LiSi Glass Ceramic Press Ceramic Hybrid Ceramic Composite Resin e.max

Titanium Metal Stainless Steel Co-Cr (Cobalt Chromium) Ni-Cr (Nickel-Chromium) Wax

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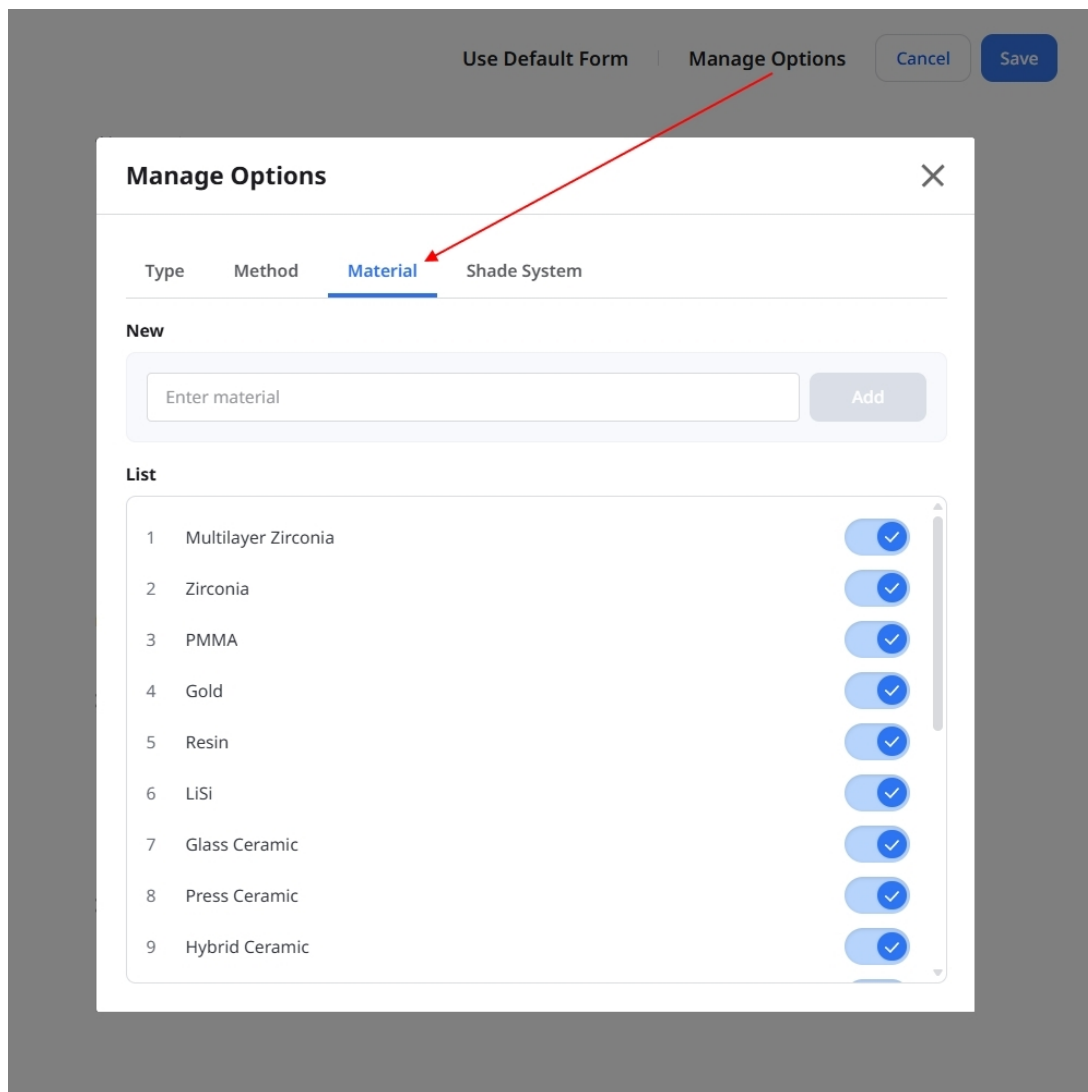
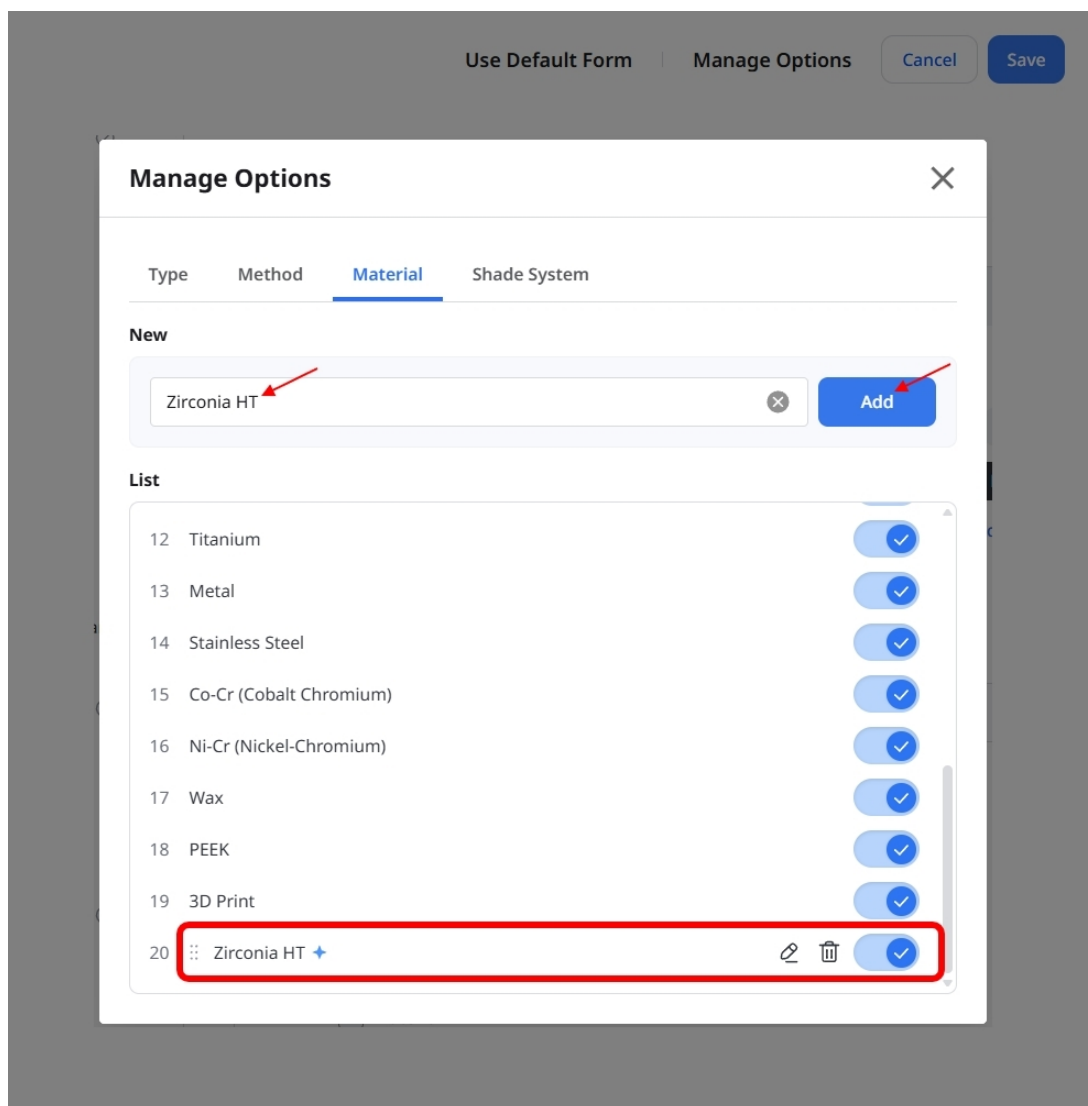
The screenshot shows a 'Manage Options' dialog box with a dark gray background. At the top, there are three buttons: 'Use Default Form', 'Manage Options' (highlighted with a red arrow), and 'Cancel'. To the right of 'Manage Options' is a blue 'Save' button. The dialog box itself has a title bar with 'Manage Options' and a close button (X). Inside, there are four tabs: 'Type', 'Method', 'Material' (selected and underlined), and 'Shade System'. Below the tabs, there is a 'New' section with a text input field labeled 'Enter material' and an 'Add' button. Below that is a 'List' section containing a scrollable list of materials, each with a blue toggle switch and a checkmark icon. The materials listed are: 1 Multilayer Zirconia, 2 Zirconia, 3 PMMA, 4 Gold, 5 Resin, 6 LiSi, 7 Glass Ceramic, 8 Press Ceramic, and 9 Hybrid Ceramic.

Material	Status
1 Multilayer Zirconia	✓
2 Zirconia	✓
3 PMMA	✓
4 Gold	✓
5 Resin	✓
6 LiSi	✓
7 Glass Ceramic	✓
8 Press Ceramic	✓
9 Hybrid Ceramic	✓

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The screenshot shows the 'Crown' configuration window. Under 'Allowed Teeth', 'Anterior', 'Premolar', and 'Molar' are all checked. In the 'Methods' section, 'Full-Contour' is selected among various options like 'Screw-Retained', 'Offset', and 'Cement-Retained'. Below this, a 'Do not use' checkbox is unchecked. A dropdown menu for 'Full-Contour' is open, showing a list of materials. 'Zirconia HT' is highlighted with a red arrow and a blue plus icon, indicating it is the selected material. Other materials listed include 'Multilayer Zirconia', 'Zirconia', 'PMMA', 'Gold', 'Resin', 'LiSi', 'Glass Ceramic', 'Press Ceramic', 'Hybrid Ceramic', 'Composite Resin', 'e.max', 'Titanium', 'Metal', 'Stainless Steel', 'Co-Cr (Cobalt Chromium)', 'Ni-Cr (Nickel-Chromium)', 'Wax', 'PEEK', and '3D Print'. A '+ Add Material' link is also visible. A 'Do not use' checkbox is also present at the bottom of the materials list.

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The screenshot shows the 'New' section of the 'Manage Options' dialog box. The 'Name' field contains 'Lab Shades'. The 'Shades' field contains 'L1, L2, L3'. Below the 'Shades' field, there is a note 'Separate values with commas.' and three buttons labeled 'L1', 'L2', and 'L3'. At the bottom, there is a blue 'Add' button.

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The screenshot shows the 'Crown' configuration window. Under 'Allowed Teeth', 'Anterior', 'Premolar', and 'Molar' are all checked. In the 'Methods' section, 'Full-Contour' is selected among various options like 'Screw-Retained', 'Offset', and 'Cement-Retained'. Below this, a 'Do not use' checkbox is present. The 'Full-Contour' section is expanded, showing a 'Materials' list where 'Zirconia HT' is selected. Other materials include 'Multilayer Zirconia', 'Zirconia', 'PMMA', 'Gold', 'Resin', 'LiSi', 'Glass Ceramic', 'Press Ceramic', 'Hybrid Ceramic', 'Composite Resin', 'e.max', 'Titanium', 'Metal', 'Stainless Steel', 'Co-Cr (Cobalt Chromium)', 'Ni-Cr (Nickel-Chromium)', 'Wax', 'PEEK', and '3D Print'. A '+ Add Material' link is also visible. The 'Zirconia HT' section is expanded, showing 'Shade Systems'. 'VITA Classic' (20 Shades), 'VITA 3D-Master' (82 Shades), 'Ivoclar' (25 Shades), 'Material' (5 Shades), and 'Lab Shades' (3 Shades) are listed. 'Lab Shades' is selected, and a red arrow points to the '+' icon next to it. Below 'Lab Shades' are buttons for 'L1', 'L2', and 'L3'.

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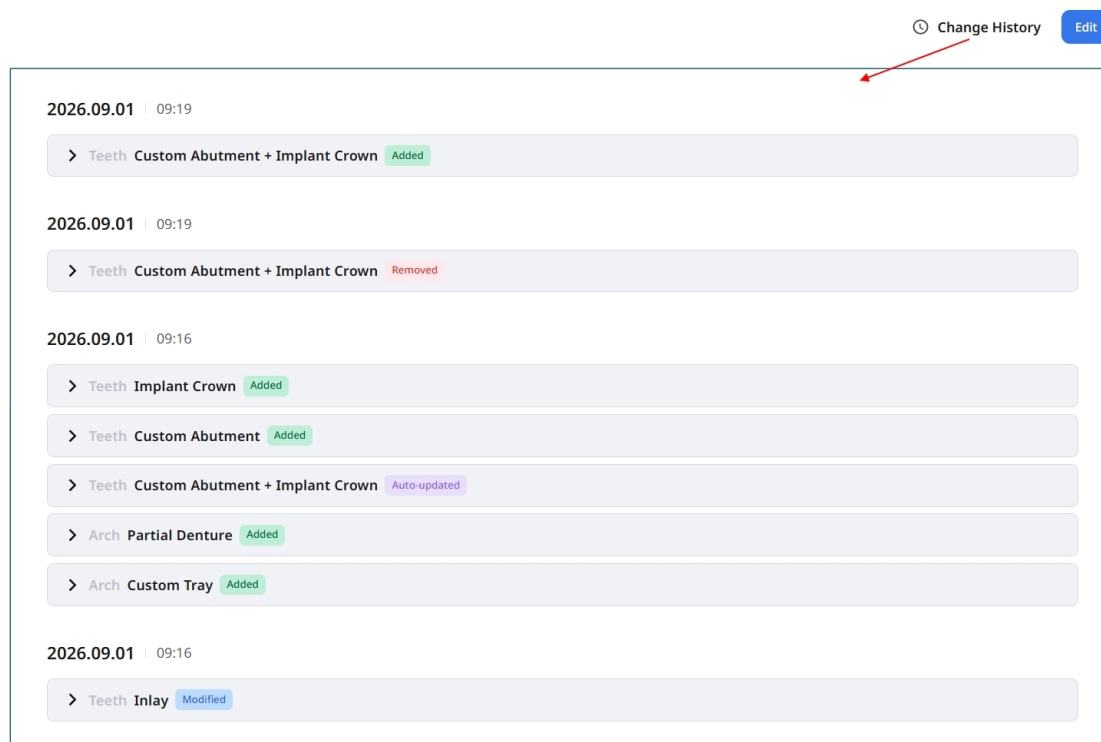
Change History Edit

2026.09.01 | 09:19

> Teeth Custom Abutment + Implant Crown Added

2026.09.01 | 09:19

> Teeth Custom Abutment + Implant Crown Removed

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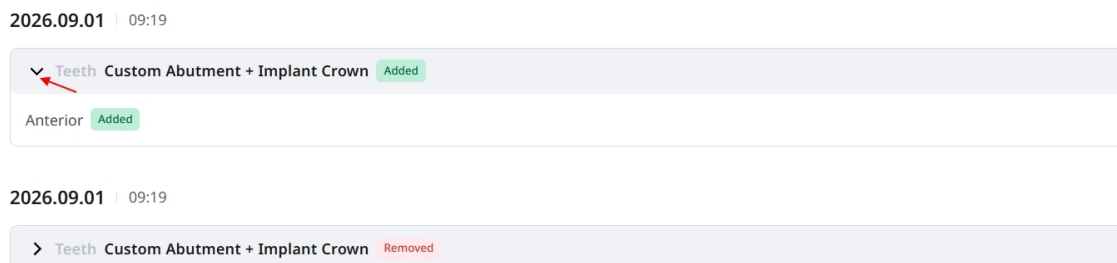
2026.09.01 | 09:19

▼ Teeth Custom Abutment + Implant Crown Added

Anterior Added

2026.09.01 | 09:19

> Teeth Custom Abutment + Implant Crown Removed

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2023-04-19

To our users,

We want to thank all Medit Link users who have supported and used our NFT service.

Unfortunately, we must inform you that this service will be terminated following the upcoming program update.

Please see the details below for more information:

- **Feature in Question**
 - Share via NFT (in Medit Link Case Talk)
- **Effective Date**
 - Starting from the update in the first half of 2023 (scheduled for April 27)

※ The update schedule is subject to change.

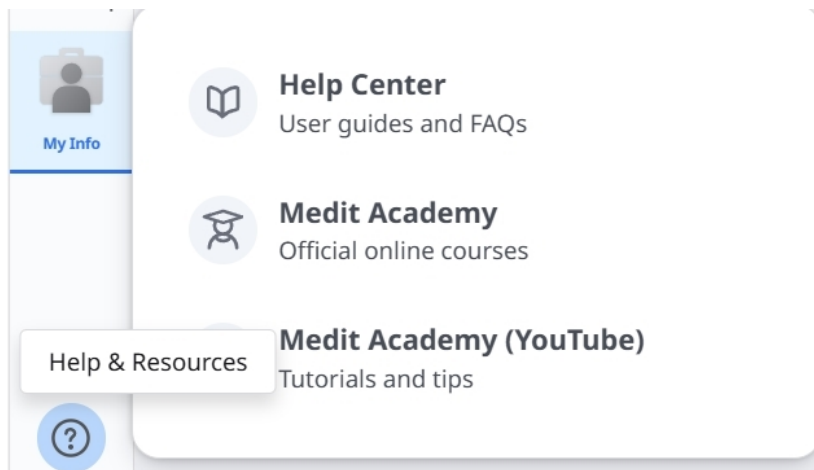
We'll continuously work to provide better service.

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The screenshot displays the Medit Link 3.4.10 Patient management interface. The top bar shows the user 'seolahjeon_' and has buttons for 'New Patient' and 'New Case'. The main area displays a list of patients with columns for Patient Name, Patient ID, Gender, Date of Birth, Number of Cases, and Registration Date. There are buttons for 'Scan' and a vertical menu icon for each patient entry. The interface also includes a sidebar with icons for Dashboard, Patient, Case Box, Order Box, Case Task, and Trash Box. The top bar shows the user 'seolahjeon_' and has buttons for 'New Patient' and 'New Case'.

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Merge all patient information to
Please select the patient you'd like to merge all data to.
After merging all data together, all information in other patient files will be deleted.

	Patient Name	Patient ID	Gender	Birth Date	Number of Cases	Date Registered	Recent Visit Date
☒	5-9_jScan_CAMOrd...				1	2020-09-18	2020-09-18
☐	test				2	2020-09-17	2020-09-17
☐	test				22	2019-07-26	2020-09-17

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Case Box 1231456 / 1231-12-04 / Female / -

All / 2018-01-01 ~ 2026-04-08 All / Contains

Form Only Need Processing Scan Completed CAD Ordered Case Completed

Total: 5

Status	Case Name	Patient Name	Form Information	Last Modified Date ↓	Scanned Date	Lab Name
☐ Form Only	1231456's Case	1231456	25-Pontic	2026-04-08 PM 4:23		
☒ Scan Completed	159's Case	1231456	16-Custom Abutment	2025-02-13 PM 4:17		
☒ Scan Completed	1231456 - clone - clone	1231456	21-Crown	2022-11-25 AM 11:27		
☒ Scan Completed	1231456 - clone	1231456	21-Crown	2022-11-25 AM 11:22		
☐ Form Only	1231456's Case	1231456	24-Crown / 25-Crown / 26-Crown / 27-Crown	2022-10-14 AM 1:19		

2 selected | Select All | Export | Status | Delete

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Total: 254

Table Options

	Status	Case Name	Patient Name	Order Date	Requested Delivery Date	Lab Name	Order ID
☐	Accepted	texture on's Case - Clone	texture on	2026-04-08 PM 5:40		MEDIT_LAB_STAGE	6012449
☒	Pending	Express Scan	Seol register	2026-04-02 PM 9:38	2026-04-09 AM 9:00	MEDIT_LAB_STAGE	6012242
☒	Case Completed	전설아의 케이스	전설아	2025-07-16 PM 5:46	2025-07-17 PM 3:00	MEDIT_LAB_STAGE	6004719
☒	Case Completed	27	1	2025-03-16 PM 3:11	2025-03-17 PM 3:00	MEDIT_LAB_STAGE	6002656
☐	Case Completed	3.3.0.136 TEST's Case	3.3.0.136 TEST	2024-04-12 PM 3:25	2024-04-13 PM 2:00	MEDIT_LAB_STAGE	5998987
☐	Case Completed	test의 케이스	Patient TEST	2024-04-09 PM 3:06	2024-04-10 PM 2:00	MEDIT_LAB_STAGE	5998955
☐	Case Completed	Mandibular Movement Case - Clone - Clone	MEDIT	2023-06-12 PM 12:12	2023-06-13 PM 2:00	MEDIT_LAB_STAGE	5996680
☐	Case Completed	Release New_2's Case	Release New_2	2023-04-25 PM 6:47	2023-04-25 PM 2:00	MEDIT_LAB_STAGE	5996180
☐	Case Completed	cX 문의's Case - Clone	cX 문의	2023-05-16 PM 5:26	2023-05-17 PM 2:00	MEDIT_LAB_STAGE	5996393
☐	Case Completed	cX 문의's Case	cX 문의	2023-05-15 PM 2:39	2023-05-16 PM 2:00	MEDIT_LAB_STAGE	5996362
☐	Case Completed	310.187's Case	310.187	2023-05-15 PM 1:33	2023-05-16 PM 2:00	MEDIT_LAB_STAGE	5996360
☐	Case Completed	310.187 test's Case - Clone			2023-05-03 PM 2:00	MEDIT_LAB_STAGE	5996306

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All

All / Contains

	Status	Case Name	Patient Name	Lab Name	Deletion Date ↓	Deleted By	Size
☐	Scan Completed	Jira_인플란트 46, 47's Case - Clone - Clone	인플란트 46, 47		2023-04-26 오전 9:12	wonhoon choi	66,654 KB
☐	Form Only	G11598492304180810429	test test		2023-04-18 오후 5:15	wonhoon choi	0 KB
☐	Form Only	G11598492304180806313	test test		2023-04-18 오후 5:15	wonhoon choi	0 KB
☐	Form Only	G11598492304180811349	tes test		2023-04-18 오후 5:15	wonhoon choi	0 KB
☐	Form Only	G11598492304180754083	Glidewell_fastdesign		2023-04-18 오후 5:00	wonhoon choi	0 KB
☐	Form Only	G11598492304180733332	Glidewell_fastdesign		2023-04-18 오후 5:00	wonhoon choi	0 KB
☐	Form Only	G11598492304180648146	test test		2023-04-18 오후 5:00	wonhoon choi	0 KB
☐	Case Completed	Multi Occlusion Normal Case's Case - Clone	Multi Occlusion Normal Case		2023-04-18 오전 10:50	wonhoon choi	23,652 KB
☐	CAD	Multi Occlusion Normal Case's Case - Clone - Clone	Multi Occlusion Normal Case		2023-04-18 오전 10:50	wonhoon choi	23,652 KB
☐	Form Only	Occlusion Test's Case	Occlusion Test		2023-04-17 오후 6:31	wonhoon choi	0 KB
☐	Form Only	Test Occlusion's Case	Test Occlusion		2023-04-17 오후 6:31	wonhoon choi	0 KB
☐	Form Only	Occlusion Test's Case	Occlusion Test		2023-04-17 오후 6:31	wonhoon choi	0 KB
☐	Form Only	Occlusion Test's Case	Occlusion Test		2023-04-17 오후 6:26	wonhoon choi	0 KB

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☒ No Date
Data Type

All
▼

☐ Case Name	Patient Name	Scanned Date ↓	Size	Status
☐ BLATRIX Charline's Case	BLATRIX Charline	2025-01-06 PM 6:28	729.30 MB	Upload Required
☐ BLATRIX Charline's Case	BLATRIX Charline	2025-01-06 PM 6:28	729.30 MB	Upload Required
☐ BLATRIX Charline's Case - 복제	BLATRIX Charline	2025-01-06 PM 6:28	733.35 MB	Upload Required
☐ BLATRIX Charline's Case - Clone	BLATRIX Charline	2025-01-06 PM 6:28	734.01 MB	Upload Required
☐ test의 케이스	Patient TEST	2023-04-26 PM 7:58	1.25 GB	Upload Required
☐ i700 HD's Case	i700 HD	2023-04-24 PM 5:23	911.24 MB	Upload Required
☐ Mandibular Movement	MEDIT	2023-04-24 PM 1:59	1.88 GB	Upload Required
☐ Mandibular Movement Case	MEDIT	2023-04-24 PM 1:59	1.75 GB	Upload Required
☐ Mandibular Movement Case - Clo...	MEDIT	2023-04-24 PM 1:59	1.75 GB	Upload Required
☐ 서라_Mandibular Movement_TEST	Mandibular Movement_export_test	2023-04-10 PM 3:39	34.73 MB	Upload Required
☐ 서라_Mandibular Movement_TEST...	Mandibular Movement_export_test	2023-04-10 PM 3:39	731.43 MB	Upload Required

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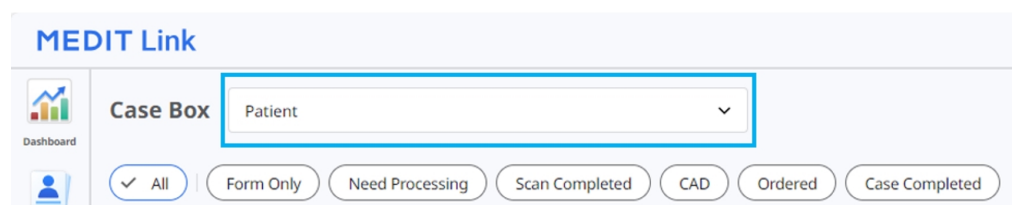
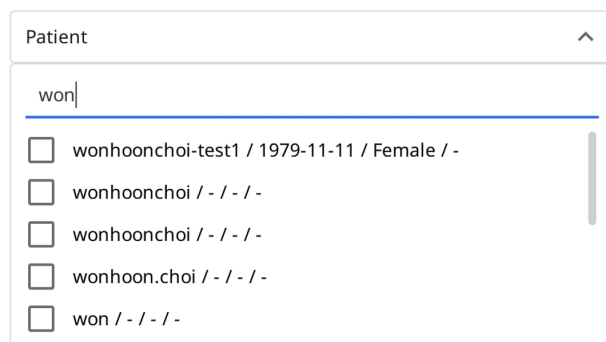
Case Box Patient

Dashboard

✓ All Form Only Need Processing Scan Completed CAD Ordered Case Completed

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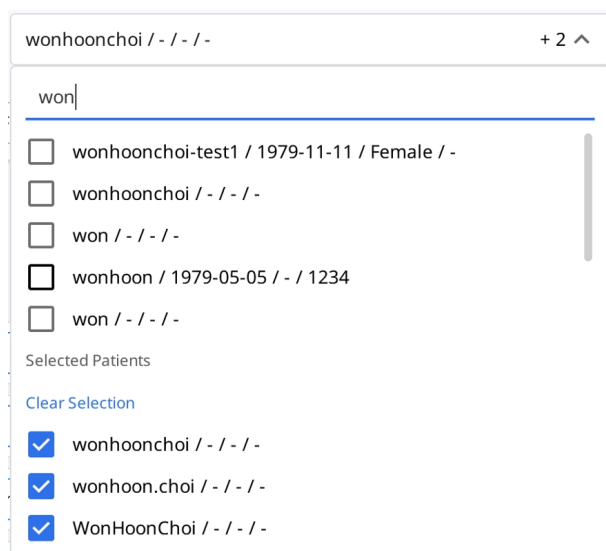
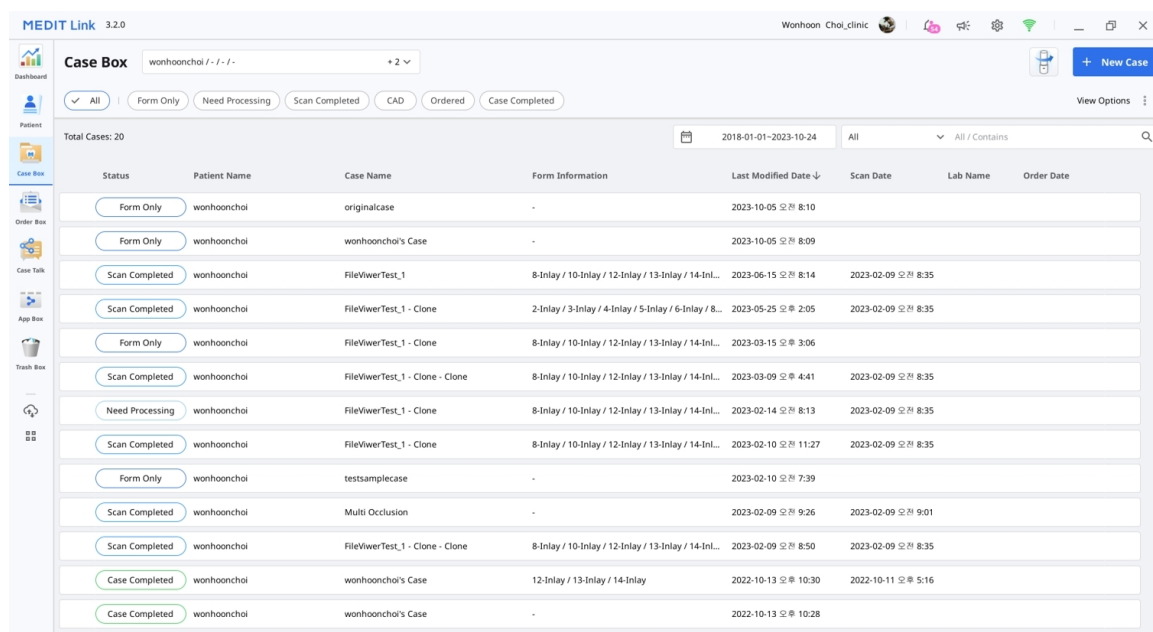
Status	Patient Name	Case Name	Form Information	Last Modified Date	Scan Date	Lab Name	Order Date
Form Only	wonhoonchoi	originalcase	-	2023-10-05 오전 8:10			
Form Only	wonhoonchoi	wonhoonchoi's Case	-	2023-10-05 오전 8:09			
Scan Completed	wonhoonchoi	FileViverTest_1	8-Inlay / 10-Inlay / 12-Inlay / 13-Inlay / 14-Inl...	2023-06-15 오전 8:14	2023-02-09 오전 8:35		
Scan Completed	wonhoonchoi	FileViverTest_1 - Clone	2-Inlay / 3-Inlay / 4-Inlay / 5-Inlay / 6-Inlay / 8...	2023-05-25 오후 2:05	2023-02-09 오전 8:35		
Form Only	wonhoonchoi	FileViverTest_1 - Clone	8-Inlay / 10-Inlay / 12-Inlay / 13-Inlay / 14-Inl...	2023-03-15 오후 3:06			
Scan Completed	wonhoonchoi	FileViverTest_1 - Clone - Clone	8-Inlay / 10-Inlay / 12-Inlay / 13-Inlay / 14-Inl...	2023-03-09 오후 4:41	2023-02-09 오전 8:35		
Need Processing	wonhoonchoi	FileViverTest_1 - Clone	8-Inlay / 10-Inlay / 12-Inlay / 13-Inlay / 14-Inl...	2023-02-14 오전 8:13	2023-02-09 오전 8:35		
Scan Completed	wonhoonchoi	FileViverTest_1 - Clone	8-Inlay / 10-Inlay / 12-Inlay / 13-Inlay / 14-Inl...	2023-02-10 오전 11:27	2023-02-09 오전 8:35		
Form Only	wonhoonchoi	testsamplcase	-	2023-02-10 오전 7:39			
Scan Completed	wonhoonchoi	Multi Occlusion	-	2023-02-09 오전 9:26	2023-02-09 오전 9:01		
Scan Completed	wonhoonchoi	FileViverTest_1 - Clone - Clone	8-Inlay / 10-Inlay / 12-Inlay / 13-Inlay / 14-Inl...	2023-02-09 오전 8:50	2023-02-09 오전 8:35		
Case Completed	wonhoonchoi	wonhoonchoi's Case	12-Inlay / 13-Inlay / 14-Inlay	2022-10-13 오후 10:30	2022-10-11 오후 5:16		
Case Completed	wonhoonchoi	wonhoonchoi's Case	-	2022-10-13 오후 10:28			

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Patient Name*

Case Name*

Creating patient information involves consent from the patient.
Patient information can be shared with 3rd parties while placing an order.

Partner

Select partner

Assignee

Select assignees

You can select up to 5 assignees.

Tags

You can add the entered tag with Enter or Spacebar.

You can add up to 10 tags of 100 characters or less.

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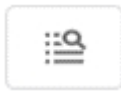
Patient List

Select a patient in the list to set the patient of the case. You can add or manage the patient information in Medit Link.

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qa_lab_01 Uses Lab Form

Full Services Dental Esthetics Crown & Bridge

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The screenshot shows the 'Form' tab of a dental software interface. The main workspace displays a dental arch diagram with teeth numbered 1 through 48. A red arrow points to the 'Form' tab in the top navigation bar. Another red arrow points to the 'Teeth' tab at the bottom of the arch diagram. On the right side, there is a sidebar with icons for Patient, Case Box, Order Box, Case Talk, App Box, and Trash Box. Below the arch diagram, there is an 'Order' section with a table of dental procedures and their details.

Type	Method	Material	Shade	Favorites
Inlay	Onlay	Coping	Denture Coping	
Crown	Pontic	Post & Core	Veneer	
Implant Crown	Custom Abutment	Surgical Guide	Offset Substructure	
Custom Abutment + Impla...	Custom Abutment + Coping	Diagnostic Wax-Up	Partial Denture	

Order List

No.	Type	Method	Material	Shade	
26	Crown	Full-Contour	Zirconia	A1	

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	Medit Scan for Clinics	Medit Scan for Labs
歯	0	0

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Share through Case Talk

Share with Anyone

Create a shared link and share it with anyone.

Expiration Date

2022-04-27

You can extend the period up to 90 days from today.

PIN Code

☐ PIN Code

Enter PIN Code

Set the PIN code using six digits for privacy protection.

☒ Share the patient's name

☒ Grant Download Permission [Premium Feature](#)

If you check this option, all users who have access to the shared link can download the files.

The premium feature is provided to all users for a limited time.
Feature availability is subject to change depending on the membership plan.

☒ Allow comments

If you check this option, you will be able to communicate through the shared link.

Comment

Please enter your comment.

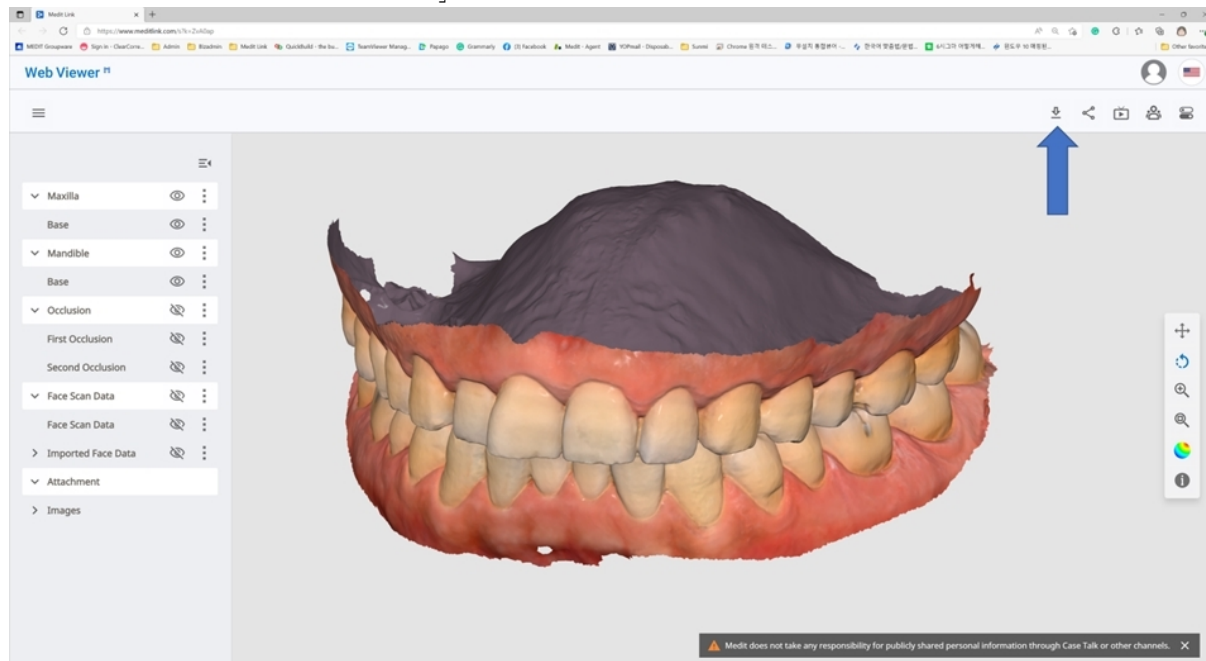
0/1024

By pressing the "Get Link" button, you agree that you understand personal information can be collected and provided to a third party (or shared in any public domain) ONLY with the consent of the person whose images on which you have worked. You also agree the consent of that person was truthfully obtained, and Medit is not responsible for any problems which may arise there after.

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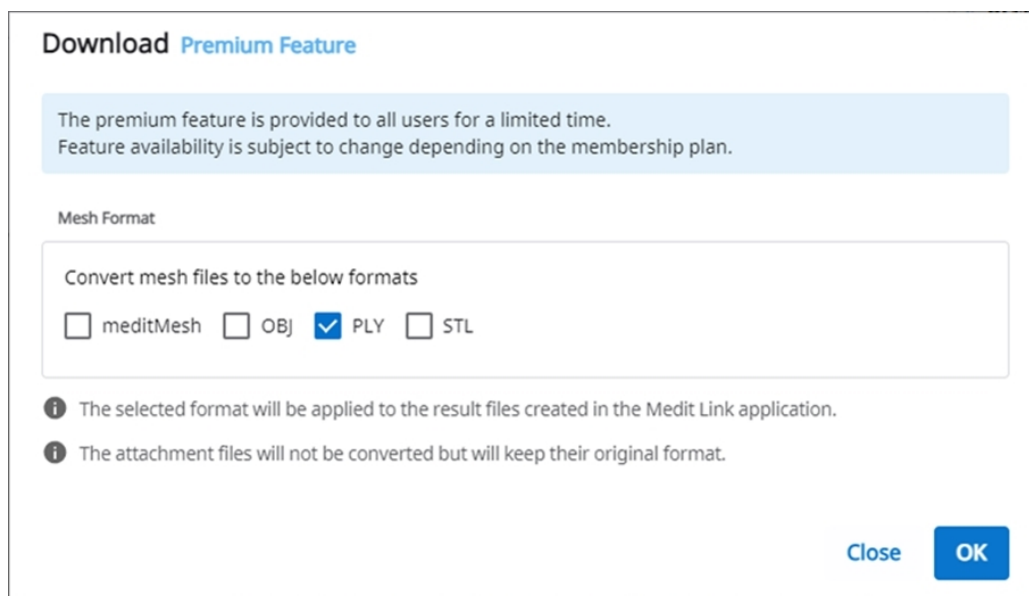
Uses Lab Form

San Diego, US

[92131] Residence Inn San Diego Rancho Bernardo

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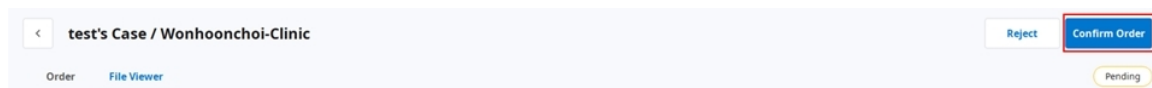
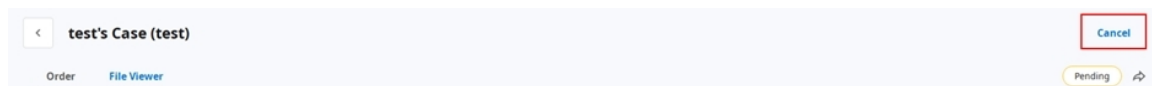
The screenshot shows a web interface for managing dental cases. The main header is 'test1's Case (test1)' with a 'Cancel Order' button. Below the header, there are two tabs: 'Order' and 'File Viewer'. The 'Order' tab is active, displaying a form with the following fields:

- Order ID: 6004858
- Assignee: (dropdown menu)
- Patient Name: test1
- Partner Name: Elysee Dental Test
- Scanned Date: (empty)
- Order Date: 2025-07-23 PM 5:17
- Requested Delivery Date: 2025-07-30
- Expected Delivery Date: 2025-07-30
- PM 05:00
- Memo: (empty text area)
- Tags: test

On the right side, there is a table titled 'Order List' and 'Shipment'.

No.	Type	Method	Material	Shad
15	Crown	Anatomic	Zirconia	A4

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Shipment Registration

Select a delivery method and enter the required information to register.
The information you enter is saved in the Shipment tab and can be checked at any time.

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Case Name
mandibular movement - Clone - Clone

Patient Name
Test Occlusion

☒ Share patient information

Requested Delivery Date
2023-04-13

Delivery Time
오전 8:15 AM 8:00 PM 8:00

Memo

Tags
You can add the entered tag with Enter or Spacebar.

☒ Include all occlusion relationships in the order
☒ Include Tags In The Order

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Tags
tag

Order List Shipment

No.	Type	Method	Material	Shad
25	Crown	Anatomic	Multi-Layered Zirconia	A3.5
26	Crown	Anatomic	Multi-Layered Zirconia	A3.5

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The screenshot shows a calendar for March 2022. The date 17 is selected. Below the calendar, there is a 'Delivery Time' section with a clock icon and a time input field showing '11:45 AM'. A slider below the input field ranges from 'AM 8:00' to 'PM 8:00'. Below the time section is a 'Comment' text area. At the bottom, there are 'Cancel' and 'Apply' buttons.

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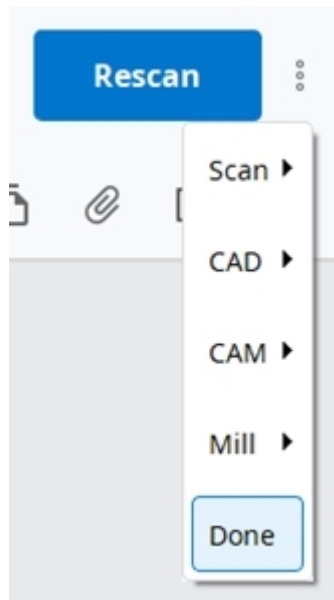
MEDIT Link 3.2.0

Work Box

Dashboard
In Box
Work Box
Case Talk
App Box
Trash Box

✓ All | Form Only | Need Processing | Scan Completed | CAD | CAM

Export Status Selected Cases: 4

Done

Status	Case Name	Patient Name	Form Information
✓ Form Only	tsstt's Case	tsstt	-
✓ Form Only	test's Case	test	-
✓ CAD	Maxilla Crown	Scan for Clinics	14-Crown
✓ CAD	Maxilla Crown	Scan for Clinics	14-Crown

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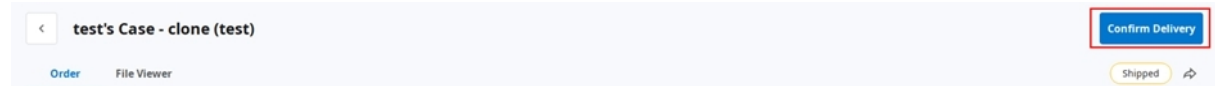
Status	Case Name	Patient Name
☐ Scan Completed	0218's Case	0205 test2mn(456)
2024-12-23		
☐ Ordered	ORG_[Clinics][3.3.6] 01 Normal - Complete	ORG_[Clinics][3.3.6] 01 Normal
☒ Scan Completed	ORG_[Clinics][3.3.6] 01 Normal - Complete	ORG_[Clinics][3.3.6] 01 Normal
☒ Scan Completed	ORG_[Clinics][3.3.6] 01 Normal - Complete	ORG_[Clinics][3.3.6] 01 Normal
☒ Scan Completed	ORG_[Clinics][3.3.6] 01 Normal - Complete	ORG_[Clinics][3.3.6] 01 Normal

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Patient Name*

Case Name*

 Creating patient information involves consent from the patient.
Patient information can be shared with 3rd parties while placing an order.

Tag

Doctorname × Specialmaterial × SpecialComment ×

 You can add up to 10 tags of 100 characters or less.

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Memo

Tags

Doctorname × Specialmaterial × SpecialComment ×

● You can add up to 10 tags of 100 characters or less.

☒ Include all occlusion relationships in the order

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Meditrial Srl
Meditrial Europe Ltd