MEDIT Link

ME-UG-702H
Revision 13 (29.09.2026)
SW version 3.5.0

Table of contents

Medit Link

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Administrator

Dentist, Manager

Have the authority to create a new dental clinic account and have full access to Medit Link.

- ✓ Scan/Order/Manage orders
- ✓ Manage members

Register

Member

Dental Hygienist, Dental Staff

Belong to the dental clinic their admin created and have limited access to Medit Link after the admin's approval.

- ✓ Scan/Order/Manage orders

Register

Lab

Administrator

Laboratory Director, Manager

Have the authority to create a new dental lab account and have full access to Medit Link.

- ✓ Scan/Manage ordered cases
- ✓ Manage members

Register

Member

Dental Technician, Laboratory Staff

Belong to the dental lab their admin created and have limited access to Medit Link after the admin's approval.

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Email Address

[Check](#)

 Please ensure your email address is active. Medit Link will send a verification email.

 Registering with a disposable and/or fake e-mail address can make you vulnerable to security issues and make it difficult to receive information from Medit Link.

Password

 Password must be 8 to 16 characters long, including a combination of at least three letters, numbers, and special characters. (Passwords are accepted in English language only.)

Confirm Password

 Please type in the same password as above.

Select Country/Region

☐ I agree to all of the Terms and Conditions.

☐ Personal Information Collection Agreement (Mandatory)

[View](#)

☐ Terms & Conditions (Mandatory)

[View](#)

☐ Agree to the Marketing and Promotion Agreement (optional)

[View](#)

 Check this option to receive marketing emails, which include product updates, software news, and promotional events.

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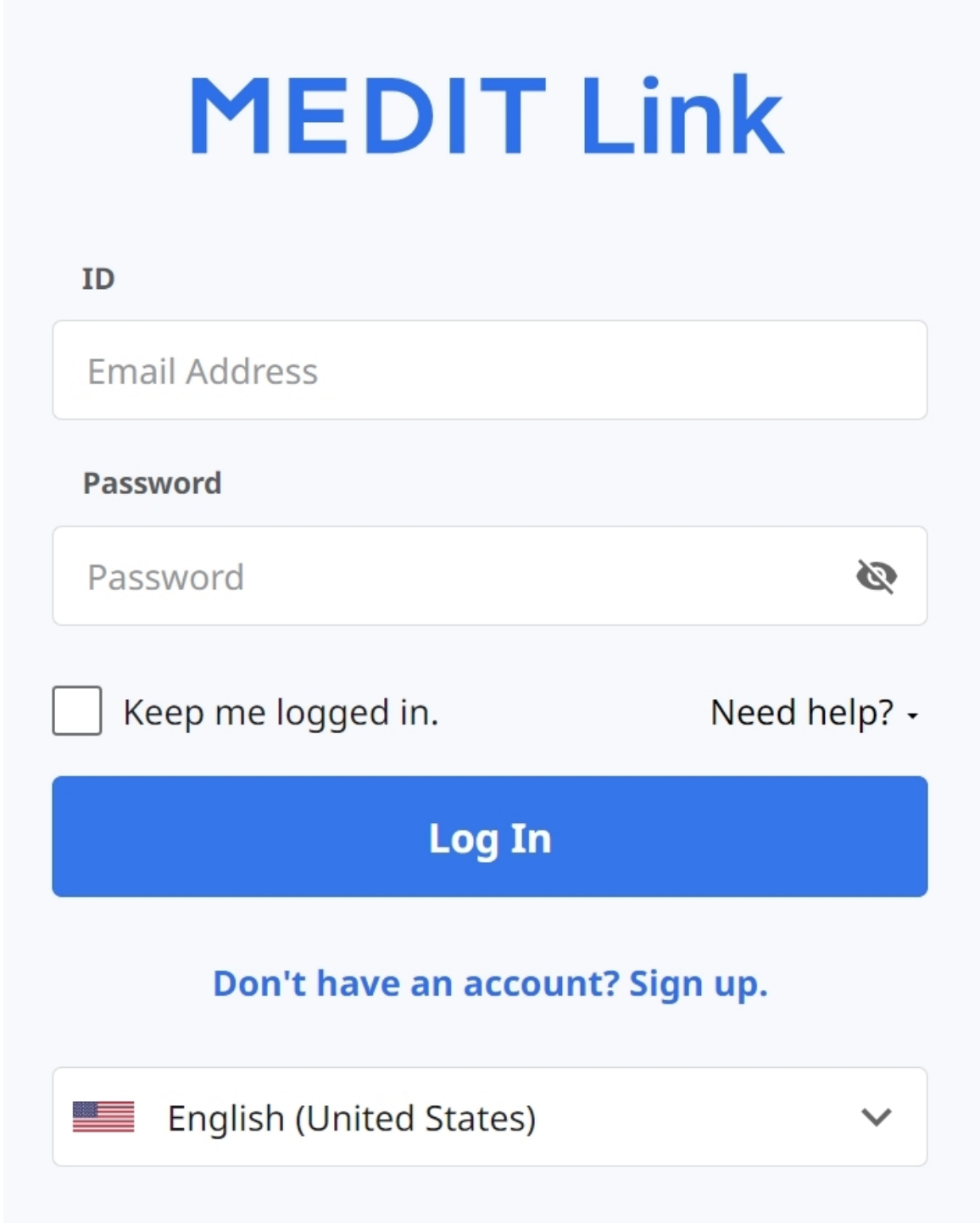
The image shows a login form for 'MEDIT Link'. At the top, the text 'MEDIT Link' is displayed in a large, bold, blue font. Below this, there are two input fields. The first is labeled 'ID' and contains the placeholder text 'Email Address'. The second is labeled 'Password' and contains the placeholder text 'Password', with a small eye icon to its right. Below the password field, there is a checkbox labeled 'Keep me logged in.' and a link 'Need help?' with a downward arrow. A large blue button labeled 'Log In' is positioned below these elements. At the bottom, there is a link 'Don't have an account? Sign up.' and a language selection dropdown menu showing 'English (United States)' with a flag icon and a downward arrow.

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
Welcome, stage-test-wh!

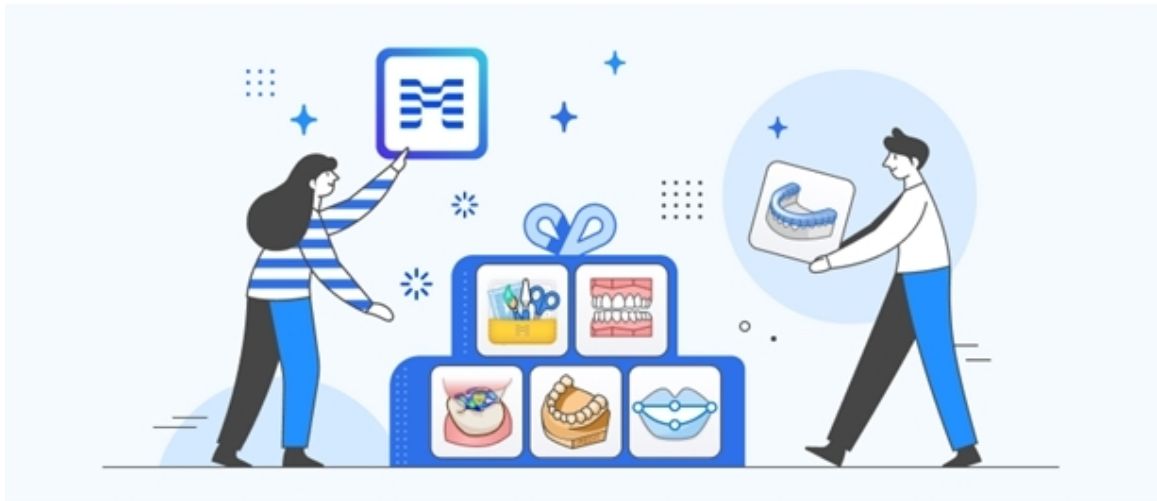Utilize Various Apps

Download Medit Link to your desktop to utilize scanning programs and various apps. All data generated is synced to the cloud and is also available on the web.

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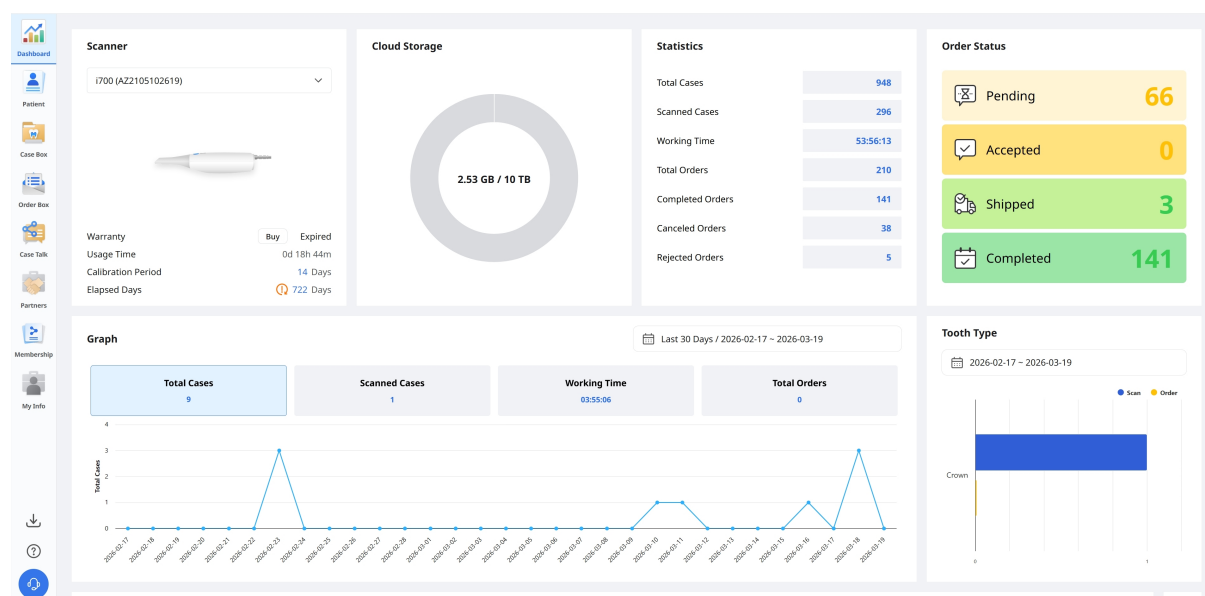
☒ Medit Link account

☐ Another email address


Enter another email address to receive the report email.

 The report will be sent on the 2nd day of every month.

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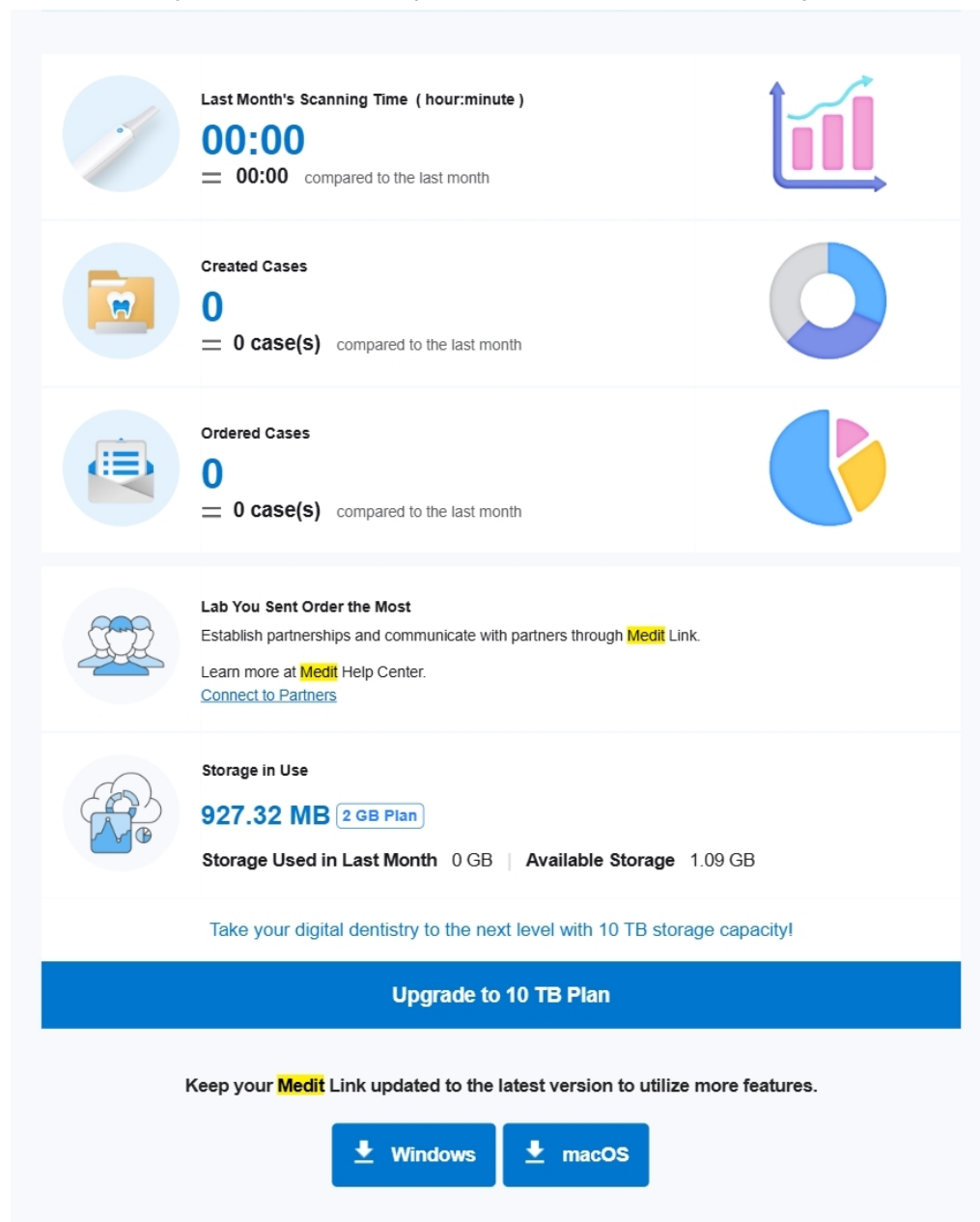This is an automatically generated email. Please do not reply to this email.
Please contact us at insights@medit.com for more information.
[Unsubscribe](#). If you no longer wish to receive this email, click "Unsubscribe." You can change the email subscription settings in [Medit Link Dashboard](#).
Copyright © Medit Corp. All rights reserved.

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
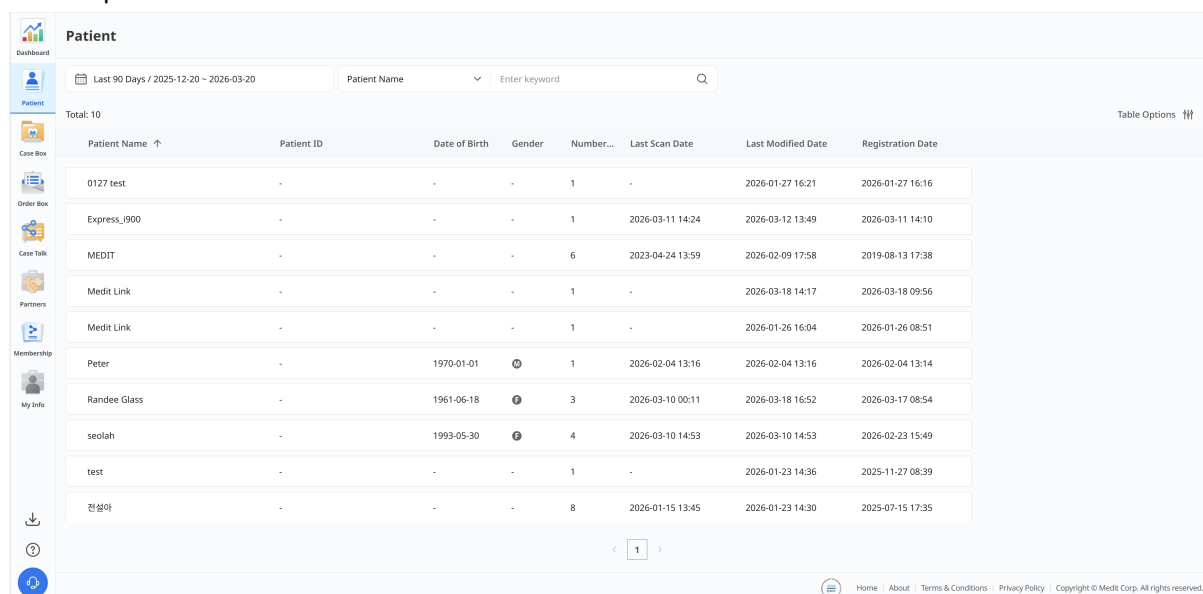
The screenshot shows a web application interface for managing patients. On the left is a sidebar with icons for Dashboard, Patient, Case Box, Order Box, Case Talk, Partners, Membership, and My Info. The main area is titled 'Patient' and includes a search bar with a date range 'Last 90 Days / 2025-12-20 ~ 2026-03-20' and a dropdown for 'Patient Name'. Below the search bar, it says 'Total: 10'. A table lists 10 patients with the following columns: Patient Name, Patient ID, Date of Birth, Gender, Number..., Last Scan Date, Last Modified Date, and Registration Date. The table is paginated, showing page 1 of 1.

Patient Name	Patient ID	Date of Birth	Gender	Number...	Last Scan Date	Last Modified Date	Registration Date
0127 test	-	-	-	1	-	2026-01-27 16:21	2026-01-27 16:16
Express_900	-	-	-	1	2026-03-11 14:24	2026-03-12 13:49	2026-03-11 14:10
MEDIT	-	-	-	6	2023-04-24 13:59	2026-02-09 17:58	2019-08-13 17:38
Medit Link	-	-	-	1	-	2026-03-18 14:17	2026-03-18 09:56
Medit Link	-	-	-	1	-	2026-01-26 16:04	2026-01-26 08:51
Peter	-	1970-01-01	♂	1	2026-02-04 13:16	2026-02-04 13:16	2026-02-04 13:14
Randee Glass	-	1961-06-18	♀	3	2026-03-10 00:11	2026-03-18 16:52	2026-03-17 08:54
seolah	-	1993-05-30	♀	4	2026-03-10 14:53	2026-03-10 14:53	2026-02-23 15:49
test	-	-	-	1	-	2026-01-23 14:36	2025-11-27 08:39
전설희	-	-	-	8	2026-01-15 13:45	2026-01-23 14:30	2025-07-15 17:35

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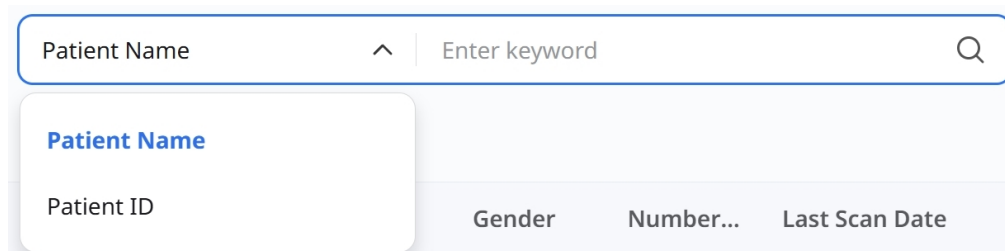The screenshot shows a search bar with the text 'Patient Name' and a dropdown menu. The dropdown menu has two options: 'Patient Name' (highlighted in blue) and 'Patient ID'. Below the dropdown menu, there are three columns: 'Gender', 'Number...', and 'Last Scan Date'.

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The screenshot shows the Medit Link Case Box interface. At the top, there's a header with the Medit Link logo, a user profile (seolahjeon), and a language dropdown (English). Below the header, the 'Case Box' section is active, showing a filter for 'Patient'. A date range 'Last 90 Days / 2025-12-20 ~ 2026-03-20' and a search bar 'Case Name' are visible. A filter bar includes 'All' (selected), 'Form Only', 'Need Processing', 'Scan Completed', 'CAD', 'Ordered', and 'Case Completed'. The table shows 18 total cases. The table columns are: Status, Case Name, Patient Name, Form Information, Last Modified..., Scanned Date, Order Date, Lab Name, and Assignee. The table lists various cases with their respective statuses and details.

Status	Case Name	Patient Name	Form Information	Last Modified...	Scanned Date	Order Date	Lab Name	Assignee
Scan Completed	Ortho Scan - Clone - Clo...	Randee Glass		2026-03-18 16:52	2026-03-10 00:11	-	-	Share
Scan Completed	Ortho Scan - Clone	Randee Glass		2026-03-18 14:19	2026-03-10 00:11	-	-	Share
Form Only	Medit Link의 케이스	Medit Link		2026-03-18 14:17	-	-	-	Share
Scan Completed	Ortho Scan	Randee Glass		2026-03-17 08:57	2026-03-10 00:11	-	-	Share
Need Processing	Express_i900의 케이스	Express_i900	23	2026-03-12 13:49	2026-03-11 14:24	-	-	Share
Form Only	Express Scan	seolah		2026-03-10 14:53	2026-03-10 14:53	-	-	Share
Form Only	Express Scan	seolah		2026-02-24 08:33	2026-02-24 08:33	-	-	Share
Need Processing	Express Scan	seolah		2026-02-23 16:38	2026-02-23 16:37	-	-	Share
Scan Completed	Express Scan	seolah	Maxilla	2026-02-23 15:50	2026-02-23 15:50	-	-	Share

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The screenshot shows the Medit Link Case Box interface with the 'Patient' filter dropdown menu open. The dropdown menu is highlighted with a red box. The filter bar includes 'All' (selected), 'Form Only', 'Need Processing', 'Scan Completed', 'CAD', 'Ordered', and 'Case Completed'.

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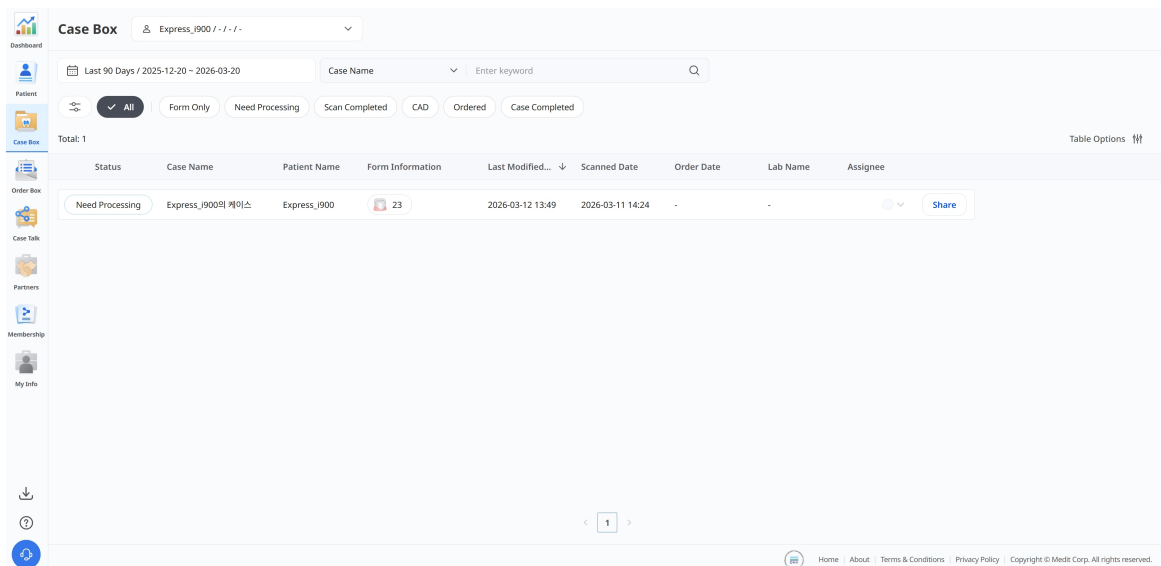
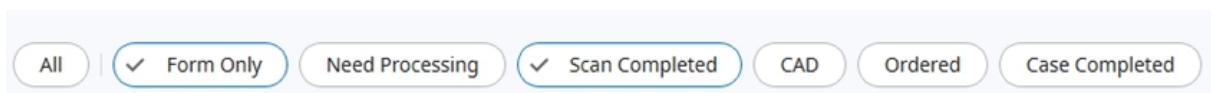
Case Box (Klinik)	Allt/Form Only/Need Processing/Scan Completed/CAD/Ordered/Case Completed
Work Box (Labb)	Allt/Form Only/Need Processing/Scan Completed/CAD/Klar

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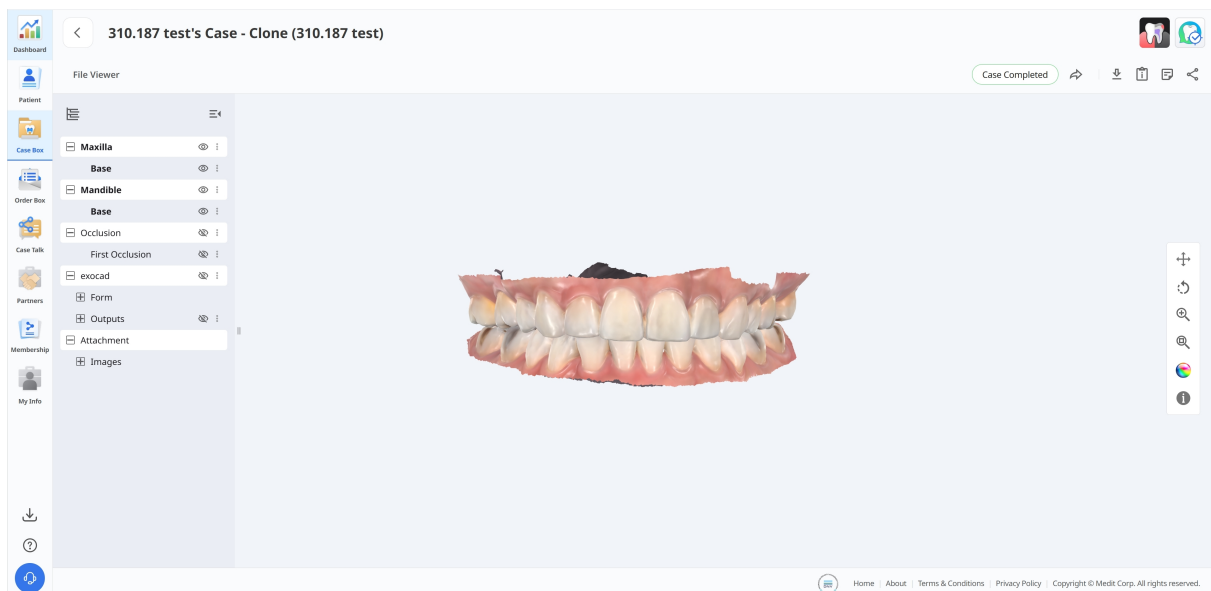
2

1

Home | About | Terms & Conditions | Privacy Policy | Copyright © Medit Corp. All rights reserved.

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
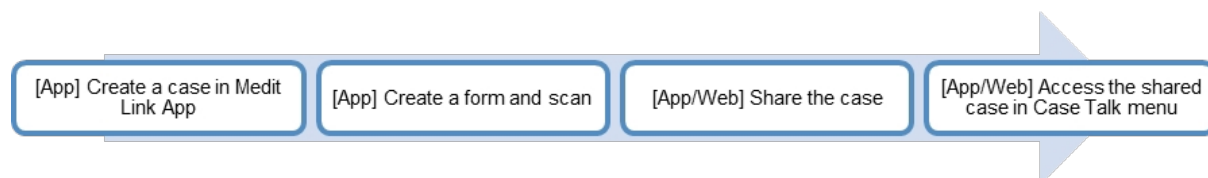
The screenshot shows the Medit Link Case Box interface. It features a sidebar with navigation options like Dashboard, Patient, Case Box, Order Box, Case Talk, Partners, Membership, and My Info. The main area displays a table of cases with columns for Status, Case Name, Patient Name, Form Information, Last Modified Date, Scanned Date, Order Date, and Lab Name. Each row has a 'Share' button on the right. The table lists several cases, including 'test's Case' and 'wonhoonchoi-test1's Case - Clone - Clone'.

Status	Case Name	Patient Name	Form Information	Last Modified Date	Scanned Date	Order Date	Lab Name	Share
Form Only	test's Case	test		2023-11-03 12:50	-	-	-	Share
Scan Completed	wonhoonchoi-test1's Case - Clone - Clone	wonhoonchoi-test1	23, 24, 25, 26	2023-11-02 19:37	2023-11-02 09:01	-	-	Share
Scan Completed	wonhoonchoi-test1's Case - Clone	wonhoonchoi-test1	23, 24, 25, 26	2023-11-02 19:21	2023-11-02 09:01	-	-	Share
Ordered	wonhoonchoi-test1's Case	wonhoonchoi-test1	23, 24, 25, 26	2023-11-02 19:21	2023-11-02 09:01	2023-11-02 ...	Wonhoonchoi_lab2	Share
Scan Completed	G11598492311021004395	test test (12345)	21	2023-11-02 19:09	-	-	-	Share
Scan Completed	test's Case	test		2023-11-01 13:32	-	-	-	Share
Scan Completed	G11598492311010319123	test test (test)	21	2023-11-01 12:20	-	-	-	Share
Scan Completed	test's Case	test		2023-10-26 17:25	2023-10-26 17:24	-	-	Share
Scan Completed	wonhoonchoi-test1's Case	wonhoonchoi-test1		2023-10-26 17:22	2023-10-26 17:21	-	-	Share
Ordered	test's Case	test	14, 15, 16, 17	2023-10-24 15:45	-	2023-10-24 ...	Wonhoonchoi_lab2	Share

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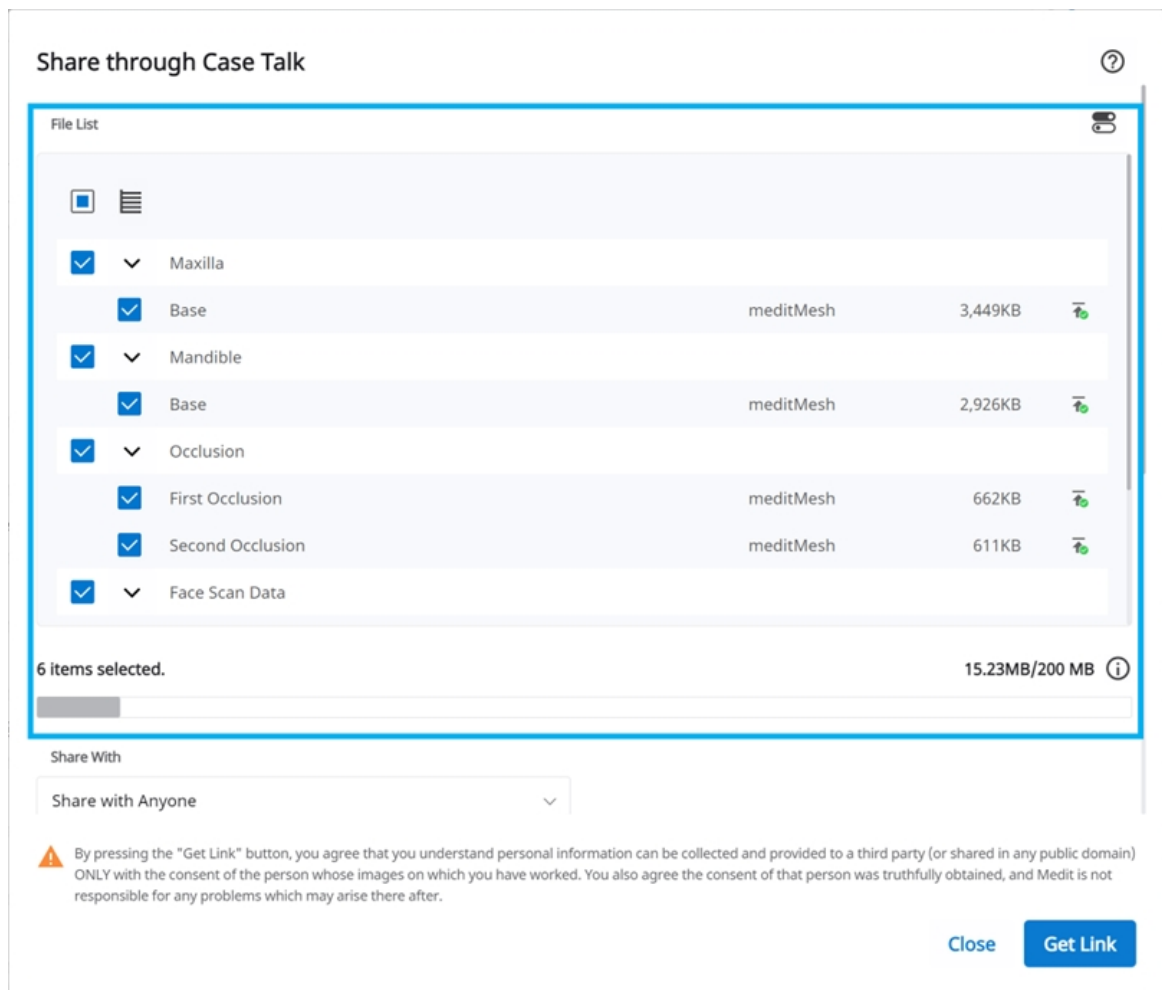Share With

Share with specific people

Create a shared link and share it with specific people.

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☒ Share the patient's name

☒ Grant Download Permission Premium Feature

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The premium feature is provided to all users for a limited time. Feature availability is subject to change depending on the membership plan.

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Invite Partner

Invite your desired partner to join Medit Link. When the lab accepts your invitation, a partnership is automatically formed.

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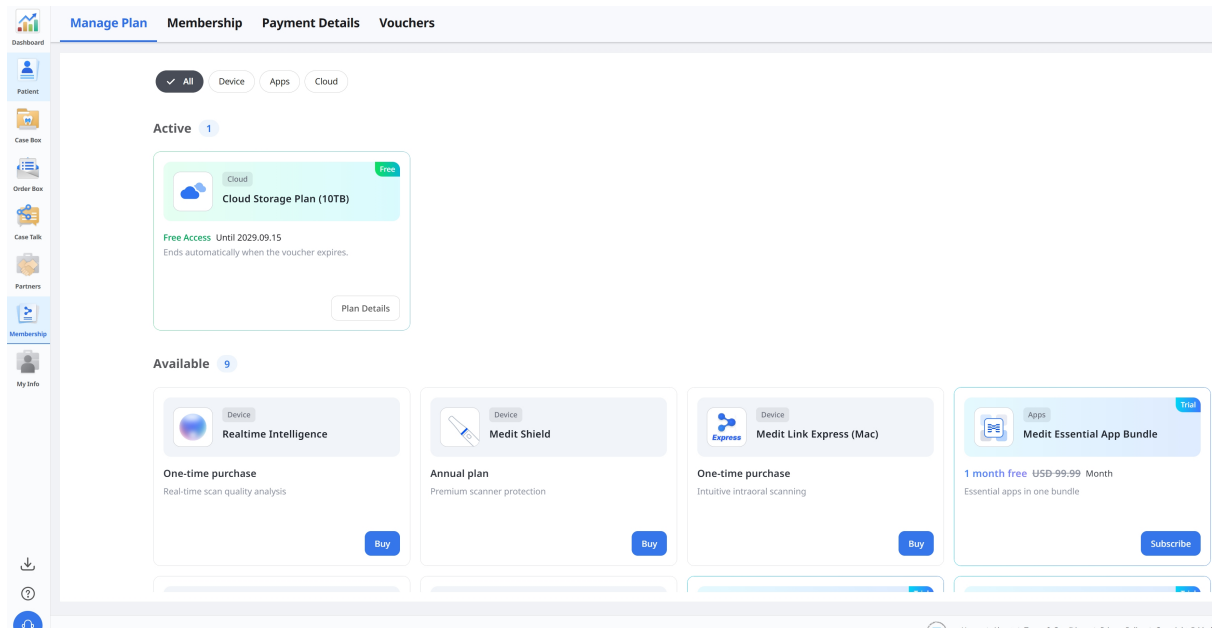
Free Access

Your free access is active.

Product	Cloud Storage Plan (10TB)
Status	Free
Expiration Date	2029.09.15
Regular Price	USD 9.99/Month

⚠ Important

- You can't subscribe to this product while the voucher is active.
- Free access ends automatically when the voucher expires. A paid subscription is required to continue.

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The screenshot displays a user interface for managing subscriptions. At the top, there are four tabs: 'Manage Plan' (selected), 'Membership', 'Payment Details', and 'Vouchers'. Below the tabs, there are three filter buttons: 'All' (selected), 'Device', and 'Apps'. The 'Active' section shows '0' active subscriptions with a message: 'No active subscriptions. Browse our services below to get started.' The 'Available' section shows '8' available services. A red arrow points from the 'Available' count to a specific service card: 'Cloud Storage Plan (10TB)'. This card is labeled 'Cloud' and 'Trial'. It offers '1 month free' (with the original price 'USD 9.99' crossed out) and 'Secure 10 TB cloud storage'. A blue 'Subscribe' button is at the bottom right of the card.

Manage Plan Membership Payment Details Vouchers

✓ All Device Apps Cloud

Active 0

No active subscriptions.
Browse our services below to get started.

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1 month free

~~USD 9.99~~ Month

- ✓ Includes 10 TB of cloud storage
- ✓ Securely store large cases and other data
- ✓ Streamline your Medit Link data management in the cloud
- ✓ Store design files without separate local backups

- Actual charges may vary depending on your bank or card issuer's fees and policies.
- For more information about Cloud Storage Plan (10TB), see [here](#).

Trial Cloud Storage: First Month Free

1 month free >

☐ I agree to automatic renewal and recurring charges.

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←

Subscribe to Medit Cloud Storage 10TB Plan

\$0.00
Then \$9.99 per day after coupon expires

Medit Cloud Storage 10TB Plan for Clinics Medit Cloud storage plan Billed daily	\$9.99
Subtotal	\$9.99
1 Month Free for New Subscribers 100% off for a day	-\$9.99
Tax	\$0.00
Total due today	\$0.00

Pay with card

Email

yhgam0829@gmail.com

Card information

1234 1234 1234 1234

MM / YY

CVC

Cardholder name

Full name on card

Country or region

South Korea

Subscribe

By confirming your subscription, you allow Medit USA Inc to charge you for future payments in accordance with their terms. You can always cancel your subscription.

Powered by | [Terms](#) | [Privacy](#)

Subscribe

×

Your Cloud Storage Plan (10TB) subscription is now active.
Enjoy full access to all features.

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
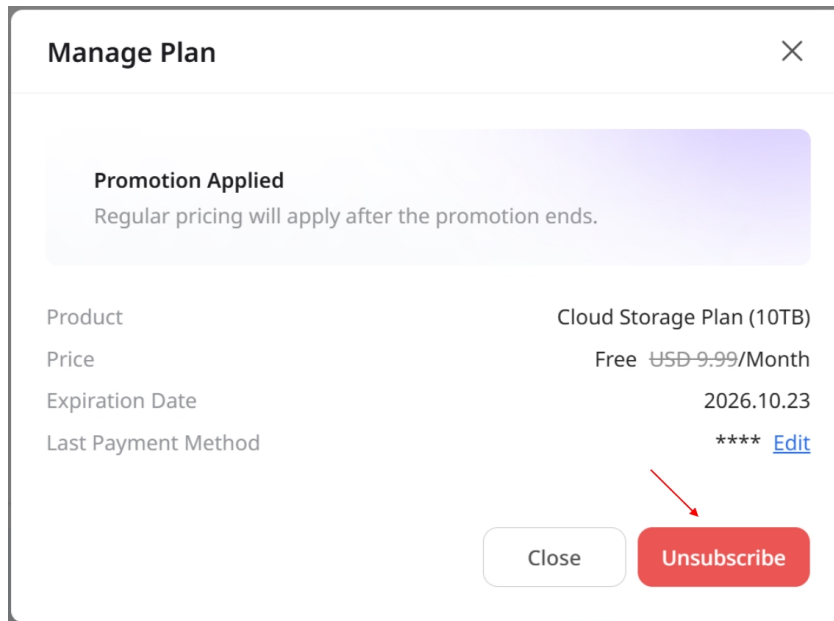
Manage Plan [Close]

Promotion Applied
Regular pricing will apply after the promotion ends.

Product	Cloud Storage Plan (10TB)
Price	Free USD 9.99 /Month
Expiration Date	2026.10.23
Last Payment Method	**** Edit

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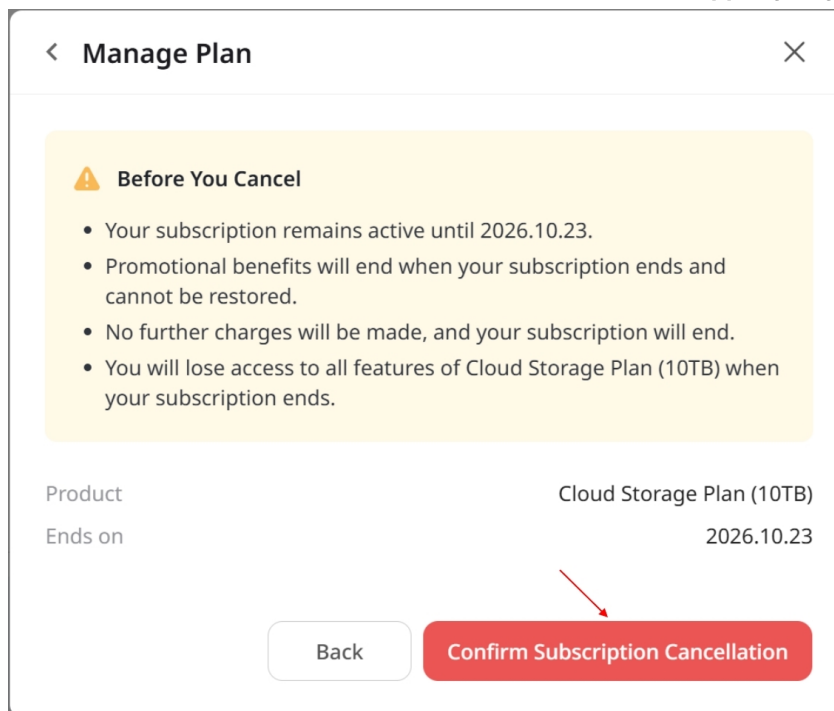Manage Plan [Close]


⚠ Before You Cancel

- Your subscription remains active until 2026.10.23.
- Promotional benefits will end when your subscription ends and cannot be restored.
- No further charges will be made, and your subscription will end.
- You will lose access to all features of Cloud Storage Plan (10TB) when your subscription ends.

Product	Cloud Storage Plan (10TB)
Ends on	2026.10.23

[Back] [Confirm Subscription Cancellation]

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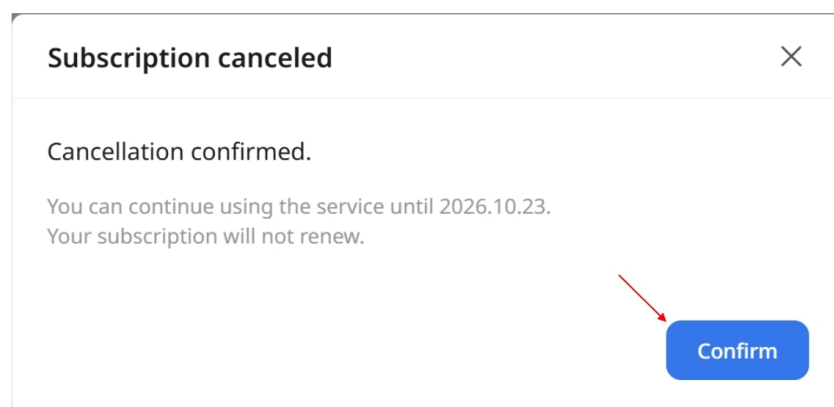Subscription canceled [Close]

Cancellation confirmed.

You can continue using the service until 2026.10.23.
Your subscription will not renew.

[Confirm]

Membership

The screenshot displays the 'Membership' section of a web application. The top navigation bar includes 'Manage Plan', 'Membership' (active), 'Payment Details', and 'Vouchers'. The left sidebar contains icons for Dashboard, Patient, Case Box, Order Box, Case Table, Partners, Membership (active), and My Info. The main content area is titled 'Medit Cloud Storage Plan' and shows the current plan as '2 GB Plan' for 'Free'. A promotional banner offers to 'Unlock 10 TB today for just USD 9.99/month!' with a 'Subscribe' button. Below this, sections for '2 GB Plan - Standard Benefits' and '10 TB Plan - Exclusive Benefits' both state 'No coupons available.' The right sidebar, 'Storage Overview & Cleanup', shows a date range 'All / 2018-01-01 - 2026-09-22', 'My Cloud Storage' at '0.00 KB / 2.00 GB' with an 'Upgrade Storage' link, and a 'Result Files & Attachments' section at '0.00 KB' with a 'Free Up Space' button. The footer contains links for Home, About, Terms & Conditions, Privacy Policy, and Copyright © Medit Corp. All rights reserved.

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2 GB Plan
Free

Try 10 TB free for the first month!

Only \$9.99/month afterward, billed automatically.
Cancel anytime via your subscription page.

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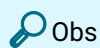
Apps

2026-06-24 - 2026-09-22

Payment Date	Product	Payment Type	Duration	Details	Amount of Payment	Status	
2026-09-21	Cloud Storage Plan (10TB)	Subscription	2026-09-21 - 2026-10-21	Monthly payment	USD 0.00	Payment Completed	
2026-09-21	Cloud Storage Plan (10TB)	Subscription	2026-09-21 - 2026-10-21	Monthly payment	USD 0.00	Payment Completed	
2026-09-21	Cloud Storage Plan (10TB)	Subscription	2026-09-21 - 2026-10-21	Monthly payment	USD 0.00	Payment Completed	
2026-09-21	Cloud Storage Plan (10TB)	Subscription	2026-09-21 - 2026-10-21	Monthly payment	USD 0.00	Payment Completed	
2026-09-21	Cloud Storage Plan (10TB)	Subscription	2026-09-21 - 2026-10-21	Monthly payment	USD 0.00	Payment Completed	

< 1 >

Home

About

Terms & Conditions


Privacy Policy

Copyright © Medix Corp. All rights reserved.

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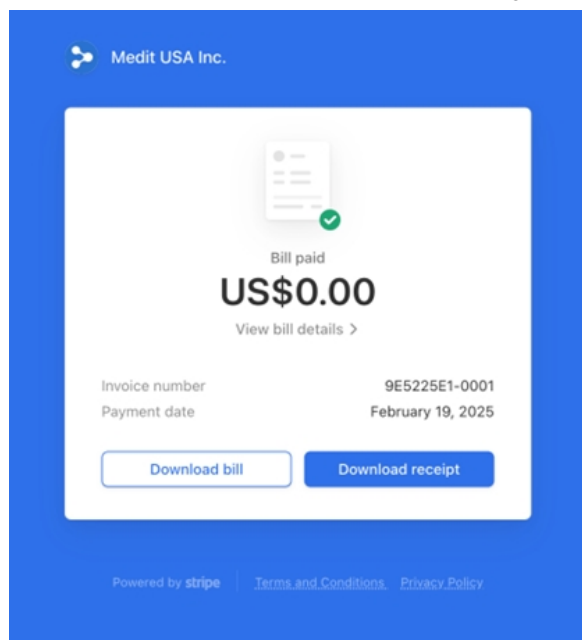ITEM

FEBRUARY 19 - FEBRUARY 20, 2025	
Medit Cloud Storage 10TB Plan for Clinics	US\$0.60
Quantity 1 piece	
subtotal	US\$0.60
1 Month Free for New Subscribers(100% 할인)	-US\$0.60
Payment amount	US\$0.00
Amount paid	US\$0.00
Remaining amount	US\$0.00

Do you have any questions? [Contact Medit USA Inc.](#)

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Manage Plan Membership Payment Details **Vouchers**

Register Code

Enter voucher code Register

Vouchers 13

All

Product	Promotion	Code Number	Registration Date	Valid Until	From	To	Status
Cloud Storage Plan	Basic EJ_Basic_Cloud Storage Plan (10TB)	2QJWR0FP	2026-09-21	2026-12-10	2026-10-21	2026-11-21	Used
Cloud Storage Plan	Basic EJ_Basic_Cloud Storage Plan (10TB)	QQ3QXJ3T	2026-09-21	2026-12-10	-	-	Registered
Cloud Storage Plan	Basic EJ_Basic_Cloud Storage Plan (10TB)	RSP4GZOI	2026-09-21	2026-12-10	-	-	Registered
Cloud Storage Plan	Basic EJ_Basic_Cloud Storage Plan (10TB)	JYAYH1XH	2026-09-21	2026-12-10	-	-	Registered
Cloud Storage Plan	Basic EJ_Basic_Cloud Storage Plan (10TB)	RSP4GZOI	2026-09-21	2026-12-10	-	-	Registered
Cloud Storage Plan	Redeem EJ_Reedeem_Cloud Storage Plan (10TB)	ZGYUDFRY	2026-09-21	2100-01-01	-	-	Registered
Cloud Storage Plan	Redeem EJ_Reedeem_Cloud Storage Plan (10TB)	1MGYSVHF	2026-09-21	2100-01-01	2026-09-21	2026-10-21	Used
Cloud Storage Plan	Redeem EJ_Reedeem_Expired1	3FWHQA	2026-09-21	2026-09-21	-	-	Expired
Cloud Storage Plan	Redeem EJ_Reedeem_Expired1	COUR8QKR	2026-09-21	2026-09-21	-	-	Expired
Cloud Storage Plan	Redeem EJ_Reedeem_Terminated1	PZPJX8	2026-09-21	2026-12-21	2026-09-21	2026-10-21	Used
Cloud Storage Plan	Redeem SQA_Reedeem_Cloud Storage Plan_36M	GFZKBWJ5	2026-09-15	2100-01-01	2026-09-15	2029-09-15	Activated
Cloud Storage Plan	Redeem SQA_Reedeem_Cloud Storage Plan_36M	VE00ZB46	2026-09-14	2100-01-01	2026-09-14	2029-09-14	Used

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Vouchers 2

All

Product	Promotion	Code Number
Cloud Storage Plan	Redeem SQA_Reedeem_Cloud Storage Plan_36M	GFZ!

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Voucher CodeZG'
ProductCloud Storage Plan
Free Access30 Days

⚠ Important

- Once activated, the voucher cannot be canceled.
- You can't subscribe to this product while the voucher is active.

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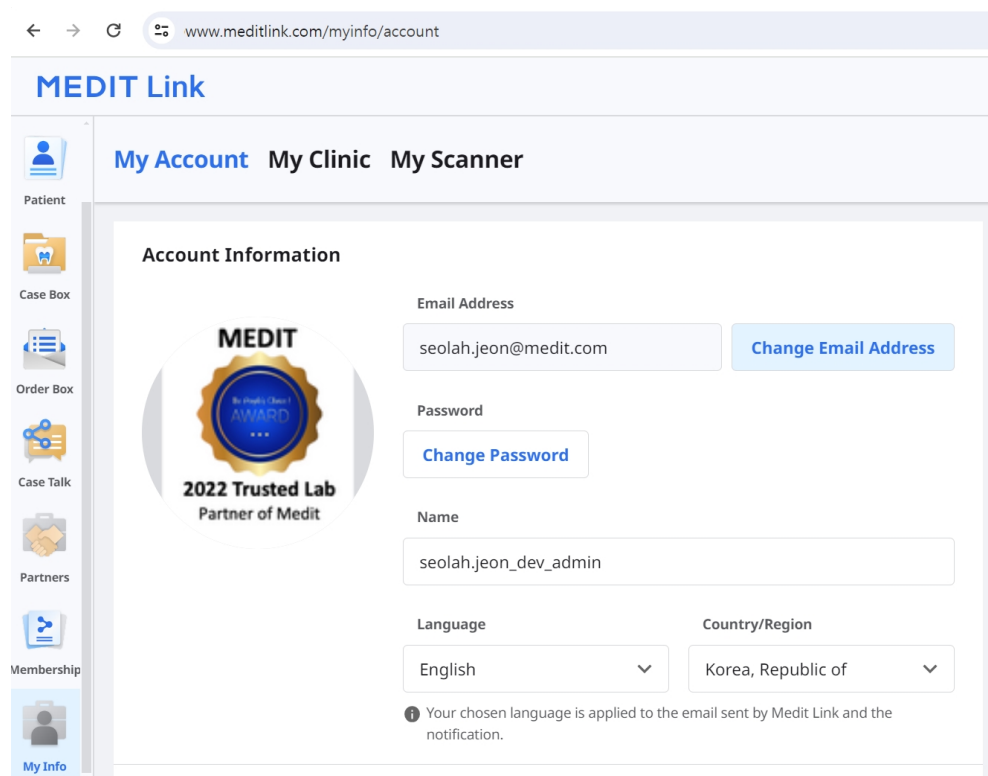
The screenshot shows a web browser window with the URL www.meditlink.com/myinfo/account. The page is titled "MEDIT Link" and has a navigation bar with "My Account", "My Clinic", and "My Scanner". A left sidebar contains icons for "Patient", "Case Box", "Order Box", "Case Talk", "Partners", "Membership", and "My Info". The main content area is titled "Account Information" and displays the following details:

- Email Address:** seolah.jeon@medit.com (with a "Change Email Address" button)
- Password:** (with a "Change Password" button)
- Name:** seolah.jeon_dev_admin
- Language:** English (dropdown menu)
- Country/Region:** Korea, Republic of (dropdown menu)

Below the dropdowns, a note states: "Your chosen language is applied to the email sent by MEDIT Link and the notification." On the left, there is a circular badge that reads "MEDIT 2022 Trusted Lab Partner of MEDIT".

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Edit

Cover and Logo

Looking for digitalized dental laboratories?

Find the most ordered Labs around your clinic with 'Featured Part

Services

Select

Specialties (choose up to 3)

Select

About Us

e.g., Welcome to OOO Dental Clinic.
We are a licensed dental clinic situated in California specializing mainly in aesthetic restorations and implantology.

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The scanner information is automatically registered once your scanner is connected to a Medit scan program.

Home

About

Terms & Conditions

Privacy Policy

Copyright © Medit Corp. All rights reserved.

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The image shows two instances of the 'Manage Options' dialog box. The left instance has the 'Teeth' tab selected, and the right instance has the 'Arch' tab selected. Both dialogs show a list of dental procedures with status toggles.

Teeth	Arch
1 Inlay	1 Denture
2 Onlay	2 Partial Denture
3 Coping	3 Custom Tray
4 Denture Coping	4 Orthodontic Appliance
5 Crown	5 Bite Splint
6 Pontic	6 Clear Aligner
7 Post & Core	7 Night Guard
8 Veneer	8 Tooth Positioner
9 Implant Crown	9 Study Model

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The image shows the 'Manage Options' dialog box with a confirmation modal open. The modal has a title 'Edit' and a message: 'The updated order form will be available to clinics immediately.' There are 'Cancel' and 'Save' buttons at the bottom of the modal. A red arrow points from the 'Save' button in the modal to the 'Save' button in the top right corner of the 'Manage Options' dialog.

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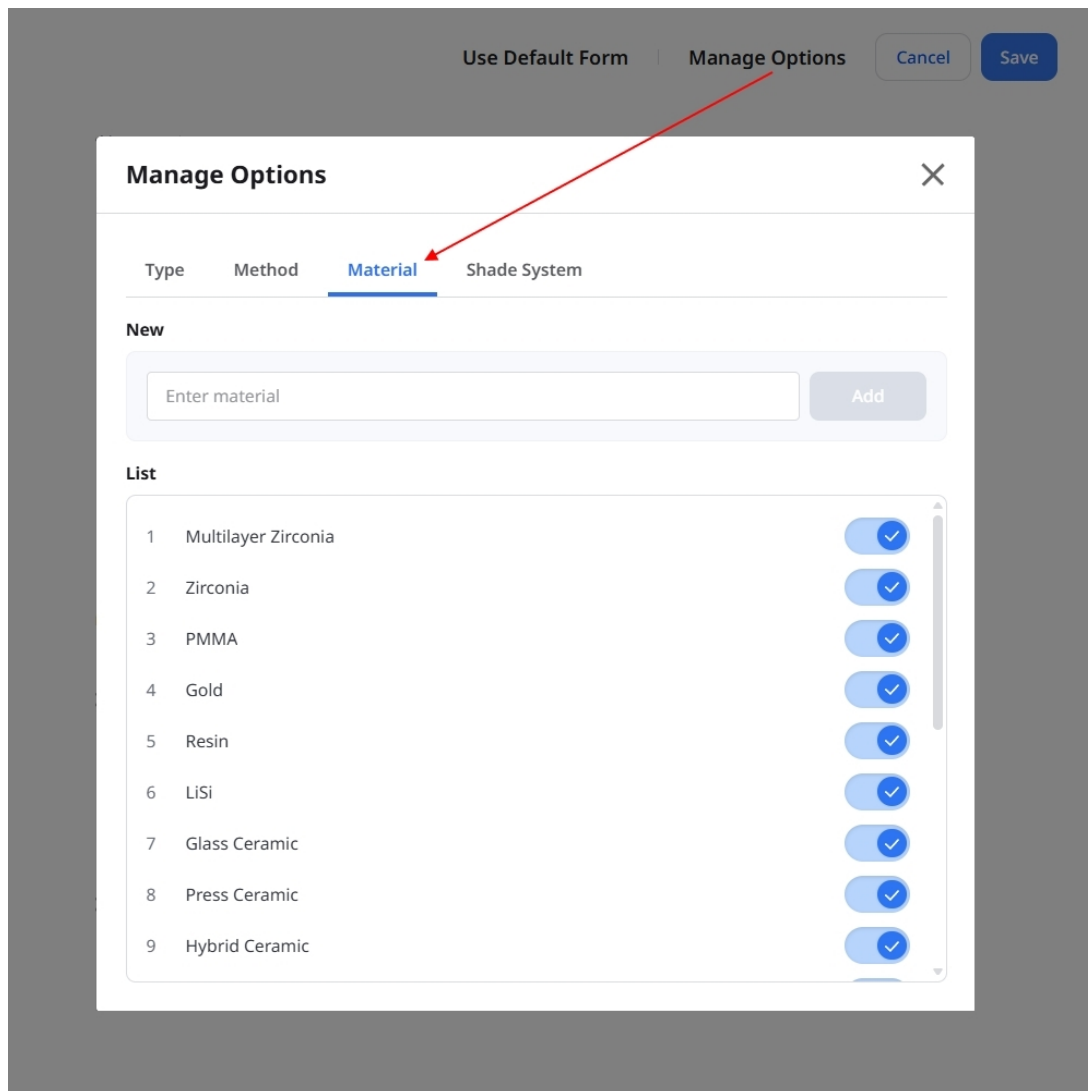
Zirconia HT

List

12	Titanium	☒
13	Metal	☒
14	Stainless Steel	☒
15	Co-Cr (Cobalt Chromium)	☒
16	Ni-Cr (Nickel-Chromium)	☒
17	Wax	☒
18	PEEK	☒
19	3D Print	☒
20	Zirconia HT edit delete	☒

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
The screenshot shows a software interface for configuring a crown. The title is "Crown". Under "Allowed Teeth", there are three checked options: Anterior, Premolar, and Molar. Under "Methods", there is a search bar and a grid of buttons. The "Full-Contour" button is highlighted in dark grey. Below the methods, there is a "Do not use" checkbox which is unchecked. A dropdown menu is open for "Full-Contour 1", showing a list of materials. In this list, "Zirconia HT" is highlighted with a blue plus sign, and a red arrow points to it from the "Co-Cr (Cobalt Chromium)" button. Other materials in the list include Multilayer Zirconia, Zirconia, PMMA, Gold, Resin, LiSi, Glass Ceramic, Press Ceramic, Hybrid Ceramic, Composite Resin, e.max, Titanium, Metal, Stainless Steel, Ni-Cr (Nickel-Chromium), Wax, PEEK, and 3D Print. There is also a "+ Add Material" link and a "Hide" button at the bottom of the materials list.

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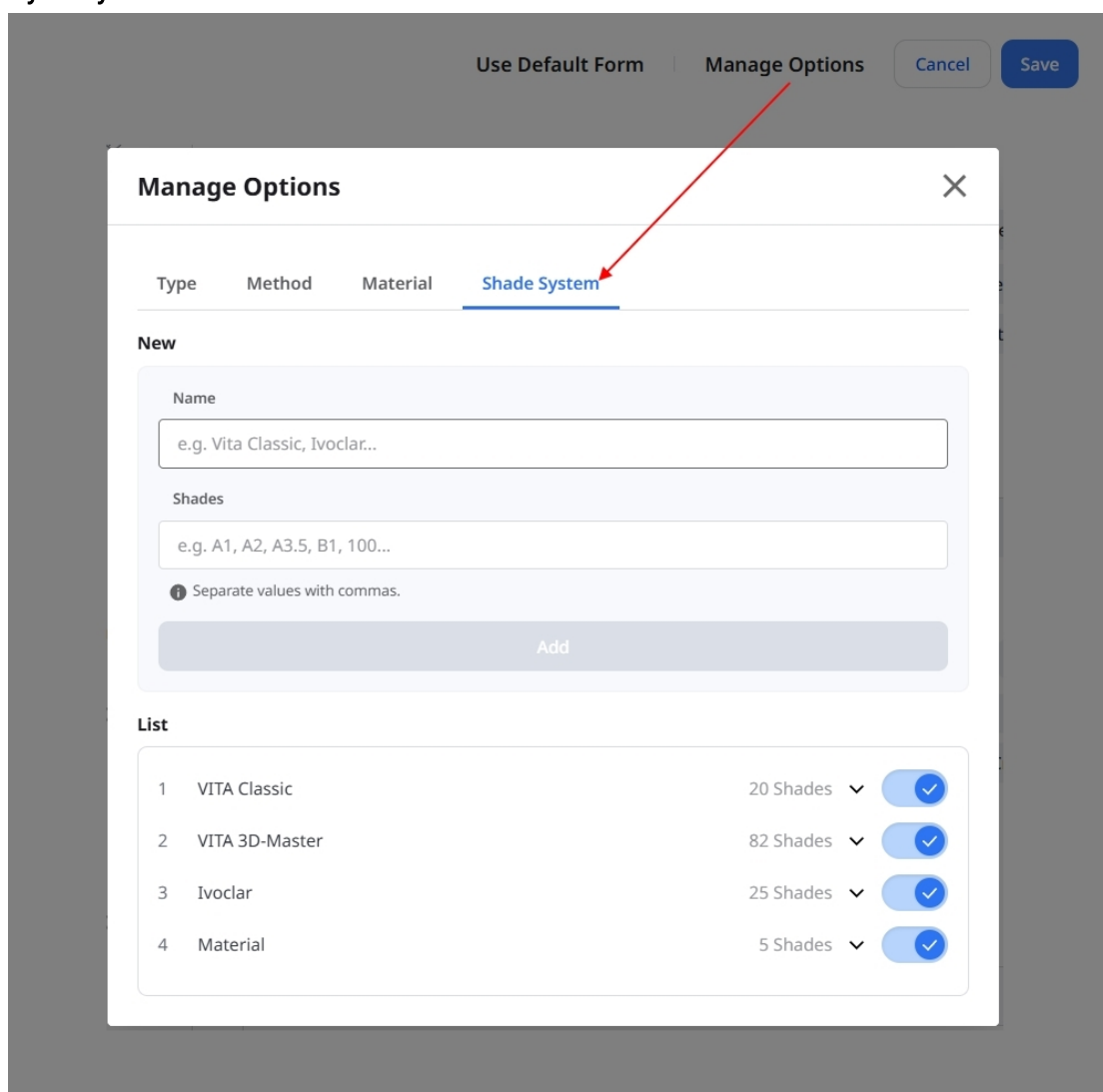The screenshot shows the 'Manage Options' dialog box. At the top, there are tabs: 'Use Default Form', 'Manage Options' (selected), 'Cancel', and 'Save'. Below the tabs, the 'Shade System' tab is active. The dialog is divided into two sections: 'New' and 'List'.

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
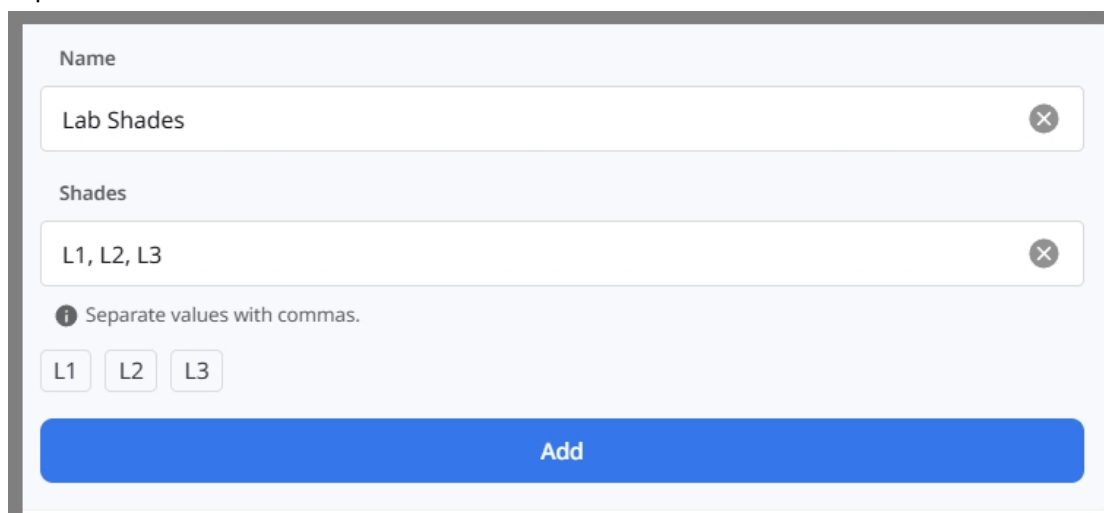
The screenshot shows the 'New' section of the 'Manage Options' dialog box. The 'Name' field contains 'Lab Shades'. The 'Shades' field contains 'L1, L2, L3'. Below the 'Shades' field, there is a hint: *i* Separate values with commas. Below the hint, there are three buttons: 'L1', 'L2', and 'L3'. At the bottom, there is a blue 'Add' button.

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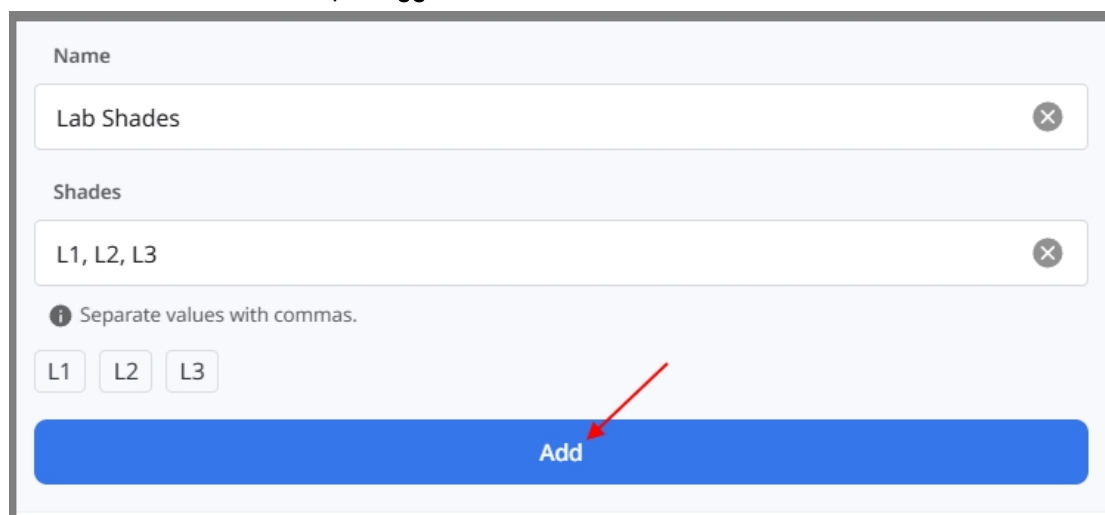The form is titled "Name" and "Shades". The "Name" field contains "Lab Shades". The "Shades" field contains "L1, L2, L3". Below the "Shades" field, there is a hint: "i Separate values with commas." and three buttons labeled "L1", "L2", and "L3". At the bottom of the form is a blue button labeled "Add", which is pointed to by a red arrow.

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The screenshot shows a web interface for configuring a crown. The main section is titled 'Crown'. Under 'Allowed Teeth', 'Anterior', 'Premolar', and 'Molar' are all checked. Under 'Methods', 'Full-Contour' is selected among several options. Below this, a search bar and a list of materials are shown. 'Zirconia HT' is selected, and a red arrow points to the '+' icon next to it. Under 'Shade Systems', 'Lab Shades' is selected, and a red arrow points to the '+' icon next to it. The 'Lab Shades' system shows three options: L1, L2, and L3.

Crown

Allowed Teeth

☒ Anterior ☒ Premolar ☒ Molar

Methods

Search

Screw-Retained **Full-Contour** Offset Cement-Retained

SCRP Temporary Cut-back Pontic Press Preform Screw-Retained Layered

Screw-Retained PFZ Screw-Retained PFG Screw-Retained PFM Cement-Retained Layered

▼ More

☐ Do not use

▼ **Full-Contour** 1

Materials

Search

Multilayer Zirconia Zirconia PMMA Gold Resin

LiSi Glass Ceramic Press Ceramic Hybrid Ceramic Composite Resin e.max

Titanium Metal Stainless Steel Co-Cr (Cobalt Chromium) Ni-Cr (Nickel-Chromium) Wax

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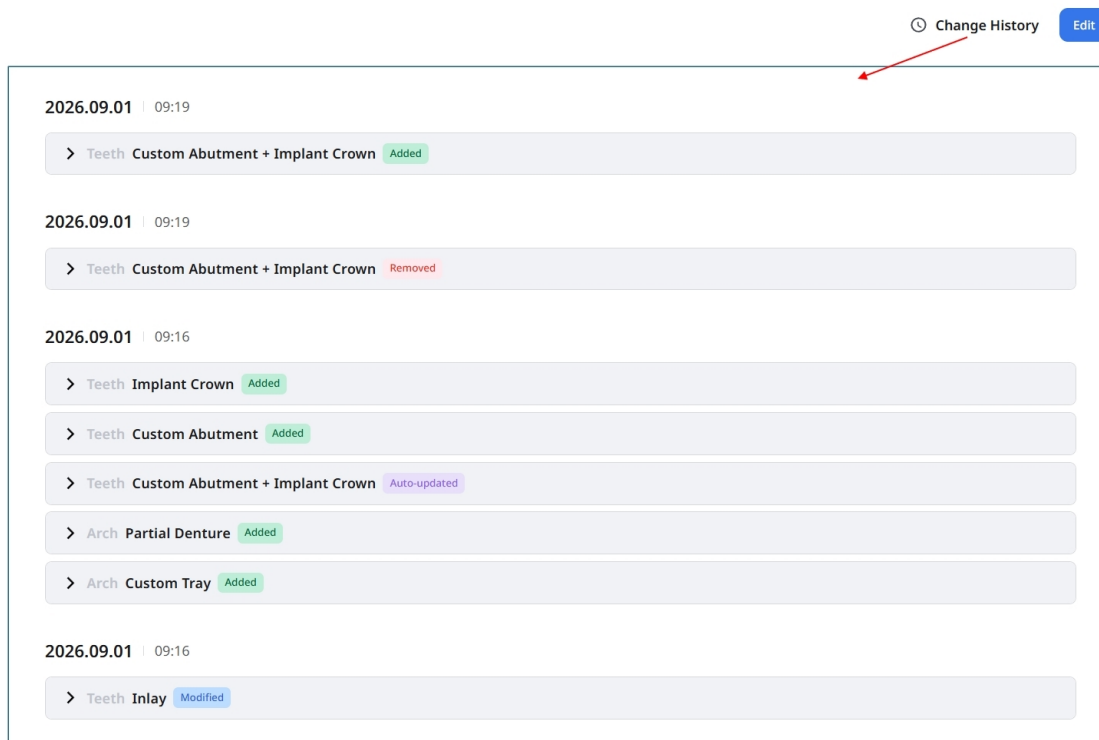
2026.09.01 | 09:19

- > Teeth Custom Abutment + Implant Crown Added

2026.09.01 | 09:19

- > Teeth Custom Abutment + Implant Crown Removed


2026.09.01 | 09:16

- > Teeth Implant Crown Added
- > Teeth Custom Abutment Added
- > Teeth Custom Abutment + Implant Crown Auto-updated
- > Arch Partial Denture Added
- > Arch Custom Tray Added

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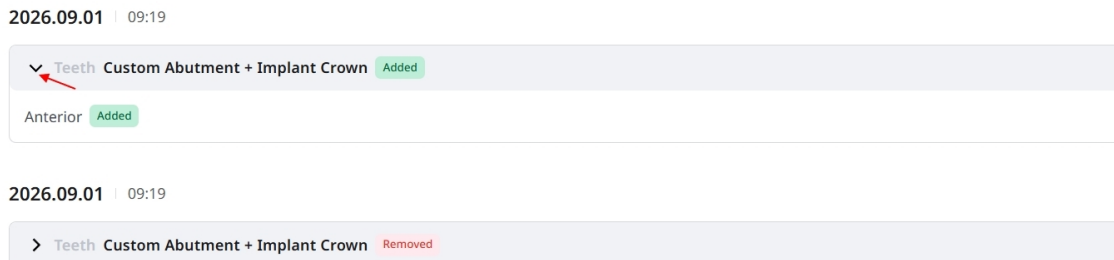2026.09.01 | 09:19

- ▼ Teeth Custom Abutment + Implant Crown Added
 - Anterior Added

2026.09.01 | 09:19

- > Teeth Custom Abutment + Implant Crown Removed

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Notice
NFT Service Termination Notice
2023-04-19

To our users,

We want to thank all Medit Link users who have supported and used our NFT service.

Unfortunately, we must inform you that this service will be terminated following the upcoming program update.

Please see the details below for more information:

- **Feature in Question**
 - Share via NFT (in Medit Link Case Talk)
- **Effective Date**
 - Starting from the update in the first half of 2023 (scheduled for April 27)

※ The update schedule is subject to change.

We'll continuously work to provide better service.

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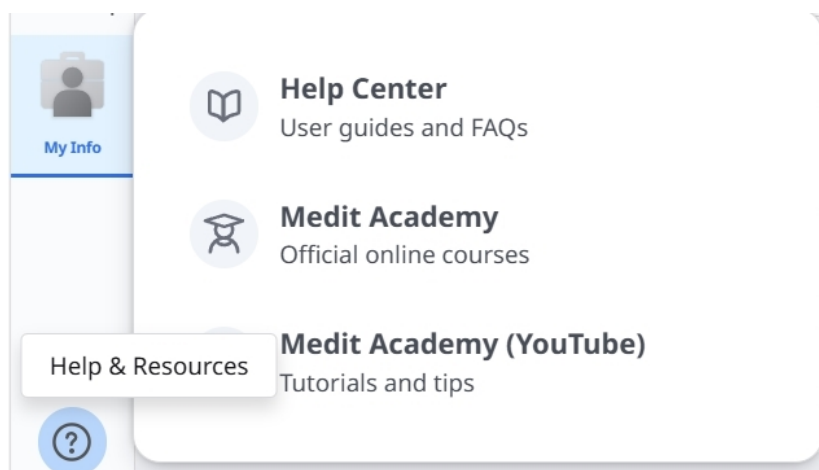
Översikt

The screenshot shows the Medit Link app interface. At the top, there's a header with 'MEDIT Link Beta 3.4.10' and a user profile. Below this is a sidebar with icons for Dashboard, Patient, Case Box, Order Box, Case Talk, and Trash Box. The main area displays a table of patient records. Red numbers 1, 2, 3, and 4 are overlaid on the interface to indicate specific areas of interest.

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Register Offline Data

Register

Finds and registers offline data on this device that is not registered in Medit Link.

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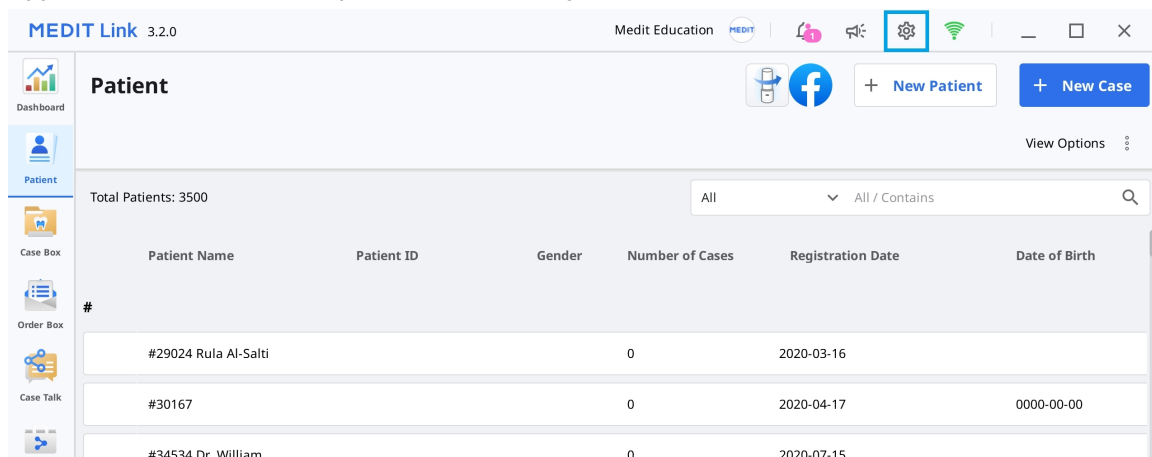
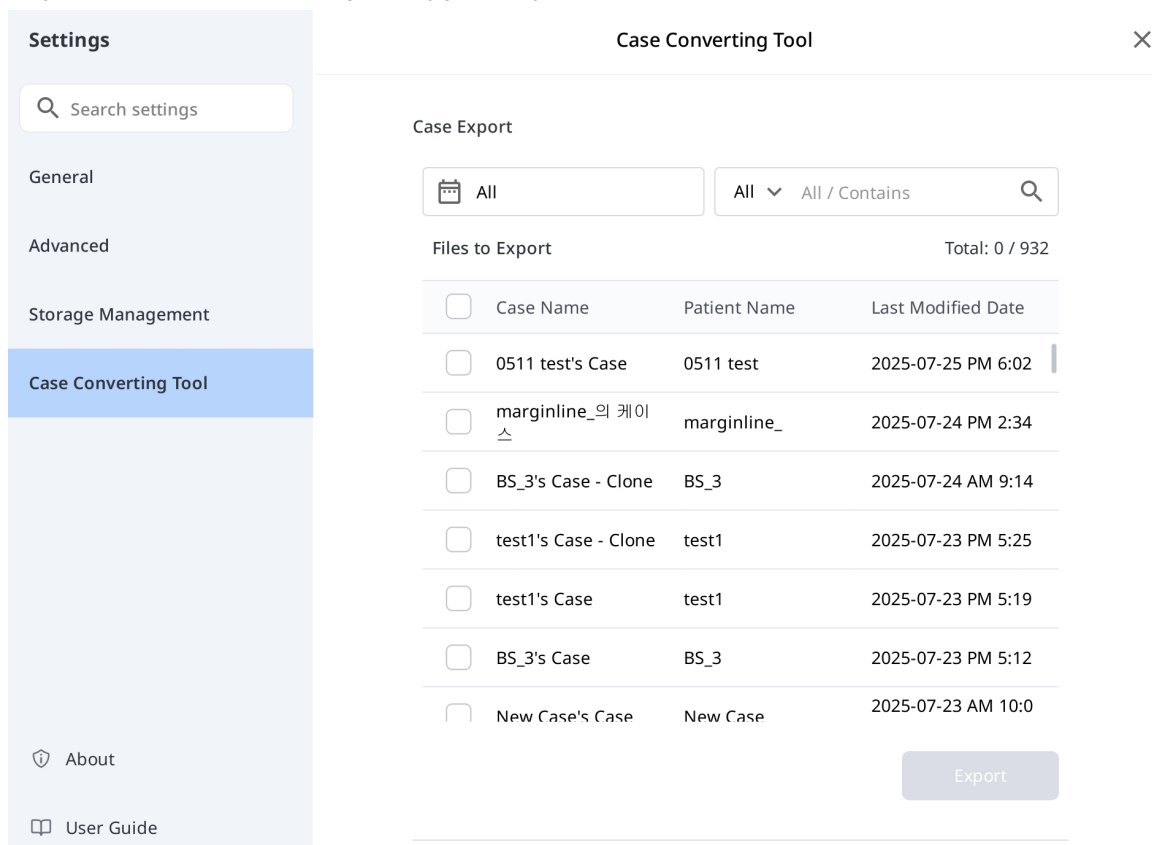
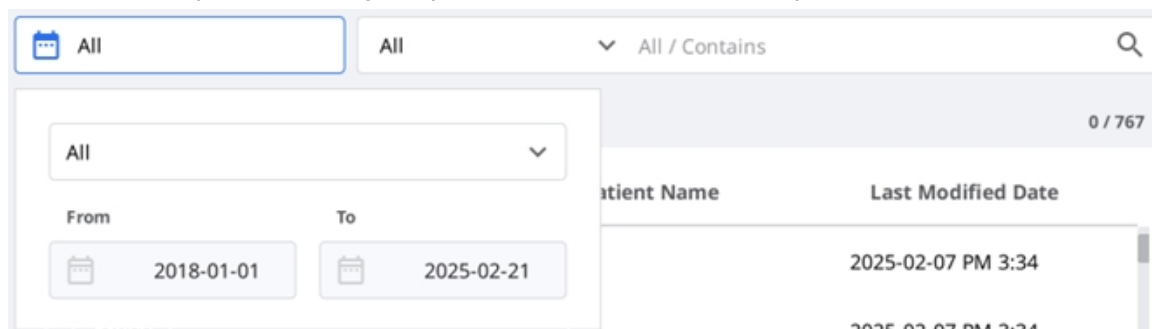
Question

Exporting successfully completed.
Do you want to open the folder containing the exported file?

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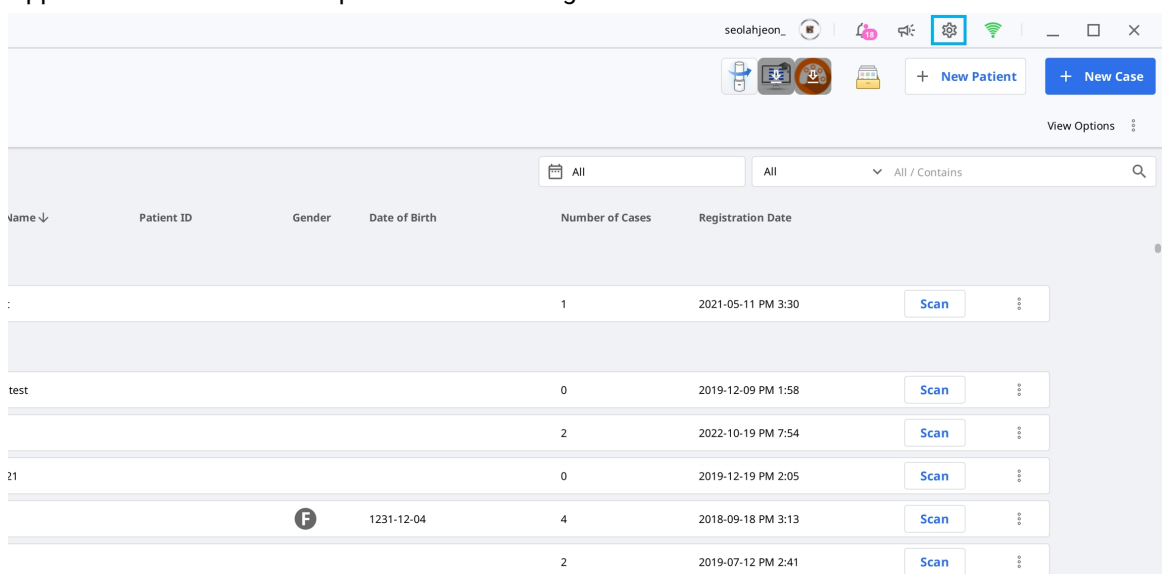
Case Import

A new case is created for every imported case. It contains case data, patient information, and scan file. However, it does not contain order information.

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 Meditrial Srl
Via Po 9 00198, Rome Italy

 Meditrial Europe Ltd
Bahnhofstrasse 23 6300 Zug, Switzerland

 Meditrial UK Ltd
19 The Circle, Queen Elizabeth Street, London, SE1 2JE, United Kingdom

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Merge all patient information to

Please select the patient you'd like to merge all data to.
After merging all data together, all information in other patient files will be deleted.

	Patient Name	Patient ID	Gender	Birth Date	Number of Cases	Date Registered	Recent Visit Date
☒	5-9_iScan_CAMOrd...				1	2020-09-18	2020-09-18
☐	test				2	2020-09-17	2020-09-17
☐	test				22	2019-07-26	2020-09-17

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Filter: All / 2018-01-01 ~ 2026-04-08 | All / Contains | All / Form Only | Form Only | Need Processing | Scan Completed | CAD | Ordered | Case Completed

Total: 5

Status	Case Name	Patient Name	Form Information	Last Modified Date	Scanned Date	Lab Name
☐ Form Only	1231456's Case	1231456	25-Pontic	2026-04-08 PM 4:23		
☒ Scan Completed	159's Case	1231456	16-Custom Abutment	2025-02-13 PM 4:17		
☒ Scan Completed	1231456 - clone - clone	1231456	21-Crown	2022-11-25 AM 11:27		
☒ Scan Completed	1231456 - clone	1231456	21-Crown	2022-11-25 AM 11:22		
☐ Form Only	1231456's Case	1231456	24-Crown / 25-Crown / 26-Crown / 27-Crown	2022-10-14 AM 1:19		

2 selected | Select All | Export | Status | Delete

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Order Box/In Box

Order Box

Dashboard | Patient | Case Box | Order Box | Case Talk | Trash Box

Filter: All / 2018-01-01 ~ 2026-04-08 | All / Contains

Status: All | Pending | Canceled | Rejected | Accepted | Shipped | Case Completed

Total: 254

Status	Case Name	Patient Name	Order Date	Requested Delivery Date	Lab Name	Order ID
Accepted	texture on's Case - Clone	texture on	2026-04-08 PM 5:40		MEDIT_LAB_STAGE	6012449
Pending	Express Scan	Seol register	2026-04-02 PM 9:38	2026-04-09 AM 9:00	MEDIT_LAB_STAGE	6012242
Case Completed	전설이의 케이스	전설이	2025-07-16 PM 5:46	2025-07-17 PM 3:00	MEDIT_LAB_STAGE	6004719
Case Completed	27	1	2025-03-16 PM 3:11	2025-03-17 PM 3:00	MEDIT_LAB_STAGE	6002656
Case Completed	3.3.0.136 TEST's Case	3.3.0.136 TEST	2024-04-12 PM 3:25	2024-04-13 PM 2:00	MEDIT_LAB_STAGE	5998987
Case Completed	test의 케이스	Patient TEST	2024-04-09 PM 3:06	2024-04-10 PM 2:00	MEDIT_LAB_STAGE	5998955
Case Completed	Mandibular Movement Case - Clone - Clone	MEDIT	2023-06-12 PM 12:12	2023-06-13 PM 2:00	MEDIT_LAB_STAGE	5996680
Case Completed	Release New_2's Case	Release New_2	2023-04-25 PM 6:47	2023-04-25 PM 2:00	MEDIT_LAB_STAGE	5996180
Case Completed	cX 문의's Case - Clone	cX 문의	2023-05-16 PM 5:26	2023-05-17 PM 2:00	MEDIT_LAB_STAGE	5996393
Case Completed	cX 문의's Case	cX 문의	2023-05-15 PM 2:39	2023-05-16 PM 2:00	MEDIT_LAB_STAGE	5996362
Case Completed	310.187's Case	310.187	2023-05-15 PM 1:33	2023-05-16 PM 2:00	MEDIT_LAB_STAGE	5996360
Case Completed	310.187 test's Case - Clone			2023-05-03 PM 2:00	MEDIT_LAB_STAGE	5996306

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Expected Delivery Date: 2025-07-22 PM 12:00

Order List:

No.	Type	Method	Material	Shade
23	Crown	Anatomic	Multi-Layered Zi...	None

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Status	Case Name	Patient Name	Lab Name	Deletion Date ↓	Deleted By	Size
☐ Scan Completed	jira_임플란트 46, 47's Case - Clone - Clone	임플란트 46, 47		2023-04-26 오전 9:12	wonhoon choi	66,654 KB
☐ Form Only	G11598492304180810429	test test		2023-04-18 오후 5:15	wonhoon choi	0 KB
☐ Form Only	G11598492304180806313	test test		2023-04-18 오후 5:15	wonhoon choi	0 KB
☐ Form Only	G11598492304180811349	tes test		2023-04-18 오후 5:15	wonhoon choi	0 KB
☐ Form Only	G11598492304180754083	Gildewell_fastdesign		2023-04-18 오후 5:00	wonhoon choi	0 KB
☐ Form Only	G11598492304180733332	Gildewell_fastdesign		2023-04-18 오후 5:00	wonhoon choi	0 KB
☐ Form Only	G11598492304180648146	test test		2023-04-18 오후 5:00	wonhoon choi	0 KB
☐ Case Completed	Multi Occlusion Normal Case's Case - Clone	Multi Occlusion Normal Case		2023-04-18 오전 10:50	wonhoon choi	23,652 KB
☐ CAD	Multi Occlusion Normal Case's Case - Clone - Clone	Multi Occlusion Normal Case		2023-04-18 오전 10:50	wonhoon choi	23,652 KB
☐ Form Only	Occlusion Test's Case	Occlusion Test		2023-04-17 오후 6:31	wonhoon choi	0 KB
☐ Form Only	Test Occlusion 's Case	Test Occlusion		2023-04-17 오후 6:31	wonhoon choi	0 KB
☐ Form Only	Occlusion Test's Case	Occlusion Test		2023-04-17 오후 6:31	wonhoon choi	0 KB
☐ Form Only	Occlusion Test's Case	Occlusion Test		2023-04-17 오후 6:26	wonhoon choi	0 KB

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Total Cases: 210

☒ No Date
 Data Type

☐ Case Name	Patient Name	Scanned Date ↓	Size	Status
☐ BLATRIX Charline's Case	BLATRIX Charline	2025-01-06 PM 6:28	729.30 MB	Upload Required
☐ BLATRIX Charline's Case	BLATRIX Charline	2025-01-06 PM 6:28	729.30 MB	Upload Required
☐ BLATRIX Charline's Case - 복제	BLATRIX Charline	2025-01-06 PM 6:28	733.35 MB	Upload Required
☐ BLATRIX Charline's Case - Clone	BLATRIX Charline	2025-01-06 PM 6:28	734.01 MB	Upload Required
☐ test의 케이스	Patient TEST	2023-04-26 PM 7:58	1.25 GB	Upload Required
☐ i700 HD's Case	i700 HD	2023-04-24 PM 5:23	911.24 MB	Upload Required
☐ Mandibular Movement	MEDIT	2023-04-24 PM 1:59	1.88 GB	Upload Required
☐ Mandibular Movement Case	MEDIT	2023-04-24 PM 1:59	1.75 GB	Upload Required
☐ Mandibular Movement Case - Clo...	MEDIT	2023-04-24 PM 1:59	1.75 GB	Upload Required
☐ 서라_Mandibular Movement_TEST	Mandibular Movement_export_test	2023-04-10 PM 3:39	34.73 MB	Upload Required
☐ 서라_Mandibular Movement_TEST...	Mandibular Movement_export_test	2023-04-10 PM 3:39	731.43 MB	Upload Required

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
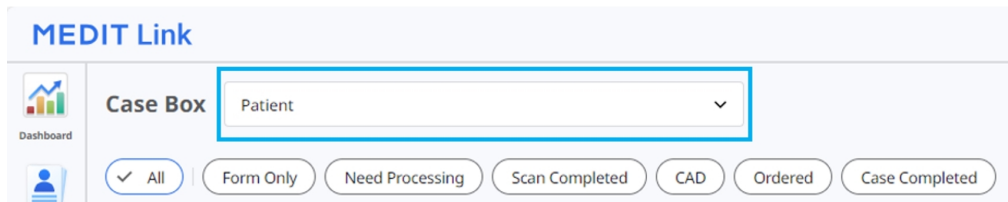
Case Box Patient

Dashboard

All Form Only Need Processing Scan Completed CAD Ordered Case Completed

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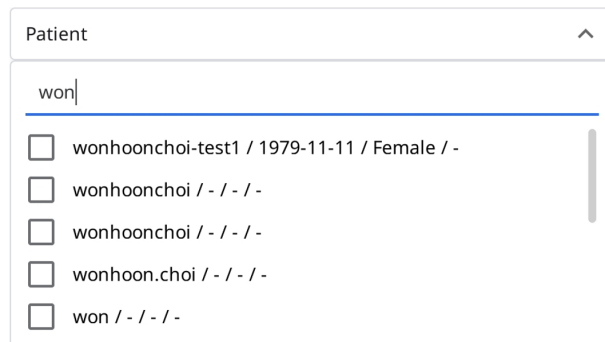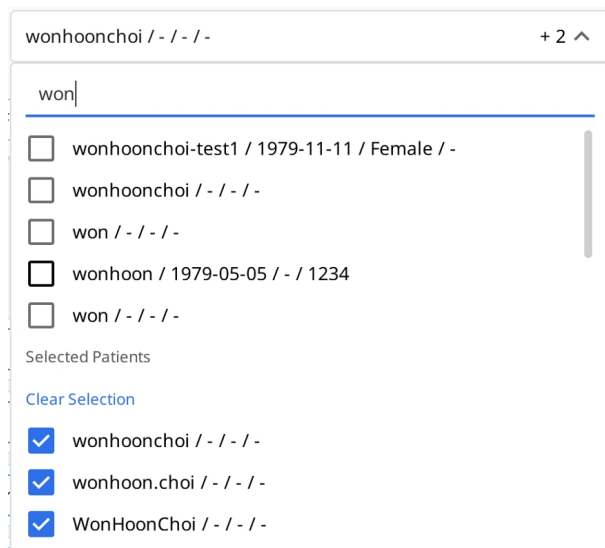
Status	Patient Name	Case Name	Form Information	Last Modified Date	Scan Date	Lab Name	Order Date
Form Only	wonhoonchoi	originalcase	-	2023-10-05 오전 8:10			
Form Only	wonhoonchoi	wonhoonchoi's Case	-	2023-10-05 오전 8:09			
Scan Completed	wonhoonchoi	FileVierTest_1	8-Inlay / 10-Inlay / 12-Inlay / 13-Inlay / 14-Inl...	2023-06-15 오전 8:14	2023-02-09 오전 8:35		
Scan Completed	wonhoonchoi	FileVierTest_1 - Clone	2-Inlay / 3-Inlay / 4-Inlay / 5-Inlay / 6-Inlay / 8...	2023-05-25 오후 2:05	2023-02-09 오전 8:35		
Form Only	wonhoonchoi	FileVierTest_1 - Clone	8-Inlay / 10-Inlay / 12-Inlay / 13-Inlay / 14-Inl...	2023-03-15 오후 3:06			
Scan Completed	wonhoonchoi	FileVierTest_1 - Clone - Clone	8-Inlay / 10-Inlay / 12-Inlay / 13-Inlay / 14-Inl...	2023-03-09 오후 4:41	2023-02-09 오전 8:35		
Need Processing	wonhoonchoi	FileVierTest_1 - Clone	8-Inlay / 10-Inlay / 12-Inlay / 13-Inlay / 14-Inl...	2023-02-14 오전 8:13	2023-02-09 오전 8:35		
Scan Completed	wonhoonchoi	FileVierTest_1 - Clone	8-Inlay / 10-Inlay / 12-Inlay / 13-Inlay / 14-Inl...	2023-02-10 오전 11:27	2023-02-09 오전 8:35		
Form Only	wonhoonchoi	testsamplecase	-	2023-02-10 오전 7:39			
Scan Completed	wonhoonchoi	Multi Occlusion	-	2023-02-09 오전 9:26	2023-02-09 오전 9:01		
Scan Completed	wonhoonchoi	FileVierTest_1 - Clone - Clone	8-Inlay / 10-Inlay / 12-Inlay / 13-Inlay / 14-Inl...	2023-02-09 오전 8:50	2023-02-09 오전 8:35		
Case Completed	wonhoonchoi	wonhoonchoi's Case	12-Inlay / 13-Inlay / 14-Inlay	2022-10-13 오후 10:30	2022-10-11 오후 5:16		
Case Completed	wonhoonchoi	wonhoonchoi's Case	-	2022-10-13 오후 10:28			

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Patient Name*

Case Name*

 Creating patient information involves consent from the patient.
Patient information can be shared with 3rd parties while placing an order.

Partner

Select partner

Assignee

Select assignees

 You can select up to 5 assignees.

Tags

You can add the entered tag with Enter or Spacebar.

 You can add up to 10 tags of 100 characters or less.

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Patient List

Select a patient in the list to set the patient of the case. You can add or manage the patient information in Medit Link.

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
qa_lab_01 Uses Lab Form

Full Services Dental Esthetics Crown & Bridge

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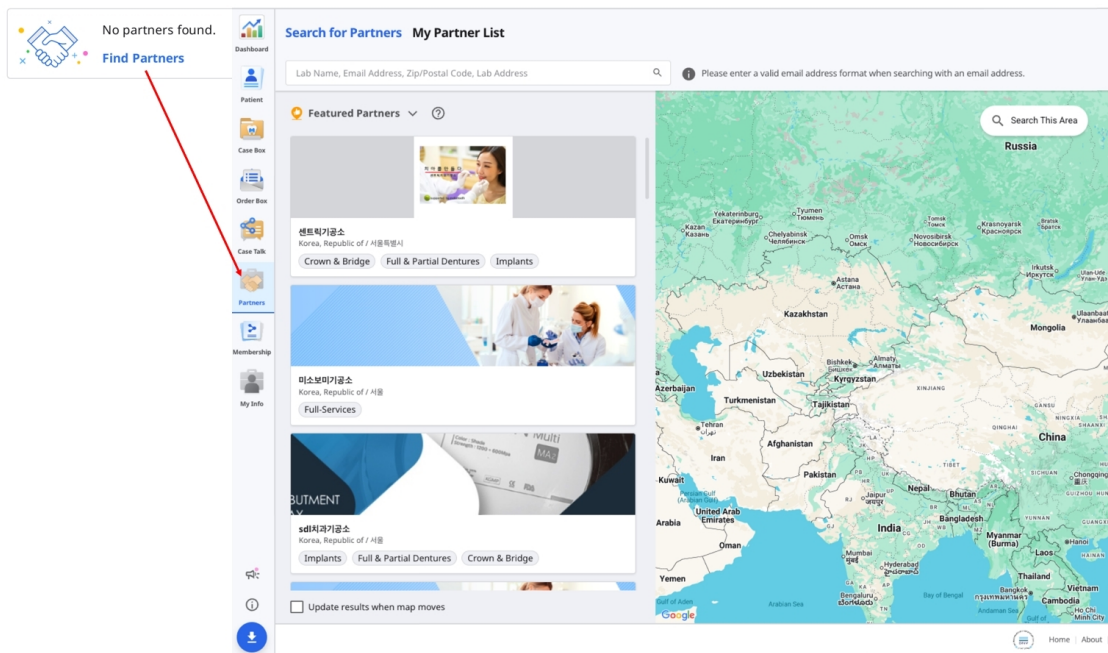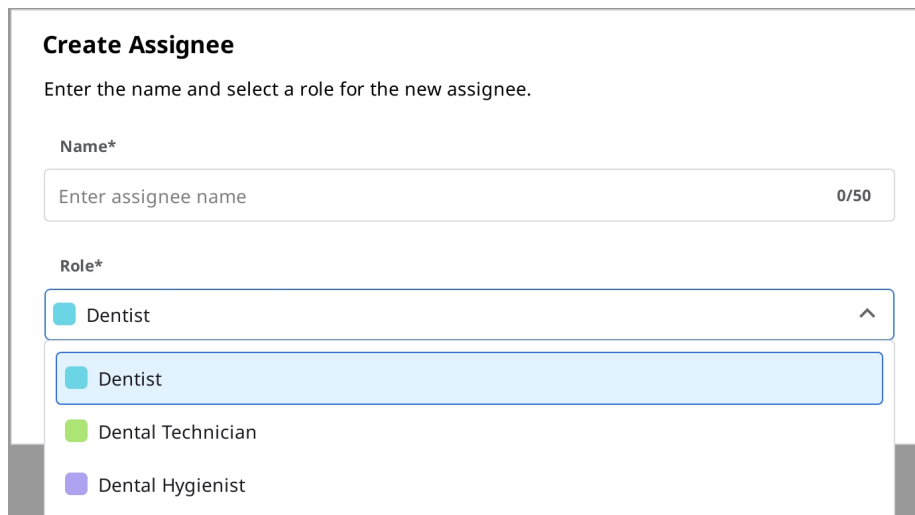
Create Assignee

Enter the name and select a role for the new assignee.

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2025-07-30

Tags

You can add the entered tag with Enter or Spacebar.

ⓘ You can add up to 10 tags of 100 characters or less.

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Options

Export Files To
C:/Users/Wonhoon Choi Laptop/Desktop/ExportTest/3/

Folder Name
Case Creation Date - Case Name

File Name
Case Creation Date - Case Name

Type
Medit

Mesh Format
Convert mesh files to the below formats

☐ meditMesh ☐ OBJ ☒ PLY ☐ STL

☐ Do not convert the attached files

Image Format
☐ Convert images to the below format

☐ JPG ☐ JPEG ☒ PNG ☐ BMP

☐ Compress Files
☐ Combine Individual Mesh
☐ Export All Occlusal Relationship
☐ Save Case Information as PDF

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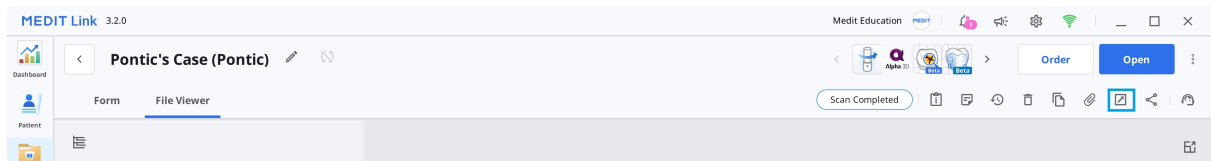
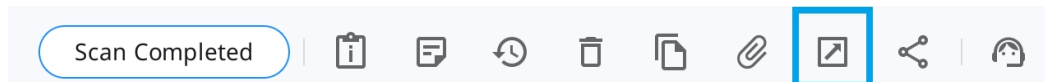
Options **File List**

Export Files To
C:/Users/Medit-Wonhoonchoi/Desktop/test/

Folder Name
Multi Occlusion - Scan - Before Scanning

File Name
Multi Occlusion - Scan - Before Scanning

CAD Type
Maestro


Mesh Format
Convert mesh files to the below formats
☐ meditMesh ☐ OBJ ☒ PLY ☐ STL
☒ Do not convert the attached files

Image Format
☐ Convert images to the below format
☐ JPG ☐ JPEG ☒ PNG ☐ BMP

☐ Compress Files
☐ Combine Individual Mesh
☐ Export All Occlusal Relationship
☐ Save Case Information as PDF [Preview](#)

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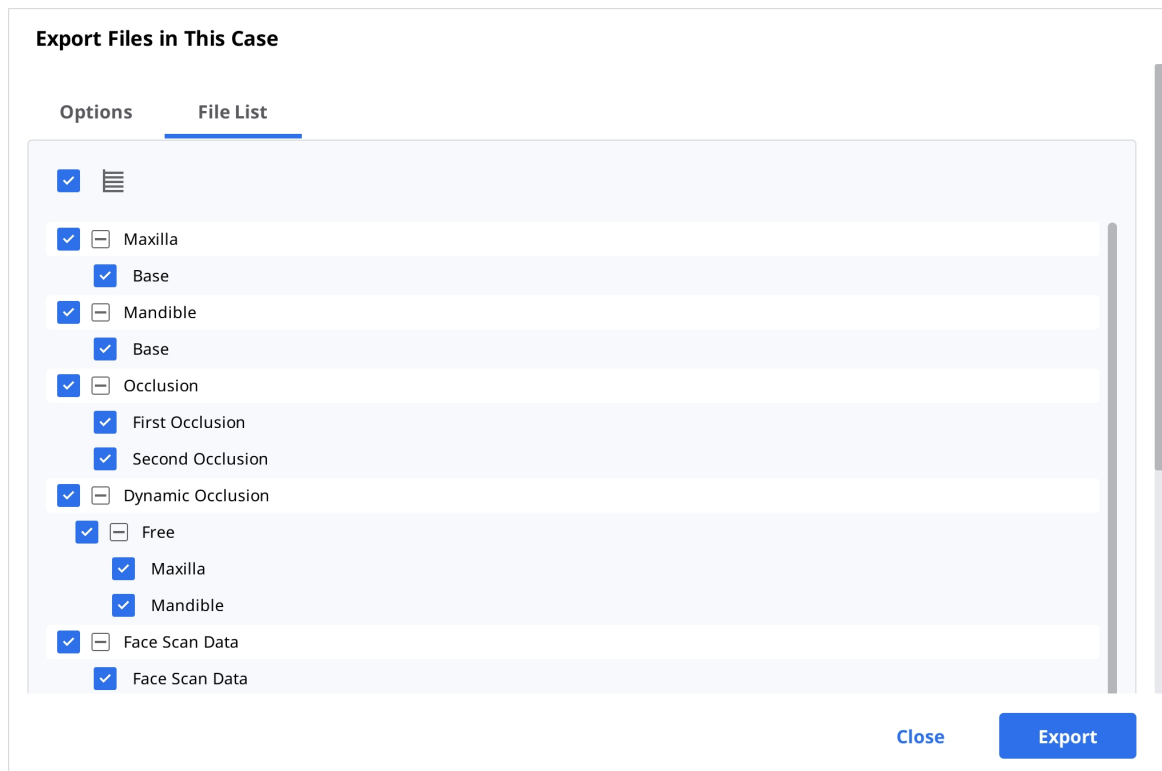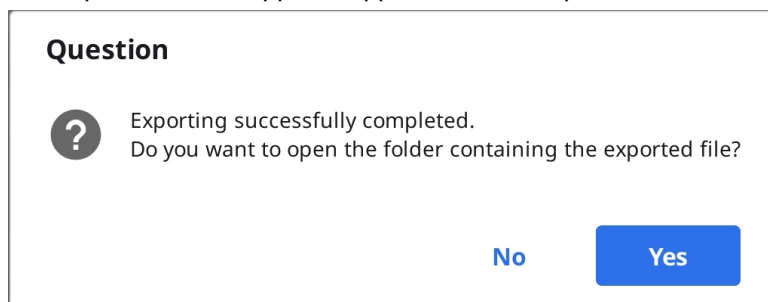
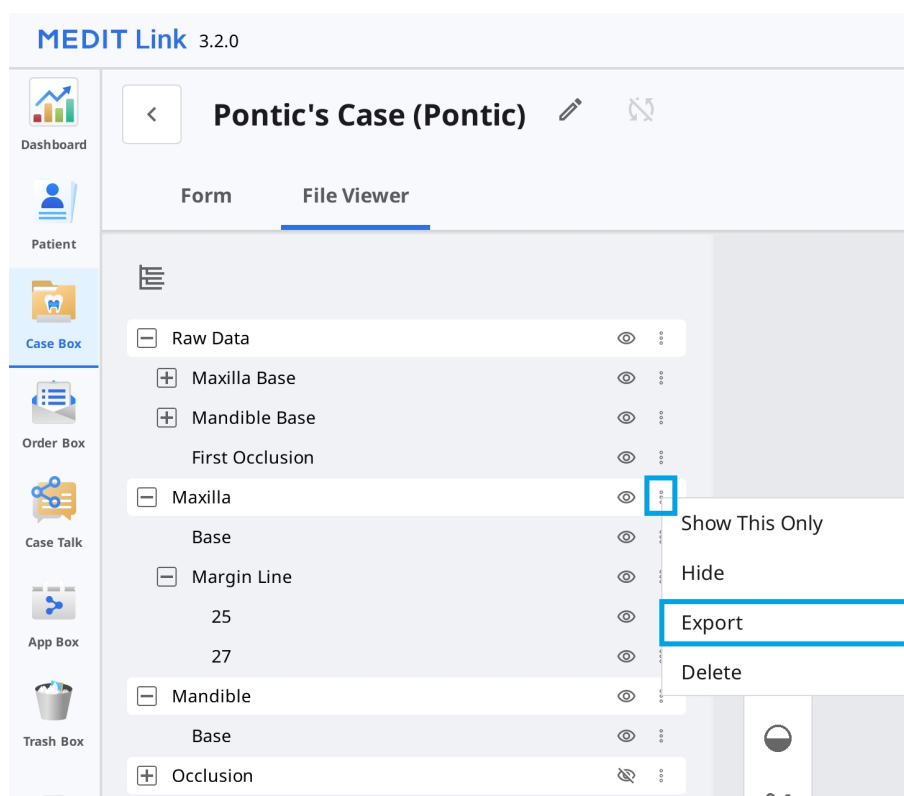
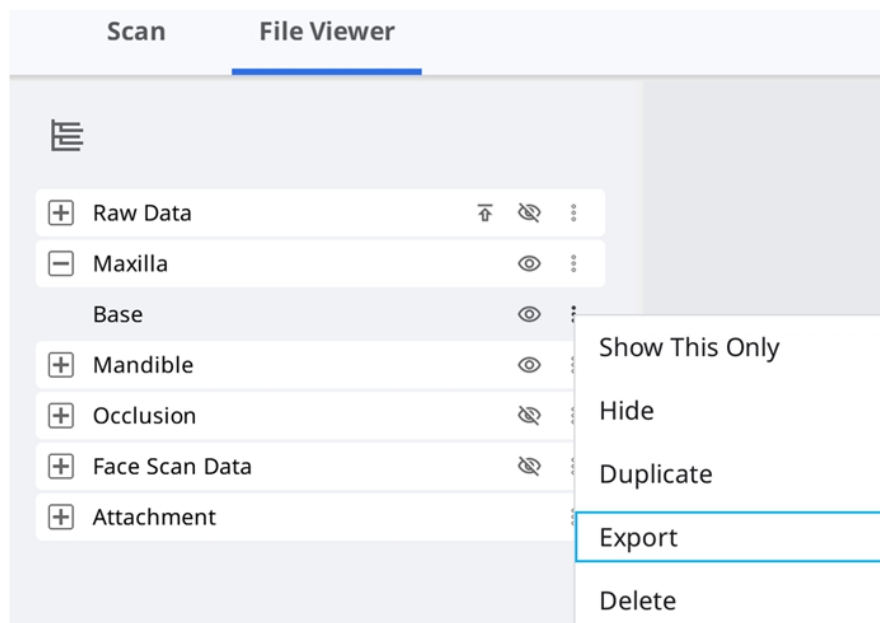
C:/Users/Medit-Wonhoonchoi/Desktop/test/

Folder Name

Multi Occlusion - Scan - Before Scanning

File Name

Multi Occlusion - Scan - Before Scanning

CAD Type

Maestro

Mesh Format

Convert mesh files to the below formats

☐ meditMesh

☐ OBJ

☒ PLY

☐ STL

☒ Do not convert the attached files

Image Format

☐ Convert images to the below format

☐ JPG

☐ JPEG

☒ PNG

☐ BMP

☐ Compress Files

☐ Combine Individual Mesh

☐ Export All Occlusal Relationship

☐ Save Case Information as PDF [Preview](#)

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Options **File List**

☐ Maxilla
 ☐ Base
 ☒ Mandible
 ☒ Base

☐ Occlusion
 ☐ First Occlusion
 ☐ Second Occlusion

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Question

? Exporting successfully completed.
Do you want to open the folder containing the exported file?

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Case Box Patient

All / 2018-01-01 - 2026-09-23 All / Contains

☒ All
 ☐ Form Only
 ☐ Need Processing
 ☐ Scan Completed
 ☐ CAD
 ☐ Ordered
 ☐ Case Completed

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Share through Case Talk

Share with Anyone

Create a shared link and share it with anyone.

Expiration Date

2022-04-27

You can extend the period up to 90 days from today.

PIN Code

☐ PIN Code

Enter PIN Code

Set the PIN code using six digits for privacy protection.

☒ Share the patient's name

☒ Grant Download Permission [Premium Feature](#)

If you check this option, all users who have access to the shared link can download the files.

The premium feature is provided to all users for a limited time.
Feature availability is subject to change depending on the membership plan.

☒ Allow comments

If you check this option, you will be able to communicate through the shared link.

Comment

Please enter your comment.


0/1024

By pressing the "Get Link" button, you agree that you understand personal information can be collected and provided to a third party (or shared in any public domain) ONLY with the consent of the person whose images on which you have worked. You also agree the consent of that person was truthfully obtained, and Medit is not responsible for any problems which may arise there after.

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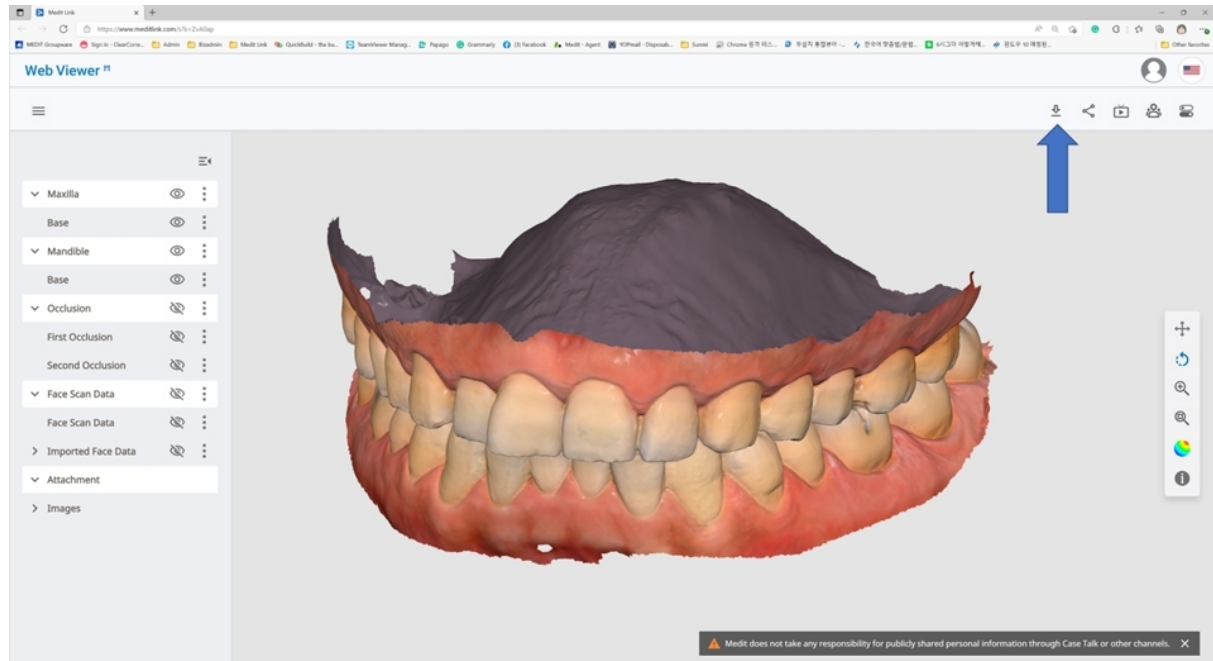Download Premium Feature

The premium feature is provided to all users for a limited time.
Feature availability is subject to change depending on the membership plan.

Mesh Format

Convert mesh files to the below formats

☐ meditMesh ☐ OBJ ☒ PLY ☐ STL

i The selected format will be applied to the result files created in the Medit Link application.

i The attachment files will not be converted but will keep their original format.

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Uses Lab Form

San Diego, US

[92131] Residence Inn San Diego Rancho Bernardo

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
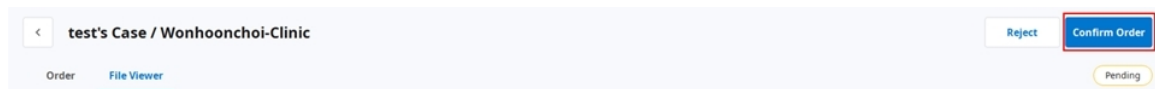
The screenshot shows a header bar with a back arrow, the text 'test's Case / Wonhoonchoi-Clinic', and two buttons: 'Reject' and 'Confirm Order'. The 'Confirm Order' button is highlighted with a red rectangular box. Below the header, there are tabs for 'Order' and 'File Viewer', with 'File Viewer' being the active tab. In the bottom right corner, there is a 'Pending' status indicator.

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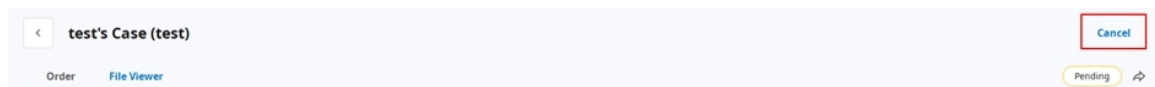
The screenshot shows a header bar with a back arrow, the text 'test's Case (test)', and a 'Cancel' button highlighted with a red rectangular box. Below the header, there are tabs for 'Order' and 'File Viewer', with 'File Viewer' being the active tab. In the bottom right corner, there is a 'Pending' status indicator.

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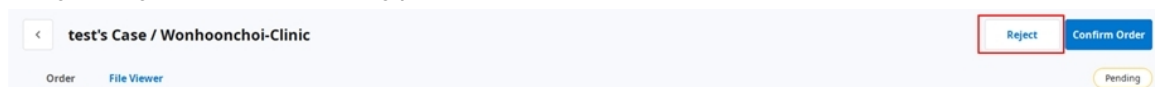The screenshot shows a header bar with a back arrow, the text 'test's Case / Wonhoonchoi-Clinic', and two buttons: 'Reject' and 'Confirm Order'. The 'Reject' button is highlighted with a red rectangular box. Below the header, there are tabs for 'Order' and 'File Viewer', with 'File Viewer' being the active tab. In the bottom right corner, there is a 'Pending' status indicator.

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The screenshot shows a web interface for managing a case. The title is 'test1's Case (test1)'. There are tabs for 'Order' and 'File Viewer'. The 'Order' tab is active, showing a form with fields for Order ID (6004858), Assignee (a dropdown menu), Patient Name (test1), Partner Name (Elysee Dental Test), Scanned Date, Order Date (2025-07-23 PM 5:17), Requested Delivery Date (2025-07-30 PM 05:00), and Expected Delivery Date (2025-07-30 PM 05:00). There is also a Memo field and a Tags field with the tag 'test'. On the right, there is a table titled 'Order List' and 'Shipment'.

No.	Type	Method	Material	Shad
15	Crown	Anatomic	Zirconia	A4

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Shipment Registration

Select a delivery method and enter the required information to register.
The information you enter is saved in the Shipment tab and can be checked at any time.

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Case Name
mandibular movement - Clone - Clone

Patient Name
Test Occlusion

☒ Share patient information

Requested Delivery Date
2023-04-13

Delivery Time
오전 8:15 AM 8:00 PM 8:00

Memo

Tags
You can add the entered tag with Enter or Spacebar.

☒ Include all occlusion relationships in the order
☒ Include Tags In The Order

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Tags
tag

Order List Shipment

No.	Type	Method	Material	Shad
25	Crown	Anatomic	Multi-Layered Zirconia	A3.5
26	Crown	Anatomic	Multi-Layered Zirconia	A3.5

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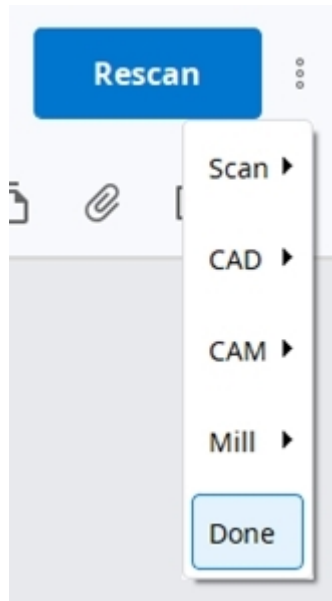
MEDIT Link 3.2.0

Work Box

Dashboard | In Box | **Work Box** | Case Talk | App Box | Trash Box

✓ All | Form Only | Need Processing | Scan Completed | CAD | CAM

[-] Export ▾ Status ▾ [X] Selected Cases: 4

Done

Status	Case Name	Patient Name	Form Information
✓ Form Only	tsstt's Case	tsstt	-
✓ Form Only	test's Case	test	-
✓ CAD	Maxilla Crown	Scan for Clinics	14-Crown
✓ CAD	Maxilla Crown	Scan for Clinics	14-Crown

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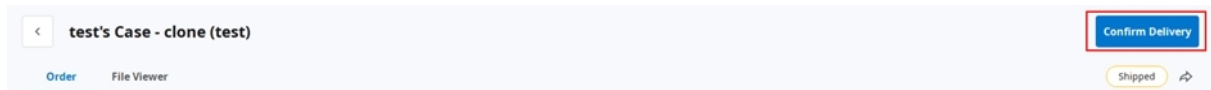
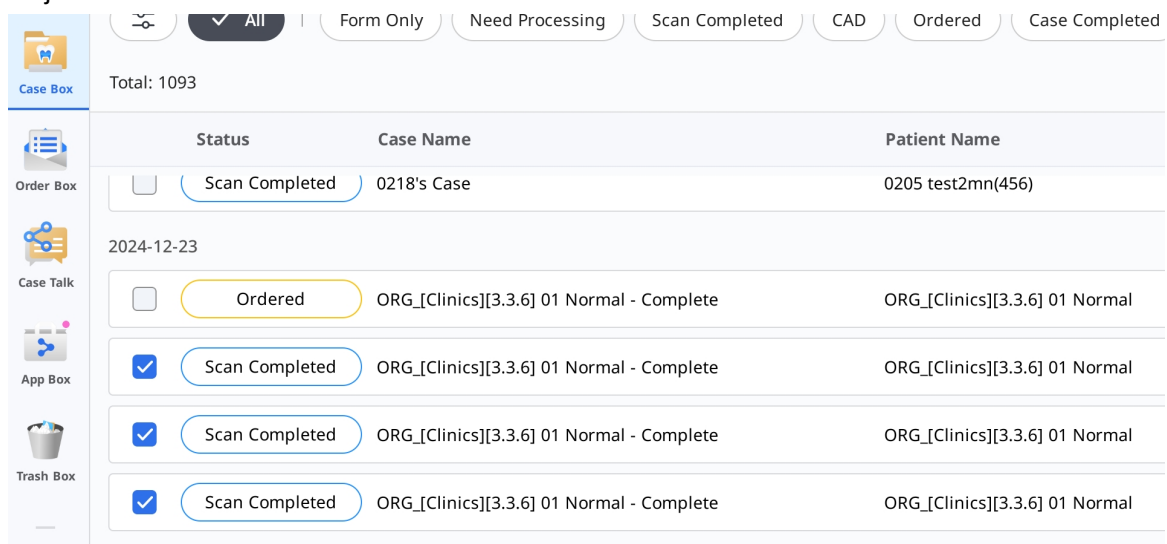
Status	Case Name	Patient Name
Scan Completed	0218's Case	0205 test2mn(456)
2024-12-23		
Ordered	ORG_[Clinics][3.3.6] 01 Normal - Complete	ORG_[Clinics][3.3.6] 01 Normal
Scan Completed	ORG_[Clinics][3.3.6] 01 Normal - Complete	ORG_[Clinics][3.3.6] 01 Normal
Scan Completed	ORG_[Clinics][3.3.6] 01 Normal - Complete	ORG_[Clinics][3.3.6] 01 Normal
Scan Completed	ORG_[Clinics][3.3.6] 01 Normal - Complete	ORG_[Clinics][3.3.6] 01 Normal

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